

### **REMARKS**

Applicant is requesting amendment of the specification, including the abstract, under 37 C.F.R. § 1.121(b)(3). Accordingly, submitted herewith in compliance with 37 C.F.R. §§ 1.125(b) and 1.125(c) is a substitute specification, as well as a marked-up version of the substitute specification showing all amendments made relative to the specification as originally filed. The substitute specification submitted herewith includes the abstract (but not claims, which are not amended). The substitute specification includes no new matter.


The specification has been amended to make reference to the SEQ ID numbers of the numerous sequences disclosed therein, thereby bringing it into compliance with 37 C.F.R. §1.821(d). In addition, paragraph [0001] of the specification has been amended to include a reference to the complete application number of international application PCT/US03/27308, an application previously unambiguously identified by title, filing date, and attorney docket number in the application as originally filed. Minor typographical errors in the specification have also been corrected as shown in the marked-up copy provided herewith. Minor errors in paragraph numbering have also been corrected. Thus, no new matter has been added by these amendments.

In the event the U.S. Patent and Trademark Office determines that an extension and/or other relief is required, applicant petitions for any required relief including extensions of time and authorizes the Commissioner to charge the cost of such petitions and/or other fees due in connection with the filing of this document to **Deposit Account No. 03-1952** referencing docket no. 543312000420.

Respectfully submitted,

Dated: May 23, 2005

By:

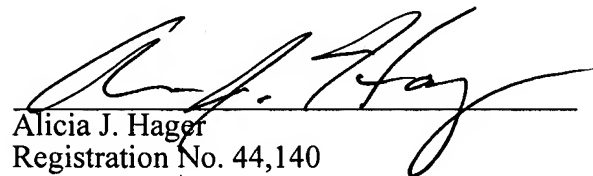  
Alicia J. Hager  
Registration No. 44,140

Morrison & Foerster LLP  
755 Page Mill Road  
Palo Alto, CA 94304  
Telephone: (650) 813-4296  
Facsimile: (650) 494-0792

## **METHODS FOR DETECTION OF GENETIC DISORDERS**

### **CROSS-REFERENCE TO RELATED APPLICATIONS**

**[0001]** This application is a continuation-in-part of PCT/US03/06198, filed February 28, 2003, which claims benefit under 35 U.S.C. §119(e) of U.S. Provisional Patent Application No. 60/378,354, filed May 8, 2002, and which is a continuation-in-part of U.S. Patent Application No. 10/093,618, filed March 11, 2002, which claims benefit under 35 U.S.C. §119(e) of U.S. Provisional Patent Application No. 60/360,232, filed March 1, 2002; this application is also a continuation of PCT/US03/\_\_\_\_\_, ~~Method for Detection of Genetic Disorders, 27308~~, filed August 29, 2003 (~~Attorney Docket No. 54331-20004.41~~); this application further is a continuation-in-part of U.S. Patent Application No. 10/376,770, filed February 28, 2003, which claims benefit under 35 U.S.C. §119(e) of U.S. Provisional Patent Application No. 60/378,354, filed May 8, 2002, and which is a continuation-in-part of U.S. Patent Application No. 10/093,618, filed March 11, 2002, which claims benefit under 35 U.S.C. §119(e) of U.S. Provisional Patent Application No. 60/360,232, filed March 1, 2002. The contents of these applications are incorporated by reference in their entirety.

### **BACKGROUND OF THE INVENTION**

#### **FIELD OF THE INVENTION**

**[0002]** The present invention is directed to a method for the detection of genetic disorders including chromosomal abnormalities and mutations. The present invention provides a rapid, non-invasive method for determining the sequence of DNA from a fetus. The method is especially useful for detection of chromosomal abnormalities in a fetus including translocations, transversions, monosomies, trisomies, and other ~~aneuploidies~~aneuploidies, deletions, additions, amplifications, translocations and rearrangements.

#### **BACKGROUND ART**

**[0003]** Chromosomal abnormalities are responsible for a significant portion of genetic defects in liveborn humans. The nucleus of a human cell contains forty-six (46) chromosomes, which contain the genetic instructions, and determine the operations of the cell. Half of the forty-six chromosomes originate from each parent. Except for the sex chromosomes, which are quite different from each other in normal males, the chromosomes from the mother and the chromosomes from the father make a matched set. The

pairs were combined when the egg was fertilized by the sperm. Occasionally, an error occurs in either the formation or combination of chromosomes, and the fertilized egg is formed with too many or too few chromosomes, or with chromosomes that are mixed in some way. Because each chromosome contains many genes, chromosomal abnormalities are likely to cause serious birth defects, affecting many body systems and often including developmental disability (e.g., mental retardation).

**[0004]** Cells mistakenly can rejoin broken ends of chromosomes, both spontaneously and after exposure to chemical compounds, carcinogens, and irradiation. When rejoining occurs within a chromosome, a chromosome segment between the two breakpoints becomes inverted and is classified as an inversion. With inversions, there is no loss of genetic material; however, inversions can cause disruption of a critical gene, or create a fusion gene that induces a disease related condition.

**[0005]** In a reciprocal translocation, two non-homologous chromosomes break and exchange fragments. In this scenario, two abnormal chromosomes result: each consists of a part derived from the other chromosome and lacks a part of itself. If the translocation is of a balanced type, the individual will display no abnormal phenotypes. However, during germ-cell formation in the translocation-bearing individuals, the proper distribution of chromosomes in the egg or sperm occasionally fails, resulting in miscarriage, malformation, or mental retardation of the offspring.

**[0006]** In a Robertsonian translocation, the centromeres of two acrocentric (a chromosome with a non-centrally located centromere) chromosomes fuse to generate one large metacentric chromosome. The karyotype of an individual with a centric fusion has one less than the normal diploid number of chromosomes.

**[0007]** Errors that generate too many or too few chromosomes can also lead to disease phenotypes. For example, a missing copy of chromosome X (monosomy X) results in Turner's Syndrome, while an additional copy of chromosome 21 results in Down's Syndrome. Other diseases such as Edward's Syndrome, and Patau Syndrome are caused by an additional copy of chromosome 18, and chromosome 13, respectively.

**[0008]** One of the most common chromosome abnormalities is known as Down syndrome. The estimated incidence of Down's syndrome is between 1 in 1,000 to 1 in 1,100 live births. Each year approximately 3,000 to 5,000 children are born in the U.S. with this chromosomal disorder. The vast majority of children with Down syndrome (approximately 95 percent) have an extra chromosome 21. Most often, the extra chromosome originates from the mother. However, in about 3-4 percent of people with Down syndrome, a translocation between chromosome 21 and either 14 or 22 is responsible for the genetic abnormality. Finally, another chromosome problem, called mosaicism, is noted in about 1 percent



of individuals with Down's syndrome. In this case, some cells have 47 chromosomes and others have 46 chromosomes. ~~Mosaicism~~Mosaicism is thought to be the result of an error in cell division soon after conception.

**[0009]** Chromosomal abnormalities are ~~congenital~~congenital, and therefore, prenatal diagnosis can be used to determine the health and condition of an unborn fetus. Without knowledge gained by prenatal diagnosis, there could be an untoward outcome for the fetus or the mother or both. Congenital anomalies account for 20 to 25% of perinatal deaths. Specifically, prenatal diagnosis is helpful for managing the remaining term of the pregnancy, planning for possible complications with the birth process, preparing for problems that can occur in the newborn infant, and finding conditions that may affect future pregnancies.

**[0010]** There are a variety of non-invasive and invasive techniques available for prenatal diagnosis including ultrasonography, amniocentesis, chorionic villus sampling (CVS), fetal blood cells in maternal blood, maternal serum alpha-fetoprotein, maternal serum beta-HCG, and maternal serum estriol. However, the techniques that are non-invasive are less specific, and the techniques with high specificity and high sensitivity are highly invasive. Furthermore, most techniques can be applied only during specific time periods during pregnancy for greatest utility.

### **Ultrasonography**

**[0011]** This is a harmless, non-invasive procedure. High frequency sound waves are used to generate visible images from the pattern of the echoes made by different tissues and organs, including the fetus in the amniotic cavity. The developing embryo can be visualized at about 6 weeks of gestation. The major internal organs and extremities can be assessed to determine if any are abnormal at about 16 to 20 weeks gestation.

**[0012]** An ultrasound examination can be useful to determine the size and position of the fetus, the amount of amniotic fluid, and the appearance of fetal anatomy; however, there are limitations to this procedure. Subtle abnormalities, such as Down syndrome, where the morphologic abnormalities are often not marked, but only subtle, may not be detected at all.

### **Amniocentesis**

**[0013]** This is a highly invasive procedure in which a needle is passed through the mother's lower abdomen into the amniotic cavity inside the uterus. This procedure can be performed at about 14 weeks gestation. For prenatal diagnosis, most amniocenteses are performed between 14 and 20 weeks gestation. However, an ultrasound examination is performed, prior to amniocentesis, to determine gestational age, position of the fetus and placenta, and determine if enough amniotic fluid is present. Within the amniotic

fluid are fetal cells (mostly derived from fetal skin) which can be grown in culture for chromosomal, biochemical, and molecular biologic analyses.

**[0014]** Large chromosomal abnormalities, such as extra or missing chromosomes or chromosome fragments, can be detected by karyotyping, which involves the identification and analysis of all 46 chromosomes from a cell and arranges them in their matched pairs, based on subtle differences in size and structure. In this systematic display, abnormalities in chromosome number and structure are apparent. This procedure typically takes 7-10 days for completion.

**[0015]** While amniocentesis can be used to provide direct genetic information, risks are associated with the procedure including fetal loss and maternal Rh sensitization. The increased risk for fetal mortality following amniocentesis is about 0.5% above what would normally be expected. Rh negative mothers can be treated with RhoGam.

### **Chorionic Villus Sampling (CVS)**

**[0016]** In this procedure, a catheter is passed via the vagina through the cervix and into the uterus to the developing placenta with ultrasound guidance. The introduction of the catheter allows cells from the placental chorionic villi to be obtained and analyzed by a variety of techniques, including chromosome analysis to determine the karyotype of the fetus. The cells can also be cultured for biochemical or molecular biologic analysis. Typically, CVS is performed between 9.5 and 12.5 weeks gestation.

**[0017]** CVS has the disadvantage of being an invasive procedure, and it has a low but significant rate of morbidity for the fetus; this loss rate is about 0.5 to 1 % higher than for women undergoing amniocentesis. Rarely, CVS can be associated with limb defects in the fetus. Also, the possibility of maternal Rh sensitization is present. Furthermore, there is also the possibility that maternal blood cells in the developing placenta will be sampled instead of fetal cells and confound chromosome analysis.

### **Maternal Serum Alpha-Fetoprotein (MSAFP)**

**[0018]** The developing fetus has two major blood proteins--albumin and alpha-fetoprotein (AFP). The mother typically has only albumin in her blood, and thus, the MSAFP test can be utilized to determine the levels of AFP from the fetus. Ordinarily, only a small amount of AFP gains access to the amniotic fluid and crosses the placenta to mother's blood. However, if the fetus has a neural tube defect, then more AFP escapes into the amniotic fluid. Neural tube defects include anencephaly (failure of closure at the cranial end of the neural tube) and spina bifida (failure of closure at the caudal end of the neural tube). The incidence of such defects is about 1 to 2 births per 1000 in the United States. Also, if

there are defects in the fetal abdominal wall, the AFP from the fetus will end up in maternal blood in higher amounts.

**[0019]** The amount of MSAFP increases with gestational age, and thus for the MSAFP test to provide accurate results, the gestational age must be known with certainty. Also, the race of the mother and presence of gestational diabetes can influence the level of MSAFP that is to be considered normal. The MSAFP is typically reported as multiples of the mean (MoM). The greater the MoM, the more likely a defect is present. The MSAFP test has the greatest sensitivity between 16 and 18 weeks gestation, but can be used between 15 and 22 weeks gestation. The MSAFP tends to be lower when Down's Syndrome or other chromosomal abnormalities is present.

**[0020]** While the MSAFP test is non-invasive, the MSAFP is not 100% specific. MSAFP can be elevated for a variety of reasons that are not related to fetal neural tube or abdominal wall defects. The most common cause for an elevated MSAFP is a wrong estimation of the gestational age of the fetus. Therefore, results from an MSAFP test are never considered definitive and conclusive.

#### **Maternal Serum Beta-HCG**

**[0021]** Beginning at about a week following conception and implantation of the developing embryo into the uterus, the trophoblast will produce detectable beta-HCG (the beta subunit of human chorionic gonadotropin), which can be used to diagnose pregnancy. The beta-HCG also can be quantified in maternal serum, and this can be useful early in pregnancy when threatened abortion or ectopic pregnancy is suspected, because the amount of beta-HCG will be lower than normal.

**[0022]** In the middle to late second trimester, the beta-HCG can be used in conjunction with the MSAFP to screen for chromosomal abnormalities, in particular for Down syndrome. An elevated beta-HCG coupled with a decreased MSAFP suggests Down syndrome. High levels of HCG suggest trophoblastic disease (molar pregnancy). The absence of a fetus on ultrasonography along with an elevated HCG suggests a hydatidiform mole.

#### **Maternal Serum Estriol**

**[0023]** The amount of estriol in maternal serum is dependent upon a viable fetus, a properly functioning placenta, and maternal well-being. Dehydroepiandrosterone (DHEA) is made by the fetal adrenal glands, and is metabolized in the placenta to estriol. The estriol enters the maternal circulation and is excreted by the maternal kidney in urine or by the maternal liver in the bile. Normal levels of estriol, measured in the third trimester, will give an indication of general well-being of the fetus. If the

estriol level drops, then the fetus is threatened and an immediate delivery may be necessary. Estriol tends to be lower when Down syndrome is present and when there is adrenal hypoplasia with anencephaly.

### **The Triple Screen Test**

**[0024]** The triple screen test comprises analysis of maternal serum alpha-feto-protein (MSAFP), human chorionic gonadotrophin (hCG), and unconjugated estriol (uE3). The blood test is usually performed 16-18 weeks after the last menstrual period. While the triple screen test is non-invasive, abnormal test results are not indicative of a birth defect. Rather, the test only indicates an increased risk and suggests that further testing is needed. For example, 100 out of 1,000 women will have an abnormal result from the triple screen test. However, only 2-3 of the 100 women will have a fetus with a birth defect. This high incidence of false positives causes tremendous stress and unnecessary anxiety to the expectant mother.

### **Fetal Cells Isolated From Maternal Blood**

**[0025]** The presence of fetal nucleated cells in maternal blood makes it possible to use these cells for noninvasive prenatal diagnosis (Walknowska, et al., Lancet 1:1119-1122, 1969; Lo et al., Lancet 2:1363-65, 1989; Lo et al., Blood 88:4390-95, 1996). The fetal cells can be sorted and analyzed by a variety of techniques to look for particular DNA sequences (Bianchi et al., Am. J. Hum. Genet. 61:822-29, (1997); Bianchi et al., PNAS 93:705-08, (1996)). Fluorescence in-situ hybridization (FISH) is one technique that can be applied to identify particular chromosomes of the fetal cells recovered from maternal blood and diagnose aneuploid conditions such as trisomies and monosomy X. Also, it has been reported that the number of fetal cells in maternal blood increases in aneuploid pregnancies.

**[0026]** The method of FISH uses DNA probes labeled with colored fluorescent tags that allow detection of specific chromosomes or genes under a microscope. Using FISH, subtle genetic abnormalities that cannot be detected by standard karyotyping are readily identifiable. This procedure typically takes 24-48 hours to complete. Additionally, using a panel of multi-colored DNA FISH probes, abnormal chromosome copy numbers can be seen.

**[0027]** While improvements have been made for the isolation and enrichment of fetal cells, it is still difficult to get many fetal blood cells. There may not be enough to reliably determine anomalies of the fetal karyotype or assay for other abnormalities. Furthermore, most techniques are time consuming, require high-inputs of labor, and are difficult to implement for a high throughput fashion.

### **Fetal DNA From Maternal Blood**

**[0028]** Fetal DNA has been detected and quantitated in maternal plasma and serum (Lo et al., Lancet 350:485-487 (1997); Lo et al., Am. J. hum. Genet. 62:768-775 (1998)). Multiple fetal cell types occur in the maternal circulation, including fetal granulocytes, lymphocytes, nucleated red blood cells, and trophoblast cells (Pertl, and Bianchi, Obstetrics and Gynecology 98: 483-490 (2001)). Fetal DNA can be detected in the serum at the seventh week of gestation, and increases with the term of the pregnancy. The fetal DNA present in the maternal serum and plasma is comparable to the concentration of DNA obtained from fetal cell isolation protocols.

**[0029]** Circulating fetal DNA has been used to determine the sex of the fetus (Lo et al., Am. J. hum. Genet. 62:768-775 (1998)). Also, fetal rhesus D genotype has been detected using fetal DNA. However, the diagnostic and clinical applications of circulating fetal DNA is limited to genes that are present in the fetus but not in the mother (Pertl and Bianchi, Obstetrics and Gynecology 98: 483-490 (2001)). Thus, a need still exists for a non-invasive method that can determine the sequence of fetal DNA and provide definitive diagnosis of chromosomal abnormalities in a fetus.

### **BRIEF SUMMARY OF THE INVENTION**

**[0030]** The invention is directed to a method for detection of genetic disorders including mutations and chromosomal abnormalities. In some embodiments, the present invention is used to detect mutations, and chromosomal abnormalities including but not limited to translocation, transversion, monosomy, trisomy, and other aneuploidies, deletion, addition, amplification, fragment, translocation, and rearrangement. Numerous abnormalities can be detected simultaneously. The present invention also provides a non-invasive method to determine the sequence of fetal DNA from a sample of a pregnant female. The present invention can be used to detect any alternation in gene sequence as compared to the wild type sequence including but not limited to point mutation, reading frame shift, transition, transversion, addition, insertion, deletion, addition-deletion, frame-shift, missense, reverse mutation, and microsatellite alteration. The present invention also provides a method for isolating free nucleic acid from a sample containing nucleic acid. The present invention also provides ~~compositions~~compositions and kits.

**[0031]** In one aspect, the invention is directed to methods for detecting chromosomal abnormalities. In one embodiment, the present invention is directed to a method for detecting chromosomal abnormalities, said method comprising quantitating the relative amount of the alleles at a heterozygous locus of interest, where the heterozygous locus of interest was previously identified by determining the

sequence of alleles at a locus of interest from template DNA, wherein said relative amount is expressed as a ratio, and wherein said ratio indicates the presence or absence of a chromosomal abnormality.

**[0032]** In some embodiments, determining the sequence includes using a method that is allele specific PCR, mass spectrometry, hybridization, primer extension, fluorescence resonance energy transfer (FRET), sequencing, Sanger dideoxy sequencing, DNA microarray, GeneCHIP arrays, HuSNP arrays, CodeLink Arrays, BeadArray Technology, MassARRAY, MassEXTEND, SNP-IT, TaqMan, InvaderStrand Assay, southern blot, slot blot, dot blot, or MALDI-TOF mass spectrometry.

**[0033]** In some embodiments, template DNA is obtained from human, non-human, mammal, reptile, cattle, cat, dog, goat, swine, pig, monkey, ape, gorilla, bull, cow, bear, horse, sheep, poultry, mouse, rat, fish, dolphin, whale, or shark. In an embodiment, the template DNA is obtained from a human source. In a preferred embodiment, the template DNA is obtained from a pregnant human female. In some embodiments, the template DNA is obtained from a sample that is a cell, fetal cell, tissue, blood, serum, plasma, saliva, urine, tear, vaginal secretion, sweat, umbilical cord blood, chorionic villi, amniotic fluid, embryonic tissue, an embryo, a two-celled embryo, a four-celled embryo, an eight celled embryo, a 16-celled embryo, a 32-celled embryo, a 64-celled embryo, a 128-celled embryo, a 256-celled embryo, a 512-celled embryo, a 1024-celled embryo, lymph fluid, cerebrospinal fluid, mucosa secretion, peritoneal fluid, ascitic fluid, fecal matter, or body exudates. In these embodiments, the sample may be mixed with an agent that inhibits cell lysis to inhibit the lysis of cells, if cells are present, where the agent is a membrane stabilizer, a cross-linker, or a cell lysis inhibitor. In some of these embodiments, agent is a cell lysis inhibitor, and may be glutaraldehyde, derivatives of glutaraldehyde, formaldehyde, formalin, or derivatives of formaldehyde. In some of these embodiments the sample is blood and in one embodiment the sample is blood from a pregnant female, e.g., a human female. In the latter embodiment, the fetus may be at a gestational age selected from the group consisting of: 0-4, 4-8, 8-12, 12-16, 16-20, 20-24, 24-28, 28-32, 32-36, 36-40, 40-44, 44-48, 48-52, or more than 52 weeks. In some of these embodiments, the template DNA may be obtained from plasma or from serum from the blood. In these embodiments, the template DNA may include a mixture of maternal DNA and fetal DNA, and in one embodiment, prior to determining the sequence of alleles of a locus of interest from template DNA, maternal DNA is sequenced to identify a homozygous locus of interest, and the homozygous locus of interest is the locus of interest analyzed in the template DNA.[[.]] In another embodiment, maternal DNA is sequenced to identify a heterozygous locus of interest, and the heterozygous locus of interest is the locus of interest analyzed in the template DNA.

**[0034]** In embodiments, alleles of multiple loci of interest are sequenced and their relative amounts quantitated and expressed as a ratio. In one embodiment, the sequence of alleles of one to tens to hundreds to thousands of loci of interest on a single chromosome on template DNA is determined. In another embodiment, the sequence of alleles of one to tens to hundreds to thousands of loci of interest on multiple chromosomes is determined.

**[0035]** In an embodiment, the locus of interest is suspected of containing a single nucleotide polymorphism or mutation. The method can be used for determining sequences of multiple loci of interest concurrently. The template DNA can comprise multiple loci from a single chromosome. The template DNA can comprise multiple loci from different chromosomes. The loci of interest on template DNA can be amplified in one reaction. Alternatively, each of the loci of interest on template DNA can be amplified in a separate reaction. The amplified DNA can be pooled together prior to digestion of the amplified DNA. Each of the labeled DNA containing a locus of interest can be separated prior to determining the sequence of the locus of interest. In one embodiment, at least one of the loci of interest is suspected of containing a single nucleotide polymorphism or a mutation.

**[0036]** There is no limitation as to the chromosomes that can be compared. The ratio for the alleles at a heterozygous locus of interest on any chromosome can be compared to the ratio for the alleles at a heterozygous locus of interest on any other chromosome. In another embodiment, the ratio of alleles at a heterozygous locus of interest on a chromosome is compared to the ratio of alleles at a heterozygous locus of interest on two, three, four or more than four chromosomes. In another embodiment, the ratio of alleles at multiple loci of interest on a chromosome is compared to the ratio of alleles at multiple loci of interest on two, three, four, or more than four chromosomes. In embodiments, the ratio for alleles at heterozygous loci of interest on a chromosome are summed and compared to the ratio for alleles at heterozygous loci of interest on a different chromosome, where a difference in ratios indicates the presence of a chromosomal abnormality. In some of these embodiments, the chromosomes that are compared are human chromosomes such as chromosome 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X, or Y. In one of the latter embodiments, the ratio for the alleles at heterozygous loci of interest of chromosomes 13, 18, and 21 are compared. In another embodiment, the sequence of one to tens to hundreds to thousands of loci of interest on the template DNA obtained from a sample of a pregnant female is determined. In one embodiment, the loci of interest are on one chromosome. In another embodiment, the loci of interest are on multiple chromosomes.

**[0037]** In some embodiments, determining the sequence of the alleles comprises amplifying alleles of a locus of interest on a template DNA using a first and a second primer, where the second primer

contains a recognition site for a restriction enzyme such that digestion with the restriction enzyme generates a 5' overhang containing the locus of interest; digesting the amplified DNA with the restriction enzyme that recognizes the recognition site on the second primer; incorporating a nucleotide into the digested DNA by using the 5' overhang containing the locus of interest as a template; and determining the sequence of the alleles of the locus of interest by determining the sequence of the DNA into which the nucleotide was incorporated. In one embodiment, determination of the sequence of the locus of interest in comprises detecting a nucleotide.

**[0038]** In other embodiments, determining the sequence of alleles comprises amplifying alleles of a locus of interest on a template DNA using a first and second primers, where the second primer contains a recognition site for a restriction enzyme such that digestion with the restriction enzyme generates a 5' overhang containing the locus of interest; digesting the amplified DNA with the restriction enzyme that recognizes the recognition site on the second primer; incorporating nucleotides into the digested DNA of (b), where a nucleotide that terminates elongation, and is complementary to the locus of interest of an allele, is incorporated into the 5' overhang of said allele, and a nucleotide complementary to the locus of interest of a different allele is incorporated into the 5' overhang of said different allele, and said terminating nucleotide, which is complementary to a nucleotide in the 5' overhang of said different allele, is incorporated into the 5' overhang of said different allele; and determining the sequence of the alleles of a locus of interest by determining the sequence of the DNA into which the complementary nucleotides have been incorporated. In one embodiment, determination of the sequence of the locus of interest comprises detecting a nucleotide.

**[0039]** The incorporation of a nucleotide may be accomplished by a DNA polymerase, including but not limited to E. coli DNA polymerase, Klenow fragment of E. coli DNA polymerase I, T7 DNA polymerase, T4 DNA polymerase, T5 DNA polymerase, Klenow class polymerases, Taq polymerase, bacteriophage 29, REDTaq™ Genomic DNA polymerase, Pfu DNA polymerase, Vent DNA polymerase or sequenase. Incorporation of a nucleotide may include incorporation of a labeled nucleotide, or labeled and unlabeled nucleotides. One nucleotide, two nucleotides, three nucleotides, four nucleotides, five nucleotides, or more than five nucleotides can be incorporated. A combination of labeled and unlabeled nucleotides can be incorporated. The labeled nucleotide may be a dideoxynucleotide triphosphate (also referred to as “dideoxy”) or deoxynucleotide triphosphate (also referred to as “deoxy”). The unlabeled nucleotide may be a dideoxynucleotide triphosphate or deoxynucleotide triphosphate. Labeled nucleotides may be labeled with a radioactive molecule, fluorescent molecule, antibody, antibody fragment, hapten, carbohydrate, biotin, derivative of biotin, phosphorescent moiety, luminescent moiety,



electrochemiluminescent moiety, chromatic moiety, and moiety having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity. In one embodiment, the labeled nucleotide is labeled with a fluorescent molecule. The incorporation of a fluorescent labeled nucleotide may further comprise using a mixture of fluorescent and unlabeled nucleotides.

**[0040]** In one embodiment, the determination of the sequence of the locus of interest comprises detecting the incorporated nucleotide. The detection method includes but is not limited to gel electrophoresis, capillary electrophoresis, microchannel electrophoresis, polyacrylamide gel electrophoresis, fluorescence detection, fluorescence polarization, DNA sequencing, Sanger dideoxy sequencing, ELISA, mass spectrometry, time of flight mass spectrometry, quadrupole mass spectrometry, magnetic sector mass spectrometry, electric sector mass spectrometry, fluorometry, infrared spectrometry, ultraviolet spectrometry, potentiostatic amperometry, DNA hybridization, DNA microarray, GeneChip arrays, HuSNP arrays, BeadArrays, MassExtend, SNP-IT, TaqMan assay, Invader assay, MassCleave, southern blot, slot blot, or dot blot.

**[0041]** In embodiments, first and second primers contain a portion of a restriction enzyme recognition site that contains a variable nucleotide, where the full restriction enzyme recognition site is generated after amplification. In some embodiments, the 3' region of said primers can contain mismatches with the template DNA, and digestion with said restriction enzyme generates a 5' overhang containing the locus of interest. In some embodiments, the restriction enzyme recognition site is for a restriction enzyme that includes but is not limited to BsaJ I, BssK I, Dde I, EcoN I, Fnu4H I, Hinf I, or ScrF I. In some embodiments, the restriction enzyme cuts DNA at a distance from the recognition site. In some of these embodiments, the recognition site is for a Type IIS restriction enzyme. In some of these embodiments, the Type IIS restriction enzyme includes but is not limited to Alw I, Alw26 I, Bbs I, Bbv I, BceA I, Bmr I, Bsa I, Bst71 I, BsmA I, BsmB I, BsmF I, BspM I, Ear I, Fau I, Fok I, Hga I, Ple I, Sap I, SSfaN I, or SthI32 I.

**[0042]** In some embodiments, the recognition site for restriction enzymes includes but is not limited to BsaJ I (5' C<sup>1</sup>CNNGG 3'), BssK I (5' <sup>1</sup>CCNGG 3'), Dde I (5' C<sup>1</sup>TNAG 3'), EcoN I (5' CCTNN<sup>1</sup>NNNAGG 3' (SEQ ID NO: 7)), Fnu4H I (5' GC<sup>1</sup>NGC 3'), Hinf I (5' G<sup>1</sup>ANTC 3'), PflF I (5' GACN<sup>1</sup>NNGTC 3'), Sau96 I (5' G<sup>1</sup>GNCC 3'), ScrF I (5' CC<sup>1</sup>NGG 3'), TthI 11 I (5' GACN<sup>1</sup>NNGTC 3'), and more preferably Fnu4H I and EcoN I, is generated after amplification.

**[0043]** The first and/or second primer can contain a tag at the 5' terminus. In some embodiments, the first primer contains a tag at the 5' terminus. The tag can be used to separate the amplified DNA from the template DNA. The tag can be used to separate the amplified DNA containing the labeled nucleotide

from the amplified DNA that does not contain the labeled nucleotide. The tag can be any chemical moiety including but not limited to radioisotope, fluorescent reporter molecule, chemiluminescent reporter molecule, antibody, antibody fragment, hapten, biotin, derivative of biotin, photobiotin, iminobiotin, digoxigenin, avidin, enzyme, acridinium, sugar, enzyme, apoenzyme, homopolymeric oligonucleotide, hormone, ferromagnetic moiety, paramagnetic moiety, diamagnetic moiety, phosphorescent moiety, luminescent moiety, electrochemiluminescent moiety, chromatic moiety, moiety having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity, or combinations thereof. In some embodiments, the tag is biotin. The biotin tag is used to separate amplified DNA from the template DNA using a streptavidin matrix. The streptavidin matrix may be coated on wells of a microtiter plate.

**[0044]** In some embodiments, the annealing length of the second primer is selected from the group consisting of 35-30, 30-25, 25-20, 20-15, 15, 14, 13, 12, 11, 10, 9, 8, 7, 6, 5, 4, and less than 4 bases.

**[0045]** In embodiments, the method of amplification includes but is not limited to polymerase chain reaction, self-sustained sequence reaction, ligase chain reaction, rapid amplification of cDNA ends, polymerase chain reaction and ligase chain reaction, Q-beta phage amplification, strand displacement amplification, or splice overlap extension polymerase chain reaction. In some embodiments, the method of amplification is PCR. In some of these embodiments, an annealing temperature for cycle 1 of PCR is about the melting temperature of the portion of the 3' region of the second primer that anneals to the template DNA. In some of the latter embodiments, an annealing temperature for cycle 2 of PCR is about the melting temperature of the portion of the 3' region of the first primer that anneals to the template DNA. In some of the latter embodiments, an annealing temperature for the remaining cycles of PCR is at about the melting temperature of the entire second primer.

**[0046]** In another aspect, the invention provides methods of determining the sequence of a locus of interest from a sample comprising free fetal DNA, where an agent that inhibits cell lysis has been added to the sample to inhibit lysis of cells, if cells are present, where the agent is a membrane stabilizer, a cross-linker, or a cell lysis inhibitor.

**[0047]** In some embodiments, the agent is a cell lysis inhibitor, and in some of these embodiments, the cell lysis inhibitor includes but is not limited to glutaraldehyde, derivatives of glutaraldehyde, formaldehyde, derivatives of formaldehyde, or formalin. In embodiments, the sample includes but is not limited to tissue, cell, blood, serum, plasma, urine, or vaginal secretion. In some embodiments, the sample is blood. In some of these embodiments, the template DNA is isolated from the serum, in other embodiments the template DNA is isolated from plasma. In some embodiments, the sample contains free

maternal template DNA and free fetal template DNA. In some embodiments, prior to determining the sequence, template DNA was isolated. In some embodiments, prior to determining the sequence of the locus of interest on fetal DNA, the sequence of the locus of interest on maternal template DNA was determined. In some embodiments, prior to determining the sequence of the locus of interest on fetal DNA, the sequence of the locus of interest on paternal template DNA was determined. In some embodiments, the locus of interest is a single nucleotide polymorphism. In other embodiments, the locus of interest is a mutation. In some embodiments, the sequence of multiple loci of interest is determined. In some of these embodiments, the multiple loci of interest are on multiple chromosomes.

**[0048]** In some embodiments, the sequence is determined by: (a) amplifying a locus of interest on a template DNA using a first and second primers, where the second primer contains a recognition site for a restriction enzyme such that digestion with the restriction enzyme generates a 5' overhang containing the locus of interest; (b) digesting the amplified DNA with the restriction enzyme that recognizes the recognition site on the second primer; (c) incorporating a nucleotide into the digested DNA of (b) by using the 5' overhang containing the locus of interest as a template; and (d) determining the sequence of the locus of interest by determining the sequence of the DNA of (c).

**[0049]** In other embodiments, the sequence is determined by: (a) amplifying alleles of a locus of interest on a template DNA using a first and second primers, where the second primer contains a recognition site for a restriction enzyme such that digestion with the restriction enzyme generates a 5' overhang containing the locus of interest; (b) digesting the amplified DNA with the restriction enzyme that recognizes the recognition site on the second primer; (c) incorporating nucleotides into the digested DNA of (b), where a nucleotide that terminates elongation, and is complementary to the locus of interest of an allele, is incorporated into the 5' overhang of said allele, and a nucleotide complementary to the locus of interest of a different allele is incorporated into the 5' overhang of said different allele, and the terminating nucleotide, which is complementary to a nucleotide in the 5' overhang of said different allele, is incorporated into the 5' overhang of said different allele; and (d) determining the sequence of the alleles of a locus of interest by determining the sequence of the DNA of (c).

**[0050]** In some embodiments, the restriction enzyme cuts DNA at a distance from the recognition site. In some of these embodiments, the recognition site includes but is for a Type IIS restriction enzyme, for example Alw I, Alw26 I, Bbs I, Bbv I, BceA I, Bmr I, Bsa I, Bst71 I, BsmA I, BsmB I, BsmF I, BspM I, Ear I, Fau I, Fok I, Hga I, Ple I, Sap I, SSfaN I, or Sthi32 I.

**[0051]** In some embodiments, the method of amplification maybe, for example, polymerase chain reaction, self-sustained sequence reaction, ligase chain reaction, rapid amplification of cDNA ends,

polymerase chain reaction and ligase chain reaction, Q-beta phage amplification, strand displacement amplification, or splice overlap extension polymerase chain reaction. In some embodiments, the method of amplification is by PCR. In some of these embodiments, an annealing temperature for cycle 1 of PCR is about the melting temperature of the portion of the 3' region of the second primer that anneals to the template DNA. In some of the latter embodiments, an annealing temperature for cycle 2 of PCR is about the melting temperature of the portion of the 3' region of the first primer that anneals to the template DNA. In some of the latter embodiments, an annealing temperature for the remaining cycles of PCR is at about the melting temperature of the entire second primer.

**[0052]** In some embodiments, the sequence of a locus of interest was determined using allele specific PCR, mass spectrometry, hybridization, primer extension, fluorescence polarization, fluorescence resonance energy transfer (FRET), fluorescence detection, sequencing, Sanger dideoxy sequencing, DNA microarray, southern blot, slot blot, dot blot, or MALDI-TOF mass spectrometry.

**[0053]** In some embodiments, the sequence of a locus of interest is determined by (1) amplification of the locus of interest; (2) hybridization of amplified loci to GeneCHIP array (3) washing GeneCHIP array; (4) staining the GeneCHIP array with detectable reagents; and (5) scanning GeneCHIP array. In some of these embodiments, the amplification method in (1) is polymerase chain reaction, self-sustained sequence reaction, ligase chain reaction, rapid amplification of cDNA ends, polymerase chain reaction and ligase chain reaction, Q-beta phage amplification, strand displacement amplification, or splice overlap extension polymerase chain reaction. In some embodiments, the method of amplification is by PCR. In some embodiments, the staining method comprises streptavidin phycoerythrin and biotinylated anti-streptavidin. In some embodiments, an agent that inhibits cell lysis has been added to the sample to inhibit the lysis of cells, if present, where the agent is membrane stabilizer, cross-linker, or cell lysis inhibitor. In some embodiments, the agent is a cell lysis inhibitor. In some of these embodiments, the cell lysis inhibitor is formalin at a percentage selected from the group consisting of: 0.0001-0.03%, 0.03-0.05%, 0.05-0.08%, 0.08-0.1%, 0.1-0.3%, 0.3-0.5%, 0.5-0.7%, 0.7-0.9%, 0.9-1.2%, 1.2-1.5%, 1.5-2%, or 2-3%. In some embodiments, the concentration of formalin in the sample is 0.1%.

**[0054]** In some embodiments, the sequence of a locus of interest is determined by (1) amplification of the locus of interest; (2) amplicon fragmentation; (3) hybridization of fragmented amplicons to CodeLink Arrays; (4) extension reaction to incorporate a nucleotide; and (5) detection of incorporated nucleotides. In some of these embodiments, the amplification method is polymerase chain reaction, self-sustained sequence reaction, ligase chain reaction, rapid amplification of cDNA ends, polymerase chain reaction and ligase chain reaction, Q-beta phage amplification, strand displacement amplification, or

splice overlap extension polymerase chain reaction. In some embodiments, the method of amplification is by PCR. In some embodiments, the amplicon fragmentation is by exonuclease digestion. In some embodiments, the incorporated nucleotide is a dideoxynucleotide or deoxynucleotide. In some embodiments, the incorporated nucleotide is labeled with a radioactive molecule, fluorescent molecule, antibody, antibody fragment, hapten, carbohydrate, biotin, derivative of biotin, phosphorescent moiety, luminescent moiety, electrochemiluminescent moiety, chromatic moiety, and moiety having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity. In some embodiments, the labeled nucleotide is labeled with a fluorescent molecule. In some embodiments, an agent that inhibits cell lysis has been added to the sample to inhibit the lysis of cells, if present, where the agent is membrane stabilizer, cross-linker, or cell lysis inhibitor. In some embodiments, the agent is a cell lysis inhibitor. In some of these embodiments, the cell lysis inhibitor is formalin at a percentage selected from the group consisting of: 0.0001-0.03%, 0.03-0.05%, 0.05-0.08%, 0.08-0.1%, 0.1-0.3%, 0.3-0.5%, 0.5-0.7%, 0.7-0.9%, 0.9-1.2%, 1.2-1.5%, 1.5-2%, or 2-3%. In some embodiments, the concentration of formalin in the sample is 0.1%.

**[0055]** In some embodiments, the sequence of a locus of interest is determined by using BeadArray Technology. In some embodiments, an agent that inhibits cell lysis has been added to the sample to inhibit the lysis of cells, if present, where the agent is membrane stabilizer, cross-linker, or cell lysis inhibitor. In some embodiments, the agent is a cell lysis inhibitor. In some of these embodiments, the cell lysis inhibitor is formalin at a percentage selected from the group consisting of: 0.0001-0.03%, 0.03-0.05%, 0.05-0.08%, 0.08-0.1%, 0.1-0.3%, 0.3-0.5%, 0.5-0.7%, 0.7-0.9%, 0.9-1.2%, 1.2-1.5%, 1.5-2%, or 2-3%. In some embodiments, the concentration of formalin in the sample is 0.1%.

**[0056]** In some embodiments, the sequence of a locus of interest is determined by (1) amplification of the locus of interest; (2) dephosphorylation of the unused reagents in (1); (3) in vitro transcription reaction of the products of (2); (4) RNase A cleavage of the products of (3); (5) mixing the products of (4) with CleanResin; (6) transfer products of (5) to SpectroCHIP; and (7) analysis of the SpectroCHIP. In some of these embodiments, the amplification method is polymerase chain reaction, self-sustained sequence reaction, ligase chain reaction, rapid amplification of cDNA ends, polymerase chain reaction and ligase chain reaction, Q-beta phage amplification, strand displacement amplification, or splice overlap extension polymerase chain reaction. In some embodiments, the method of amplification is by PCR. In some embodiments, the dephosphorylation reaction is catalyzed by shrimp alkaline phosphatase. In some embodiments, an agent that inhibits cell lysis has been added to the sample to inhibit the lysis of cells, if present, where the agent is membrane stabilizer, cross-linker, or cell lysis inhibitor. In some

embodiments, the agent is a cell lysis inhibitor. In some of these embodiments, the cell lysis inhibitor is formalin at a percentage selected from the group consisting of: 0.0001-0.03%, 0.03-0.05%, 0.05-0.08%, 0.08-0.1%, 0.1-0.3%, 0.3-0.5%, 0.5-0.7%, 0.7-0.9%, 0.9-1.2%, 1.2-1.5%, 1.5-2%, or 2-3%. In some embodiments, the concentration of formalin in the sample is 0.1%.

**[0057]** In some embodiments, the sequence of a locus of interest is determined by (1) amplification of a locus of interest; (2) dephosphorylation of the unused reagents in (1); (3) hybridization of a primer to the locus of interest; (4) incorporation of a nucleotide; (5) mixing the products of (4) with CleanResin; (6) transfer products of (5) to SpectroCHIP; and (7) analysis of the SpectroCHIP. In some of these embodiments, the amplification method is polymerase chain reaction, self-sustained sequence reaction, ligase chain reaction, rapid amplification of cDNA ends, polymerase chain reaction and ligase chain reaction, Q-beta phage amplification, strand displacement amplification, or splice overlap extension polymerase chain reaction. In some embodiments, the method of amplification is by PCR. In some embodiments, the dephosphorylation reaction is catalyzed by shrimp alkaline phosphatase. In some embodiments, hybridization of primer is adjacent to the locus of interest. In some embodiments, the incorporated nucleotide is a dideoxynucleotide or deoxynucleotide. In some embodiments, the incorporated nucleotide is labeled with radioactive molecule, fluorescent molecule, antibody, antibody fragment, hapten, carbohydrate, biotin, derivative of biotin, phosphorescent moiety, luminescent moiety, electrochemiluminescent moiety, chromatic moiety, and moiety having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity. In some embodiments, the labeled nucleotide is labeled with a fluorescent molecule. In some embodiments, an agent that inhibits cell lysis has been added to the sample to inhibit the lysis of cells, if present, where the agent is membrane stabilizer, cross-linker, or cell lysis inhibitor. In some embodiments, the agent is a cell lysis inhibitor. In some of these embodiments, the cell lysis inhibitor is formalin at a percentage selected from the group consisting of: 0.0001-0.03%, 0.03-0.05%, 0.05-0.08%, 0.08-0.1%, 0.1-0.3%, 0.3-0.5%, 0.5-0.7%, 0.7-0.9%, 0.9-1.2%, 1.2-1.5%, 1.5-2%, or 2-3%. In some embodiments, the concentration of formalin in the sample is 0.1%.

**[0058]** In some embodiments, the sequence of a locus of interest is determined by (1) amplification of the locus of interest; (2) exonuclease treatment of the products of (1); (3) single stranded DNA of (2) is annealed to an oligonucleotide; (4) incorporation of a nucleotide using the annealed template and primer of (3); (5) detection of the incorporated nucleotide. In some embodiments, the amplification method is by polymerase chain reaction, self-sustained sequence reaction, ligase chain reaction, rapid amplification of cDNA ends, polymerase chain reaction and ligase chain reaction, Q-beta phage

amplification, strand displacement amplification, or splice overlap extension polymerase chain reaction. In some embodiments, the method of amplification is by PCR. In some embodiments, the primer hybridizes adjacent to the locus of interest. In some embodiment, the incorporated nucleotide is a dideoxynucleotide or deoxynucleotide. In some embodiments, the incorporation reaction comprises two terminating nucleotides and two non-terminating nucleotides. In some embodiments, the incorporated nucleotide is labeled with radioactive molecule, fluorescent molecule, antibody, antibody fragment, hapten, carbohydrate, biotin, derivative of biotin, phosphorescent moiety, luminescent moiety, electrochemiluminescent moiety, chromatic moiety, and moiety having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity. In some embodiments, the terminating nucleotides are labeled with radioactive molecule, fluorescent molecule, antibody, antibody fragment, hapten, carbohydrate, biotin, derivative of biotin, phosphorescent moiety, luminescent moiety, electrochemiluminescent moiety, chromatic moiety, and moiety having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity. In some embodiments, the labeled nucleotide is labeled with a fluorescent molecule. In some embodiments, the terminating nucleotides are labeled with a fluorescent molecule. In some embodiments, an agent that inhibits cell lysis has been added to the sample to inhibit the lysis of cells, if present, where the agent is membrane stabilizer, cross-linker, or cell lysis inhibitor. In some embodiments, the agent is a cell lysis inhibitor. In some of these embodiments, the cell lysis inhibitor is formalin at a percentage selected from the group consisting of: 0.0001-0.03%, 0.03-0.05%, 0.05-0.08%, 0.08-0.1%, 0.1-0.3%, 0.3-0.5%, 0.5-0.7%, 0.7-0.9%, 0.9-1.2%, 1.2-1.5%, 1.5-2%, or 2-3%. In some embodiments, the concentration of formalin in the sample is 0.1%.

**[0059]** In some embodiments, the sequence of a locus of interest is determined by (1) amplification of the locus of interest, wherein the amplification reaction comprises a forward primer, a reverse primer, and a probe that anneals to the locus of interest, which is within the region of the amplicon; and (2) detection of the PCR products, wherein the amount of PCR product is used to determine the presence or absence of a specific genetic sequence. In some embodiments, the amplification is by PCR. In some embodiments, the probe contains a reporter dye at the 5' end and the 3' end contains a quenching dye. In some embodiments, the PCR products are detected using the ABI 7700 Sequence Detection System. In some embodiments, an agent that inhibits cell lysis has been added to the sample to inhibit the lysis of cells, if present, where the agent is membrane stabilizer, cross-linker, or cell lysis inhibitor. In some embodiments, the agent is a cell lysis inhibitor. In some of these embodiments, the cell lysis inhibitor is formalin at a percentage selected from the group consisting of: 0.0001-0.03%, 0.03-0.05%, 0.05-0.08%,

0.08-0.1%, 0.1-0.3%, 0.3-0.5%, 0.5-0.7%, 0.7-0.9%, 0.9-1.2%, 1.2-1.5%, 1.5-2%, or 2-3%. In some embodiments, the concentration of formalin in the sample is 0.1%.

**[0060]** In another aspect, the invention provides methods for determining the sequence of a locus of interest in a sample containing fetal DNA.

**[0061]** In some embodiments, the method for determining the sequence includes (a) amplifying a locus of interest on a template DNA using a first and second primers, where the second primer contains a recognition site for a restriction enzyme such that digestion with the restriction enzyme generates a 5' overhang containing the locus of interest; (b) digesting the amplified DNA with the restriction enzyme that recognizes the recognition site on the second primer; (c) incorporating a nucleotide into the digested DNA of (b) by using the 5' overhang containing the locus of interest as a template; and (d) determining the sequence of the locus of interest by determining the sequence of the DNA of (c).

**[0062]** In other embodiments, the method for determining the sequence includes (a) amplifying alleles of a locus of interest on a template DNA using a first and second primers, where the second primer contains a recognition site for a restriction enzyme such that digestion with the restriction enzyme generates a 5' overhang containing the locus of interest; (b) digesting the amplified DNA with the restriction enzyme that recognizes the recognition site on the second primer; (c) incorporating nucleotides into the digested DNA of (b), where a nucleotide that terminates elongation, and is complementary to the locus of interest of an allele, is incorporated into the 5' overhang of said allele, and a nucleotide complementary to the locus of interest of a different allele is incorporated into the 5' overhang of said different allele, and said terminating nucleotide, which is complementary to a nucleotide in the 5' overhang of said different allele, is incorporated into the 5' overhang of said different allele; and (d) determining the sequence of the alleles of a locus of interest by determining the sequence of the DNA of (c).

**[0063]** In embodiments, the sample is cell, tissue, blood, serum, plasma, saliva, urine, tears, vaginal secretion, sweat, umbilical cord blood, chorionic villi, amniotic fluid, embryonic tissue, embryo, a two-celled embryo, a four-celled embryo, an eight-celled embryo, a 16-celled embryo, a 32-celled embryo, a 64-celled embryo, a 128-celled embryo, a 256-celled embryo, a 512-celled embryo, a 1024-celled embryo, lymph fluid, cerebrospinal fluid, mucosa secretion, peritoneal fluid, ascitic fluid, fecal matter, or body exudates.

**[0064]** In another aspect, the invention provides methods for preparing a sample for analysis that include isolating free nucleic acid from a sample that contains nucleic acid, where an agent that inhibits cell lysis has been added to the sample to inhibit lysis of cells, if cells are present, where the agent is



membrane stabilizer, cross-linker, or cell lysis inhibitor. In this aspect, the portion of the sample that is to be analyzed is the free nucleic acid, not the cellular portion. In an embodiment, the present invention provides a method for isolating nucleic acid said method comprising (a) obtaining a sample containing nucleic acid; (b) adding a cell lysis inhibitor, cell membrane stabilizer, or cross-linker to the sample of (a); and (c) isolating nucleic acid. In an embodiment, the method is used for isolating free nucleic acid. In an embodiment, the method is used for isolating free fetal nucleic acid. In another embodiment, the present invention provides a method for isolating free fetal nucleic acid said method comprising (a) obtaining a sample containing nucleic acid; (b) adding a cell lysis inhibitor, cell membrane stabilizer, or cross-linker to the sample of (a); (c) isolating the plasma from the blood sample, wherein the plasma is isolated by centrifuging the blood sample; and (d) removing the supernatant, which contains the plasma, using procedures to minimize disruption of the "buffy-coat."

**[0065]** In some embodiments, the agent is cell lysis inhibitor, and in some of these embodiments, the cell lysis inhibitor is glutaraldehyde, derivatives of glutaraldehyde, formaldehyde, formalin, and derivatives of formaldehyde, crosslinkers, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA, DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS or compounds listed in Table XXIII.[[.]]

**[0066]** In some embodiments the cell lysis inhibitor is formalin. In some of these embodiments, the final concentration of formalin in the sample is 0.0001-0.03%, 0.03-0.05%, 0.05-0.08%, 0.08-0.1%, 0.1-0.3%, 0.3-0.5%, 0.5-0.7%, 0.7-0.9%, 0.9-1.2%, 1.2-1.5%, 1.5-2%, or 2-3%. In one embodiment, the final concentration of formalin in the sample is 0.1%.

**[0067]** An agent that stabilizes cell membranes may be added to the sample including but not limited to aldehydes, urea formaldehyde, phenol formaldehyde, DMAE (dimethylaminoethanol), cholesterol, cholesterol derivatives, high concentrations of magnesium, vitamin E, and vitamin E derivatives, calcium, calcium gluconate, taurine, niacin, hydroxylamine derivatives, bimoclomol, sucrose, astaxanthin, glucose, amitriptyline, isomer A hopane tetral phenylacetate, isomer B hopane tetral phenylacetate, citicoline, inositol, vitamin B, vitamin B complex, cholesterol hemisuccinate, sorbitol, calcium, coenzyme Q, ubiquinone, vitamin K, vitamin K complex, menaquinone, zonegran, zinc, ginkgo biloba extract, diphenylhydantoin, perftoran, polyvinylpyrrolidone, phosphatidylserine, tegretol, PABA, disodium

cromglycate, nedocromil sodium, phenytoin, zinc citrate, mexilit, dilantin, sodium hyaluronate, or polaxamer 188.

**[0068]** In another embodiment, an agent that prevents DNA destruction is added to the sample including but not limited to DNase inhibitors, zinc chloride, ethylenediaminetetraacetic acid, guanidine-HCl, guanidine isothiocyanate, N-lauroylsarcosine, and Na-dodecylsulphate.

**[0069]** In some embodiments, the sample is obtained from human, non-human, mammal, reptile, cattle, cat, dog, goat, swine, pig, monkey, ape, gorilla, bull, cow, bear, horse, sheep, poultry, mouse, rat, fish, dolphin, whale, or shark. In some of these embodiments, the sample is obtained from a human source.

**[0070]** In some embodiments, the sample containing nucleic acid is obtained from any nucleic acid containing source including but not limited to a cell, fetal cell, tissue, blood, serum, plasma, saliva, urine, tear, vaginal secretion, breast fluid, breast milk, sweat, umbilical cord blood, chorionic villi, amniotic fluid, embryonic tissue, embryo, a two-celled embryo, a four-celled embryo, an eight-celled embryo, a 16-celled embryo, a 32-celled embryo, a 64-celled embryo, a 128-celled embryo, a 256-celled embryo, a 512-celled embryo, a 1024-celled embryo lymph fluid, cerebrospinal fluid, mucosa secretion, peritoneal fluid, ascitic fluid, fecal matter, or body exudates. In some of these embodiments, the sample is blood.

**[0071]** In embodiments the sample is from a pregnant female. In an embodiment, the sample is obtained from a pregnant human female. In an embodiment, the sample is blood obtained from a pregnant female and, e.g., the nucleic acid is isolated from plasma obtained from blood of a pregnant female; the plasma is generated using procedures designed to minimize the amount of maternal cell lysis. In some of these embodiments, the blood is obtained from a human pregnant female when the fetus is at a gestational age of 0-4, 4-8, 8-12, 12-16, 16-20, 20-24, 24-28, 28-32, 32-36, 36-40, 40-44, 44-48, 48-52, or more than 52 weeks. In some of these embodiments, the sample is obtained from plasma from the blood.

**[0072]** In some embodiments, the isolation of nucleic acid includes a centrifugation step; e.g., in some embodiments free nucleic acid is isolated from plasma obtained from blood, for example from a pregnant female. In some embodiments, the centrifugation step is performed with the centrifuge braking power set to zero (the centrifuge comes to a stop by natural deceleration). In some embodiments, the centrifugation step is performed at a speed of 0-50 rpm, 50-100 rpm, 100-200 rpm, 200-300 rpm, 300-400 rpm, 400-500 rpm, 500-600 rpm, 600-700 rpm, 700-800 rpm, 800-900 rpm, 900-1000 rpm, 1000-2000 rpm, 2000-3000 rpm, 3000-4000 rpm, 4000-5000 rpm, 5000-6000 rpm, 6000-7000 rpm, 7000-8000 rpm, or greater than 8000 rpm. In one embodiment, the blood, e.g., from the pregnant female, is

centrifuged at a speed less than 4000 rpm. In another embodiment, the acceleration power of the centrifuge is not used.

**[0073]** In another aspect, the invention provides a method for detecting a chromosomal abnormality by (a) determining the sequence of alleles of a locus of interest from template DNA, and (b) quantitating the relative amount of the alleles at a heterozygous locus of interest that was identified from the locus of interest of (a), wherein said relative amount is expressed as a ratio, and wherein said ratio indicates the presence or absence of a chromosomal abnormality.

**[0074]** In yet another aspect, the invention provides compositions.

**[0075]** In one embodiment, the invention provides a composition containing fetal DNA and maternal DNA, where the percentage of free fetal DNA in the total free DNA of the composition is about 15-16% fetal DNA, about 16-17% fetal DNA, about 17-18% fetal DNA, about 18-19% fetal DNA, about 19-20% fetal DNA, about 20-21% fetal DNA, about 21-22% fetal DNA, about 22-23% fetal DNA, about 23-24% fetal DNA, about 24-25% fetal DNA, about 25-35% fetal DNA, about 35-45% fetal DNA, about 45-55% fetal DNA, about 55-65% fetal DNA, about 65-75% fetal DNA, about 75-85% fetal DNA, about 85-90% fetal DNA, about 90-91% fetal DNA, about 91-92% fetal DNA, about 92-93% fetal DNA, about 93-94% fetal DNA, about 94-95% fetal DNA, about 95-96% fetal DNA, about 96-97% fetal DNA, about 97-98% fetal DNA, about 98-99% fetal DNA, or about 99-99.7% fetal DNA.

**[0076]** In another embodiment, the invention provides a composition containing fetal DNA and maternal DNA, where the percentage of free fetal DNA in the total free DNA of the composition is about 15-16% fetal DNA, about 16-17% fetal DNA, about 17-18% fetal DNA, about 18-19% fetal DNA, about 19-20% fetal DNA, about 20-21% fetal DNA, about 21-22% fetal DNA, about 22-23% fetal DNA, about 23-24% fetal DNA, about 24-25% fetal DNA, about 25-35% fetal DNA, about 35-45% fetal DNA, about 45-55% fetal DNA, about 55-65% fetal DNA, about 65-75% fetal DNA, about 75-85% fetal DNA, about 85-90% fetal DNA, about 90-91% fetal DNA, about 91-92% fetal DNA, about 92-93% fetal DNA, about 93-94% fetal DNA, or about 94-95% fetal DNA.

**[0077]** In yet another aspect, the invention provides a prenatal diagnostic method including analyzing a composition comprising fetal DNA and maternal DNA, where the percentage of free fetal DNA in the total free DNA of the composition is about 15-16% fetal DNA, about 16-17% fetal DNA, about 17-18% fetal DNA, about 18-19% fetal DNA, about 19-20% fetal DNA, about 20-21% fetal DNA, about 21-22% fetal DNA, about 22-23% fetal DNA, about 23-24% fetal DNA, about 24-25% fetal DNA, about 25-35% fetal DNA, about 35-45% fetal DNA, about 45-55% fetal DNA, about 55-65% fetal DNA, about 65-75% fetal DNA, about 75-85% fetal DNA, about 85-90% fetal DNA, about 90-91% fetal DNA,

about 91-92% fetal DNA, about 92-93% fetal DNA, about 93-94% fetal DNA, or about 94-95% fetal DNA.

**[0078]** In still yet another aspect, the invention provides a kit for use in any of the methods of the invention, where the kit contains a set of primers used in the method, where the second primer contains a sequence that generates a recognition site for a restriction enzyme such that digestion with the restriction enzyme generates a 5' overhang containing the locus of interest, and a set of instructions.

#### **BRIEF DESCRIPTION OF THE FIGURES**

**[0079]** FIG. 1A. ~~A Schematic~~ A schematic diagram depicting a double stranded DNA molecule. A pair of primers, depicted as bent arrows, flank the locus of interest, depicted as a triangle symbol at base N14. The locus of interest can be a single nucleotide polymorphism, point mutation, insertion, deletion, translocation, etc. Each primer contains a restriction enzyme recognition site about 10 by from the 5' terminus depicted as region "a" in the first primer and as region "d" in the second primer. Restriction recognition site "a" can be for any type of restriction enzyme but recognition site "d" is for a restriction enzyme, which cuts "n" nucleotides away from its recognition site and leaves a 5' overhang and a recessed 3' end. Examples of such enzymes include but are not limited to BceAI and BsmF I. The 5' overhang serves as a template for incorporation of a nucleotide into the 3' recessed end.

**[0080]** The first primer is shown modified with biotin at the 5' end to aid in purification. The sequence of the 3' end of the primers is such that the primers anneal at a desired distance upstream and downstream of the locus of interest. The second primer anneals close to the locus of interest; the annealing site, which is depicted as region "c," is designed such that the 3' end of the second primer anneals one base away from the locus of interest. The second primer can anneal any distance from the locus of interest provided that digestion with the restriction enzyme, which recognizes the region "d" on this primer, generates a 5' overhang that contains the locus of interest. The first primer annealing site, which is depicted as region "b," is about 20 bases.

**[0081]** FIG. 1B. A schematic diagram depicting the annealing and extension steps of the first cycle of amplification by PCR. The first cycle of amplification is performed at about the melting temperature of the 3' region, which anneals to the template DNA, of the second primer, depicted as region "c," and is 13 base pairs in this example. At this temperature, both the first and second primers anneal to their respective complementary strands and begin extension, depicted by dotted lines. In this first cycle, the second primer extends and copies the region b where the first primer can anneal in the next cycle.

**[0082]** FIG. 1C. A schematic diagram depicting the annealing and extension steps following denaturation in the second cycle of amplification of PCR. The second cycle of amplification is performed at a higher annealing temperature (TM2), which is about the melting temperature of the 20 by of the 3' region of the first primer that anneals to the template DNA, depicted as region "b." Therefore at TM2, the first primer, which contains region b' which is complementary to region b, can bind to the DNA that was copied in the first cycle of the reaction. However, at TM2 the second primer cannot anneal to the original template DNA or to DNA that was copied in the first cycle of the reaction because the annealing temperature is too high. The second primer can anneal to 13 bases in the original template DNA but TM2 is calculated at about the melting temperature of 20 bases.

**[0083]** FIG. 1D. A schematic diagram depicting the annealing and extension reactions after denaturation during the third cycle of amplification. In this cycle, the annealing temperature, TM3, is about the melting temperature of the entire second primer, including regions "c" and "d." The length of regions "c" + "d" is about 27-33 by long, and thus TM3 is significantly higher than TM1 and TM2. At this higher TM the second primer, which contain regions c' and d', anneals to the copied DNA generated in cycle 2.

**[0084]** FIG. 1E. A schematic diagram depicting the annealing and extension reactions for the remaining cycles of amplification. The annealing temperature for the remaining cycles is TM3, which is about the melting temperature of the entire second primer. At TM3, the second primer binds to templates that contain regions c' and d' and the first primer binds to templates that contain regions a' and b. By raising the annealing temperature successively in each cycle for the first three cycles, from TM1, TM2, and TM3, nonspecific amplification is significantly reduced.

**[0085]** FIG. 1F. A schematic diagram depicting the amplified locus of interest bound to a solid matrix.

**[0086]** FIG. 1G. A schematic diagram depicting the bound, amplified DNA after digestion with restriction enzyme "d." The "downstream" end is released into the supernatant, and can be removed by washing with any suitable buffer. The upstream end containing the locus of interest remains bound to the solid matrix.

**[0087]** FIG. 1H. A schematic diagram depicting the bound amplified DNA, after "filling in" with a labeled ddNTP. A DNA polymerase is used to "fill in" the base (N'14) that is complementary to the locus of interest (N14). In this example, only ddNTPs are present in this reaction, such that only the locus of interest or SNP of interest is filled in.

[0088] FIG. 1I. A schematic diagram depicting the labeled, bound DNA after digestion with restriction enzyme “a.” The labeled DNA is released into the supernatant, which can be collected to identify the base that was incorporated.

[0089] FIG. 2. A schematic diagram depicting double stranded DNA templates with  $n$  number of loci of interest and  $n$  number of primer pairs,  $x_1, y_1$  to  $x_n, y_n$ , specifically annealed such that a primer flanks each locus of interest. The first primers are biotinylated at the 5' end, depicted by •, and contain a restriction enzyme recognition site, “a”, which can be any type of restriction enzyme. The second primers contain a restriction enzyme recognition site, “d,” where “d” is a recognition site for a restriction enzyme that cuts “ $n$ ” nucleotides away from its recognition site, and generates a 5' overhang containing the locus of interest and a recessed 3' end. The second primers anneal adjacent to the respective loci of interest. The exact position of the restriction enzyme site “d” in the second primers is designed such that digesting the PCR product of each locus of interest with restriction enzyme “d” generates a 5' overhang containing the locus of interest and a 3' recessed end. The annealing sites of the first primers are about 20 bases long and are selected such that each successive first primer is further away from its respective second primer. For example, if at locus 1 the 3' ends of the first and second primers are  $Z$  base pairs apart, then at locus 2, the 3' ends of the first and second primers are  $Z + K$  base pairs apart, where  $K = 1, 2, 3$  or more than three bases. Primers for locus  $N$  are  $Z_{N-1} + K$  base pairs apart. The purpose of making each successive first primer further apart from their respective second primers is such that the “filled in” restriction fragments (generated after amplification, purification, digestion and labeling as described in FIGS. 1B-1I) differ in size and can be resolved, for example by electrophoresis, to allow detection of each individual locus of interest.

[0090] FIG. 3A-3C: PCR amplification of SNPs using multiple annealing temperatures. A sample containing genomic DNA templates from thirty-six human volunteers was analyzed for the following four SNPs: SNP HC21S00340 (lane 1), identification number as assigned in the Human Chromosome 21 cSNP Database, located on chromosome 21; SNP TSC 0095512 (lane 2), located on chromosome 1, SNP TSC 0214366 (lane 3), located on chromosome 1; and SNP TSC 0087315 (lane 4), located on chromosome 1. Each SNP was amplified by PCR using three different annealing temperature protocols, herein referred to as the low stringency annealing temperature; medium stringency annealing temperature; and high stringency annealing temperature. Regardless of the annealing temperature protocol, each SNP was amplified for 40 cycles of PCR. The denaturation step for each PCR reaction was performed for 30 seconds at 95°C. 3A. Photograph of a gel demonstrating PCR amplification of the 4 different SNPs using the low stringency annealing temperature protocol. 3B. Photograph of a gel demonstrating PCR

amplification of the 4 different SNPs using medium stringency annealing temperature protocol. 3C. Photograph of a gel demonstrating PCR amplification of the 4 different SNPs using the high stringency annealing temperature protocol.

[0091] FIG. 4A. (From top to bottom: SEQ ID NOS: 17, 667, 668, 18.) A depiction of the DNA sequence of SNP HC21S00027, as assigned by the Human Chromosome 21 cSNP database, located on chromosome 21. A first primer and a second primer are indicated above and below, respectively, the sequence of HC21S00027. The first primer is biotinylated and contains the restriction enzyme recognition site for EcoRI. The second primer contains the restriction enzyme recognition site for BsmF I and contains 13 bases that anneal to the DNA sequence. The SNP is indicated by R (A/G) and r (T/C) (complementary to R).

[0092] FIG. 4B. (From top to bottom: SEQ ID NOS: 17, 667, 668, 19.) A depiction of the DNA sequence of SNP HC21S00027, as assigned by the Human Chromosome 21 cSNP database, located on chromosome 21. A first primer and a second primer are indicated above and below, respectively, the sequence of HC21S00027. The first primer is biotinylated and contains the restriction enzyme recognition site for EcoRI. The second primer contains the restriction enzyme recognition site for BceA I and has 13 bases that anneal to the DNA sequence. The SNP is indicated by R (A/G) and r (T/C) (complementary to R).

[0093] FIG. 4C. (From top to bottom: SEQ ID NOS: 11, 669, 670, 20.) A depiction of the DNA sequence of SNP TSC0095512 from chromosome 1. The first primer and the second primer are indicated above and below, respectively, the sequence of TSC0095512. The first primer is biotinylated and contains the restriction enzyme recognition site for EcoRI. The second primer contains the restriction enzyme recognition site for BsmF I and has 13 bases that anneal to the DNA sequence. The SNP is indicated by S (G/C) and s (C/G) (complementary to S).

[0094] FIG. 4D. (From top to bottom: SEQ ID NOS: 11, 669, 670, 12.) A depiction of the DNA sequence of SNP TSC0095512 from chromosome 1. The first primer and the second primer are indicated above and below, respectively, the sequence of TSC0095512. The first primer is biotinylated and contains the restriction enzyme recognition site for EcoRI. The second primer contains the restriction enzyme recognition site for BceA I and has 13 bases that anneal to the DNA sequence. The SNP is indicated by S (G/C) and s (C/G) (complementary to S).

[0095] FIGS. 5A-5D. (FIG. 5A: SEQ ID NOS: 671 (top) and 672 (bottom); FIG. 5B: SEQ ID NOS: 673 (top) and 674 (bottom); FIG. 5C: SEQ ID NOS: 675 (top) and 676 (bottom); FIG. 5D: SEQ ID NOS: 677 (top) and 678 (bottom)). A schematic diagram depicting the nucleotide sequences of SNP

HC21S00027 (FIGS. 5A and 5B) and SNP TSC0095512 (FIGS. 5C and 5D) after amplification with the primers described in FIGS. 4A-4D. Restriction sites in the primer sequence are indicated in bold.

**[0096]** FIGS. 6A-6D. A schematic diagram depicting the nucleotide sequences of each amplified SNP after digestion with the appropriate Type IIS restriction enzyme. FIGS. 6A and 6B (SEQ ID NOS: 679 (upper left), 680 (upper right), 681 (lower left) and 682 (lower right)) and 6C and 6D (SEQ ID NOS: 679 (upper left), 683 (upper right), 684 (lower left) and 685 (lower right)) depict fragments of SNP HC21S00027 digested with the Type IIS restriction enzymes BsmF I and BceA I, respectively. FIGS. 6C and 6D (SEQ ID NOS: 686 (upper left), 687 (upper right), 688 (lower left) and 689 (lower right)) and 6E and 6F (SEQ ID NOS: 686 (upper left), 690 (upper right), 691 (lower left) and 692 (lower right)) depict fragments of SNP TSC0095512 digested with the Type IIS restriction enzymes BsmF I and BceA I, respectively.

**[0097]** FIGS. 7A-7D. A schematic diagram depicting the incorporation of a fluorescently labeled nucleotide using the 5' overhang of the digested SNP site as a template to "fill in" the 3' recessed end. FIGS. 7A and 7B (SEQ ID NOS: 693 (top) and 694 (bottom)) and 7C and 7D (SEQ ID NOS: 693 (top) and 695 (bottom)) depict the digested SNP HC21S00027 locus with an incorporated labeled ddNTP (\*R<sup>dd</sup> = fluorescent dideoxy nucleotide). FIGS. 7C and 7D (SEQ ID NOS: 696 (top) and 697 (bottom)) and 7E and 7F (SEQ ID NOS: 696 (top) and 698 (bottom)) depict the digested SNP TSC0095512 locus with an incorporated labeled ddNTP (\*S<sup>dd</sup> = fluorescent dideoxy nucleotide). The use of ddNTPs ensures that the 3' recessed end is extended by one nucleotide, which is complementary to the nucleotide of interest or SNP site present in the 5' overhang.

**[0098]** FIG. 7E. (From top to bottom: SEQ ID NOS: 693, 694, 699, 694, 700, 694, 701, 694.) A schematic diagram depicting the incorporation of dNTPs and a ddNTP into the 5' overhang containing the SNP site. SNP HC21S00007 was digested with BsmF I, which generates a four base 5' overhang. The use of a mixture of dNTPs and ddNTPs allows the 3' recessed end to be extended one nucleotide (a ddNTP is incorporated first); two nucleotides (a dNTP is incorporated followed by a ddNTP); three nucleotides (two dNTPs are incorporated, followed by a ddNTP); or four nucleotides (three dNTPs are incorporated, followed by a ddNTP). All four products can be separated by size, and the incorporated nucleotide detected (\*R<sup>dd</sup> = fluorescent dideoxy nucleotide). Detection of the first nucleotide, which corresponds to the SNP or locus site, and the next three nucleotides provides an additional level of quality assurance. The SNP is indicated by R (A/G) and r (T/C) (complementary to R).

**[0099]** FIGS. 8A-8D. Release of the "filled in" SNP from the solid support matrix, i.e. streptavidin coated well. SNP HC21S00027 is shown in FIGS. 8A and 8B (SEQ ID NOS: 702 (upper left), 703



(upper right), 704 (lower left) and 705 (lower right)) and 8B (SEQ ID NOS: 702 (upper left), 703 (upper right), 704 (lower left) and 706 (lower right)), while SNP TSC0095512 is shown in FIGS. 8C and 8D (SEQ ID NOS: 707 (upper left), 708 (upper right), 709 (lower left) and 710 (lower right)) and 8D (SEQ ID NOS: 707 (upper left), 708 (upper right), 709 (lower left) and 711 (lower right)). The “filled in” SNP is free in solution, and can be detected.

**[0100]** FIG. 9A. Sequence analysis of SNP HC21S00027 digested with BceAI. Four “fill in” reactions are shown; each reaction contained one fluorescently labeled nucleotide, ddGTP, ddATP, ddTTP, or ddCTP, and unlabeled ddNTPs. The 5’ overhang generated by digestion with BceA I and the expected nucleotides at this SNP site are indicated.

**[0101]** FIG. 9B. Sequence analysis of SNP TSC0095512. SNP TSC0095512 was amplified with a second primer that contained the recognition site for BceA I, and in a separate reaction, with a second primer that contained the recognition site for BsmF I. Four fill in reactions are shown for each PCR product; each reaction contained one fluorescently labeled nucleotide, ddGTP, ddATP, ddTTP, or ddCTP, and unlabeled ddNTPs. The 5’ overhang generated by digestion with BceA I and with BsmF I and the expected nucleotides are indicated.

**[0102]** FIG. 9C. Sequence analysis of SNP TSC0264580 after amplification with a second primer that contained the recognition site for BsmF I. Four fill in reactions are shown; each reaction contained one fluorescently labeled nucleotide, which was ddGTP, ddATP, ddTTP, or ddCTP and unlabeled ddNTPs. Two different 5’ overhangs are depicted: one represents the DNA molecules that were cut 11 nucleotides away on the sense strand and 15 nucleotides away on the antisense strand and the other represents the DNA molecules that were cut 10 nucleotides away on the sense strand and 14 nucleotides away on the antisense strand. The expected nucleotides also are indicated.

**[0103]** FIG. 9D. Sequence analysis of SNP HC21S00027 amplified with a second primer that contained the recognition site for BsmF I. A mixture of labeled ddNTPs and unlabeled dNTPs was used to fill in the 5’ overhang generated by digestion with BsmF I. Two different 5’ overhangs are depicted: one represents the DNA molecules that were cut 11 nucleotides away on the sense strand and 15 nucleotides away on the antisense strand and the other represents the DNA molecules that were cut 10 nucleotides away on the sense strand and 14 nucleotides away on the antisense strand. The nucleotide upstream from the SNP, the nucleotide at the SNP site (the sample contained DNA templates from 36 individuals; both nucleotides would be expected to be represented in the sample), and the three nucleotides downstream of the SNP are indicated.

**[0104]** FIG. 10. Sequence analysis of multiple SNPs. SNPs HC21S00131, and HC21S00027, which are located on chromosome 21, and SNPs TSC0087315, SNP TSC0214366, SNP TSC0413944, and SNP TSC0095512, which are on chromosome 1, were amplified in separate PCR reactions with second primers that contained a recognition site for BsmF I. The primers were designed so that each amplified locus of interest was of a different size. After amplification, the reactions were pooled into a single sample, and all subsequent steps of the method performed (as described for FIGS. 1F-1I) on that sample. Each SNP and the nucleotide found at each SNP are indicated.

**[0105]** FIG. 11A-11B. Quantification of the percentage of fetal DNA in maternal blood. Blood was obtained from a pregnant human female with informed consent. DNA was isolated and serial dilutions were made to determine the percentage of fetal DNA present in the sample. The SRY gene, which is located on chromosome Y, was used to detect fetal DNA. The cystic fibrosis gene, which is located on chromosome 7, was used to detect both maternal and fetal DNA. 11 A. Amplification of the SRY gene and the cystic fibrosis gene using a DNA template isolated from a blood sample that was treated with EDTA. 11B. Amplification of the SRY gene and the cystic fibrosis gene using a DNA template that was isolated from a blood sample that was treated with formalin and EDTA.

**[0106]** FIG. 12. Genetic analysis of an individual previously genotyped with Trisomy 21 (Down's Syndrome). Blood was collected, with informed consent, from an individual who had previously been genotyped with trisomy 21. DNA was isolated and two SNPs on chromosome 21 and two SNPs on chromosome 13 were genotyped. As shown in the photograph of the gel, the SNPs at chromosome 21 show disproportionate ratios of the two nucleotides. Visual inspection of the gel demonstrates that one nucleotide of the two nucleotides at the SNP sites analyzed for chromosome 21 is of greater intensity, suggesting it is not present in a 50:50 ratio. However, visual inspection of the gel suggests that the nucleotides at the heterozygous SNP sites analyzed on chromosome 13 are present in the expected 50:50 ratio.

**[0107]** FIG. 13. Sequence determination of both alleles of SNPs TSC0837969, TSC0034767, TSC1130902, TSC0597888, TSC0195492, TSC0607185 using one fluorescently labeled nucleotide. Labeled ddGTP was used in the presence of unlabeled dATP, dCTP, dTTP to fill-in the overhang generated by digestion with BsmF I. The nucleotide preceding the variable site on the strand that was filled-in was not guanine, and the nucleotide after the variable site on the strand that was filled in was not guanine. The nucleotide two bases after the variable site on the strand that was filled-in was guanine. Alleles that contain guanine at variable site are filled in with labeled ddGTP. Alleles that do not contain

guanine are filled in with unlabeled dATP, dCTP, or dTTP, and the polymerase continues to incorporate nucleotides until labeled ddGTP is filled in at position 3 complementary to the overhang.

**[0108]** FIG. 14. Identification of SNPs with alleles that are variable within the population. The sequences of both alleles of seven SNPs located on chromosome 13 were determined using a template DNA comprised of DNA obtained from two hundred and forty five individuals. Labeled ddGTP was used in the presence of unlabeled dATP, dCTP, dTTP to fill-in the overhang generated by digestion with BsmF I. The nucleotide preceding the variable site on the strand that was filled-in was not guanine, and the nucleotide after the variable site on the strand that was filled in was not guanine. The nucleotide two bases after the variable site on the strand that was filled-in was guanine. Alleles that contain guanine at variable site are filled in with labeled ddGTP. Alleles that do not contain guanine are filled in with unlabeled dATP, dCTP, or dTTP, and the polymerase continues to incorporate nucleotides until labeled ddGTP is filled in at position 3 complementary to the overhang.

**[0109]** FIG. 15. Determination of the ratio for one allele to the other allele at heterozygous SNPs. The observed nucleotides for SNP TSC0607185 are cytosine (referred to as allele 1) and thymidine (referred to as allele 2) on the sense strand. The ratio of allele 2 to allele 1 was calculated using template DNA isolated from five individuals. The ratio of allele 2 to allele 1 (allele 2 / allele 1) was consistently 1:1.

**[0110]** The observed nucleotides for SNP TSC1130902 are guanine (referred to as allele 1) and adenine (referred to as allele 2) on the sense strand. The ratio of allele 2 to allele 1 was calculated using template DNA isolated from five individuals. The ratio of allele 2 to allele 1 (allele 2 / allele 1) was consistently 75:25.

**[0111]** FIG. 16. The percentage of allele 2 to allele 1 at SNP TSC0108992 remains linear when calculated on template DNA containing an extra copy of chromosome 21. SNP TSC0108992 was amplified using template DNA from four individuals, and two separate fill-in reactions (labeled as A and B) were performed for each PCR reaction (labeled 1 through 4). The calculated percentage of allele 2 to allele 1 on template DNA from normal individuals was 0.47. The deviation from the theoretically predicted percentage of 0.50 remained linear on template DNA isolated from an individual with Down's syndrome.

**[0112]** FIG. 17A. Analysis of a SNP located on chromosome 21 from template DNA isolated from an individual with a normal genetic karyotype. SNP TSC0108992 was amplified using the methods described herein, and after digestion with the type IIS restriction enzyme BsmF I, the 5' overhang was filled in using labeled ddTTP, and unlabeled dATP, dCTP, and dGTP. Three separate PCR reactions

were performed, and each PCR reaction was split into two samples. The percentage of allele 2 at the SNP site ( $\text{allele 2} / (\text{allele 2} + \text{allele 1})$ ) was calculated, which resulted in mean of 0.50.

**[0113]** FIG. 17B. Analysis of a SNP located on chromosome 21 from template DNA isolated from an individual with a trisomy 21 genetic karyotype. SNP TSC0108992 was amplified using the methods described herein, and after digestion with the type IIS restriction enzyme BsmF I, the 5' overhang was filled in using labeled ddTTP, and unlabeled dATP, dCTP, and dGTP. Three separate PCR reactions were performed, and each PCR reaction was split into two samples. The percentage of allele 2 at the SNP site ( $\text{allele 2} / (\text{allele 2} + \text{allele 1})$ ) was calculated, which resulted in mean of 0.30.

**[0114]** FIG. 17C. Analysis of a SNP located on chromosome 21 from a mixture comprised of template DNA from an individual with Trisomy 21, and template DNA from an individual with a normal genetic karyotype in a ratio of 3:1 (Trisomy 21: Normal). SNP TSC0108992 was amplified from the mixture of template DNA using the methods described herein, and after digestion with the type IIS restriction enzyme BsmF I, the 5' overhang was filled in using labeled ddTTP, and unlabeled dATP, dCTP, and dGTP. Three separate PCR reactions were performed, and each PCR reaction was split into two samples. The percentage of allele 2 at the SNP site ( $\text{allele 2} / (\text{allele 2} + \text{allele 1})$ ) was calculated, which resulted in mean of 0.319.

**[0115]** FIG. 17D. Analysis of a SNP located on chromosome 21 from a mixture comprised of template DNA from an individual with Trisomy 21, and template DNA from an individual with a normal genetic karyotype in a ratio of 1:1 (Trisomy 21: Normal). SNP TSC0108992 was amplified from the mixture of template DNA using the methods described herein, and after digestion with the type IIS restriction enzyme BsmF I, the 5' overhang was filled in using labeled ddTTP, and unlabeled dATP, dCTP, and dGTP. Three separate PCR reactions were performed, and each PCR reaction was split into two samples. The percentage of allele 2 at the SNP site ( $\text{allele 2} / (\text{allele 2} + \text{allele 1})$ ) was calculated, which resulted in mean of 0.352.

**[0116]** FIG. 17E. Analysis of a SNP located on chromosome 21 from a mixture comprised of template DNA from an individual with Trisomy 21, and template DNA from an individual with a normal genetic karyotype in a ratio of 1:2.3 (Trisomy 21: Normal). SNP TSC0108992 was amplified from the mixture of template DNA using the methods described herein, and after digestion with the type IIS restriction enzyme BsmF I, the 5' overhang was filled in using labeled ddTTP, and unlabeled dATP, dCTP, and dGTP. Three separate PCR reactions were performed, and each PCR reaction was split into two samples. The percentage of allele 2 at the SNP site ( $\text{allele 2} / (\text{allele 2} + \text{allele 1})$ ) was calculated, which resulted in mean of 0.382.

**[0117]** FIG. 17F. Analysis of a SNP located on chromosome 21 from a mixture comprised of template DNA from an individual with Trisomy 21, and template DNA from an individual with a normal genetic karyotype in a ratio of 1:4 (Trisomy 21: Normal). SNP TSC0108992 was amplified from the mixture of template DNA using the methods described herein, and after digestion with the type IIS restriction enzyme BsmF I, the 5' overhang was filled in using labeled ddTTP, and unlabeled dATP, dCTP, and dGTP. Three separate PCR reactions were performed, and each PCR reaction was split into two samples. The percentage of allele 2 at the SNP site ( $\text{allele 2} / (\text{allele 2} + \text{allele 1})$ ) was calculated, which resulted in mean of 0.397.

**[0118]** FIG. 18A. Agarose gel analysis of nine (9) SNPs amplified from template DNA. Each of the nine SNPs were amplified from genomic DNA using the methods described herein. Lane 1 corresponds to SNP TSC0397235, lane 2 corresponds to TSC0470003, lane 3 corresponds to TSC1649726, lane 4 corresponds to TSC1261039, lane 5 corresponds to TSC0310507, lane 6 corresponds to TSC1650432, lane 7 corresponds to TSC1335008, lane 8 corresponds to TSC0128307, and lane 9 corresponds to TSC0259757.

**[0119]** FIG. 18B. The original template DNA was amplified using 12 base primers that annealed to various regions on chromosome 13. One hundred different primer sets were used to amplify regions throughout chromosome 13. For each of the nine SNPs, a primer that annealed approximately 130 bases from the locus of interest and 130 bases downstream of the locus of interest were used. This amplification reaction, which contained a total of 100 different primer sets, was used to amplify the regions containing the loci of interest. The resulting PCR product was used in a subsequent PCR reaction, wherein each of the nine SNPs were individually amplified using a first primer and a second primer, wherein the second primer contained the binding site for the type IIS restriction enzyme BsmF I. SNPs were loaded in the same order as FIG. 18A.

**[0120]** FIG. 19A. Quantification of the percentage of allele 2 to allele 1 for SNP TSC047003 on original template DNA (IA) and multiplexed template DNA (M1-M3), wherein the DNA was first amplified using 12 base primers that annealed 150 bases upstream and downstream of the loci of interest. Then, three separate PCR reactions were performed on the multiplexed template DNA, using a first and second primer.

**[0121]** FIG. 19B. Quantification of the percentage of allele 2 to allele 1 for SNP TSC1261039 on original template DNA (IA) and multiplexed template DNA (M1-M3), wherein the DNA was first amplified using 12 base primers that annealed 150 bases upstream and downstream of the loci of interest.

Then, three separate PCR reactions were performed on the multiplexed template DNA, using a first and second primer.

**[0122]** FIG. 19C. Quantification of the percentage of allele 2 to allele 1 for SNP TSC310507 on original template DNA (IA) and multiplexed template DNA (M1-M3), wherein the DNA was first amplified using 12 base primers that annealed 150 bases upstream and downstream of the loci of interest. Then, three separate PCR reactions were performed on the multiplexed template DNA, using a first and second primer.

**[0123]** FIG. 19D. Quantification of the percentage of allele 2 to allele 1 for SNP TSC1335008 on original template DNA (IA) and multiplexed template DNA (M1-M3), wherein the DNA was first amplified using 12 base primers that annealed 150 bases upstream and downstream of the loci of interest. Then, three separate PCR reactions were performed on the multiplexed template DNA, using a first and second primer.

**[0124]** FIG. 20. Detection of fetal DNA from plasma DNA isolated from a pregnant female. Four SNPs wherein the maternal DNA was homozygous were analyzed on the plasma DNA. The maternal DNA was homozygous for adenine at TSC0838335 (lane 1), while the plasma DNA displayed a heterozygous pattern (lane 2). The guanine allele represented the fetal DNA, which was clearly distinguished from the maternal signal. Both the maternal DNA and the plasma DNA were homozygous for adenine at TSC0418134 (lanes 3 and 4). The maternal DNA was homozygous for guanine at TSC0129188 (lane 5), while the plasma DNA displayed a heterozygous pattern (lane 6). The adenine allele represented the fetal DNA. Both the maternal DNA and the plasma DNA were homozygous for adenine at TSC0501389 (lanes 7 and 8).

#### **DETAILED DESCRIPTION OF THE INVENTION**

**[0125]** The present invention provides a method for detecting genetic disorders, including but not limited to mutations, insertions, deletions, and chromosomal abnormalities, and is especially useful for the detection of genetic disorders of a fetus. The method is especially useful for detection of a translocation, addition, amplification, transversion, inversion, aneuploidy, polyploidy, monosomy, trisomy, trisomy 21, trisomy 13, trisomy 14, trisomy 15, trisomy 16, trisomy 18, trisomy 22, triploidy, tetraploidy, and sex chromosome abnormalities including but not limited to XO, XXY, XYY, and XXX. The method also provides a non-invasive technique for determining the sequence of fetal DNA and identifying mutations within the fetal DNA.

**[0126]** The invention is directed to a method for detecting chromosomal abnormalities, the method comprising: (a) determining the sequence of alleles of a locus of interest on a template DNA; and (b) quantitating a ratio for the alleles at a heterozygous locus of interest that was identified from the locus of interest of (a), wherein said ratio indicates the presence or absence of a chromosomal abnormality.

**[0127]** In another embodiment, the present invention provides a non-invasive method for determining the sequence of a locus of interest on fetal DNA, said method comprising: (a) obtaining a sample from a pregnant female; (b) adding a cell lysis inhibitor, cell membrane stabilizer or cross-linker to the sample of (a); (c) obtaining template DNA from the sample of (b), wherein said template DNA comprises fetal DNA and maternal DNA; and (d) determining the sequence of a locus of interest on template DNA.

**[0128]** In another ~~embodiment, the~~ embodiment, the present invention is directed to a method for isolating DNA, said method comprising (a) obtaining a sample containing nucleic acid; (b) adding a cell lysis inhibitor, cell membrane stabilizer or cross-linker to sample of (a); and (c) isolating the DNA.

**[0129]** In another ~~embodiment, the~~ embodiment, the present invention is directed to a method for isolating free DNA, said method comprising (a) obtaining a sample containing nucleic acid; (b) adding a cell lysis inhibitor, cell membrane stabilizer or cross-linker to sample of (a); and (c) isolating the DNA.

**[0130]** In another embodiment, the present invention is directed to a method for isolating free DNA from a sample containing nucleic acid to which a cell lysis inhibitor, cell membrane stabilizer or cross-linker has been added, said method ~~comprising~~ comprising isolating the DNA.

**[0131]** In another ~~embodiment~~ embodiment, the present invention is directed to a method for isolating free fetal DNA, said method comprising (a) obtaining a sample containing nucleic acid; (b) adding a cell lysis inhibitor, cell membrane stabilizer or cross-linker to sample of (a); and (c) isolating the DNA. In another embodiment, the DNA is isolated using any technique suitable in the art including but not limited to cesium chloride gradients, gradients, sucrose gradients, glucose gradients, centrifugation protocols, boiling, Qiagen purification systems, QIA DNA blood purification kit, HiSpeed Plasmid Maxi Kit, QIAfilter plasmid kit, Promega DNA purification systems, MangeSil Paramagnetic Particle based systems, Wizard SV technology, Wizard Genomic DNA purification kit, Amersham purification systems, GFX Genomic Blood DNA purification kit, Invitrogen Life Technologies Purification Systems, CONCERT purification system, Mo Bio Laboratories purification systems, UltraClean BloodSpin Kits, and UltraClean Blood DNA Kit.

**[0132]** In another ~~embodiment~~ embodiment, the present invention is directed to a method for isolating free fetal DNA from a sample containing nucleic acid to which a cell lysis inhibitor, cell

membrane stabilizer or cross-linker has been added, said method comprising isolating the DNA. In a preferred embodiment, the free fetal DNA is isolated from plasma or serum obtained from the blood of a pregnant female.

**[0133]** In another embodiment, the DNA is isolated using techniques and/or protocols that substantially reduce the amount of maternal DNA in the sample including but not limited to centrifuging the samples, with the braking power for the centrifuge set to zero (the brake on the centrifuge is not used), transferring the supernatant to a new tube with minimal or no disturbance of the “buffy-coat,” and transferring only a portion of the supernatant to a new tube. In a preferred embodiment, both acceleration power and braking power for the centrifuge are set to zero.

**[0134]** In another embodiment, the DNA is isolated using techniques and/or protocols that substantially reduce the amount of maternal DNA in the sample including but not limited to centrifuging the samples, with the acceleration power for the centrifuge set to zero, transferring the supernatant to a new tube with minimal or no disturbance of the “buffy-coat,” and transferring only a portion of the supernatant to a new tube.

**[0135]** In another embodiment, the “buffy-coat” is removed from the tube prior to removal of the supernatant using any applicable method including but not limited to using a syringe or needle to withdraw the “buffy-coat.”

**[0136]** In another embodiment, the braking power for the centrifuge is set at a percentage including but not limited to 1-5%, 5-10%, 10-20%, 20-30%, 30-40%, 40-50%, 50-60%, 60-70%, 70-80%, 80-90%, 90-95%, 95-99% of maximum braking power.

**[0137]** In another embodiment, the acceleration power for the centrifuge is set at a percentage including but not limited to 1-5%, 5-10%, 10-20%, 20-30%, 30-40%, 40-50%, 50-60%, 60-70%, 70-80%, 80-90%, 90-95%, 95-99% of maximum acceleration power.

**[0138]** In another embodiment, the present invention is directed to a composition comprising free fetal DNA and free maternal DNA, wherein the composition comprises a relationship of free fetal DNA to free maternal DNA including but not limited to at least about 15% free fetal DNA, at least about 20% free fetal DNA, at least about 30% free fetal DNA, at least about 40% free fetal DNA, at least about 50% free fetal DNA, at least about 60% free fetal DNA, at least about 70% free fetal DNA, at least about 80% free fetal DNA, at least about 90% free fetal DNA, at least about 91% free fetal DNA, at least about 92% free fetal DNA, at least about 93% free fetal DNA, at least about 94% free fetal DNA, at least about 95% free fetal DNA, at least about 96% free fetal DNA, at least about 97% free fetal DNA, at least about 98% free fetal DNA, at least about 99% free fetal DNA, and at least about 99.5% free fetal DNA.



**[0139]** In another embodiment, the present invention is directed to a method of using a composition comprising free fetal DNA and free maternal DNA for prenatal diagnostics, wherein the composition comprises a relationship of free fetal DNA to free maternal DNA including but not limited to at least about 15% free fetal DNA, at least about 20% free fetal DNA, at least about 30% free fetal DNA, at least about 40% free fetal DNA, at least about 50% free fetal DNA, at least about 60% free fetal DNA, at least about 70% free fetal DNA, at least about 80% free fetal DNA, at least about 90% free fetal DNA, at least about 91% free fetal DNA, at least about 92% free fetal DNA, at least about 93% free fetal DNA, at least about 94% free fetal DNA, at least about 95% free fetal DNA, at least about 96% free fetal DNA, at least about 97% free fetal DNA, at least about 98% free fetal DNA, at least about 99% free fetal DNA, and at least about 99.5% free fetal DNA.

**[0140]** In another embodiment, the present invention is directed to a composition comprising free fetal DNA and free maternal DNA, wherein the composition comprises a relationship of free fetal DNA to free maternal DNA including but not limited to about 13-15% free fetal DNA, about 15-16% free fetal DNA, about 16-17% free fetal DNA, about 17-18% free fetal DNA, about 18-19% free fetal DNA, about 19-20% free fetal DNA, about 20-21% free fetal DNA, about 21-22% free fetal DNA, about 22-23% free fetal DNA, about 23-24% free fetal DNA, about 24-25% free fetal DNA, about 25-35% free fetal DNA, about 35-45% free fetal DNA, about 45-55% free fetal DNA, about 55-65% free fetal DNA, about 65-75% free fetal DNA, about 75-85% free fetal DNA, about 85-90% free fetal DNA, about 90-91% free fetal DNA, about 91-92% free fetal DNA, about 92-93% free fetal DNA, about 93-94% free fetal DNA, about 94-95% free fetal DNA, about 95-96% free fetal DNA, about 96-97% free fetal DNA, about 97-98% free fetal DNA, about 98-99% free fetal DNA, and about 99-99.7% free fetal DNA.

**[0141]** In another embodiment, the present invention is directed to a method of using a composition comprising free fetal DNA and free maternal DNA for prenatal diagnostics, wherein the composition comprises a relationship of free fetal DNA to free maternal DNA including but not limited to about 13-15% free fetal DNA, about 15-16% free fetal DNA, about 16-17% free fetal DNA, about 17-18% free fetal DNA, about 18-19% free fetal DNA, about 19-20% free fetal DNA, about 20-21% free fetal DNA, about 21-22% free fetal DNA, about 22-23% free fetal DNA, about 23-24% free fetal DNA, about 24-25% free fetal DNA, about 25-35% free fetal DNA, about 35-45% free fetal DNA, about 45-55% free fetal DNA, about 55-65% free fetal DNA, about 65-75% free fetal DNA, about 75-85% free fetal DNA, about 85-90% free fetal DNA, about 90-91% free fetal DNA, about 91-92% free fetal DNA, about 92-93% free fetal DNA, about 93-94% free fetal DNA, about 94-95% free fetal DNA, about 95-96% free

fetal DNA, about 96-97% free fetal DNA, about 97-98% free fetal DNA, about 98-99% free fetal DNA, or about 99-99.7% free fetal DNA.

**[0142]** In another embodiment, the present invention is directed to a composition comprising free fetal DNA and free maternal DNA, wherein the composition comprises a relationship of free fetal DNA to free maternal DNA including but not limited a maximum of 13%-15% free fetal DNA, a maximum of 15-18% free fetal DNA, a maximum of 18-20% free fetal DNA, a maximum of 20-40% free fetal DNA, a maximum of 40-50% free fetal DNA, a maximum of 50-60% free fetal DNA, a maximum of 60-70% free fetal DNA, a maximum of 70-80% free fetal DNA, a maximum of 80-90% free fetal DNA, a ~~maximum~~maximum of 90-92% free fetal DNA, a ~~maximum~~maximum of 92-94% free fetal DNA, a maximum of 94-95% free fetal DNA, a ~~maximum~~maximum of 95-96% free fetal DNA, a ~~maximum~~maximum of 96-97% free fetal DNA, a ~~maximum~~maximum of 97-98% free fetal DNA, a maximum of 98-99% free fetal DNA, a ~~maximum~~maximum of 99-99.5% free fetal DNA, and a ~~maximum~~maximum of 99.5-99.9% free fetal DNA.

**[0143]** In another embodiment, the present invention is directed to a method of using a composition comprising free fetal DNA and free maternal DNA for prenatal diagnostics, wherein the composition comprises a relationship of free fetal DNA to free maternal DNA including but not limited a maximum of 13%-15% free fetal DNA, a maximum of 15-18% free fetal DNA, a maximum of 18-20% free fetal DNA, a maximum of 20-40% free fetal DNA, a maximum of 40-50% free fetal DNA, a maximum of 50-60% free fetal DNA, a maximum of 60-70% free fetal DNA, a maximum of 70-80% free fetal DNA, a maximum of 80-90% free fetal DNA, a maximum of 90-92% free fetal DNA, a maximum of 92-94% free fetal DNA, a maximum of 94-95% free fetal DNA, a maximum of 95-96% free fetal DNA, a maximum of 96-97% free fetal DNA, a maximum of 97-98% free fetal DNA, a maximum of 98-99% free fetal DNA, a maximum of 99-99.5% free fetal DNA, and a maximum of 99.5-99.9% free fetal DNA.

### **DNA Template**

**[0144]** By a “locus of interest” is intended a selected region of nucleic acid that is within a larger region of nucleic acid. A locus of interest can include but is not limited to 1-100, 1-50, 1-20, or 1-10 nucleotides, preferably 1-6, 1-5, 1-4, 1-3, 1-2, or 1 nucleotide(s).

**[0145]** As used herein, an “allele” is one of several alternate forms of a gene or non-coding regions of DNA that occupy the same position on a chromosome. The term allele can be used to describe DNA from any organism including but not limited to bacteria, viruses, fungi, protozoa, molds, yeasts, plants, humans, non-humans, animals, and archeabacteria.

**[0146]** For example, bacteria typically have one large strand of DNA. The term allele with respect to bacterial DNA refers to the form of a gene found in one cell as compared to the form of the same gene in a different bacterial cell of the same species.

**[0147]** Alleles can have the identical sequence or can vary by a single nucleotide or more than one nucleotide. With regard to organisms that have two copies of each chromosome, if both chromosomes have the same allele, the condition is referred to as homozygous. If the alleles at the two chromosomes are different, the condition is referred to as heterozygous. For example, if the locus of interest is SNP X on chromosome 1, and the maternal chromosome contains an adenine at SNP X (A allele) and the paternal chromosome contains a guanine at SNP X (G allele), the individual is heterozygous at SNP X.

**[0148]** As used herein, sequence means the identity of one nucleotide or more than one contiguous nucleotides in a polynucleotide. In the case of a single nucleotide, e.g., a SNP, “sequence” and “identity” are used interchangeably herein.

**[0149]** The term “chromosomal abnormality” refers to a deviation between the structure of the subject chromosome and a normal homologous chromosome. The term “normal” refers to the predominate karyotype or banding pattern found in healthy individuals of a particular species. A chromosomal abnormality can be numerical or structural, and includes but is not limited to aneuploidy, polyploidy, inversion, a trisomy, a monosomy, duplication, deletion, deletion of a part of a chromosome, addition, addition of a part of chromosome, insertion, a fragment of a chromosome, a region of a chromosome, chromosomal rearrangement, and translocation. A chromosomal abnormality can be correlated with presence of a pathological condition or with a predisposition to develop a pathological condition. As defined herein, a single nucleotide polymorphism (“SNP”) is not a chromosomal abnormality.

**[0150]** As used herein, incorporation of a nucleotide by a polymerase is referred to as an elongation reaction or a fill-in reaction interchangeably.

**[0151]** As used herein with respect to individuals, “mutant alleles” refers to variant alleles that are associated with a disease state.

**[0152]** The term “template” refers to any nucleic acid molecule that can be used for amplification in the invention. RNA or DNA that is not naturally double stranded can be made into double stranded DNA so as to be used as template DNA. Any double stranded DNA or preparation containing multiple, different double stranded DNA molecules can be used as template DNA to amplify a locus or loci of interest contained in the template DNA.

**[0153]** The template DNA can be obtained from any source including but not limited to humans, non-humans, mammals, reptiles, cattle, cats, dogs, goats, swine, pigs, monkeys, apes, gorillas, bulls, cows, bears, horses, sheep, poultry, mice, rats, fish, dolphins, whales, and sharks.

**[0154]** The template DNA can be from any appropriate sample including but not limited to, nucleic acid-containing samples of tissue, bodily fluid (for example, blood, serum, plasma, saliva, urine, tears, peritoneal fluid, ascitic fluid, vaginal secretion, breast fluid, breast milk, lymph fluid, cerebrospinal fluid or mucosa secretion), umbilical cord blood, chorionic villi, amniotic fluid, an embryo, a two-celled embryo, a four-celled embryo, an eight-celled embryo, a 16-celled embryo, a 32-celled embryo, a 64-celled embryo, a 128-celled embryo, a 256-celled embryo, a 512-celled embryo, a 1024-celled embryo, embryonic tissues, lymph fluid, cerebrospinal fluid, mucosa secretion, or other body exudate, fecal matter, an individual cell or extract of the such sources that contain the nucleic acid of the same, and subcellular structures such as mitochondria, using protocols well established within the art.

**[0155]** In one embodiment, the template DNA can be obtained from a sample of a pregnant female.

**[0156]** In another embodiment, the template DNA can be obtained from an embryo. In a preferred embodiment, the template DNA can be obtained from a single-cell of an embryo.

**[0157]** In one embodiment, the template DNA is fetal DNA. Fetal DNA can be obtained from sources including but not limited to maternal blood, maternal serum, maternal plasma, fetal cells, umbilical cord blood, chorionic villi, amniotic fluid, urine, saliva, cells or tissues.

**[0158]** In another embodiment, a cell lysis inhibitor is added to the sample including but not limited to formaldehyde, formaldehyde derivatives, formalin, glutaraldehyde, glutaraldehyde derivatives, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA, DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS or compounds listed in Table XXIII. In another embodiment, two, three, four, five or more than five cell lysis inhibitors can be added to the sample. In a preferred embodiment, formalin is present in the sample at a percentage including but not limited to 0.0001-0.03%, 0.03-0.05%, 0.05-0.08%, 0.08-0.1%, 0.1-0.3%, 0.3-0.5%, 0.5-0.7%, 0.7-0.9%, 0.9-1.2%, 1.2-1.5%, 1.5-2%, 2-3%, 3-5%, and greater than 5%. In another embodiment, any combination of cross-linker, cell membrane stabilizer, or cell lysis inhibitor can be added to the sample including but not limited to a cross-linker and a cell membrane stabilizer, a cross-linker and a cell lysis inhibitor, and a cell membrane stabilizer and a cell lysis inhibitor. More than one cross-linker can be used with more than one

cell membrane stabilizer. More than one cross-linker can be used with more than one cell lysis inhibitor. More than one cell membrane stabilizer can be used with more than cell lysis inhibitor.

**[0159]** In another embodiment, the cell lysis inhibitor is added to the sample such that lysis is less than about 10% of the cells. In a preferred embodiment, the cell lysis inhibitor is added to the sample such that lysis is less than about 5% of the cells. In a most preferred embodiment, the cell lysis inhibitor is added to the sample such that lysis is less than about 1% of the cells.

**[0160]** In another embodiment, a cell membrane stabilizer is added to the sample such that lysis is less than about 10% of the cells. In a preferred embodiment, the cell membrane stabilizer is added to the sample such that lysis is less than about 5% of the cells. In a most preferred embodiment, the cell membrane stabilizer is added to the sample such that lysis is less than about 1% of the cells.

**[0161]** In another embodiment, a cross-linker is added to the sample such that lysis is less than about 10% of the cells. In a preferred embodiment, the cross-linker is added to the sample such that lysis is less than about 5% of the cells. In a most preferred embodiment, the cross-linker is added to the sample such that lysis is less than about 1% of the cells.

**[0162]** In another embodiment, the cell lysis inhibitor, cross-linker or cell membrane stabilizer is added to the sample in an applicable time period including but not limited to 1-10 seconds, 10-30 seconds, 30-60 seconds, 1-5 minutes, 5-10 minutes, 10-20 minutes, 20-30 minutes, 30-40 minutes, 40-50 minutes, 60-90 minutes, 90-180 minutes or greater than 180 minutes after collection of the sample. In another embodiment, the cell lysis inhibitor, cross-linker, or cell membrane stabilizer is present in the apparatus to which the sample is collected including but not limited to a glass tube, a plastic tube, a circular container, an eppendorf tube, an IV bag, or any other appropriate collection device. In another embodiment, after the addition of the cell lysis inhibitor, cell membrane stabilizer, or cross-linker, the sample is left at about room temperature for the period of time to allow the reagent to function, including but not limited to 1-5, 5-10, 10-20, 20-40, 40-60, 60-90, 90-120, 120-150, 150-180, 180-240, 240-300 or greater than 300 minutes.

**[0163]** In another embodiment, the template DNA contains both maternal DNA and fetal DNA. In a preferred embodiment, template DNA is obtained from blood of a pregnant female. Blood is collected using any standard technique for blood-drawing including but not limited to venipuncture. For example, blood can be drawn from a vein from the inside of the elbow or the back of the hand. Blood samples can be collected from a pregnant female at any time during fetal gestation. For example, blood samples can be collected from human females at 1-4, 4-8, 8-12, 12-16, 16-20, 20-24, 24-28, 28-32, 32-36, 36-40, or 40-44 weeks of fetal gestation, and preferably between 8-28 weeks of fetal gestation.

**[0164]** The blood sample is centrifuged to separate the plasma from the maternal cells. The plasma and maternal cell fractions are transferred to separate tubes and re-centrifuged. The plasma fraction contains cell-free fetal DNA and maternal DNA. Any standard DNA isolation technique can be used to isolate the fetal DNA and the maternal DNA including but not limited to QIAamp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183).

**[0165]** In a preferred embodiment, blood can be collected into an apparatus containing a magnesium chelator including but not limited to EDTA, and is stored at 4°C. Optionally, a calcium chelator, including but not limited to EGTA, can be added.

**[0166]** In another embodiment, a cell lysis inhibitor is added to the maternal blood including but not limited to formaldehyde, formaldehyde derivatives, formalin, glutaraldehyde, glutaraldehyde derivatives, a protein cross-linker, a nucleic acid cross-linker, a protein and nucleic acid cross-linker, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA, DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS, or compounds listed in Table XXIII.

**[0167]** In another embodiment, an agent that stabilizes cell membranes may be added to the maternal blood samples to reduce maternal cell lysis including but not limited to aldehydes, urea formaldehyde, phenol formaldehyde, DMAE (dimethylaminoethanol), cholesterol, cholesterol derivatives, high concentrations of magnesium, vitamin E, and vitamin E derivatives, calcium, calcium gluconate, taurine, niacin, hydroxylamine derivatives, bimocromol, sucrose, astaxanthin, glucose, amitriptyline, isomer A hopane tetral phenylacetate, isomer B hopane tetral phenylacetate, citicoline, inositol, vitamin B, vitamin B complex, cholesterol hemisuccinate, sorbitol, calcium, coenzyme Q, ubiquinone, vitamin K, vitamin K complex, menaquinone, zonegran, zinc, ginkgo biloba extract, diphenylhydantoin, perftoran, polyvinylpyrrolidone, phosphatidylserine, tegretol, PABA, disodium cromoglycate, nedocromil sodium, phenytoin, zinc citrate, mexitil, dilantin, sodium hyaluronate, or polaxamer 188.

**[0168]** In another embodiment, the template DNA is obtained from the plasma or serum of the blood of the pregnant female. The percentage of fetal DNA in maternal plasma is between 0.39-11.9% (*Pertl, and Bianchi, Obstetrics and Gynecology* 98: 483-490 (2001)). The majority of the DNA in the plasma sample is maternal, which makes using the DNA for genotyping the fetus difficult. However, methods that increase the percentage of fetal DNA in the maternal plasma allow the sequence of the fetal DNA to be determined, and allow for the detection of genetic disorders including mutations, insertions, deletions,

and chromosomal abnormalities. The addition of cell lysis inhibitors, cell membrane stabilizers or cross-linkers to the maternal blood sample can increase the relative percentage of fetal DNA. While lysis of both maternal and fetal cells is inhibited, the vast majority of cells are maternal, and thus by reducing the lysis of maternal cells, there is a relative increase in the percentage of free fetal DNA. See Example 4.

**[0169]** In another embodiment, any blood drawing technique, method, protocol, or equipment that reduce the amount of cell lysis can be used, including but not limited to a large bore needle, a shorter length needle, a needle coating that increases laminar flow, e.g., teflon, a modification of the bevel of the needle to increase laminar flow, or techniques that reduce the rate of blood flow. The fetal cells likely are destroyed in the maternal blood by the mother's immune system. However, it is likely that a large portion of the maternal cell lysis occurs as a result of the blood draw or processing of the blood sample. Thus, methods that prevent or reduce cell lysis will reduce the amount of maternal DNA in the sample, and increase the relative percentage of free fetal DNA.

**[0170]** In another embodiment, an agent that preserves or stabilizes the structural integrity of cells can be used to reduce the amount of cell lysis.

**[0171]** In another embodiment, any protocol that reduces the amount of free maternal DNA in the maternal blood can be used prior to obtaining the sample. In another embodiment, prior to obtaining the sample, the pregnant female rests without physical activity for a period of time including but not limited to 0-5, 5-10, 10-15, 15-20, 20-25, 25-30, 30-35, 35-40, 40-45, 45-50, 50-55, 55-60, 60-120, 120-180, 180-240, 240-300, 300-360, 360-420, 420-480, 480-540, 540-600, 600-660, 660-720, 720-780, 780-840, 840-900, 900-1200, 1200-1500, 1500-1800, 1800-2100, 2100-2400, 2400-2700, 2700-3000, 3000-3300, 3300-3600, 3600-3900, 3900-4200, 4200-4500, and greater than 4500 minutes. In another embodiment, the sample is obtained from the pregnant female after her body has reached a relaxed state. The period of rest prior to obtaining the sample may reduce the amount of maternal nucleic acid in the sample. In another embodiment, the sample is obtained from the pregnant female in the a.m., including but not limited to 4-5 am, 5-6 am, 6-7 am, 7-8 am, 8-9 am, 9-10 am, 10-11 am, and 11-12 am.

**[0172]** In another embodiment, the sample is obtained from the pregnant female after she has slept for a period of time including but not limited to 0-1, 1-2, 2-3, 3-4, 4-5, 5-6, 6-7, 7-8, 8-9, 9-10, 10-11, 11-12, or greater than 12 hours.

**[0173]** In another embodiment, prior to obtaining the sample, the pregnant female exercises for a period of time followed by a period of rest. In another embodiment, the period of exercise includes but is not limited to 0-15, 15-30, 30-45, 45-60, 60-120, 120-240, or greater than 240 minutes.

[0174] In another embodiment, agents that prevent the destruction of DNA, including but not limited to a DNase inhibitor, zinc chloride, ethylenediaminetetraacetic acid, guanidine-HCl, guanidine isothiocyanate, N-lauroylsarcosine, and Na-dodecylsulphate, can be added to the blood sample.

[0175] In another embodiment, fetal DNA is obtained from a fetal cell, wherein said fetal cell can be isolated from sources including but not limited to maternal blood, umbilical cord blood, chorionic villi, amniotic fluid, embryonic tissues and mucous obtained from the cervix or vagina of the mother.

[0176] In a preferred embodiment, fetal cells are isolated from maternal peripheral blood. An antibody specific for fetal cells can be used to purify the fetal cells from the maternal serum (Mueller et al., *Lancet* 336: 197-200 (1990); Ganshirt-Ahlert et al., *Am. J. Obstet. Gynecol.* 166: 1350-1355 (1992)). Flow cytometry techniques can also be used to enrich fetal cells (Herzenberg et al., *PNAS* 76: 1453-1455 (1979); Bianchi et al., *PNAS* 87: 3279-3283 (1990); Bruch et al., *Prenatal Diagnosis* 11: 787-798 (1991)). U.S. Pat. No. 5,432,054 also describes a technique for separation of fetal nucleated red blood cells, using a tube having a wide top and a narrow, capillary bottom made of polyethylene. Centrifugation using a variable speed program results in a stacking of red blood cells in the capillary based on the density of the molecules. The density fraction containing low density red blood cells, including fetal red blood cells, is recovered and then differentially hemolyzed to preferentially destroy maternal red blood cells. A density gradient in a hypertonic medium is used to separate red blood cells, now enriched in the fetal red blood cells from lymphocytes and ruptured maternal cells. The use of a hypertonic solution shrinks the red blood cells, which increases their density, and facilitate purification from the more dense lymphocytes. After the fetal cells have been isolated, fetal DNA can be purified using standard techniques in the art.

[0177] The nucleic acid that is to be analyzed can be any nucleic acid, e.g., genomic, plasmid, cosmid, yeast artificial chromosomes, artificial or man-made DNA, including unique DNA sequences, and also DNA that has been reverse transcribed from an RNA sample, such as cDNA. The sequence of RNA can be determined according to the invention if it is capable of being made into a double stranded DNA form to be used as template DNA.

[0178] The terms "primer" and "oligonucleotide primer" are interchangeable when used to discuss an oligonucleotide that anneals to a template and can be used to prime the synthesis of a copy of that template.

[0179] "Amplified" DNA is DNA that has been "copied" once or multiple times, e.g. by polymerase chain reaction. When a large amount of DNA is available to assay, such that a sufficient number of copies of the locus of interest are already present in the sample to be assayed, it may not be necessary to



“amplify” the DNA of the locus of interest into an even larger number of replicate copies. Rather, simply “copying” the template DNA once using a set of appropriate primers, which may contain hairpin structures that allow the restriction enzyme recognition sites to be double stranded, can suffice.

**[0180]** “Copy” as in “copied DNA” refers to DNA that has been copied once, or DNA that has been amplified into more than one copy.

**[0181]** In one embodiment, the nucleic acid is amplified directly in the original sample containing the source of nucleic acid. It is not essential that the nucleic acid be extracted, purified or isolated; it only needs to be provided in a form that is capable of being amplified. Hybridization of the nucleic acid template with primer, prior to amplification, is not required. For example, amplification can be performed in a cell or sample lysate using standard protocols well known in the art. DNA that is on a solid support, in a fixed biological preparation, or otherwise in a composition that contains non-DNA substances and that can be amplified without first being extracted from the solid support or fixed preparation or non-DNA substances in the composition can be used directly, without further purification, as long as the DNA can anneal with appropriate primers, and be copied, especially amplified, and the copied or amplified products can be recovered and utilized as described herein.

**[0182]** In a preferred embodiment, the nucleic acid is extracted, purified or isolated from non-nucleic acid materials that are in the original sample using methods known in the art prior to amplification.

**[0183]** In another embodiment, the nucleic acid is extracted, purified or isolated from the original sample containing the source of nucleic acid and prior to amplification, the nucleic acid is fragmented using any number of methods well known in the art including but not limited to enzymatic digestion, manual shearing, or sonication. For example, the DNA can be digested with one or more restriction enzymes that have a recognition site, and especially an eight base or six base pair recognition site, which is not present in the loci of interest. Typically, DNA can be fragmented to any desired length, including 50, 100, 250, 500, 1,000, 5,000, 10,000, 50,000 and 100,000 base pairs long. In another embodiment, the DNA is fragmented to an average length of about 1000 to 2000 base pairs. However, it is not necessary that the DNA be fragmented.

**[0184]** Fragments of DNA that contain the loci of interest can be purified from the fragmented DNA before amplification. Such fragments can be purified by using primers that will be used in the amplification (see “Primer Design” section below) as hooks to retrieve the loci of interest, based on the ability of such primers to anneal to the loci of interest. In a preferred embodiment, tag-modified primers are used, such as e.g. biotinylated primers.

**[0185]** By purifying the DNA fragments containing the loci of interest, the specificity of the amplification reaction can be improved. This will minimize amplification of nonspecific regions of the template DNA. Purification of the DNA fragments can also allow multiplex PCR (Polymerase Chain Reaction) or amplification of multiple loci of interest with improved specificity.

**[0186]** The loci of interest that are to be sequenced can be selected based upon sequence alone. In humans, over 1.42 million single nucleotide polymorphisms (SNPs) have been described (*Nature* 409:928-933 (2001); The SNP Consortium LTD). On the average, there is one SNP every 1.9 kb of human genome. However, the distance between loci of interest need not be considered when selecting the loci of interest to be sequenced according to the invention. If more than one locus of interest on genomic DNA is being analyzed, the selected loci of interest can be on the same chromosome or on different chromosomes.

**[0187]** In a preferred embodiment, the selected loci of interest can be clustered to a particular region on a chromosome. Multiple loci of interest can be located within a region of DNA such that even with any breakage or fragmentation of the DNA, the multiple loci of interest remain linked. For example, if the DNA is obtained and by natural forces is broken into fragments of 5 Kb, multiple loci of interest can be selected within the 5 Kb regions. This allows each fragment, as measured by the loci of interest within that fragment, to serve as an experimental unit, and will reduce any possible experimental noise of comparing loci of interest on multiple chromosomes.

**[0188]** The loci of interest on a chromosome can be any distance from each other including but not limited to 10-50, 50-100, 100-150, 150-200, 200-250, 250-500, 500-750, 750-1000, 1000-1500, 1500-2000, 2000-2500, 2500-3000, 3000-3500, 3500-4000, 4000-4500, 4500-5000, 5000-10,000 and greater than 10,000 base pairs.

**[0189]** In a preferred embodiment, the length of sequence that is amplified is preferably different for each locus of interest so that the loci of interest can be separated by size.

**[0190]** In fact, it is an advantage of the invention that primers that copy an entire gene sequence need not be utilized. Rather, the copied locus of interest is preferably only a small part of the total gene or a small part of a non-coding region of DNA. There is no advantage to sequencing the entire gene as this can increase cost and delay results. Sequencing only the desired bases or loci of interest maximizes the overall efficiency of the method because it allows for the sequence of the maximum number of loci of interest to be determined in the fastest amount of time and with minimal cost.

**[0191]** Because a large number of sequences can be analyzed together, the method of the invention is especially amenable to the large-scale screening of a number of loci of interest.

[0192] Any number of loci of interest can be analyzed and processed, especially at the same time, using the method of the invention. The sample(s) can be analyzed to determine the sequence at one locus of interest or at multiple loci of interest at the same time. The loci of interest can be present on a single chromosome or on multiple chromosomes.

[0193] Alternatively, 2, 3, 4, 5, 6, 7, 8, 9, 10-20, 20-25, 25-30, 30-35, 35-40, 40-45, 45-50, 50-100, 100-250, 250-500, 500-1,000, 1,000-2,000, 2,000-3,000, 3,000-5,000, 5,000-10,000, 10,000-50,000 or more than 50,000 loci of interest can be analyzed at the same time when a global genetic screening is desired. Such a global genetic screening might be desired when using the method of the invention to provide a genetic fingerprint to identify an individual or for SNP genotyping.

[0194] The locus of interest to be copied can be within a coding sequence or outside of a coding sequence. Preferably, one or more loci of interest that are to be copied are within a gene. In a preferred embodiment, the template DNA that is copied is a locus or loci of interest that is within a genomic coding sequence, either intron or exon. In a highly preferred embodiment, exon DNA sequences are copied. The loci of interest can be sites where mutations are known to cause disease or predispose to a disease state. The loci of interest can be sites of single nucleotide polymorphisms. Alternatively, the loci of interest that are to be copied can be outside of the coding sequence, for example, in a transcriptional regulatory region, and especially a promoter, enhancer, or repressor sequence.

#### **Method for Determining the Sequence of a Locus of Interest**

[0195] Any method that provides information on the sequence of a nucleic acid can be used including but not limited to allele specific PCR, PCR, gel electrophoresis, ELISA, mass spectrometry, MALDI-TOF mass spectrometry hybridization, primer extension, fluorescence detection, fluorescence resonance energy transfer (FRET), fluorescence polarization, DNA sequencing, Sanger dideoxy sequencing, DNA sequencing gels, capillary electrophoresis on an automated DNA sequencing machine, microchannel electrophoresis, microarray, southern blot, slot blot, dot blot, single primer linear nucleic acid amplification, as described in U.S. Patent No. 6,251,639, SNP-IT, GeneChips, HuSNP, BeadArray, TaqMan assay, Invader assay, MassExtend, or MassCleave™ (hMC) method.

[0196] The preferred method of determining the sequence has previously been described in U.S. Application No. 10/093,618, filed on March 11, 2002, hereby incorporated by reference in its entirety.

#### **I. Primer Design**

[0197] Published sequences, including consensus sequences, can be used to design or select primers for use in amplification of template DNA. The selection of sequences to be used for the construction of

primers that flank a locus of interest can be made by examination of the sequence of the loci of interest, or immediately thereto. The recently published sequence of the human genome provides a source of useful consensus sequence information from which to design primers to flank a desired human gene locus of interest.

**[0198]** By “flanking” a locus of interest is meant that the sequences of the primers are such that at least a portion of the 3’ region of one primer is complementary to the antisense strand of the template DNA and upstream from the locus of interest site (forward primer), and at least a portion of the 3’ region of the other primer is complementary to the sense strand of the template DNA and downstream of the locus of interest (reverse primer). A “primer pair” is intended a pair of forward and reverse primers. Both primers of a primer pair anneal in a manner that allows extension of the primers, such that the extension results in amplifying the template DNA in the region of the locus of interest.

**[0199]** Primers can be prepared by a variety of methods including but not limited to cloning of appropriate sequences and direct chemical synthesis using methods well known in the art (Narang *et al.*, *Methods Enzymol.* 68:90 (1979); Brown *et al.*, *Methods Enzymol.* 68:109 (1979)). Primers can also be obtained from commercial sources such as Operon Technologies, Amersham Pharmacia Biotech, Sigma, and Life Technologies. The primers can have an identical melting temperature. The lengths of the primers can be extended or shortened at the 5’ end or the 3’ end to produce primers with desired melting temperatures. In a preferred embodiment, one of the primers of the prime pair is longer than the other primer. In a preferred embodiment, the 3’ annealing lengths of the primers, within a primer pair, differ. Also, the annealing position of each primer pair can be designed such that the sequence and length of the primer pairs yield the desired melting temperature. The simplest equation for determining the melting temperature of primers smaller than 25 base pairs is the Wallace Rule ( $T_d = 2(A+T) + 4(G+C)$ ). Computer programs can also be used to design primers, including but not limited to Array Designer Software (Arrayit Inc.), Oligonucleotide Probe Sequence Design Software for Genetic Analysis (Olympus Optical Co.), NetPrimer, and DNAsis from Hitachi Software Engineering. The  $T_m$  (melting or annealing temperature) of each primer is calculated using software programs such as Net Primer (free web based program at <http://premierbiosoft.com/netprimer/netprlaunch/netprlaunch.html>; internet address as of April 17, 2002).

**[0200]** In another embodiment, the annealing temperature of the primers can be recalculated and increased after any cycle of amplification, including but not limited to cycle 1, 2, 3, 4, 5, cycles 6-10, cycles 10-15, cycles 15-20, cycles 20-25, cycles 25-30, cycles 30-35, or cycles 35-40. After the initial cycles of amplification, the 5’ half of the primers is incorporated into the products from each loci of

interest, thus the  $T_M$  can be recalculated based on both the sequences of the 5' half and the 3' half of each primer.

**[0201]** For example, in FIG. 1B, the first cycle of amplification is performed at about the melting temperature of the 3' region, which anneals to the template DNA, of the second primer (region "c"), which is 13 bases. After the first cycle, the annealing temperature can be raised to  $T_{M2}$ , which is about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer, which is depicted as region "b." The second primer cannot bind to the original template DNA because it only anneals to 13 bases in the original DNA template, and  $T_{M2}$  is about the melting temperature of approximately 20 bases, which is the 3' annealing region of the first primer (FIG. 1C). However, the first primer can bind to the DNA that was copied in the first cycle of the reaction. In the third cycle, the annealing temperature is raised to  $T_{M3}$ , which is about the melting temperature of the entire sequence of the second primer, which is depicted as regions "c" and "d." The DNA template produced from the second cycle of PCR contains both regions c' and d', and therefore, the second primer can anneal and extend at  $T_{M3}$  (FIG. 1D). The remaining cycles are performed at  $T_{M3}$ . The entire sequence of the first primer (a + b') can anneal to the template from the third cycle of PCR, and extend (FIG. 1E). Increasing the annealing temperature will decrease non-specific binding and increase the specificity of the reaction, which is especially useful if amplifying a locus of interest from human genomic DNA, which is about  $3 \times 10^9$  base pairs long.

**[0202]** As used herein, the term "about" with regard to annealing temperatures is used to encompass temperatures within 10 degrees celcius of the stated temperatures.

**[0203]** In one embodiment, one primer pair is used for each locus of interest. However, multiple primer pairs can be used for each locus of interest.

**[0204]** In one embodiment, primers are designed such that one or both primers of the primer pair contain sequence in the 5' region for one or more restriction endonucleases (restriction enzyme).

**[0205]** As used herein, with regard to the position at which restriction enzymes digest DNA, the "sense" strand is the strand reading 5' to 3' in the direction in which the restriction enzyme cuts. For example, BsmF I recognizes the following sequences:

5' GGGAC(N) <sub>10</sub> 3' (SEQ ID NO: 1)	5' (N) <sub>14</sub> GTCCC 3' (SEQ ID NO: 2)
3' CCCTG(N) <sub>14</sub> 5' (SEQ ID NO: 2)	3'(N) <sub>10</sub> CAGGG 5' (SEQ ID NO: 1)

**[0206]** The sense strand is the strand containing the "GGGAC" sequence as it reads 5' to 3' in the direction that the restriction enzyme cuts.

**[0207]** As used herein, with regard to the position at which restriction enzymes digest DNA, the “antisense” strand is the strand reading 3’ to 5’ in the direction in which the restriction enzyme cuts.

**[0208]** In another embodiment, one of the primers in a primer pair is designed such that it contains a restriction enzyme recognition site for a restriction enzyme that cuts “n” nucleotides away from the recognition site, and produces a recessed 3’ end and a 5’ overhang that contains the locus of interest (herein referred to as a “second primer”). “N” is a distance from the recognition site to the site of the cut by the restriction enzyme. In other words, the second primer of a primer pair contains a recognition site for a restriction enzyme that does not cut DNA at the recognition site but cuts “n” nucleotides away from the recognition site. For example, if the recognition sequence is for the restriction enzyme BceA I, the enzyme will cut ten (10) nucleotides from the recognition site on the sense strand, and twelve (12) nucleotides away from the recognition site on the antisense strand.

**[0209]** The 3’ region and preferably, the 3’ half, of the primers is designed to anneal to a sequence that flanks the loci of interest (FIG. 1A). The second primer can anneal any distance from the locus of interest provided that digestion with the restriction enzyme that recognizes the restriction enzyme recognition site on this primer generates a 5’ overhang that contains the locus of interest. The 5’ overhangs can be of any size, including but not limited to 1, 2, 3, 4, 5, 6, 7, 8, and more than 8 bases.

**[0210]** In a preferred embodiment, the 3’ end of the primer that anneals closer to the locus of interest (second primer) can anneal 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, or more than 14 bases from the locus of interest or at the locus of interest.

**[0211]** In a preferred embodiment, the second primer is designed to anneal closer to the locus of interest than the other primer of a primer pair (the other primer is herein referred to as a “first primer”). The second primer can be a forward or reverse primer and the first primer can be a reverse or forward primer, respectively. Whether the first or second primer should be the forward or reverse primer can be determined by which design will provide better sequencing results.

**[0212]** For example, the primer that anneals closer to the locus of interest can contain a recognition site for the restriction enzyme BsmF I, which cuts ten (10) nucleotides from the recognition site on the sense strand, and fourteen (14) nucleotides from the recognition site on the antisense strand. In this case, the primer can be designed so that the restriction enzyme recognition site is 13 bases, 12 bases, 10 bases or 11 bases from the locus of interest. If the recognition site is 13 bases from the locus of interest, digestion with BsmF I will generate a 5’ overhang (RXXX), wherein the locus of interest (R) is the first nucleotide in the overhang (reading 3’ to 5’), and X is any nucleotide. If the recognition site is 12 bases from the locus of interest, digestion with BsmF I will generate a 5’ overhang (XRXX), wherein the locus

of interest (R) is the second nucleotide in the overhang (reading 3' to 5'). If the recognition site is 11 bases from the locus of interest, digestion with BsmF I will generate a 5' overhang (XXRX), wherein the locus of interest (R) is the third nucleotide in the overhang (reading 3' to 5'). The distance between the restriction enzyme recognition site and the locus of interest should be designed so that digestion with the restriction enzyme generates a 5' overhang, which contains the locus of interest. The effective distance between the recognition site and the locus of interest will vary depending on the choice of restriction enzyme.

**[0213]** In another embodiment, the primer that anneals closer to the locus of interest site, relative to the other primer, can be designed so that the restriction enzyme that generates the 5' overhang, which contains the locus of interest, will see the same sequence at the cut site, independent of the nucleotide at the locus of interest site. For example, if the primer that anneals closer to the locus of interest is designed so that the recognition site for the restriction enzyme BsmF I (5' GGGAC 3') is thirteen bases from the locus of interest, the restriction enzyme will cut the antisense strand one base from the locus of interest. The nucleotide at the locus of interest is adjacent to the cut site, and may vary from DNA molecule to DNA molecule. If it is desired that the nucleotides adjacent to the cut site be identical, the primer can be designed so that the restriction enzyme recognition site for BsmF I is twelve bases away from the locus of interest site. Digestion with BsmF I will generate a 5' overhang, wherein the locus of interest site is in the second position of the overhang (reading 3' to 5') and is no longer adjacent to the cut site. Designing the primer so that the restriction enzyme recognition site is twelve (12) bases from the locus of interest site allows the nucleotides adjacent to the cut site to be the same, independent of the nucleotide at the locus of interest. Also, primers that have been designed so that the restriction enzyme recognition site, BsmF I, is eleven (11) or ten (10) bases from the locus of interest site will allow the nucleotides adjacent to the cut site to be the same, independent of the nucleotide at the locus of interest. Similar strategies of primer design can be employed with other restriction enzymes so that the nucleotides adjacent to the cut site will be the same, independent of the nucleotide at the loci of interest.

**[0214]** The 3' end of the first primer (either the forward or the reverse) can be designed to anneal at a chosen distance from the locus of interest. Preferably, for example, this distance is between 1-10, 10-25, 25-50, 50-75, 75-100, 100-150, 150-200, 200-250, 250-300, 300-350, 350-400, 400-450, 450-500, 500-550, 550-600, 600-650, 650-700, 700-750, 750-800, 800-850, 850-900, 900-950, 950-1000 and greater than 1000 bases away from the locus of interest. The annealing sites of the first primers are chosen such that each successive upstream primer is further and further away from its respective downstream primer.

**[0215]** For example, if at locus of interest 1 the 3' ends of the first and second primers are Z bases apart, then at locus of interest 2, the 3' ends of the upstream and downstream primers are Z + K bases apart, where K = 1, 2, 3, 4, 5-10, 10-20, 20-30, 30-40, 40-50, 50-60, 60-70, 70-80, 80-90, 90-100, 100-200, 200-300, 300-400, 400-500, 500-600, 600-700, 700-800, 800-900, 900-1000, or greater than 1000 bases (FIG 2). The purpose of making the first primers further and further apart from their respective second primers is so that the PCR products of all the loci of interest differ in size and can be separated, e.g., on a sequencing gel. This allows for multiplexing by pooling the PCR products in later steps.

**[0216]** In one embodiment, the 5' region of the first or second primer can have a recognition site for any type of restriction enzyme. In a preferred embodiment, the 5' region of the first and/or second primer has at least one restriction enzyme recognition site that is different from the restriction enzyme recognition site that is used to generate the 5' overhang, which contains the locus of interest.

**[0217]** In one embodiment, the 5' region of the first primer can have a recognition site for any type of restriction enzyme. In a preferred embodiment, the first primer has at least one restriction enzyme recognition site that is different from the restriction enzyme recognition site in the second primer. In another preferred embodiment, the first primer anneals further away from the locus of interest than the second primer.

**[0218]** In a preferred embodiment, the second primer contains a restriction enzyme recognition sequence for a Type IIS restriction enzyme including but not limited to BceA I and BsmF I, which produce a two base 5' overhang and a four base 5' overhang, respectively. Restriction enzymes that are Type IIS are preferred because they recognize asymmetric base sequences (not palindromic like the orthodox Type II enzymes). Type IIS restriction enzymes cleave DNA at a specified position that is outside of the recognition site, typically up to 20 base pairs outside of the recognition site. These properties make Type IIS restriction enzymes, and the recognition sites thereof, especially useful in the method of the invention. Preferably, the Type IIS restriction enzymes used in this method leave a 5' overhang and a recessed 3'.

**[0219]** A wide variety of Type IIS restriction enzymes are known and such enzymes have been isolated from bacteria, phage, archaeobacteria and viruses of eukaryotic algae and are commercially available (Promega, Madison WI; New England Biolabs, Beverly, MA; Szybalski W. et al., Gene 100:13-26, 1991). Examples of Type IIS restriction enzymes that would be useful in the method of the invention include, but are not limited to enzymes such as those listed in Table I.



TABLE I

Enzyme-Source	Recognition/Cleavage Site	Supplier
Alw I - <i>Acinetobacter lwoffii</i>	GGATC(4/5)	NE Biolabs
Alw26 I - <i>Acinetobacter lwoffii</i>	GTCTC(1/5)	Promega
Bbs I - <i>Bacillus laterosporus</i>	GAAGAC(2/6)	NE Biolabs
Bbv I - <i>Bacillus brevis</i>	GCAGC(8/12)	NE Biolabs
BceA I - <i>Bacillus cereus</i> 1315	IACGGC(12/14)	NE Biolabs
Bmr I - <i>Bacillus megaterium</i>	CTGGG(5/4)	NE Biolabs
Bsa I - <i>Bacillus stearothermophilus</i> 6-55	GGTCTC(1/5)	NE Biolabs
Bst71 I - <i>Bacillus stearothermophilus</i> 71	GCAGC(8/12)	Promega
BsmA I - <i>Bacillus stearothermophilus</i> A664	GTCTC(1/5)	NE Biolabs
BsmB I - <i>Bacillus stearothermophilus</i> B61	CGTCTC(1/5)	NE Biolabs
BsmF I - <i>Bacillus stearothermophilus</i> F	GGGAC(10/14)	NE Biolabs
BspM I - <i>Bacillus species</i> M	ACCTGC(4/8)	NE Biolabs
Ear I - <i>Enterobacter aerogenes</i>	CTCTTC(1/4)	NE Biolabs
Fau I - <i>Flavobacterium aquatile</i>	CCCGC(4/6)	NE Biolabs
Fok I - <i>Flavobacterium okeonokoites</i>	GGATG(9/13)	NE Biolabs
Hga I - <i>Haemophilus gallinarum</i>	GACGC(5/10)	NE Biolabs
Ple I - <i>Pseudomonas lemoignei</i>	GAGTC(4/5)	NE Biolabs
Sap I - <i>Saccharopolyspora species</i>	GCTCTTC(1/4)	NE Biolabs
SfaN I - <i>Streptococcus faecalis</i> ND547	GCATC(5/9)	NE Biolabs
Sth132 I - <i>Streptococcus thermophilus</i> ST132	CCCG(4/8)	No commercial supplier (Gene 195:201-206 (1997))

[0220] In one embodiment, a primer pair has sequence at the 5' region of each of the primers that provides a restriction enzyme recognition site that is unique for one restriction enzyme.

[0221] In another embodiment, a primer pair has sequence at the 5' region of each of the primers that provide a restriction site that is recognized by more than one restriction enzyme, and especially for more than one Type IIS restriction enzyme. For example, certain consensus sequences can be recognized by more than one enzyme. For example, BsgI, Eco571 and BpmI all recognize the consensus (G/C)TGnAG and cleave 16 by away on the antisense strand and 14 by away on the sense strand. A primer that provides such a consensus sequence would result in a product that has a site that can be recognized by any of the restriction enzymes BsgI, Eco571 and BpmI.

[0222] Other restriction enzymes that cut DNA at a distance from the recognition site, and produce a recessed 3' end and a 5' overhang include Type III restriction enzymes.

[0223] For example, the restriction enzyme EcoP15I recognizes the sequence 5' CAGCAG 3' and cleaves 25 bases downstream on the sense strand and 27 bases on the antisense strand. It will be further appreciated by a person of ordinary skill in the art that new restriction enzymes are continually being discovered and can readily be adopted for use in the subject invention.

**[0224]** In another embodiment, the second primer can contain a portion of the recognition sequence for a restriction enzyme, wherein the full recognition site for the restriction enzyme is generated upon amplification of the template DNA such that digestion with the restriction enzyme generates a 5' overhang containing the locus of interest. For example, the recognition site for BsmF I is 5' GGGACN<sub>10</sub> 3' (SEQ ID NO: 1). The 3' region, which anneals to the template DNA, of the second primer can end with the nucleotides "GGG," which do not have to be complementary with the template DNA. If the 3' annealing region is about 10-20 bases, even if the last three bases do not anneal, the primer will extend and, generate a BsmF I site.

Second primer: 5' GGAAATTCCATGATGCGTGGG→ (SEQ ID NO: 3)

Template DNA 3' CCTTTAAGGTACTACGCAN<sub>1</sub>N<sub>2</sub>N<sub>3</sub>TG 5' (SEQ ID NO: 27)

5' GGAAATTCCATGATGCCTN<sub>1</sub>N<sub>2</sub>N<sub>3</sub>AC 3' (SEQ ID NO: 4)

**[0225]** The second primer can be designed to anneal to the template DNA, wherein the next two bases of the template DNA are thymidine and guanine, such that an adenosine and cytosine are incorporated into the primer forming a recognition site for BsmF I, 5' GGGACN<sub>10</sub> 3' (SEQ ID NO: 1). The second primer can be designed to anneal in such a manner that digestion with BsmF I generates a 5' overhang containing the locus of interest.

**[0226]** In another embodiment, the second primer can contain an entire or full recognition site for a restriction enzyme or a portion of a recognition site, which generates a full recognition site upon primer-dependent replication of the template DNA such that digestion with a restriction enzyme that cuts at the recognition site and generates a 5' overhang that contains the locus of interest. For example, the restriction enzyme BsaJ I binds the following recognition site: 5' C<sup>1</sup>CN<sub>1</sub>N<sub>2</sub>GG 3'. The second primer can be designed such that the 3' region, which anneals to the template DNA of the primer ends with "CC", the SNP of interest is represented by "N<sub>1</sub>", and the template sequence downstream of the SNP is "N<sub>2</sub>GG."

Second primer: 5' GGAAATTCCATGATGCGTACC→ (SEQ ID NO: 5)

Template DNA 3' CCTTTAAGGTACTACGCATGGN<sub>1</sub>N<sub>2</sub>CC 5' (SEQ ID NO: 28)

5' GGAAATTCCATGATGCCTACCN<sub>1</sub>N<sub>2</sub>GG 3' (SEQ ID NO: 6)

**[0227]** After digestion with BsaJ I, a 5' overhang of the following sequence would be generated:

5' C 3'

3' GGN<sub>1</sub>N<sub>2</sub>CC 5'

**[0228]** If the nucleotide guanine is not reported at the locus of interest, the 3' recessed end can be filled in with unlabeled cytosine, which is complementary to the first nucleotide in the overhang. After removing the excess cytosine, labeled ddNTPs can be used to fill in the next nucleotide, N<sub>1</sub>, which

represents the locus of interest. Other restriction enzymes can be used including but not limited to BssK I (5' 'CCNGG 3'), Dde I (5' C<sup>†</sup>TNAG 3'), EcoN I (5' CCTNN<sup>†</sup>NNNAGG 3' (SEQ ID NO: 7)), Fnu4H I (5' GC<sup>†</sup>NGC 3'), Hinf I (5' G<sup>†</sup>ANTC 3') PflF I (5' GACN<sup>†</sup>NNGTC 3'), Sau96 I (5' G<sup>†</sup>GNCC 3'), ScrF I (5' CC<sup>†</sup>NGG 3'), and Tth1 11 I (5' GACN<sup>†</sup>NNGTC 3').

**[0229]** It is not necessary that the 3' region, which anneals to the template DNA, of the second primer be 100% complementary to the template DNA. For example, the last 1, 2, or 3 nucleotides of the 3' end of the second primer can be mismatches with the template DNA. The region of the primer that anneals to the template DNA will target the primer, and allow the primer to extend. Even if the last two nucleotides are not complementary to the template DNA, the primer will extend and generate a restriction enzyme recognition site. For example, the last two nucleotides in the second primer are "CC." The second primer anneals to the template DNA, and allows extension even if "CC" is not complementary to the nucleotides Na, and Nb, on the template DNA.

Second primer:	5' GGAAATTCATGATGCGTACC→ (SEQ ID NO: 5)
Template DNA	3' CCTTTAAGGTACTACGCATN <sub>a</sub> N <sub>b</sub> N <sub>1</sub> N <sub>2</sub> CC 5' (SEQ ID NO: 29)
	5' GGAAATTCATGATGCCTAN <sub>a</sub> N <sub>b</sub> N <sub>1</sub> N <sub>2</sub> GG 3' (SEQ ID NO: 8)

**[0230]** After digestion with BsaJ I, a 5' overhang of the following sequence would be generated:


5' C	3'
3' GGN <sub>1</sub> N <sub>2</sub> CC 5'	

**[0231]** If the nucleotide guanine is not reported at the locus of interest, the 5' overhang can be filled in with unlabeled cytosine. The excess cytosine can be rinsed away, and filled in with labeled ddNTPs. The first nucleotide incorporated (N<sub>1</sub>') corresponds to the locus of interest. If guanine is reported at the locus of interest, the loci of interest can be filled in with unlabeled cytosine and a nucleotide downstream of the locus of interest can be detected. For example, assume N<sub>2</sub> is adenine. If the locus of interest is guanine, unlabeled cytosine can be used in the fill in reaction. After removing the cytosine, a fill in reaction with labeled thymidine can be used. The labeled thymidine will be incorporated only if the locus of interest was a guanine. Thus, the sequence of the locus of interest can be determined by detecting a nucleotide downstream of the locus of interest.


**[0232]** In another embodiment, the first and second primers contain a portion of a recognition sequence for a restriction enzyme, wherein the full recognition site for the restriction enzyme is generated upon amplification of the template DNA such that digestion with the restriction enzyme generates a 5' overhang containing the locus of interest. The recognition site for any restriction enzyme that contains one or more than one variable nucleotide can be generated including but not limited to the restriction

enzymes BssK I (5' <sup>1</sup>CCNGG 3'), Dde I (5' C<sup>1</sup>TNAG 3'), Econ I (5' CCTNN<sup>1</sup>NNNAGG 3' (SEQ ID NO: 7)), Fnu4H I (5' GC<sup>1</sup>NGC 3'), Hinf I (5' G<sup>1</sup>ANTC 3'), PflI I (5' GACN<sup>1</sup>NNGTC 3'), Sau96 I (5' G<sup>1</sup>GNCC 3'), ScrF I (5' CC<sup>1</sup>NGG 3'), and Tth111 I (5' GACN<sup>1</sup>NNGTC 3').

**[0233]** In a preferred embodiment, the 3' regions of the first and second primers contain the partial sequence for a restriction enzyme, wherein the partial sequence contains 1, 2, 3, 4 or more than 4 mismatches with the template DNA; these mismatches create the restriction enzyme recognition site. The number of mismatches that can be tolerated at the 3' end depends on the length of the primer. For example, if the locus of interest is represented by N<sub>1</sub>, a first primer can be designed to be complementary to the template DNA, depicted below as region "a." The 3' region of the first primer ends with "CC," which is not complementary to the template DNA. The second primer is designed to be complementary to the template DNA, which is depicted below as region "b". The 3' region of the second primer ends with "CC," which is not complementary to the template DNA.

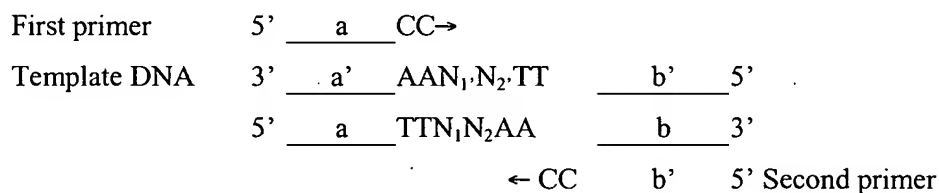


**[0234]** After one round of amplification the following products would be generated:

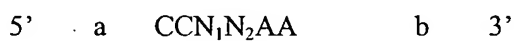

and




**[0235]** In cycle two, the primers can anneal to the templates that were generated from the first cycle of PCR:

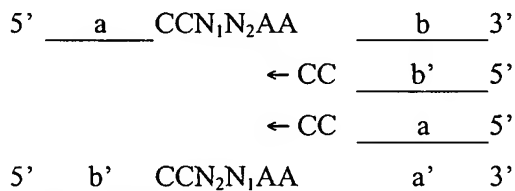

**[0236]** After cycle two of PCR, the following products would be generated:



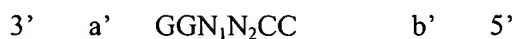

**[0237]**    The restriction enzyme recognition site for BsaJ I is generated, and after digestion with BsaJ I, a 5' overhang containing the locus of interest is created. The locus of interest can be detected as described in detail below.

**[0238]**    In another embodiment, a primer pair has sequence at the 5' region of each of the primers that provides two or more restriction sites that are recognized by two or more restriction enzymes.

**[0239]**    In a most preferred embodiment, a primer pair has different restriction enzyme recognition sites at the 5' regions, especially 5' ends, such that a different restriction enzyme is required to cleave away any undesired sequences. For example, the first primer for locus of interest "A" can contain sequence recognized by a restriction enzyme, "X," which can be any type of restriction enzyme, and the second primer for locus of interest "A," which anneals closer to the locus of interest, can contain sequence for a restriction enzyme, "Y," which is a Type IIS restriction enzyme that cuts "n" nucleotides away and leaves a 5' overhang and a recessed 3' end. The 5' overhang contains the locus of interest. After binding the amplified DNA to streptavidin coated wells, one can digest with enzyme "Y," rinse, then fill in with labeled nucleotides and rinse, and then digest with restriction enzyme "X," which will release the DNA fragment containing the locus of interest from the solid matrix. The locus of interest can be analyzed by detecting the labeled nucleotide that was "filled in" at the locus of interest, e.g. SNP site.

**[0240]**    In another embodiment, the second primers for the different loci of interest that are being amplified according to the invention contain recognition sequence in the 5' regions for the same restriction enzyme and likewise all the first primers also contain the same restriction enzyme recognition site, which is a different enzyme from the enzyme that recognizes the second primers.

**[0241]**    In another embodiment, the second primers for the multiple loci of interest that are being amplified according to the invention contain restriction enzyme recognition sequences in the 5' regions for different restriction enzymes.

**[0242]**    In another embodiment, the first primers for the multiple loci of interest that are being amplified according to the invention contain restriction enzyme recognition sequences in the 5' regions for different restriction enzymes. Multiple restriction enzyme sequences provide an opportunity to influence the order in which pooled loci of interest are released from the solid support. For example, if 50 loci of interest are amplified, the first primers can have a tag at the extreme 5' end to aid in purification and a restriction enzyme recognition site, and the second primers can contain a recognition site for a type

IIS restriction enzyme. For example, several of the first primers can have a restriction enzyme recognition site for EcoR I, other first primers can have a recognition site for Pst I, and still other first primers can have a recognition site for BamH I. After amplification, the loci of interest can be bound to a solid support with the aid of the tag on the first primers. By performing the restriction digests one restriction enzyme at a time, one can serially release the amplified loci of interest. If the first digest is performed with EcoR I, the loci of interest amplified with the first primers containing the recognition site for EcoR I will be released, and collected while the other loci of interest remain bound to the solid support. The amplified loci of interest can be selectively released from the solid support by digesting with one restriction enzyme at a time. The use of different restriction enzyme recognition sites in the first primers allows a larger number of loci of interest to be amplified in a single reaction tube.

**[0243]** In a preferred embodiment, any region 5' of the restriction enzyme digestion site of each primer can be modified with a functional group that provides for fragment manipulation, processing, identification, and/or purification. Examples of such functional groups, or tags, include but are not limited to biotin, derivatives of biotin, carbohydrates, haptens, dyes, radioactive molecules, antibodies, and fragments of antibodies, peptides, and immunogenic molecules.

**[0244]** In another embodiment, the template DNA can be replicated once, without being amplified beyond a single round of replication. This is useful when there is a large amount of the DNA available for analysis such that a large number of copies of the loci of interest are already present in the sample, and further copies are not needed. In this embodiment, the primers are preferably designed to contain a "hairpin" structure in the 5' region, such that the sequence doubles back and anneals to a sequence internal to itself in a complementary manner. When the template DNA is replicated only once, the DNA sequence comprising the recognition site would be single-stranded if not for the "hairpin" structure. However, in the presence of the hairpin structure, that region is effectively double stranded, thus providing a double stranded substrate for activity by restriction enzymes.

**[0245]** To the extent that the reaction conditions are compatible, all the primer pairs to analyze a locus or loci of interest of DNA can be mixed together for use in the method of the invention. In a preferred embodiment, all primer pairs are mixed with the template DNA in a single reaction vessel. Such a reaction vessel can be, for example, a reaction tube, or a well of a microtiter plate.

**[0246]** Alternatively, to avoid competition for nucleotides and to minimize primer dimers and difficulties with annealing temperatures for primers, each locus of interest or small groups of loci of interest can be amplified in separate reaction tubes or wells, and the products later pooled if desired. For example, the separate reactions can be pooled into a single reaction vessel before digestion with the

restriction enzyme that generates a 5' overhang, which contains the locus of interest or SNP site, and a 3' recessed end. Preferably, the primers of each primer pair are provided in equimolar amounts. Also, especially preferably, each of the different primer pairs is provided in equimolar amounts relative to the other pairs that are being used.

**[0247]** In another embodiment, combinations of primer pairs that allow efficient amplification of their respective loci of interest can be used (see e.g. FIG. 2). Such combinations can be determined prior to use in the method of the invention. Multi-well plates and PCR machines can be used to select primer pairs that work efficiently with one another. For example, gradient PCR machines, such as the Eppendorf Mastercycler® gradient PCR machine, can be used to select the optimal annealing temperature for each primer pair. Primer pairs that have similar properties can be used together in a single reaction tube.

**[0248]** In another embodiment, a multi-sample container including but not limited to a 96-well or more plate can be used to amplify a single locus of interest with the same primer pairs from multiple template DNA samples with optimal PCR conditions for that locus of interest. Alternatively, a separate multi-sample container can be used for amplification of each locus of interest and the products for each template DNA sample later pooled. For example, gene A from 96 different DNA samples can be amplified in microtiter plate 1, gene B from 96 different DNA samples can be amplified in microtiter plate 2, etc., and then the amplification products can be pooled.

**[0249]** The result of amplifying multiple loci of interest is a preparation that contains representative PCR products having the sequence of each locus of interest. For example, if DNA from only one individual is used as the template DNA and if hundreds of disease-related loci of interest were amplified from the template DNA, the amplified DNA would be a mixture of small, PCR products from each of the loci of interest. Such a preparation could be further analyzed at that time to determine the sequence at each locus of interest or at only some loci of interest. Additionally, the preparation could be stored in a manner that preserves the DNA and can be analyzed at a later time. Information contained in the amplified DNA can be revealed by any suitable method including but not limited to fluorescence detection, sequencing, gel electrophoresis, and mass spectrometry (see "Detection of Incorporated Nucleotide" section below).

## **II. Amplification of Loci of Interest**

**[0250]** The template DNA can be amplified using any suitable method known in the art including but not limited to PCR (polymerase chain reaction), 3SR (self-sustained sequence reaction), LCR (ligase chain reaction), RACE-PCR (rapid amplification of cDNA ends), PLCR (a combination of polymerase

chain reaction and ligase chain reaction), Q-beta phage amplification (Shah *et al.*, *J. Medical Micro.* 33: 1435-41 (1995)), SDA (strand displacement amplification), SOE-PCR (splice overlap extension PCR), and the like. These methods can be used to design variations of the releasable primer mediated cyclic amplification reaction explicitly described in this application. In the most preferred embodiment, the template DNA is amplified using PCR (PCR: A Practical Approach, M. J. McPherson, et al., IRL Press (1991); PCR Protocols: A Guide to Methods and Applications, Innis, et al., Academic Press (1990); and PCR Technology: Principals and Applications of DNA Amplification, H. A. Erlich, Stockton Press (1989)). PCR is also described in numerous U.S. patents, including U.S. Pat. Nos. 4,683,195; 4,683,202; 4,800,159; 4,965,188; 4,889,818; 5,075,216; 5,079,352; 5,104,792; 5,023,171; 5,091,310; and 5,066,584.

**[0251]** The components of a typical PCR reaction include but are not limited to a template DNA, primers, a reaction buffer (dependent on choice of polymerase), dNTPs (dATP, dTTP, dGTP, and dCTP) and a DNA polymerase. Suitable PCR primers can be designed and prepared as discussed above (see "Primer Design" section above). Briefly, the reaction is heated to 95°C for 2 min. to separate the strands of the template DNA, the reaction is cooled to an appropriate temperature (determined by calculating the annealing temperature of designed primers) to allow primers to anneal to the template DNA, and heated to 72°C for two minutes to allow extension.

**[0252]** In a preferred embodiment, the annealing temperature is increased in each of the first three cycles of amplification to reduce non-specific amplification. See also Example 1, below. The  $TM_1$  of the first cycle of PCR is about the melting temperature of the 3' region of the second primer that anneals to the template DNA. The annealing temperature can be raised in cycles 2-10, preferably in cycle 2, to  $TM_2$ , which is about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer. If the annealing temperature is raised in cycle 2, the annealing temperature remains about the same until the next increase in annealing temperature. Finally, in any cycle subsequent to the cycle in which the annealing temperature was increased to  $TM_2$ , preferably cycle 3, the annealing temperature is raised to  $TM_3$ , which is about the melting temperature of the entire second primer. After the third cycle, the annealing temperature for the remaining cycles can be at about  $TM_3$  or can be further increased. In this example, the annealing temperature is increased in cycles 2 and 3. However, the annealing temperature can be increased from a low annealing temperature in cycle 1 to a high annealing temperature in cycle 2 without any further increases in temperature or the annealing temperature can progressively change from a low annealing temperature to a high annealing temperature in any number of incremental steps. For example, the annealing temperature can be changed in cycles 2, 3, 4, 5, 6, etc.



**[0253]** After annealing, the temperature in each cycle is increased to an “extension” temperature to allow the primers to “extend” and then following extension the temperature in each cycle is increased to the denaturization temperature. For PCR products less than 500 base pairs in size, one can eliminate the extension step in each cycle and just have denaturization and annealing steps. A typical PCR reaction consists of 25-45 cycles of denaturation, annealing and extension as described above. However, as previously noted, one cycle of amplification (one copy) can be sufficient for practicing the invention.

**[0254]** In another embodiment, multiple sets of primers wherein a primer set comprises a forward primer and a reverser primer, can be used to amplify the template DNA for 1-5, 5-10, 10-15, 15-20 or more than 20 cycles, and then the amplified product is further amplified in a reaction with a single primer set or a subset of the multiple primer sets. In a preferred embodiment, a low concentration of each primer set is used to minimize primer-dimer formation. A low concentration of starting DNA can be amplified using multiple primer sets. Any number of primer sets can be used in the first amplification reaction including but not limiting to 1-10, 10-20, 20-30, 30-40, 40-50, 50-60, 60-70, 70-80, 80-90, 90-100, 100-150, 150-200, 200-250, 250-300, 300-350, 350-400, 400-450, 450-500, 500-1000, and greater than 1000. In another embodiment, the amplified product is amplified in a second reaction with a single primer set. In another embodiment, the amplified product is further amplified with a subset of the multiple primer pairs including but not limited to 2-10, 10-20, 20-30, 30-40, 40-50, 50-60, 60-70, 70-80, 80-90, 90-100, 100-150, 150-200, 200-250, and more than 250.

**[0255]** The multiple primer sets will amplify the loci of interest, such that a minimal amount of template DNA is not limiting for the number of loci that can be detected. For example, if template DNA is isolated from a single cell or the template DNA is obtained from a pregnant female, which comprises both maternal template DNA and fetal template DNA, low concentrations of each primer set can be used in a first amplification reaction to amplify the loci of interest. The low concentration of primers reduces the formation of primer-dimer and increases the probability that the primers will anneal to the template DNA and allow the polymerase to extend. The optimal number of cycles performed with the multiple primer sets is determined by the concentration of the primers. Following the first amplification reaction, additional primers can be added to further amplify the loci of interest. Additional amounts of each primer set can be added and further amplified in a single reaction. Alternatively, the amplified product can be further amplified using a single primer set in each reaction or a subset of the multiple primers sets. For example, if 150 primer sets were used in the first amplification reaction, subsets of 10 primer sets can be used to further amplify the product from the first reaction.

[0256] Any DNA polymerase that catalyzes primer extension can be used including but not limited to E. coli DNA polymerase, Klenow fragment of E. coli DNA polymerase 1, T7 DNA polymerase, T4 DNA polymerase, Taq polymerase, Pfu DNA polymerase, Vent DNA polymerase, bacteriophage 29, REDTaq™ Genomic DNA polymerase, or sequenase. Preferably, a thermostable DNA polymerase is used. A “hot start” PCR can also be performed wherein the reaction is heated to 95°C for two minutes prior to addition of the polymerase or the polymerase can be kept inactive until the first heating step in cycle 1. “Hot start” PCR can be used to minimize nonspecific amplification. Any number of PCR cycles can be used to amplify the DNA, including but not limited to 2, 5, 10, 15, 20, 25, 30, 35, 40, or 45 cycles. In a most preferred embodiment, the number of PCR cycles performed is such that equimolar amounts of each loci of interest are produced.

### III. Purification of Amplified DNA

[0257] Purification of the amplified DNA is not necessary for practicing the invention. However, in one embodiment, if purification is preferred, the 5' end of the primer (first or second primer) can be modified with a tag that facilitates purification of the PCR products. In a preferred embodiment, the first primer is modified with a tag that facilitates purification of the PCR products. The modification is preferably the same for all primers, although different modifications can be used if it is desired to separate the PCR products into different groups.

[0258] The tag can be any chemical moiety including but not limited to a radioisotope, fluorescent reporter molecule, chemiluminescent reporter molecule, antibody, antibody fragment, hapten, biotin, derivative of biotin, photobiotin, iminobiotin, digoxigenin, avidin, enzyme, acridinium, sugar, enzyme, apoenzyme, homopolymeric oligonucleotide, hormone, ferromagnetic moiety, paramagnetic moiety, diamagnetic moiety, phosphorescent moiety, luminescent moiety, electrochemiluminescent moiety, chromatic moiety, moiety having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity, or combinations thereof.

[0259] As one example, the 5' ends of the primers can be biotinylated (Kandpal *et al.*, *Nucleic Acids Res.* 18:1789-1795 (1990); Kaneoka *et al.*, *Biotechniques* 10:30-34 (1991); Green *et al.*, *Nucleic Acids Res.* 18:6163-6164 (1990)). The biotin provides an affinity tag that can be used to purify the copied DNA from the genomic DNA or any other DNA molecules that are not of interest. Biotinylated molecules can be purified using a streptavidin coated matrix as shown in FIG. 1F, including but not limited to Streptawell, transparent, High-Bind plates from Roche Molecular Biochemicals (catalog number 1 645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog).

**[0260]** The PCR product of each locus of interest is placed into separate wells of a Streptavidin coated plate. Alternatively, the PCR products of the loci of interest can be pooled and placed into a streptavidin coated matrix, including but not limited to the Streptawell, transparent, High-Bind plates from Roche Molecular Biochemicals (catalog number 1 645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog).

**[0261]** The amplified DNA can also be separated from the template DNA using non-affinity methods known in the art, for example, by polyacrylamide gel electrophoresis using standard protocols.

#### **IV. Digestion of Amplified DNA**

**[0262]** The amplified DNA can be digested with a restriction enzyme that recognizes a sequence that had been provided on the first or second primer using standard protocols known within the art (FIGS. 6A-6D). Restriction enzyme digestions are performed using standard protocols well known within the art. The enzyme used depends on the restriction recognition site generated with the first or second primer. See "Primer Design" section, above, for details on restriction recognition sites generated on primers.

**[0263]** Type IIS restriction enzymes are extremely useful in that they cut approximately 10-20 base pairs outside of the recognition site. Preferably, the Type IIS restriction enzymes used are those that generate a 5' overhang and a recessed 3' end, including but not limited to BceA I and BsmF I (see e.g. Table 1). In a most preferred embodiment, the second primer (either forward or reverse) contains a restriction enzyme recognition sequence for BsmF I or BceA I. The Type IIS restriction enzyme BsmF I recognizes the nucleic acid sequence GGGAC, and cuts 14 nucleotides from the recognition site on the antisense strand and 10 nucleotides from the recognition site on the sense strand. Digestion with BsmF I generates a 5' overhang of four (4) bases.

**[0264]** For example, if the second primer is designed so that after amplification the restriction enzyme recognition site is 13 bases from the locus of interest, then after digestion, the locus of interest is the first base in the 5' overhang (reading 3' to 5'), and the recessed 3' end is one base from the locus of interest. The 3' recessed end can be filled in with a nucleotide that is complementary to the locus of interest. One base of the overhang can be filled in using dideoxynucleotides. However, 1, 2, 3, or 4 bases of the overhang can be filled in using deoxynucleotides or a mixture of dideoxynucleotides and deoxynucleotides.

**[0265]** The restriction enzyme BsmF I cuts DNA ten (10) nucleotides from the recognition site on the sense strand and fourteen (14) nucleotides from the recognition site on the antisense strand. However,

in a sequence dependent manner, the restriction enzyme BsmF I also cuts eleven (11) nucleotides from the recognition site on the sense strand and fifteen (15) nucleotides from the recognition site on the antisense strand. Thus, two populations of DNA molecules exist after digestion: DNA molecules cut at 10/14 and DNA molecules cut at 11/15. If the recognition site for BsmF I is 13 bases from the locus of interest in the amplified product, then DNA molecules cut at the 11/15 position will generate a 5' overhang that contains the locus of interest in the second position of the overhang (reading 3' to 5'). The 3' recessed end of the DNA molecules can be filled in with labeled nucleotides. For example, if labeled dideoxynucleotides are used, the 3' recessed end of the molecules cut at 11/15 would be filled in with one base, which corresponds to the base upstream from the locus of interest, and the 3' recessed end of molecules cut at 10/14 would be filled in with one base, which corresponds to the locus of interest. The DNA molecules that have been cut at the 10/14 position and the DNA molecules that have been cut at the 11/15 position can be separated by size, and the incorporated nucleotides detected. This allows detection of both the nucleotide before the locus of interest, detection of the locus of interest, and potentially the three bases after the locus of interest.

**[0266]** Alternatively, if the base upstream from the locus of interest and the locus of interest are different nucleotides, then the 3' recessed end of the molecules cut at 11/15 can be filled in with deoxynucleotide that is complementary to the upstream base. The remaining deoxynucleotide is washed away, and the locus of interest site can be filled in with either labeled deoxynucleotides, unlabeled deoxynucleotides, labeled dideoxynucleotides, or unlabeled dideoxynucleotides. After the fill in reaction, the nucleotide can be detected by any suitable method. Thus, after the first fill in reaction with dNTP, the 3' recessed end of the molecules cut at 10/14 and 11/15 is upstream from the locus of interest. The 3' recessed end can now be filled in one base, which corresponds to the locus of interest, two bases, three bases or four bases.

**[0267]** The restriction enzyme BceA I recognizes the nucleic acid sequence ACGGC and cuts 12 (twelve) nucleotides from the recognition site on the sense strand and 14 (fourteen) nucleotides from the recognition site on the antisense strand. If the distance from the recognition site for BceA I on the second primer is designed to be thirteen (13) bases from the locus of interest (see FIGS. 4A-4D), digestion with BceA I will generate a 5' overhang of two bases, which contains the locus of interest, and a recessed 3' end that is upstream from the locus of interest. The locus of interest is the first nucleotide in the 5' overhang (reading 3' to 5').

**[0268]** Alternative cutting is also seen with the restriction enzyme BceA I, although at a much lower frequency than is seen with BsmF I. The restriction enzyme BceA I can cut thirteen (13) nucleotides

from the recognition site on the sense strand and fifteen (15) nucleotides from the recognition site on the antisense strand. Thus, two populations of DNA molecules exist: DNA molecules cut at 12/14 and DNA molecules cut at 13/15. If the restriction enzyme recognition site is 13 bases from the locus of interest in the amplified product, DNA molecules cut at the 13/15 position yield a 5' overhang, which contains the locus of interest in the second position of the overhang (reading 3' to 5'). Labeled dideoxynucleotides can be used to fill in the 3' recessed end of the DNA molecules. The DNA molecules cut at 13/15 will have the base upstream from the locus of interest filled in, and the DNA molecules cut at 12/14 will have the locus of interest site filled in. The DNA molecules cut at 13/15 and those cut at 12/14 can be separated by size, and the incorporated nucleotide detected. Thus, the alternative cutting can be used to obtain additional sequence information.

**[0269]** Alternatively, if the two bases in the 5' overhang are different, the 3' recessed end of the DNA molecules, which were cut at 13/15, can be filled in with the deoxynucleotide complementary to the first base in the overhang, and excess deoxynucleotide washed away. After filling in, the 3' recessed end of the DNA molecules that were cut at 12/14 and the DNA molecules that were cut at 13/15 are upstream from the locus of interest. The 3' recessed ends can be filled with either labeled dideoxynucleotides, unlabeled dideoxynucleotides, labeled deoxynucleotides, or unlabeled deoxynucleotides.

**[0270]** If the primers provide different restriction sites for certain of the loci of interest that were copied, all the necessary restriction enzymes can be added together to digest the copied DNA simultaneously. Alternatively, the different restriction digests can be made in sequence, for example, using one restriction enzyme at a time, so that only the product that is specific for that restriction enzyme is digested.

**[0271]** Optimal restriction enzyme digestion conditions, including but not limited to the concentration of enzyme, temperature, buffer conditions, and the time of digestion can be optimized for each restriction enzyme. For example, the alternative cutting seen with the type IIS restriction enzyme BsmF I can be reduced, if desired, by performing the restriction enzyme digestion at lower temperatures including but not limited to 25-16°, 16-12°C, 12-8°C, 8-4°C, or 4-0°C.

## **V. Incorporation of Labeled Nucleotides**

**[0272]** Digestion with the restriction enzyme that recognizes the sequence on the second primer generates a recessed 3' end and a 5' overhang, which contains the locus of interest (FIG. 1G). The recessed 3' end can be filled in using the 5' overhang as a template in the presence of unlabeled or labeled nucleotides or a combination of both unlabeled and labeled nucleotides. The nucleotides can be labeled

with any type of chemical group or moiety that allows for detection including but not limited to radioactive molecules, fluorescent molecules, antibodies, antibody fragments, haptens, carbohydrates, biotin, derivatives of biotin, phosphorescent moieties, luminescent moieties, electrochemiluminescent moieties, chromatic moieties, and moieties having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity. The nucleotides can be labeled with one or more than one type of chemical group or moiety. Each nucleotide can be labeled with the same chemical group or moiety. Alternatively, each different nucleotide can be labeled with a different chemical group or moiety. The labeled nucleotides can be dNTPs, ddNTPs, or a mixture of both dNTPs and ddNTPs. The unlabeled nucleotides can be dNTPs, ddNTPs or a mixture of both dNTPs and ddNTPs.

**[0273]** Any combination of nucleotides can be used to incorporate nucleotides including but not limited to unlabeled deoxynucleotides, labeled deoxynucleotides, unlabeled dideoxynucleotides, labeled dideoxynucleotides, a mixture of labeled and unlabeled deoxynucleotides, a mixture of labeled and unlabeled dideoxynucleotides, a mixture of labeled deoxynucleotides and labeled dideoxynucleotides, a mixture of labeled deoxynucleotides and unlabeled dideoxynucleotides, a mixture of unlabeled deoxynucleotides and unlabeled dideoxynucleotides, a mixture of unlabeled deoxynucleotides and labeled dideoxynucleotides, dideoxynucleotide analogues, deoxynucleotide analogues, a mixture of dideoxynucleotide analogues and deoxynucleotide analogues, phosphorylated nucleoside analogues, 2'-deoxynucleotide-5'-triphosphate, and modified 2'-deoxynucleotide-5'-triphosphate.

**[0274]** For example, as shown in FIG. 1H, in the presence of a polymerase, the 3' recessed end can be filled in with fluorescent ddNTP using the 5' overhang as a template. The incorporated ddNTP can be detected using any suitable method including but not limited to fluorescence detection.

**[0275]** All four nucleotides can be labeled with different fluorescent groups, which will allow one reaction to be performed in the presence of all four labeled nucleotides. Alternatively, four separate "fill in" reactions can be performed for each locus of interest; each of the four reactions will contain a different labeled nucleotide (e.g. ddATP\*, ddTTP\*, ddGTP\*, or ddCTP\*, where \* indicates a labeled nucleotide). Each nucleotide can be labeled with different chemical groups or the same chemical groups. The labeled nucleotides can be dideoxynucleotides or deoxynucleotides.

**[0276]** In another embodiment, nucleotides can be labeled with fluorescent dyes including but not limited to fluorescein, pyrene, 7-methoxycoumarin, Cascade Blue.TM., Alexa Fluor 350, Alexa Fluor 430, Alexa Fluor 488, Alexa Fluor 532, Alexa Fluor 546, Alexa Fluor 568, Alexa Fluor 594, Alexa Fluor 633, Alexa Fluor 647, Alexa Fluor 660, Alexa Fluor 680, AMCA-X, dialkylaminocoumarin, Pacific Blue, Marina Blue, BODIPY 493/503, BODIPY FL-X, DTAF, Oregon Green 500, Dansyl-X, 6-FAM, Oregon Green 488,

Oregon Green 514, Rhodamine Green-X, Rhodol Green, Calcein, Eosin, ethidium bromide, NBD, TET, 2', 4', 5', 7' tetrabromosulfonefluorescein, BODIPY-R6G, BODIPY-FI BR2, BODIPY 530/550, HEX, BODIPY 558/568, BODIPY-TMR-X., PyMPO, BODIPY 564/570, TAMRA, BODIPY 576/589, Cy3, Rhodamine Red-x, BODIPY 581/591, carboxyXrhodamine, Texas Red-X, BODIPY-TR-X., Cy5, SpectrumAqua, SpectrumGreen #1, SpectrumGreen #2, SpectrumOrange, SpectrumRed, or naphthofluorescein.

**[0277]** In another embodiment, the "fill in" reaction can be performed with fluorescently labeled dNTPs, wherein the nucleotides are labeled with different fluorescent groups. The incorporated nucleotides can be detected by any suitable method including but not limited to Fluorescence Resonance Energy Transfer (FRET).

**[0278]** In another embodiment, a mixture of both labeled ddNTPs and unlabeled dNTPs can be used for filling in the recessed 3' end of the SNP or locus of interest. Preferably, the 5' overhang consists of more than one base, including but not limited to 2, 3, 4, 5, 6 or more than 6 bases. For example, if the 5' overhang consists of the sequence "XGAA," wherein X is the locus of interest, e.g. SNP, then filling in with a mixture of labeled ddNTPs and unlabeled dNTPs will produce several different DNA fragments. If a labeled ddNTP is incorporated at position "X," the reaction will terminate and a single labeled base will be incorporated. If however, an unlabeled dNTP is incorporated, the polymerase continues to incorporate other bases until a labeled ddNTP is incorporated. If the first two nucleotides incorporated are dNTPs, and the third is a ddNTP, the 3' recessed end will be extended by three bases. This DNA fragment can be separated from the other DNA fragments that were extended by 1, 2, or 4 bases by size. A mixture of labeled ddNTPs and unlabeled dNTPs will allow all bases of the overhang to be filled in, and provides additional sequence information about the locus of interest, e.g. SNP (see FIGS. 7E and 9D).

**[0279]** After incorporation of the labeled nucleotide, the amplified DNA can be digested with a restriction enzyme that recognizes the sequence provided by the first primer. For example, in FIG 1I, the amplified DNA is digested with a restriction enzyme that binds to region "a," which releases the DNA fragment containing the incorporated nucleotide from the streptavidin matrix.

**[0280]** Alternatively, one primer of each primer pair for each locus of interest can be attached to a solid support matrix including but not limited to a well of a microtiter plate. For example, streptavidin-coated microtiter plates can be used for the amplification reaction with a primer pair, wherein one primer is biotinylated. First, biotinylated primers are bound to the streptavidin-coated microtiter plates. Then, the plates are used as the reaction vessel for PCR amplification of the loci of interest. After

the amplification reaction is complete, the excess primers, salts, and template DNA can be removed by washing. The amplified DNA remains attached to the microtiter plate. The amplified DNA can be digested with a restriction enzyme that recognizes a sequence on the second primer and generates a 5' overhang, which contains the locus of interest. The digested fragments can be removed by washing. After digestion, the SNP site or locus of interest is exposed in the 5' overhang. The recessed 3' end is filled in with a labeled nucleotide, including but not limited to, fluorescent ddNTP in the presence of a polymerase. The labeled DNA can be released into the supernatant in the microtiter plate by digesting with a restriction enzyme that recognizes a sequence in the 5' region of the first primer.

**[0281]** In another embodiment, one nucleotide can be used to determine the sequence of multiple alleles of a gene. A nucleotide that terminates the elongation reaction can be used to determine the sequence of multiple alleles of a gene. At one allele, the terminating nucleotide is complementary to the locus of interest in the 5' overhang of said allele. The nucleotide is incorporated and terminates the reaction. At a different allele, the terminating nucleotide is not complementary to the locus of interest, which allows a non-terminating nucleotide to be incorporated at the locus of interest of the different allele. However, the terminating nucleotide is complementary to a nucleotide downstream from the locus of interest in the 5' overhang of said different allele. The sequence of the alleles can be determined by analyzing the patterns of incorporation of the terminating nucleotide. The terminating nucleotide can be labeled or unlabeled.

**[0282]** In a another embodiment, the terminating nucleotide is a nucleotide that terminates or hinders the elongation reaction including but not limited to a dideoxynucleotide, a dideoxynucleotide derivative, a dideoxynucleotide analog, a dideoxynucleotide homolog, a dideoxynucleotide with a sulfur chemical group, a deoxynucleotide, a deoxynucleotide derivative, a deoxynucleotide homolog, a deoxynucleotide analog, a deoxynucleotide with a sulfur chemical group, arabinoside triphosphate, a ~~arabinoside~~an arabinoside triphosphate analog, a ~~arabinoside~~an arabinoside triphosphate homolog, or an arabinoside derivative.

**[0283]** In another embodiment, a terminating nucleotide labeled with one signal generating moiety tag, including but not limited to a fluorescent dye, can be used to determine the sequence of the alleles of a locus of interest. The use of a single nucleotide labeled with one signal generating moiety tag eliminates any difficulties that can arise when using different fluorescent moieties. In addition, using one nucleotide labeled with one signal generating moiety tag to determine the sequence of alleles of a locus of interest reduces the number of reactions, and eliminates pipetting errors.



**[0284]** For example, if the second primer contains the restriction enzyme recognition site for BsmFI, digestion will generate a 5' overhang of 4 bases. The second primer can be designed such that the locus of interest is located in the first position of the overhang. A representative overhang is depicted below, where R represents the locus of interest:

	5' CAC				
	3' GTG	R	T	G	G
Overhang position		1	2	3	4

**[0285]** One nucleotide with one signal generating moiety tag can be used to determine whether the variable site is homozygous or heterozygous. For example, if the variable site is adenine (A) or guanine (G), then either adenine or guanine can be used to determine the sequence of the alleles of the locus of interest, provided that there is an adenine or guanine in the overhang at position 2, 3, or 4.

**[0286]** For example, if the nucleotide in position 2 of the overhang is thymidine, which is complementary to adenine, then labeled ddATP, unlabeled dCTP, dGTP, and dTTP can be used to determine the sequence of the alleles of the locus of interest. The ddATP can be labeled with any signal generating moiety including but not limited to a fluorescent dye. If the template DNA is homozygous for adenine, then labeled ddATP\* will be incorporated at position 1 complementary to the overhang at the alleles, and no nucleotide incorporation will be seen at position 2, 3 or 4 complementary to the overhang.

Allele 1	5' CCC	A*			
	3' GGG	T	T	G	G
Overhang position		1	2	3	4
Allele 2	5' CCC	A*			
	3' GGG	T	T	G	G
Overhang position		1	2	3	4

**[0287]** One signal will be seen corresponding to incorporation of labeled ddATP at position 1 complementary to the overhang, which indicates that the individual is homozygous for adenine at this position. This method of labeling eliminates any difficulties that may arise from using different dyes that have different quantum coefficients.

**[0288]** Homozygous guanine:

If the template DNA is homozygous for guanine, then no ddATP will be incorporated at position 1 complementary to the overhang, but ddATP will be incorporated at the first available position, which in this case is position 2 complementary to the overhang. For example, if the second position in the overhang corresponds to a thymidine, then:

Allele 1	5' CCC	<b>G</b>	<b>A*</b>		
	3' GGG	C	T	G	G
Overhang position		1	2	3	4
Allele 2	5' CCC	<b>G</b>	<b>A*</b>		
	3' GGG	C	T	G	G
Overhang position		1	2	3	4

**[0289]** One signal will be seen corresponding to incorporation of ddATP at position 2 complementary to the overhang, which indicates that the individual is homozygous for guanine. The molecules that are filled in at position 2 complementary to the overhang will have a different molecular weight than the molecules filled in at position 1 complementary to the overhang.

**[0290]** Heterozygous condition:

Allele 1	5' CCC	<b>A*</b>			
	3' GGG	T	T	G	G
Overhang position		1	2	3	4
Allele 2	5' CCC	<b>G</b>	<b>A*</b>		
	3' GGG	C	T	G	G
Overhang position		1	2	3	4

**[0291]** Two signals will be seen; the first signal corresponds to the ddATP filled in at position one complementary to the overhang and the second signal corresponds to the ddATP filled in at position 2 complementary to the overhang. The two signals can be separated based on molecular weight; allele 1 and allele 2 will be separated by a single base pair, which allows easy detection and quantitation of the signals. Molecules filled in at position one can be distinguished from molecules filled in at position two using any method that discriminates based on molecular weight including but not limited to gel electrophoresis, capillary gel electrophoresis, DNA sequencing, and mass spectrometry. It is not necessary that the nucleotide be labeled with a chemical moiety; the DNA molecules corresponding to the different alleles can be separated based on molecular weight.

**[0292]** If position 2 of the overhang is not complementary to adenine, it is possible that positions 3 or 4 may be complementary to adenine. For example, position 3 of the overhang may be complementary to the nucleotide adenine, in which case labeled ddATP may be used to determine the sequence of both alleles.

**[0293]** Homozygous for adenine:

Allele 1	5' CCC	<b>A*</b>			
----------	--------	-----------	--	--	--

	3' GGG	T	G	T	G
Overhang position		1	2	3	4
Allele 2	5' CCC	A*			
	3' GGG	T	G	T	G
Overhang position		1	2	3	4

**[0294]**     Homozygous for guanine:

Allele 1	5' CCC	G	C	A*	
	3' GGG	C	G	T	G
Overhang position		1	2	3	4
Allele 2	5' CCC	G	C	A*	
	3' GGG	C	G	T	G
Overhang position		1	2	3	4

**[0295]**     Heterozygous:

Allele 1	5' CCC	A*			
	3' GGG	T	G	T	G
Overhang position		1	2	3	4
Allele 2	5' CCC	G	C	A*	
	3' GGG	C	G	T	G
Overhang position		1	2	3	4

**[0296]**     Two signals will be seen; the first signal corresponds to the ddATP filled in at position 1 complementary to the overhang and the second signal corresponds to the ddATP filled in at position 3 complementary to the overhang. The two signals can be separated based on molecular weight; allele 1 and allele 2 will be separated by two bases, which can be detected using any method that discriminates based on molecular weight.

**[0297]**     Alternatively, if positions 2 and 3 are not complementary to adenine (*i.e* positions 2 and 3 of the overhang correspond to guanine, cytosine, or adenine) but position 4 is complementary to adenine, labeled ddATP can be used to determine the sequence of both alleles.

**[0298]**     Homozygous for adenine:

Allele 1	5' CCC	A*			
	3' GGG	T	G	G	T
Overhang position		1	2	3	4
Allele 2	5' CCC	A*			

	3' GGG	T	G	G	T
Overhang position		1	2	3	4

**[0299]** One signal will be seen that corresponds to the molecular weight of molecules filled in with ddATP at position one complementary to the overhang, which indicates that the individual is homozygous for adenine at the variable site.

**[0300]** Homozygous for guanine:

Allele 1	5' CCC	G	C	C	A*
	3' GGG	C	G	G	T
Overhang position		1	2	3	4
Allele 2	5' CCC	G	C	C	A*
	3' GGG	C	G	G	T
Overhang position		1	2	3	4

**[0301]** One signal will be seen that corresponds to the molecular weight of molecules filled in at position 4 complementary to the overhang, which indicates that the individual is homozygous for guanine.

**[0302]** Heterozygous:

Allele 1	5' CCC	A*			
	3' GGG	T	G	G	T
Overhang position		1	2	3	4
Allele 2	5' CCC	G	C	C	A*
	3' GGG	C	G	G	T
Overhang position		1	2	3	4

**[0303]** Two signals will be seen; the first signal corresponds to the ddATP filled in at position one complementary to the overhang and the second signal corresponds to the ddATP filled in at position 4 complementary to the overhang. The two signals can be separated based on molecular weight; allele 1 and allele 2 will be separated by three bases, which allows detection and quantitation of the signals. The molecules filled in at position 1 and those filled in at position 4 can be distinguished based on molecular weight.

**[0304]** As discussed above, if the variable site contains either adenine or guanine, either labeled adenine or labeled guanine can be used to determine the sequence of both alleles. If positions 2, 3, or 4 of the overhang are not complementary to adenine but one of the positions is complementary to a guanine, then labeled ddGTP can be used to determine whether the template DNA is homozygous or heterozygous

for adenine or guanine. For example, if position 3 in the overhang corresponds to a cytosine then the following signals will be expected if the template DNA is homozygous for guanine, homozygous for adenine, or heterozygous:

**[0305]**      Homozygous for guanine:

Allele 1	5' CCC	<b>G*</b>				
	3' GGG	C	T	C	T	
Overhang position		1	2	3	4	
Allele 2	5' CCC	<b>G*</b>				
	3' GGG	C	T	C	T	
Overhang position		1	2	3	4	

**[0306]**      One signal will be seen that corresponds to the molecular weight of molecules filled in with ddGTP at position one complementary to the overhang, which indicates that the individual is homozygous for guanine.

**[0307]**      Homozygous for adenine:

Allele 1	5' CCC	<b>A</b>	<b>A</b>	<b>G*</b>		
	3' GGG	T	T	C	T	
Overhang position		1	2	3	4	
Allele 2	5' CCC	<b>A</b>	<b>A</b>	<b>G*</b>		
	3' GGG	T	T	C	T	
Overhang position		1	2	3	4	

**[0308]**      One signal will be seen that corresponds to the molecular weight of molecules filled in at position 3 complementary to the overhang, which indicates that the individual is homozygous for adenine at the variable site.

**[0309]**      Heterozygous:

Allele 1	5' CCC	<b>G*</b>				
	3' GGG	C	T	C	T	
Overhang position		1	2	3	4	
Allele 2	5' CCC	<b>A</b>	<b>A</b>	<b>G*</b>		
	3' GGG	T	T	C	T	
Overhang position		1	2	3	4	

**[0310]**      Two signals will be seen; the first signal corresponds to the ddGTP filled in at position one complementary to the overhang and the second signal corresponds to the ddGTP filled in at position 3

complementary to the overhang. The two signals can be separated based on molecular weight; allele 1 and allele 2 will be separated by two bases, which allows easy detection and quantitation of the signals.

**[0311]** In another embodiment, the nucleotide labeled with a single chemical moiety, which is used to determine the sequence of alleles of interest, can be analyzed by a variety of methods including but not limited to fluorescence detection, DNA sequencing gel, capillary electrophoresis on an automated DNA sequencing machine, microchannel electrophoresis, and other methods of sequencing, mass spectrometry, time of flight mass spectrometry, quadrupole mass spectrometry, magnetic sector mass spectrometry, electric sector mass spectrometry infrared spectrometry, ultraviolet spectrometry, potentiostatic amperometry or by DNA hybridization techniques including Southern Blots, Slot Blots, Dot Blots, and DNA microarrays, wherein DNA fragments would be useful as both "probes" and "targets," ELISA, fluorimetry, Fluorescence Resonance Energy Transfer (FRET), SNP-IT, GeneChips, HuSNP, BeadArray, TaqMan assay, Invader assay, MassExtend, or MassCleave™ (hMC) method.

**[0312]** Some type IIS restriction enzymes also display alternative cutting as discussed above. For example, BsmFI will cut at 10/14 and 11/15 from the recognition site. However, the cutting patterns are not mutually exclusive; if the 11/15 cutting pattern is seen at a particular sequence, 10/14 cutting is also seen. If the restriction enzyme BsmF I cuts at 10/14 from the recognition site, the 5' overhang will be  $X_1X_2X_3X_4$ . If BsmF I cuts 11/15 from the recognition site, the 5' overhang will be  $X_0X_1X_2X_3$ . If position  $X_0$  of the overhang is complementary to the labeled nucleotide, the labeled nucleotide will be incorporated at position  $X_0$  and provides an additional level of quality assurance. It provides additional sequence information.

**[0313]** For example, if the variable site is adenine or guanine, and position 3 in the overhang is complementary to adenine, labeled ddATP can be used to determine the genotype at the variable site. If position 0 of the 11/15 overhang contains the nucleotide complementary to adenine, ddATP will be filled in and an additional signal will be seen.

**[0314]** Heterozygous:

10/14 Allele 1	5' CCA	<b>A*</b>			
	3' GGT	T	G	T	G
Overhang position		1	2	3	4
10/14 Allele 2	5' CCA	<b>G</b>	C	<b>A*</b>	
	3' GGT	C	G	T	G
Overhang position		1	2	3	4

11/15 Allele 1	5' CC	A*			
	3' GG	T	T	G	T
Overhang position		0	1	2	3
11/15 Allele 2	5' CC	A*			
	3' GG	T	C	G	T
Overhang position		0	1	2	3

**[0315]** Three signals are seen; one corresponding to the ddATP incorporated at position 0 complementary to the overhang, one corresponding to the ddATP incorporated at position 1 complementary to the overhang, and one corresponding to the ddATP incorporated at position 3 complementary to the overhang. The molecules filled in at position 0, 1, and 3 complementary to the overhang differ in molecular weight and can be separated using any technique that discriminates based on molecular weight including but not limited to gel electrophoresis, and mass spectrometry.

**[0316]** For quantitating the ratio of one allele to another allele or when determining the relative amount of a mutant DNA sequence in the presence of wild type DNA sequence, an accurate and highly sensitive method of detection must be used. The alternate cutting displayed by type IIS restriction enzymes may increase the difficulty of determining ratios of one allele to another allele because the restriction enzyme may not display the alternate cutting (11/15) pattern on the two alleles equally. For example, allele 1 may be cut at 10/14 80% of the time, and 11/15 20% of the time. However, because the two alleles may differ in sequence, allele 2 may be cut at 10/14 90% of the time, and 11/15 20% of the time.

**[0317]** For purposes of quantitation, the alternate cutting problem can be eliminated when the nucleotide at position 0 of the overhang is not complementary to the labeled nucleotide. For example, if the variable site corresponds to adenine or guanine, and position 3 of the overhang is complementary to adenine (*i.e.*, a thymidine is located at position 3 of the overhang), labeled ddATP can be used to determine the genotype of the variable site. If position 0 of the overhang generated by the 11/15 cutting properties is not complementary to adenine, (*i.e.*, position 0 of the overhang corresponds to guanine, cytosine, or adenine) no additional signal will be seen from the fragments that were cut 11/15 from the recognition site. Position 0 complementary to the overhang can be filled in with unlabeled nucleotide, eliminating any complexity seen from the alternate cutting pattern of restriction enzymes. This method provides a highly accurate method for quantitating the ratio of a variable site including but not limited to a mutation, or a single nucleotide polymorphism.

**[0318]** For instance, if SNP X can be adenine or guanine, this method of labeling allows quantitation of the alleles that correspond to adenine and the alleles that correspond to guanine, without determining if the restriction enzyme displays any differences between the alleles with regard to alternate cutting patterns.

**[0319]** Heterozygous:

10/14 Allele 1	5' CCG	A*				
	3' GGC	T	G	T	G	
Overhang position		1	2	3	4	
10/14 Allele 2	5' CCG	G	C	A*		
	3' GGC	C	G	T	G	
Overhang position		1	2	3	4	

**[0320]** The overhang generated by the alternate cutting properties of BsmF I is depicted below:

11/15 Allele 1	5' CC				
	3' GG	C	T	G	T
Overhang position		0	1	2	3
11/15 Allele 2	5' CC				
	3' GG	C	C	G	T
Overhang position		0	1	2	3

**[0321]** After filling in with labeled ddATP and unlabeled dGTP, dCTP, dTTP, the following molecules would be generated:

11/15 Allele 1	5' CC	G	A*		
	3' GG	C	T	G	T
Overhang position		0	1	2	3
11/15 Allele 2	5' CC	G	G	C	A*
	3' GG	C	C	G	T
Overhang position		0	1	2	3

**[0322]** Two signals are seen; one corresponding to the molecules filled in with ddATP at position one complementary to the overhang and one corresponding to the molecules filled in with ddATP at position 3 complementary to the overhang. Position 0 of the 11/15 overhang is filled in with unlabeled nucleotide, which eliminates any difficulty in quantitating a ratio for the nucleotide at the variable site on allele 1 and the nucleotide at the variable site on allele 2.



**[0323]** Any nucleotide can be used including adenine, adenine derivatives, adenine homologues, guanine, guanine derivatives, guanine homologues, cytosine, cytosine derivatives, cytosine homologues, thymidine, thymidine derivatives, or thymidine homologues, or any combinations of adenine, adenine derivatives, adenine homologues, guanine, guanine derivatives, guanine homologues, cytosine, cytosine derivatives, cytosine homologues, thymidine, thymidine derivatives, or thymidine homologues.

**[0324]** The nucleotide can be labeled with any chemical group or moiety, including but not limited to radioactive molecules, fluorescent molecules, antibodies, antibody fragments, haptens, carbohydrates, biotin, derivatives of biotin, phosphorescent moieties, luminescent moieties, electrochemiluminescent moieties, chromatic moieties, and moieties having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity. The nucleotide can be labeled with one or more than one type of chemical group or moiety.

**[0325]** In another embodiment, labeled and unlabeled nucleotides can be used. Any combination of deoxynucleotides and dideoxynucleotides can be used including but not limited to labeled dideoxynucleotides and labeled deoxynucleotides; labeled dideoxynucleotides and unlabeled deoxynucleotides; unlabeled dideoxynucleotides and unlabeled deoxynucleotides; and unlabeled dideoxynucleotides and labeled deoxynucleotides.

**[0326]** In another embodiment, nucleotides labeled with a chemical moiety can be used in the PCR reaction. Unlabeled nucleotides then are used to fill-in the 5' overhangs generated after digestion with the restriction enzyme. An unlabeled terminating nucleotide can be used to in the presence of unlabeled nucleotides to determine the sequence of the alleles of a locus of interest.

**[0327]** For example, if labeled dTTP was used in the PCR reaction, the following 5' overhang would be generated after digestion with BsmF I:

10/14 Allele 1	5' CT*G	A				
	3' GA C	T	G	T	G	
Overhang position		1	2	3	4	
10/14 Allele 2	5' CT*G	G	C	A		
	3' GA C	C	G	T	G	
Overhang position		1	2	3	4	

**[0328]** Unlabeled ddATP, unlabeled dCTP, unlabeled dGTP, and unlabeled dTTP can be used to fill-in the 5' overhang. Two signals will be generated; one signal corresponds to the DNA molecules filled in with unlabeled ddATP at position 1 complementary to the overhang and the second signal corresponds to DNA molecules filled in with unlabeled ddATP at position 3 complementary to the overhang. The DNA

molecules can be separated based on molecular weight and can be detected by the fluorescence of the dTTP, which was incorporated during the PCR reaction.

**[0329]** The labeled DNA loci of interest sites can be analyzed by a variety of methods including but not limited to fluorescence detection, DNA sequencing gel, capillary electrophoresis on an automated DNA sequencing machine, microchannel electrophoresis, and other methods of sequencing, mass spectrometry, time of flight mass spectrometry, quadrupole mass spectrometry, magnetic sector mass spectrometry, electric sector mass spectrometry infrared spectrometry, ultraviolet spectrometry, potentiostatic amperometry or by DNA hybridization techniques including Southern Blots, Slot Blots, Dot Blots, and DNA microarrays, wherein DNA fragments would be useful as both "probes" and "targets," ELISA, fluorimetry, Fluorescence Resonance Energy Transfer (FRET), SNP-IT, GeneChips, HuSNP, BeadArray, TaqMan assay, Invader assay, MassExtend, or MassCleave™ (hMC) method.

**[0330]** This method of labeling is extremely sensitive and allows the detection of alleles of a locus of interest that are in various ratios including but not limited to 1:1, 1:2, 1:3, 1:4, 1:5, 1:6-1:10, 1:11-1:20, 1:21-1:30, 1:31-1:40, 1:41-1:50, 1:51-1:60, 1:61-1:70, 1:71-1:80, 1:81-1:90, 1:91-1:100, 1:101-1:200, 1:250, 1:251-1:300, 1:301-1:400, 1:401-1:500, 1:501-1:600, 1:601-1:700, 1:701-1:800, 1:801-1:900, 1:901-1:1000, 1:1001-1:2000, 1:2001-1:3000, 1:3001-1:4000, 1:4001-1:5000, 1:5001-1:6000, 1:6001-1:7000, 1:7001-1:8000, 1:8001-1:9000, 1:9001-1:10,000, 1:10,001-1:20,000, 1:20,001-1:30,000, 1:30,001-1:40,000, 1:40,001-1:50,000, and greater than 1:50,000.

**[0331]** For example, this method of labeling allows one nucleotide labeled with one signal generating moiety to be used to determine the sequence of alleles at a SNP locus, or detect a mutant allele amongst a population of normal alleles, or detect an allele encoding antibiotic resistance from a bacterial cell amongst alleles from antibiotic sensitive bacteria, or detect an allele from a drug resistant virus amongst alleles from drug-sensitive virus, or detect an allele from a non-pathogenic bacterial strain amongst alleles from a pathogenic bacterial strain.

**[0332]** As shown above, a single nucleotide can be used to determine the sequence of the alleles at a particular locus of interest. This method is especially useful for determining if an individual is homozygous or heterozygous for a particular mutation or to determine the sequence of the alleles at a particular SNP site. This method of labeling eliminates any errors caused by the quantum coefficients of various dyes. It also allows the reaction to proceed in a single reaction vessel including but not limited to a well of a microtiter plate, or a single eppendorf tube.

**[0333]** This method of labeling is especially useful for the detection of multiple genetic signals in the same sample. For example, this method is useful for the detection of fetal DNA in the blood, serum, or

plasma of a pregnant female, which contains both maternal DNA and fetal DNA. The maternal DNA and fetal DNA may be present in the blood, serum or plasma at ratios such as 97:3; however, the above-described method can be used to detect the fetal DNA. This method of labeling can be used to detect two, three, four or more than four different genetic signals in the sample population

**[0334]** This method of labeling is especially useful for the detection of a mutant allele that is among a large population of wild type alleles. Furthermore, this method of labeling allows the detection of a single mutant cell in a large population of wild type cells. For example, this method of labeling can be used to detect a single cancerous cell among a large population of normal cells. Typically, cancerous cells have mutations in the DNA sequence. The mutant DNA sequence can be identified even if there is a large background of wild type DNA sequence. This method of labeling can be used to screen, detect, or diagnosis any type of cancer including but not limited to colon, renal, breast, bladder, liver, kidney, brain, lung, prostate, and cancers of the blood including leukemia.

**[0335]** This labeling method can also be used to detect pathogenic organisms, including but not limited to bacteria, fungi, viruses, protozoa, and mycobacteria. It can also be used to discriminate between pathogenic strains of microorganism and non-pathogenic strains of microorganisms including but not limited to bacteria, fungi, viruses, protozoa, and mycobacteria.

**[0336]** For example, there are several strains of *Escherichia coli* (*E. coli*), and most are non-pathogenic. However, several strains, such as *E. coli* O157 are pathogenic. There are genetic differences between non-pathogenic *E. coli* strains and pathogenic *E. coli*. The above described method of labeling can be used to detect pathogenic microorganisms in a large population of non-pathogenic organisms, which are sometimes associated with the normal flora of an individual.

## **VI. Analysis of the locus of interest**

**[0337]** The loci of interest can be analyzed by a variety of methods including but not limited to fluorescence detection, DNA sequencing gel, capillary electrophoresis on an automated DNA sequencing machine, (e.g. the ABI Prism 3100 Genetic Analyzer or the ABI Prism 3700 Genetic Analyzer), microchannel electrophoresis, and other methods of sequencing, Sanger dideoxy sequencing, mass spectrometry, time of flight mass spectrometry, quadrupole mass spectrometry, magnetic sector mass spectrometry, electric sector mass spectrometry infrared spectrometry, ultraviolet spectrometry, potentiostatic amperometry or by DNA hybridization techniques including Southern Blot, Slot Blot, Dot Blot, and DNA microarray, wherein DNA fragments would be useful as both “probes” and “targets,” ELISA, fluorimetry, fluorescence polarization, Fluorescence Resonance Energy Transfer (FRET), SNP-

IT, GeneChips, HuSNP, BeadArray, TaqMan assay, Invader assay, MassExtend, or MassCleave™ (hMC) method.

**[0338]** The loci of interest can be analyzed using gel electrophoresis followed by fluorescence detection of the incorporated nucleotide. Another method to analyze or read the loci of interest is to use a fluorescent plate reader or fluorimeter directly on the 96-well streptavidin coated plates. The plate can be placed onto a fluorescent plate reader or scanner such as the Pharmacia 9200 Typhoon to read each locus of interest.

**[0339]** Alternatively, the PCR products of the loci of interest can be pooled and after “filling in” (FIG. 10), the products can be separated by size, using any method appropriate for the same, and then analyzed using a variety of techniques including but not limited to fluorescence detection, DNA sequencing gel, capillary electrophoresis on an automated DNA sequencing machine, microchannel electrophoresis, other methods of sequencing, Sanger dideoxy sequencing, DNA hybridization techniques including Southern Blot, Slot Blot, Dot Blot, and DNA microarray, mass spectrometry, time of flight mass spectrometry, quadrupole mass spectrometry, magnetic sector mass spectrometry, electric sector mass spectrometry infrared spectrometry, ultraviolet spectrometry, potentiostatic amperometry. For example, polyacrylamide gel electrophoresis can be used to separate DNA by size and the gel can be scanned to determine the color of fluorescence in each band (using e.g., ABI 377 DNA sequencing machine or a Pharmacia Typhoon 9200).

**[0340]** In another embodiment, the sequence of the locus of interest can be determined by detecting the incorporation of a nucleotide that is 3' to the locus of interest, wherein said nucleotide is a different nucleotide from the possible nucleotides at the locus of interest. This embodiment is especially useful for the sequencing and detection of SNPs. The efficiency and rate at which DNA polymerases incorporate nucleotides varies for each nucleotide.

**[0341]** According to the data from the Human Genome Project, 99% of all SNPs are binary. The sequence of the human genome can be used to determine a nucleotide that is 3' to the SNP of interest. When a nucleotide that is 3' to the SNP site differs from the possible nucleotides at the SNP site, a nucleotide that is one or more than one base 3' to the SNP can be used to determine the sequence of the SNP site.

**[0342]** For example, suppose the sequence of SNP X on chromosome 13 is to be determined. The sequence of the human genome indicates that SNP X can either be adenosine or guanine and that a nucleotide 3' to the locus of interest is a thymidine. A primer that contains a restriction enzyme recognition site for BsmF I, which is designed to be 13 bases from the locus of interest after

amplification, is used to amplify a DNA fragment containing SNP X. Digestion with the restriction enzyme BsmF I generates a 5' overhang that contains the locus of interest, which can either be adenosine or guanine. The digestion products can be split into two "fill in" reactions: one contains dTTP, and the other reaction contains dCTP. If the locus of interest is homozygous for guanine, only the DNA molecules that were mixed with dCTP will be filled in. If the locus of interest is homozygous for adenosine, only the DNA molecules that were mixed with dTTP will be filled in. If the locus of interest is heterozygous, the DNA molecules that were mixed with dCTP will be filled in as well as the DNA molecules that were mixed with dTTP. After washing to remove the excess dNTP, the samples are filled in with labeled ddATP, which is complementary to the nucleotide (thymidine) that is 3' to the locus of interest. The DNA molecules that were filled in by the previous reaction will be filled in with labeled ddATP. If the individual is homozygous for adenosine, the DNA molecules that were mixed with dTTP subsequently will be filled in with the labeled ddATP. However, the DNA molecules that were mixed with dCTP, would not have incorporated that nucleotide, and therefore, could not incorporate the ddATP. Detection of labeled ddATP only in the molecules that were mixed with dTTP indicates that the nucleotide at SNP X on chromosome 13 is adenosine.

**[0343]** In another embodiment, large scale screening for the presence or absence of single nucleotide polymorphisms or mutations can be performed. One to tens to hundreds to thousands of loci of interest on a single chromosome or on multiple chromosomes can be amplified with primers as described above in the "Primer Design" section. The primers can be designed so that each amplified loci of interest is of a different size (FIG. 2). The multiple loci of interest can be of a DNA sample from one individual representing multiple loci of interest on a single chromosome, multiple chromosomes, multiple genes, a single gene, or any combination thereof.

**[0344]** When human data is being analyzed, the known sequence can be a specific sequence that has been determined from one individual (including e.g. the individual whose DNA is currently being analyzed), or it can be a consensus sequence such as that published as part of the human genome.

#### **Ratio of Alleles at Heterozygous Locus of Interest**

**[0345]** In one embodiment, the ratio of alleles at a heterozygous locus of interest can be calculated. The intensity of a nucleotide at the loci of interest can be quantified using any number of computer programs including but not limited to GeneScan and ImageQuant. For example, for a heterozygous SNP, there are two nucleotides, and each should be present in a 1:1 ratio. In a preferred embodiment, the ratio of multiple heterozygous SNPs can be calculated.

**[0346]** In one embodiment, the ratio for a variable nucleotide at alleles at a heterozygous locus of interest can be calculated. The intensity of each variable nucleotide present at the loci of interest can be quantified using any number of computer programs including but not limited to GeneScan and ImageQuant. For example, for a heterozygous SNP, there will be two nucleotides present, and each may be present in a 1:1 ratio. In a preferred embodiment, the ratio of multiple heterozygous SNPs can be calculated.

**[0347]** In another embodiment, the ratio of alleles at a heterozygous locus of interest on a chromosome is summed and compared to the ratio of alleles at a heterozygous locus of interest on a different chromosome. In a preferred embodiment, the ratio of alleles at multiple heterozygous loci of interest on a chromosome is summed and compared to the ratio of alleles at multiple heterozygous loci of interest on a different chromosome. The ratio obtained from SNP 1, SNP 2, SNP 3, SNP 4, etc on chromosome 1 can be summed. This ratio can then be compared to the ratio obtained from SNP A, SNP B, SNP C, SNP D, etc.

**[0348]** For example, 100 SNPs can be analyzed on chromosome 1. Of these 100 SNPs, assume 50 are heterozygous. The ratio of the alleles at heterozygous SNPs on chromosome 1 can be summed, and should give a ratio of approximately 50:50. Likewise, of 100 SNPs analyzed on chromosome 21, assume 50 are heterozygous. The ratio of alleles at heterozygous SNPs on chromosome 21 is summed. With a normal number of chromosomes, the ratio should be approximately 50:50, and thus there should be no difference between the ratio obtained from chromosome 1 and 21. However, if there is an additional copy of chromosome 21, an additional allele will be provided, and the ratio should be approximately 66:33. Thus, the ratio for nucleotides at heterozygous SNPs can be used to detect the presence or absence of chromosomal abnormalities. Any chromosomal abnormality can be detected including aneuploidy, polyploidy, inversion, a trisomy, a monosomy, duplication, deletion, deletion of a part of a chromosome, addition, addition of a part of chromosome, insertion, a fragment of a chromosome, a region of a chromosome, chromosomal rearrangement, and translocation. The method is especially useful for the detection of trisomy 13, trisomy 18, trisomy 21, XXY, and XYY.

**[0349]** The present invention provides a method to quantitate a ratio for the alleles at a heterozygous locus of interest. The loci of interest include but are not limited to single nucleotide polymorphisms, mutations. There is no need to amplify the entire sequence of a gene or to quantitate the amount of a particular gene product. The present invention does not rely on quantitative PCR.

### **Detection of Fetal Chromosomal Abnormalities**

**[0350]** As discussed above in the section entitled "DNA template," the template DNA can be obtained from a sample of a pregnant female, wherein the template DNA comprises maternal template DNA and fetal template DNA. In one embodiment, the template DNA is obtained from the blood of a pregnant female. In a preferred embodiment, the template DNA is obtained from the plasma or serum from the blood of a pregnant female.

**[0351]** In one embodiment, the template DNA from the sample from the pregnant female comprises both maternal template DNA and fetal template DNA. In another embodiment, maternal template DNA is obtained from any nucleic acid containing source including but not limited to cell, tissue, blood, serum, plasma, saliva, urine, tears, vaginal secretion, lymph fluid, cerebrospinal fluid, mucosa secretion, peritoneal fluid, ascitic fluid, fecal matter, or body exudates, and sequenced to identify homozygous or heterozygous loci of interest, which are the loci of interest analyzed on the template DNA obtained from the sample from the pregnant female.

**[0352]** In a preferred embodiment, the sequence of the alleles of multiple loci of interest on maternal template DNA is determined to identify homozygous loci of interest. In another embodiment, the sequence of the alleles of multiple loci of interest on maternal template DNA is determined to identify heterozygous loci of interest. The sequence of the alleles of multiple loci of interest on maternal template DNA can be determined in a single reaction or in multiple reactions.

**[0353]** For example, if 100 maternal loci of interest on chromosome 21 and 100 maternal loci of interest on chromosome 1 are analyzed, one would predict approximately 50 loci of interest on each chromosome to be homozygous and 50 to be heterozygous. The 50 homozygous loci of interest, or the 50 heterozygous loci of interest or the 50 homozygous and 50 heterozygous loci of interest, or any combination of the homozygous and heterozygous loci of interest on each chromosome can be analyzed using the template DNA from the sample from the pregnant female.

**[0354]** The locus of interest on the template DNA from the sample of the pregnant female is analyzed using the amplification, isolation, digestion, fill in, and detection methods described above. The same primers used to analyze the locus of interest on the maternal template DNA are used to screen the template DNA from the sample from the pregnant female. Any number of loci of interest can be analyzed on the template DNA from the sample from the pregnant female. For example, 1, 1-5, 5-10, 10-20, 20-30, 30-40, 40-50, 50-60, 60-70, 70-80, 80-90, 90-100, 100-150, 150-200, 200-250, 250-300, 300-500, 500-1000, 1000-2000, 2000-3000, 3000-4000 or more than 4000 homozygous maternal loci of interest

can be analyzed in the template DNA from the sample from the pregnant female. In a preferred embodiment, multiple loci of interest on multiple chromosomes are analyzed.

**[0355]** From the population of homozygous maternal loci of interest, there will be both heterozygous and homozygous loci of interest from the template DNA from the sample from the pregnant female; the heterozygous loci of interest can be further analyzed. At heterozygous loci of interest, the ratio of alleles can be used to determine the number of chromosomes that are present.

**[0356]** The percentage of fetal DNA present in the sample from the pregnant female can be calculated by determining the ratio of alleles at a heterozygous locus of interest on a chromosome that is not typically associated with a chromosomal abnormality. In a preferred embodiment, the ratio of alleles at multiple heterozygous loci of interest on a chromosome can be used to determine the percentage of fetal DNA. For example, chromosome 1, which is the largest chromosome in the human genome, can be used to determine the percentage of fetal DNA.

**[0357]** For example, suppose SNP X is homozygous at the maternal template DNA (A/A). At SNP X, the template DNA from the sample from the pregnant female, which can contain both fetal DNA and maternal DNA, is heterozygous (A/G). The nucleotide guanine represents the fetal DNA because at SNP X the mother is homozygous, and thus the guanine is attributed to the fetal DNA. The guanine at SNP X can be used to calculate the percentage of fetal DNA in the sample.

**[0358]** Alternatively, multiple loci of interest on two or more chromosomes can be examined to determine the percentage of fetal DNA. For example, multiple loci of interest can be examined on chromosomes 13, and 18 to determine the percentage of fetal DNA because organisms with chromosomal abnormalities at chromosome 13 and 18 are not viable.

**[0359]** Alternatively, for a male fetus, a marker on the Y chromosome can be used to determine the amount of fetal DNA present in the sample. A panel of serial dilutions can be made using the template DNA isolated from the sample from the pregnant female, and quantitative PCR analysis performed. Two PCR reactions can be performed: one PCR reaction to amplify a marker on the Y chromosome, for example SRY, and the other reaction to amplify a region on any of the autosomal chromosomes. The amount of fetal DNA can be calculated using the following formula:

**[0360]** Percent Fetal DNA: (last dilution Y chromosome detected / last dilution autosomal chromosome detected) \* 2 \* 100.

**[0361]** If at SNP A, the mother is homozygous A/A, and the fetus is heterozygous A/G, then the ratio of A:G can be used to detect chromosomal abnormalities. If the fetal DNA is fifty percent (50%) of the DNA in the maternal blood, then at SNP A where the maternal nucleotide is an adenine and the other



nucleotide is a guanine, one would expect the ratio of adenine (two adenines from the maternal template DNA and one from the fetal template DNA) to guanine (from the fetal template DNA) to be 25:75 or 0.33. However, if the fetus has a trisomy of this particular chromosome, and the additional chromosome is contributed by the mother, and thus an additional adenine nucleotide is present, then one would expect the ratio of 0.25 ( $50 \text{ (G)} / (2 \times 50 \text{ maternal A} + 2 \times 50 \text{ fetal A})$ ). Thus, there is a difference of 8% between the ratio obtained from a chromosome present in two copies, and a chromosome present in a trisomy condition. On the other hand, if the additional chromosome is contributed by the father, and thus, an additional guanine is present, then one would expect the ratio of 0.66 ( $2 \times 50 \text{ for G fetal allele} / (2 \times 50 \text{ maternal A allele} + 50 \text{ for fetal A allele})$ ).

**[0362]** However, if the fetal DNA is 40% of the DNA in the maternal blood, the expected ratio without a trisomy is 0.25 ( $40 \text{ for fetal G allele} / 2 \times 60 \text{ for maternal A allele} + 1 \times 60 \text{ for fetal A allele}$ ). If the fetus has a trisomy, and the additional chromosome is provided by the mother, the expected ratio would be 0.20 ( $40 \text{ for fetal G allele} / (2 \times 60 \text{ for maternal A allele} + 2 \times 40 \text{ for fetal A allele})$ ). A 5% difference between the ratios obtained from a chromosome present in two copies and a chromosome present in the Trisomy condition is detected.

**[0363]** In another embodiment, multiple loci of interest on multiple chromosomes can be examined. The ratios for the alleles at each heterozygous locus of interest on a chromosome can be summed and compared to the ratios for the alleles at each locus of interest on a different chromosome. The chromosomes that are compared can be of human origin, and include but are not limited to chromosomes 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X, and Y. The ratio obtained from multiple chromosomes can be compared to the ratio obtained for a single chromosome or from multiple chromosomes.

**[0364]** In one embodiment, one of the chromosomes used in the comparison can be chromosome 13, 15, 16; 18, 21, 22, X or Y. In a preferred embodiment, the ratios on chromosomes 13, 18, and 21 are compared.

**[0365]** For example, assuming 40% fetal DNA in the sample from the pregnant female, the ratio of the alleles at a heterozygous locus of interest on chromosome 1 will be 0.25 ( $40 \text{ for fetal G allele} / (2 \times 60 \text{ for maternal A allele} + 40 \text{ for fetal A allele})$ ). Likewise, the ratio of alleles at a heterozygous locus of interest on chromosome 21 will be present in a ratio of 0.25. However, in a fetus with trisomy 21 where the additional chromosome is contributed by the mother, the nucleotides at a heterozygous locus of interest on chromosome 21 will be present in a ratio of 0.20 ( $40 \text{ for fetal G allele} / (60 \times 2 \text{ for maternal A allele} + 40 \times 2 \text{ for fetal A allele})$ ). By contrast, the ratio for chromosome 1 will remain at 0.25, and thus the

5% difference in ratios will signify an additional chromosome. One to tens to hundreds to thousands of loci of interest can be analyzed.

[0366] In another embodiment, the loci of interest on the template DNA from the sample from the pregnant female can be genotyped without prior identification of the homozygous maternal loci of interest. It is not necessary to genotype the maternal template DNA prior to analysis of the template DNA containing both maternal and fetal template DNA.

[0367] The ratio of the alleles at the loci of interest can be used to determine the presence or absence of a chromosomal abnormality. The template DNA from the sample from the pregnant female contains both maternal template DNA and fetal template DNA. There are 3 possibilities at each SNP for either the maternal template DNA or the fetal template DNA: heterozygous, homozygous for allele 1, or homozygous for allele 2. The possible nucleotide ratios for a SNP that is either an adenine or a guanine are shown in Table II. The ratios presented in Table II are calculated with the fetal DNA at 50% of the DNA in the sample from the pregnant female.

**Table II. Ratios for nucleotides for a heterozygous SNP.**

Fetal SNP

Maternal SNP	A/A	G/G	A/G
A/A	100% A	N/A	75% A, 25%G
G/G	N/A	100% G	25% A, 75% G
A/G	75% A, 25%G	25% A, 75% G	50% A, 50% G

[0368] There are three nucleotide ratios: 100% of a single nucleotide, 50:50, or 75:25. These ratios will vary depending on the amount of fetal DNA present in sample from the pregnant female. However, the percentage of fetal DNA should be constant regardless of the chromosome analyzed. Therefore, if chromosomes are present in two copies, the above calculated ratios will be seen.

[0369] On the other hand, these percentages will vary when an additional chromosome is present. For example, assume that SNP X can be adenine or guanine, and that the percentage of fetal DNA in the sample from the pregnant female is 50%. Analysis of the loci of interest on chromosome 1 will provide the ratios discussed above: 100:0, 50:50, and 75:25. The possible ratios for a SNP that is A/G with an additional chromosome are provided in Table III.

**Table III: Nucleotides ratios at a SNP when an additional copy of a chromosome is present**

Maternal SNPX	Fetal SNP			
	A/A/A	G/G/G	A/G/G	A/A/G
A/A	100% A	N/A	60% A, 40%G	80% A, 20%G
G/G	N/A	100% G	20% A, 80% G	40% A, 60% G
A/G	80% A, 20% G	20% A, 80% G	40% A, 60% G	60% A, 40% G

**[0370]** The possible ratios for the alleles at a heterozygous SNP with an additional copy of a chromosome are: 0:100, 40:60, and 20:80. Two of these ratios, 40:60, and 20:80 differ from the ratios of alleles at heterozygous SNPs obtained with two copies of a chromosome. As discussed above, the ratios for the nucleotides at a heterozygous SNP depend on the amount of fetal DNA present in the sample. However, the ratios, whatever they are, will remain constant across chromosomes unless there is a chromosomal abnormality.

**[0371]** The ratio of alleles at heterozygous loci of interest on a chromosome can be compared to the ratio for alleles at heterozygous loci of interest on a different chromosome. For example, the ratio for multiple loci of interest on chromosome 1 (the ratio at SNP 1, SNP 2, SNP 3, SNP 4, etc.) can be compared to the ratio for multiple loci of interest on chromosome 21 (the ratio at SNP A, SNP B, SNP C, SNP D, etc.). Any chromosome can be compared to any other chromosome. There is no limit to the number of chromosomes that can be compared.

**[0372]** Referring back to the data in Tables II and III, the ratios for nucleotides at a heterozygous SNP on chromosome 1, which was present in two copies, were 25:75, and 50:50. On the other, the ratio for nucleotides at a heterozygous SNP on chromosome 21, which was present in three copies, were 40:60, and 20:80. The difference between these two ratios indicates a chromosomal abnormality. The ratios can be pre-calculated for the full range of varying degrees of fetal DNA present in the maternal serum. Tables II and III demonstrate that both maternal homozygous and heterozygous loci of interest can be used to detect the presence of a fetal chromosomal abnormality.

**[0373]** The above example illustrates how the ratios for nucleotides at heterozygous SNPs can be used to detect the presence of an additional chromosome. The same type of analysis can be used to detect chromosomal rearrangements, translocations, mini-chromosomes, duplications of regions of chromosomes, monosomies, deletions of regions of chromosomes, and fragments of chromosomes. The

method does not require genotyping of the mother or the father, however, it may be done to reduce the number of SNPs that need to be analyzed with the plasma sample.

[0374] The present invention does not quantitate the amount of a fetal gene product, nor is the utility of the present invention limited to the analysis of genes found on the Y chromosome. The present invention does not merely rely on the detection of a paternally inherited nucleic acid, rather, the present invention provides a method that allows the ratio of maternal to fetal alleles at loci of interest, including SNPs, to be calculated.

[0375] In another embodiment, a single allele at a locus of interest can be used to determine the presence or absence of a chromosomal abnormality and detect a genetic disorder in the fetus. In a preferred embodiment, the maternal allele at a locus of interest is used to determine the presence or absence of a chromosomal abnormality in the fetus. The biological mother can be genotyped to identify a homozygous locus of interest. Likewise, the biological father can be genotyped to identify a homozygous locus of interest. The locus of interest wherein the maternal template DNA is homozygous for one allele and the paternal template DNA is homozygous for the other allele is analyzed using the template DNA obtained from the plasma of the mother, which contains both maternal and fetal template DNA. Any number of loci of interest can be analyzed including but not limited to 1, 1-5, 5-10, 10-20, 20-30, 30-40, 40-50, 50-60, 60-70, 70-80, 80-90, 90-100, 100-150, 150-200, 200-250, 250-300, 300-500, 500-1000, 1000-2000, 2000-3000, 3000-4000, 4000-8000, 8000-16000, 16000-32000 or greater than 32000 loci of interest.

[0376] In a preferred embodiment, the signal from the maternal genome and the fetal allele, which was inherited from the mother, at the locus of interest is quantitated. For example, if the 5' overhang, which is generated after digestion with the type IIS enzyme, is filled in with a nucleotide that is fluorescently labeled, the intensity of the incorporated dye can be quantitated

Maternal Template DNA - - Homozygous for Adenine

Allele 1	5' CCG	A*				
	3' GGC	T	G	T	G	
Overhang position		1	2	3	4	
Allele 2	5' CCG	A*				
	3' GGC	T	G	T	G	
Overhang position		1	2	3	4	

Paternal Template DNA - - Homozygous for Cytosine

Allele 1	5' CCG	C*				
----------	--------	----	--	--	--	--

	3' GGC	G	G	T	G
Overhang position		1	2	3	4
Allele 2	5' CCG	C*			
	3' GGC	T	G	T	G
Overhang position		1	2	3	4

Template DNA in the plasma - - both maternal template DNA and fetal template DNA

Maternal Template DNA - - Homozygous for Adenine

Allele 1	5' CCG	A*			
	3' GGC	T	G	T	G
Overhang position		1	2	3	4
Allele 2	5' CCG	A*			
	3' GGC	T	G	T	G
Overhang position		1	2	3	4

Fetal Template DNA - -Heterozygous

Allele 1	5' CCG	A*			
	3' GGC	T	G	T	G
Overhang position		1	2	3	4
Allele 2	5' CCG	ddC			
	3' GGC	T	G	T	G
Overhang position		1	2	3	4

**[0377]** The template DNA obtained from the plasma of the pregnant female is filled in with labeled ddATP, and unlabeled ddCTP (depicted as ddC above), ddGTP, and ddTTP. The plasma DNA contains two maternal adenine alleles, and one fetal adenine allele. By filling in with labeled ddATP and unlabeled ddCTP, only the maternal allele and the fetal allele inherited from the mother are detected. The paternal allele is not detected in this manner. The fill-in reactions can be performed as described in the Examples below.

**[0378]** A single locus of interest can be analyzed or multiple loci of interest. The intensity of the maternal allele at multiple loci of interest can be quantitated. An average can be calculated for a chromosome and compared to the average obtained for a different chromosome. For example, the average intensity of the maternal allele and the fetal allele inherited from the mother at chromosome 1 can be compared to the average intensity of the maternal allele and the fetal allele inherited from the mother at

chromosomes 13, 18, or 21. In a preferred embodiment, chromosomes 13, 15, 18, 21, 22, X and Y, when applicable, are compared.

**[0379]** The signal from a locus of interest may be stronger than another locus of interest. However, there is no reason why the signal from the locus of interest on one chromosome would be stronger than the signal from the locus of interest on another chromosome. While the signal from various loci of interest may be variable, the variation should be seen across the genome. The average signal of the loci of interest should be the same when any chromosomes are compared.

**[0380]** The conditions of the PCR reaction can be optimized so that an equivalent amount of PCR product is produced. For example, the concentration of the primers, the concentration of nucleotides, and the number of cycles for each loci of interest can be optimized. In addition, the fill-in reactions can be done under conditions such that any increase in a specific allele can be detected. The fill-in reaction conditions can be optimized to detect any increase in the allele of interest including but not limited to the concentration of reagents, the time of the fill-in reaction, and the temperature of the reaction.

**[0381]** With a normal genetic karyotype, the signal at each locus of interest comprise signal from the maternal genome, and signal from the fetal allele, which was inherited from the mother. The percent of fetal DNA in the sample remains constant, regardless of the chromosome that is analyzed. For example, if at SNP X, the maternal genome is A/A, and the paternal genome is G/G, then the fetal genome will be A/G, and the fetal adenine allele will comprise a specified percentage of the signal from the adenine allele. If the percentage of fetal DNA is 20% in the maternal plasma, then the fetal adenine allele will contribute 20% of the signal for the adenine allele. The contribution of the fetal allele, which was inherited from the mother, will be constant for any locus of interest that is analyzed.

**[0382]** When there is a chromosomal abnormality, the signal from the maternal genome and the fetal allele, which was inherited from the mother, at the loci of interest will differ from the signal observed for other chromosomes. For example, with a Trisomy, the signal at the locus of interest will comprise the maternal genome and two fetal alleles, which were inherited from the mother. The signal from the loci of interest for the chromosome that is present in three copies will have the contribution of an additional fetal allele, which will alter the signal of the alleles at these loci of interest.

**[0383]** In another embodiment, a ratio can be calculated using a single allele and a standard DNA of known quantity. In a preferred embodiment, a ratio is calculated using the alleles of the maternal genome, and the fetal allele, which was inherited from the mother, and a standard DNA. The biological mother can be genotyped to identify a homozygous locus of interest. Likewise, the biological father can be genotyped to identify a homozygous locus of interest. The locus of interest wherein the maternal

template DNA is homozygous for one allele and the paternal template DNA is homozygous for the other allele is analyzed using the template DNA obtained from the plasma of the mother, which contains both maternal and fetal template DNA.

**[0384]** In a preferred embodiment, the signal from the maternal genome and the fetal allele, which was inherited from the mother, at the locus of interest is quantitated. For example, if the 5' overhang, which is generated after digestion with the type IIS enzyme, is filled in with a nucleotide that is fluorescently labeled, the intensity of the incorporated dye can be quantitated.

Template DNA in the plasma - - both maternal template DNA and fetal template DNA

Maternal Template DNA - - Homozygous for Adenine

Allele 1	5' CCG	A*			
	3' GGC	T	G	T	G
Overhang position		1	2	3	4
Allele 2	5' CCG	A*			
	3' GGC	T	G	T	G
Overhang position		1	2	3	4

Fetal Template DNA - -Heterozygous

Allele 1	5' CCG	A*			
	3' GGC	T	G	T	G
Overhang position		1	2	3	4
Allele 2	5' CCG	ddC			
	3' GGC	T	G	T	G
Overhang position		1	2	3	4

**[0385]** The template DNA obtained from the plasma of the pregnant female is filled in with labeled ddATP, and unlabeled ddCTP (depicted as ddC above), ddGTP, and ddTTP. The plasma DNA contains two maternal adenine alleles, and one fetal adenine allele. By filling in with labeled ddATP and unlabeled ddCTP, only the maternal allele and the fetal allele inherited from the mother are detected.

**[0386]** A single locus of interest or multiple loci of interest can be analyzed. For each locus of interest, a DNA molecule is designed to migrate at about the same position as the locus of interest. In a preferred embodiment, the DNA molecule is of known quantity. A ratio is calculated using the alleles of the maternal genome and the fetal allele, which was inherited from the mother, and the DNA molecule designed to migrate at about the same position as the locus of interest. For example, if the locus of interest is designed to migrate at 30 base pairs, the DNA molecule can be designed to migrate at about 30

base pairs including but not limited to 20-25, 25-30, 30-35, 35-45, and greater than 45. The alleles of the maternal genome and the fetal allele, which was inherited from the mother, and the standard DNA molecule can be analyzed in the same reaction or can be analyzed in a separate reaction. The alleles of the maternal genome and the fetal allele, which was inherited from the mother, and the standard DNA molecule can be analyzed in the same lane of a gel or can be analyzed in separate lanes of a gel. The use of standard DNA molecules of known quantity, which are designed to migrate at the same position as the loci of interest, will correct for various factors including but not limited to the intensity of the bands relative to the location on the gel.

**[0387]** The ratio of multiple loci of interest on a chromosome can be quantitated, and an average calculated. The average can be compared to the average obtained for another chromosome. The ratio is used to indicate the presence or absence of a chromosomal abnormality. Analysis of the alleles of the maternal genome and the fetal allele also allows detection of single gene or multi-gene genetic disorders.

**[0388]** Any chromosome of any organism can be analyzed using the methods of the invention. For example, in humans, chromosome 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X or Y can be analyzed using the methods of the invention. The ratio for the alleles at a heterozygous locus of interest on any chromosome can be compared to the ratio for the alleles at a heterozygous locus of interest on any other chromosome.

**[0389]** Thus, the present invention provides a non-invasive technique, which is independent of fetal cell isolation, for rapid, accurate and definitive detection of chromosome abnormalities in a fetus. The present invention also provides a non-invasive method for determining the sequence of DNA from a fetus. The present invention can be used to detect any alternation in gene sequence as compared to the wild type sequence including but not limited to point mutation, reading frame shift, transition, transversion, addition, insertion, deletion, addition-deletion, frame-shift, missense, reverse mutation, and microsatellite alteration.

#### **Detection of Fetal Chromosomal Abnormalities Using Short Tandem Repeats**

**[0390]** Short tandem repeats (STRs) are short sequences of DNA, normally of 2-5 base pairs in length, which are repeated numerous times in a head-tail manner. Tandemly repeated DNA sequences are widespread throughout the human genome, and show sufficient variability among the individuals in a population. Minisatellites have core repeats with 9-80 base pairs.

**[0391]** In another embodiment, short tandem repeats can be used to detect fetal chromosomal abnormalities. Template DNA can be obtained from a nucleic acid containing sample including but not



limited to cell, tissue, blood, serum, plasma, saliva, urine, tears, vaginal secretion, lymph fluid, cerebrospinal fluid, mucosa secretion, peritoneal fluid, ascitic fluid, fecal matter, or body exudates. In another embodiment, a cell lysis inhibitor is added to the nucleic acid containing sample. In a preferred embodiment, the template DNA is obtained from the blood of a pregnant female. In another embodiment, the template DNA is obtained from the plasma or serum from the blood of a pregnant female.

**[0392]** The template DNA obtained from the blood of the pregnant female will contain both fetal DNA and maternal DNA. The fetal DNA comprises STRs from the mother and the father. The variation in the STRs between the mother and father can be used to detect chromosomal abnormalities.

**[0393]** Primers can be designed to amplify short tandem repeats. Any method of amplification can be used including but not limited to polymerase chain reaction, self-sustained sequence reaction, ligase chain reaction, rapid amplification of cDNA ends, polymerase chain reaction and ligase chain reaction, Q-beta phage amplification, strand displacement amplification, and splice overlap extension polymerase chain reaction. In a preferred embodiment, PCR is used.

**[0394]** Any number of short tandem repeats can be analyzed including but not limited to 1-5, 5-10, 10-50, 50-100, 100-200, 200-300, 300-400, 400-500, 500-1000, and greater than 1000. The short tandem repeats can be analyzed in a single PCR reaction or in multiple PCR reactions. In a preferred embodiment, STRs from multiple chromosomes are analyzed.

**[0395]** After amplification, the PCR products can be analyzed by any number of methods including but not restricted to gel electrophoresis, and mass spectrometry. The template DNA from the pregnant female comprises STRs of maternal and paternal origin. The STRs of paternal origin represent the fetal DNA. The paternal and maternal STRs may be identical in length or the maternal and the paternal STRs may differ.

**[0396]** Heterozygous STRs are those of which the maternal and paternal differ in length. The amount of each PCR product can be quantitated for each heterozygous STR. With a normal number of chromosomes, the amount of each PCR product should be approximately equal. However, with an extra chromosome, one of the STR PCR products will be present at a greater amount.

**[0397]** For example, multiple STRs on chromosome 1 can be analyzed on the template DNA obtained from the blood of the pregnant female. Each STR, whether of maternal or paternal origin, should be present at approximately the same amount. Likewise, with two chromosome 21s, each STR should be present at approximately the same amount. However, with a trisomy 21, one of the STR PCR products, when the maternal and paternal differ in length (a heterozygous STR) should be present at a higher amount. The ratio for each heterozygous STR on one chromosome can be compared to the ratio

for each heterozygous STR on a different chromosome, wherein a difference indicates the presence or absence of a chromosomal abnormality.

### **Kits**

**[0398]** The methods of the invention are most conveniently practiced by providing the reagents used in the methods in the form of kits. A kit preferably contains one or more of the following components: written instructions for the use of the kit, appropriate buffers, salts, DNA extraction detergents, primers, nucleotides, labeled nucleotides, 5' end modification materials, and if desired, water of the appropriate purity, confined in separate containers or packages, such components allowing the user of the kit to extract the appropriate nucleic acid sample, and analyze the same according to the methods of the invention. The primers that are provided with the kit will vary, depending upon the purpose of the kit and the DNA that is desired to be tested using the kit.

**[0399]** A kit can also be designed to detect a desired or variety of single nucleotide polymorphisms, especially those associated with an undesired condition or disease. For example, one kit can comprise, among other components, a set or sets of primers to amplify one or more loci of interest associated with Huntington's disease. Another kit can comprise, among other components, a set or sets of primers for genes associated with a predisposition to develop type I or type II diabetes. Still, another kit can comprise, among other components, a set or sets of primers for genes associated with a predisposition to develop heart disease. Details of utilities for such kits are provided in the "Utilities" section below.

### **Utilities**

**[0400]** The methods of the invention can be used whenever it is desired to know the genotype of an individual. The method of the invention is especially useful for the detection of genetic disorders. The method of the invention is especially useful as a non-invasive technique for the detection of genetic disorders in a fetus. In a preferred embodiment, the method of the invention provides a method for identification of single nucleotide polymorphisms.

**[0401]** In a preferred embodiment, the method is useful for detecting chromosomal abnormalities including but not limited to trisomies, monosomies, duplications, deletions, additions, chromosomal rearrangements, translocations, and other aneuploidies. The method is especially useful for the detection of chromosomal abnormalities in a fetus.

**[0402]** In a preferred embodiment, the method of the invention provides a method for identification of the presence of a disease in a fetus, especially a genetic disease that arises as a result of the presence of a genomic sequence, or other biological condition that it is desired to identify in an individual for which it

is desired to know the same. The identification of such sequence in the fetus based on the presence of such genomic sequence can be used, for example, to determine if the fetus is a carrier or to assess if the fetus is predisposed to developing a certain genetic trait, condition or disease. The method of the invention is especially useful in prenatal genetic testing of parents and child.

[0403] Examples of diseases that can be diagnosed by this invention are listed in Table IV.

**TABLE IV**

Achondroplasia
Adrenoleukodystrophy, X-Linked
Agammaglobulinemia, X-Linked
Alagille Syndrome
Alpha-Thalassemia X-Linked Mental Retardation Syndrome
Alzheimer Disease
Alzheimer Disease, Early-Onset Familial
Amyotrophic Lateral Sclerosis Overview
Androgen Insensitivity Syndrome
Angelman Syndrome
Ataxia Overview, Hereditary
Ataxia-Telangiectasia
Becker Muscular Dystrophy also The Dystrophinopathies)
Beckwith-Wiedemann Syndrome
Beta-Thalassemia
Biotinidase Deficiency
Branchiootorenal Syndrome
BRCA1 and BRCA2 Hereditary Breast/Ovarian Cancer
Breast Cancer
CADASIL
Canavan Disease
Cancer
Charcot-Marie-Tooth Hereditary Neuropathy
Charcot-Marie-Tooth Neuropathy Type 1
Charcot-Marie-Tooth Neuropathy Type 2
Charcot-Marie-Tooth Neuropathy Type 4
Charcot-Marie-Tooth Neuropathy Type X
Cockayne Syndrome
Colon Cancer
Contractural Arachnodactyly, Congenital
Craniosynostosis Syndromes (FGFR-Related)
Cystic Fibrosis
Cystinosis
Deafness and Hereditary Hearing Loss
DRPLA (Dentatorubral-Pallidoluysian Atrophy)
DiGeorge Syndrome (also 22q11 Deletion Syndrome)
Dilated Cardiomyopathy, X-Linked

Down Syndrome (Trisomy 21)
Duchenne Muscular Dystrophy (also The Dystrophinopathies)
Dystonia, Early-Onset Primary (DYT1)
Dystrophinopathies, The
Ehlers-Danlos Syndrome, Kyphoscoliotic Form
Ehlers-Danlos Syndrome, Vascular Type
Epidermolysis Bullosa Simplex
Exostoses, Hereditary Multiple
Facioscapulohumeral Muscular Dystrophy
Factor V Leiden Thrombophilia
Familial Adenomatous Polyposis (FAP)
Familial Mediterranean Fever
Fragile X Syndrome
Friedreich Ataxia
Frontotemporal Dementia with Parkinsonism-17
Galactosemia
Gaucher Disease
Hemochromatosis, Hereditary
Hemophilia A
Hemophilia B
Hemorrhagic Telangiectasia, Hereditary
Hearing Loss and Deafness, Nonsyndromic, DFNA (Connexin 26)
Hearing Loss and Deafness, Nonsyndromic, DFNB 1 (Connexin 26)
Hereditary Spastic Paraplegia
Hermansky-Pudlak Syndrome
Hexosaminidase A Deficiency (also Tay-Sachs)
Huntington Disease
Hypochondroplasia
Ichthyosis, Congenital, Autosomal Recessive
Incontinentia Pigmenti
Kennedy Disease (also Spinal and Bulbar Muscular Atrophy)
Krabbe Disease
Leber Hereditary Optic Neuropathy
Lesch-Nyhan Syndrome Leukemias
Li-Fraumeni Syndrome
Limb-Girdle Muscular Dystrophy
Lipoprotein Lipase Deficiency, Familial
Lissencephaly
Marfan Syndrome
MELAS (Mitochondrial Encephalomyopathy, Lactic Acidosis, and Stroke-Like Episodes)
Monosomies
Multiple Endocrine Neoplasia Type 2
Multiple Exostoses, Hereditary Muscular Dystrophy, Congenital
Myotonic Dystrophy
Nephrogenic Diabetes Insipidus

Neurofibromatosis 1
Neurofibromatosis 2
Neuropathy with Liability to Pressure Palsies, Hereditary
Niemann-Pick Disease Type C
Nijmegen Breakage Syndrome Norrie Disease
Oculocutaneous Albinism Type 1
Oculopharyngeal Muscular Dystrophy
Ovarian Cancer
Pallister-Hall Syndrome
Parkin Type of Juvenile Parkinson Disease
Pelizaeus-Merzbacher Disease
Pendred Syndrome
Peutz-Jeghers Syndrome Phenylalanine Hydroxylase Deficiency
Prader-Willi Syndrome
<i>PROP 1</i> -Related Combined Pituitary Hormone Deficiency (CPHD)
Prostate Cancer
Retinitis Pigmentosa
Retinoblastoma
Rothmund-Thomson Syndrome
Smith-Lemli-Opitz Syndrome
Spastic Paraplegia, Hereditary
Spinal and Bulbar Muscular Atrophy (also Kennedy Disease)
Spinal Muscular Atrophy
Spinocerebellar Ataxia Type 1
Spinocerebellar Ataxia Type 2
Spinocerebellar Ataxia Type 3
Spinocerebellar Ataxia Type 6
Spinocerebellar Ataxia Type 7
Stickler Syndrome (Hereditary Arthroophthalmopathy)
Tay-Sachs (also GM2 Gangliosidosis)
Trisomies
Tuberous Sclerosis Complex
Usher Syndrome Type I
Usher Syndrome Type II
Velocardiofacial Syndrome (also 22q11 Deletion Syndrome)
Von Hippel-Lindau Syndrome
Williams Syndrome
Wilson Disease
X-Linked Adrenoleukodystrophy
X-Linked Agammaglobulinemia
X-Linked Dilated Cardiomyopathy (also The Dystrophinopathies)
X-Linked Hypotonic Facies Mental Retardation Syndrome

[0404] The method of the invention is useful for screening an individual at multiple loci of interest, such as tens, hundreds, or even thousands of loci of interest associated with a genetic trait or genetic

disease by sequencing the loci of interest that are associated with the trait or disease state, especially those most frequently associated with such trait or condition. The invention is useful for analyzing a particular set of diseases including but not limited to heart disease, cancer, endocrine disorders, immune disorders, neurological disorders, musculoskeletal disorders, ophthalmologic disorders, genetic abnormalities, trisomies, monosomies, transversions, translocations, skin disorders, and familial diseases.

[0405] The method of the invention can also be used to confirm or identify the relationship of a DNA of unknown sequence to a DNA of known origin or sequence, for example, for use in, maternity or paternity testing, and the like.

[0406] Having now generally described the invention, the same will become better understood by reference to certain specific examples which are included herein for purposes of illustration only and are not intended to be limiting unless ~~other-wise~~otherwise specified.

### EXAMPLES

[0407] The following examples are illustrative only and are not intended to limit the scope of the invention as defined by the claims.

#### EXAMPLE 1

[0408] DNA sequences were amplified by PCR, wherein the annealing step in cycle 1 was performed at a specified temperature, and then increased in cycle 2, and further increased in cycle 3 for the purpose of reducing non-specific amplification. The TM1 of cycle 1 of PCR was determined by calculating the melting temperature of the 3' region, which anneals to the template DNA, of the second primer. For example, in FIG. 1B, the TM1 can be about the melting temperature of region "c." The annealing temperature was raised in cycle 2, to TM2, which was about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer. For example, in FIG. 1C, the annealing temperature (TM2) corresponds to the melting temperature of region "b." In cycle 3, the annealing temperature was raised to TM3, which was about the melting temperature of the entire sequence of the second primer. For example, in FIG. 1D, the annealing temperature (TM3) corresponds to the melting temperature of region "c" + region "d". The remaining cycles of amplification were performed at TM3.

#### Preparation of Template DNA

[0409] The template DNA was prepared from a 5 ml sample of blood obtained by venipuncture from a human volunteer with informed consent. Blood was collected from 36 volunteers. Template DNA was isolated from each blood sample using QIAamp DNA Blood Midi Kit supplied by QIAGEN (Catalog

number 51183). Following isolation, the template DNA from each of the 36 volunteers was pooled for further analysis.

### Primer Design

**[0410]** The following four single nucleotide polymorphisms were analyzed: SNP HC21S00340, identification number as assigned by Human Chromosome 21 cSNP Database, (FIG. 3, lane 1) located on chromosome 21; SNP TSC 0095512 (FIG. 3, lane 2) located on chromosome 1, SNP TSC 0214366 (FIG. 3, lane 3) located on chromosome 1; and SNP TSC 0087315 (FIG. 3, lane 4) located on chromosome 1. The SNP Consortium Ltd database can be accessed at <http://snp.cshl.org/>, website address effective as of February 14, 2002.

**[0411]** SNP HC21S00340 was amplified using the following primers:

First primer:

5'TAGAATAGCACTGAATTCAGGAATACAATCATTGTCAC 3' (SEQ ID NO: 9)

Second primer:

5'ATCACGATAAACGGCCAAACTCAGGTTA3' (SEQ ID NO: 10)

**[0412]** SNP TSC0095512 was amplified using the following primers:

First primer:

5'AAGTTTAGATCAGAATTCGTGAAAGCAGAAGTTGTCTG 3' (SEQ ID NO: 11)

Second primer:

5'TCTCCAACTAACGGCTCATCGAGTAAAG 3' (SEQ ID NO: 12)

**[0413]** SNP TSC0214366 was amplified using the following primers:

First primer:

5'ATGACTAGCTATGAATTCGTTCAAGGTAGAAAATGGAA 3' (SEQ ID NO: 13)

Second primer:

5'GAGAATTAGAACGGCCCAAATCCCACTC3' (SEQ ID NO: 14)

**[0414]** SNP TSC 0087315 was amplified using the following primers:

First primer:

5'TTACAATGCATGAATTCATCTTGGTCTCTCAAAGTGC 3' (SEQ ID NO: 15)

Second primer:

5'TGGACCATAAACGGCCAAAAACTGTAAG 3' (SEQ ID NO: 16).

**[0415]** All primers were designed such that the 3' region was complementary to either the upstream or downstream sequence flanking each locus of interest and the 5' region contained a restriction enzyme

recognition site. The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI. The second primer contained the recognition site for the restriction enzyme BceA I.

### **PCR Reaction**

**[0416]** All four loci of interest were amplified from the template genomic DNA using PCR (U.S. Patent Nos. 4,683,195 and 4,683,202). The components of the PCR reaction were as follows: 40 ng of template DNA, 5  $\mu$ M first primer, 5  $\mu$ M second primer, 1 X HotStarTaq Master Mix as obtained from Qiagen (Catalog No. 203443). The HotStarTaq Master Mix contained DNA polymerase, PCR buffer, 200  $\mu$ M of each dNTP, and 1.5 mM MgCl<sub>2</sub>.

**[0417]** Amplification of each template DNA that contained the SNP of interest was performed using three different series of annealing temperatures, herein referred to as low stringency annealing temperature, medium stringency annealing temperature, and high stringency annealing temperature. Regardless of the annealing temperature protocol, each PCR reaction consisted of 40 cycles of amplification. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN. As instructed by the manufacturer, the reactions were incubated at 95°C for 15 min. prior to the first cycle of PCR. The denaturation step after each extension step was performed at 95°C for 30 sec. The annealing reaction was performed at a temperature that permitted efficient extension without any increase in temperature.

**[0418]** The low stringency annealing reaction comprised three different annealing temperatures in each of the first three cycles. The annealing temperature for the first cycle was 37°C for 30 sec.; the annealing temperature for the second cycle was 57°C for 30 sec.; the annealing temperature for the third cycle was 64°C for 30 sec. Annealing was performed at 64°C for subsequent cycles until completion.

**[0419]** As shown in the photograph of the gel (FIG. 3A), multiple bands were observed after amplification of SNP TSC 0087315 (lane 4). Amplification of SNP HC21S00340 (lane 1), SNP TSC0095512 (lane 2), and SNP TSC0214366 (lane 3) generated a single band of high intensity and one band of faint intensity, which was of higher molecular weight. When the low annealing temperature conditions were used, the correct size product was generated and this was the predominant product in each reaction.

**[0420]** The medium stringency annealing reaction comprised three different annealing temperatures in each of the first three cycles. The annealing temperature for the first cycle was 40°C for 30 seconds; the annealing temperature for the second cycle was 60°C for 30 seconds; and the annealing temperature



for the third cycle was 67°C for 30 seconds. Annealing was performed at 67°C for subsequent cycles until completion. Similar to what was observed under low stringency annealing conditions, amplification of SNP TSC0087315 (FIG. 3B, lane 4) generated multiple bands under conditions of medium stringency. Amplification of the other three SNPs (lanes 1-3) produced a single band. These results demonstrate that variable annealing temperatures can be used to cleanly amplify loci of interest from genomic DNA with a primer that has an annealing length of 13 bases.

**[0421]** The high stringency annealing reaction was comprised of three different annealing temperatures in each of the first three cycles. The annealing temperature of the first cycle was 46°C for 30 seconds; the annealing temperature of the second cycle was 65°C for 30 seconds; and the annealing temperature for the third cycle was 72°C for 30 seconds. Annealing was performed at 72°C for subsequent cycles until completion. As shown in the photograph of the gel (FIG. 3C), amplification of SNP TSC0087315 (lane 4) using the high stringency annealing temperatures generated a single band of the correct molecular weight. By raising the annealing temperatures for each of the first three cycles, non-specific amplification was eliminated. Amplification of SNP TSC0095512 (lane 2) generated a single band. SNPs HC21S00340 (lane 1), and TSC0214366 (lane 3) failed to amplify at the high stringency annealing temperatures, however, at the medium stringency annealing temperatures, these SNPs amplified as a single band. These results demonstrate that variable annealing temperatures can be used to reduce non-specific PCR products, as demonstrated for SNP TSC0087315 (FIG. 3, lane 4).

## **EXAMPLE 2**

**[0422]** SNPs on chromosomes 1 (TSC0095512), 13 (TSC0264580), and 21 (HC21S00027) were analyzed. SNP TSC0095512 was analyzed using two different sets of primers, and SNP HC21S00027 was analyzed using two types of reactions for the incorporation of nucleotides.

### **Preparation of Template DNA**

**[0423]** The template DNA was prepared from a 5 ml sample of blood obtained by venipuncture from a human volunteer with informed consent. Template DNA was isolated using the QIAmp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183). The template DNA was isolated as per instructions included in the kit. Following isolation, template DNA from thirty-six human volunteers were pooled together and cut with the restriction enzyme EcoRI. The restriction enzyme digestion was performed as per manufacturer's instructions.

## Primer Design

[0424] SNP HC21S00027 was amplified by PCR using the following primer set:

First primer:

5' ATAACCGTATGCGAATTCTATAATTTTCCTGATAAAGG 3' (SEQ ID NO: 17)

Second primer:

5' CTTAAATCAGGGGACTAGGTAAACTTCA 3' (SEQ ID NO: 18).

[0425] The first primer contained a biotin tag at the extreme 5' end, and the nucleotide sequence for the restriction enzyme EcoRI. The second primer contained the nucleotide sequence for the restriction enzyme BsmF I (FIG. 4A).

[0426] Also, SNP HC21S00027 was amplified by PCR using the same first primer but a different second primer with the following sequence:

Second primer:

5' CTTAAATCAGACGGCTAGGTAAACTTCA 3' (SEQ ID NO: 19)

[0427] This second primer contained the recognition site for the restriction enzyme BceA I (FIG. 4B).

[0428] SNP TSC0095512 was amplified by PCR using the following primers:

First primer:

5' AAGTTTAGATCAGAATTCGTGAAAGCAGAAGTTGTCTG 3' (SEQ ID NO: 11)

Second primer:

5' TCTCCAACCTAGGGACTCATCGAGTAAAG 3' (SEQ ID NO: 20).

[0429] The first primer had a biotin tag at the 5' end and contained a restriction enzyme recognition site for EcoRI. The second primer contained a restriction enzyme recognition site for BsmF I (FIG. 4C).

[0430] Also, SNP TSC0095512 was amplified using the same first primer and a different second primer with the following sequence:

Second primer:

5' TCTCCAACCTAACGGCTCATCGAGTAAAG 3' (SEQ ID NO: 12)

[0431] This second primer contained the recognition site for the restriction enzyme BceA I (FIG. 4D).

[0432] SNP TSC0264580, which is located on chromosome 13, was amplified with the following primers:

First primer:

5' AACGCCGGGCGAGAATTCAGTTTTTCAACTTGCAAGG 3' (SEQ ID NO: 21)

Second primer:

5' CTACACATATCTGGGACGTTGGCCATCC 3' (SEQ ID NO: 22).

[0433] The first primer contained a biotin tag at the extreme 5' end and had a restriction enzyme recognition site for EcoRI. The second primer contained a restriction enzyme recognition site for BsmF I.

#### PCR Reaction

[0434] All loci of interest were amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). In this example, the loci of interest were amplified in separate reaction tubes but they could also be amplified together in a single PCR reaction. For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction can be optimized for each locus of interest but in this example, 40 ng of template human genomic DNA and 5  $\mu$ M of each primer were used. Forty cycles of PCR were performed. The following PCR conditions were used:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 64°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty nine (39) times;
- (9) 72°C for 5 minutes.

[0435] In the first cycle of PCR, the annealing temperature was about the melting temperature of the 3' annealing region of the second primers, which was 37°C. The annealing temperature in the second cycle of PCR was about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer, which was 57°C. The annealing temperature in the third cycle of PCR was about the melting temperature of the entire sequence of the second primer, which was 64°C. The annealing temperature for the remaining cycles was 64°C. Escalating the annealing temperature from TM1 to TM2 to TM3 in the first three cycles of PCR greatly improves specificity. These annealing temperatures are representative, and the skilled artisan will understand the annealing temperatures for each cycle are dependent on the specific primers used.

**[0436]** The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results. The PCR products for SNP HC21S00027 and SNP TSC095512 are shown in FIGS. 5A-5D.

#### **Purification of Fragment of Interest**

**[0437]** The PCR products were separated from the genomic template DNA. Each PCR product was divided into four separate reaction wells of a Streptawell, transparent, High-Bind plate from Roche Diagnostics GmbH (catalog number 1 645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog). The first primers contained a 5' biotin tag so the PCR products bound to the Streptavidin coated wells while the genomic template DNA did not. The streptavidin binding reaction was performed using a Thermomixer (Eppendorf) at 1000 rpm for 20 min. at 37°C. Each well was aspirated to remove unbound material, and washed three times with 1X PBS, with gentle mixing (Kandpal et al., Nucl. Acids Res. 18:1789-1795 (1990); Kaneoka et al., Biotechniques 10:30-34 (1991); Green et al., Nucl. Acids Res. 18:6163-6164 (1990)).

#### **Restriction Enzyme Digestion of Isolated Fragments**

**[0438]** The purified PCR products were digested with the restriction enzyme that bound the recognition site incorporated into the PCR products from the second primer. SNP HC21S00027 (FIG. 6A and 6B) and SNP TSC0095512 (FIG. 6C and 6D) were amplified in separate reactions using two different second primers. FIG. 6A (SNP HC21S00027) and FIG. 6C (SNP TSC0095512) depict the PCR products after digestion with the restriction enzyme BsmF I (New England Biolabs catalog number R0572S). FIG. 6B (SNP HC21S00027) and FIG. 6D (SNP TSC0095512) depict the PCR products after digestion with the restriction enzyme BceA I (New England Biolabs, catalog number R0623 S). The digests were performed in the Streptawells following the instructions supplied with the restriction enzyme. SNP TSC0264580 was digested with BsmF I. After digestion with the appropriate restriction enzyme, the wells were washed three times with PBS to remove the cleaved fragments.

#### **Incorporation of Labeled Nucleotide**

**[0439]** The restriction enzyme digest described above yielded a DNA fragment with a 5' overhang, which contained the SNP site or locus of interest and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase.

**[0440]** For each SNP, four separate fill in reactions were performed; each of the four reactions contained a different fluorescently labeled dideoxynucleotide (ddATP, ddCTP, ddGTP, or ddTTP). The

following components were added to each fill in reaction: 1  $\mu$ l of a fluorescently labeled dideoxynucleotide, 0.5  $\mu$ l of unlabeled ddNTPs (40  $\mu$ M), which contained all nucleotides except the nucleotide that was fluorescently labeled, 2  $\mu$ l of 10X sequenase buffer, 0.25  $\mu$ l of Sequenase, and water as needed for a 20  $\mu$ l reaction. All of the fill in reactions were performed at 40°C for 10 min.

Non-fluorescently labeled nucleotides was purchased from Fermentas Inc. (Hanover, MD). All other labeling reagents were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565). In the presence of fluorescently labeled ddNTPs, the 3' recessed end was extended by one base, which corresponds to the SNP or locus of interest (FIG 7A-7D).

**[0441]** A mixture of labeled ddNTPs and unlabeled dNTPs also was used for the “fill in” reaction for SNP HC21S00027. The “fill in” conditions were as described above except that a mixture containing 40  $\mu$ M unlabeled dNTPs, 1  $\mu$ l fluorescently labeled ddATP, 1  $\mu$ l fluorescently labeled ddCTP, 1  $\mu$ l fluorescently labeled ddGTP, and 1  $\mu$ l ddTTP was used. The fluorescent ddNTPs were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565; Amersham did not publish the concentrations of the fluorescent nucleotides). SNP HC21S00027 was digested with the restriction enzyme BsmF I, which generated a 5' overhang of four bases. As shown in FIG. 7E, if the first nucleotide incorporated is a labeled dideoxynucleotide, the 3' recessed end is filled in by one base, allowing detection of the SNP or locus of interest. However, if the first nucleotide incorporated is a dNTP, the polymerase continues to incorporate nucleotides until a ddNTP is filled in. For example, the first two nucleotides can be filled in with dNTPs, and the third nucleotide with a ddNTP, allowing detection of the third nucleotide in the overhang. Thus, the sequence of the entire 5' overhang can be determined, which increases the information obtained from each SNP or locus of interest.

**[0442]** After labeling, each Streptawell was rinsed with 1X PBS (100  $\mu$ l) three times. The “filled in” DNA fragments were then released from the Streptawells by digestion with the restriction enzyme EcoRI, according to the manufacturer's instructions that were supplied with the enzyme (FIGS. 8A-8D).

Digestion was performed for 1 hour at 37°C with shaking at 120 rpm.

#### **Detection of the Locus of Interest**

**[0443]** After release from the streptavidin matrix, 2-3  $\mu$ l of the 10  $\mu$ l sample was loaded in a 48 well membrane tray (The Gel Company, catalog number TAM48-01). The sample in the tray was absorbed with a 48 Flow Membrane Comb (The Gel Company, catalog number AM48), and inserted into a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691).

**[0444]** The sample was electrophoresed into the gel at 3000 volts for 3 min. The membrane comb was removed, and the gel was run for 3 hours on an ABI 377 Automated Sequencing Machine. The incorporated labeled nucleotide was detected by fluorescence.

**[0445]** As shown in FIG. 9A, from a sample of thirty six (36) individuals, one of two nucleotides, either adenosine or guanine, was detected at SNP HC21S00027. These are the two nucleotides reported to exist at SNP HC21S00027 (<http://snp.cshl.org/snpsearch.shtml>).

**[0446]** One of two nucleotides, either guanine or cytosine, was detected at SNP TS00095512 (FIG. 9B). The same results were obtained whether the locus of interest was amplified with a second primer that contained a recognition site for BceA I or the second primer contained a recognition site for BsmF I.

**[0447]** As shown in FIG. 9C, one of two nucleotides was detected at SNP TSC0264580, which was either adenosine or cytosine. These are the two nucleotides reported for this SNP site (<http://snp.cshl.org/snpsearch.shtml>). In addition, a thymidine was detected one base from the locus of interest. In a sequence dependent manner, BsmF I cuts some DNA molecules at the 10/14 position and other DNA molecules, which have the same sequence, at the 11/15 position. When the restriction enzyme BsmF I cuts 11 nucleotides away on the sense strand and 15 nucleotides away on the antisense strand, the 3' recessed end is one base from the SNP site. The sequence of SNP TSC0264580 indicated that the base immediately preceding the SNP site was a thymidine. The incorporation of a labeled ddNTP into this position generated a fragment one base smaller than the fragment that was cut at the 10/14 position. Thus, the DNA molecules cut at the 11/15 position provided sequence information about the base immediately preceding the SNP site, and the DNA molecules cut at the 10/14 position provided sequence information about the SNP site.

**[0448]** SNP HC21S00027 was amplified using a second primer that contained the recognition site for BsmF I. A mixture of labeled ddNTPs and unlabeled dNTPs was used to fill in the 5' overhang generated by digestion with BsmF I. If a dNTP was incorporated, the polymerase continued to incorporate nucleotides until a ddNTP was incorporated. A population of DNA fragments, each differing by one base, was generated, which allowed the full sequence of the overhang to be determined.

**[0449]** As seen in FIG. 9D, an adenosine was detected, which was complementary to the nucleotide (a thymidine) immediately preceding the SNP or locus of interest. This nucleotide was detected because of the 11/15 cutting property of BsmF I, which is described in detail above. A guanine and an adenosine were detected at the SNP site, which are the two nucleotides reported for this SNP site (FIG. 9A). The two nucleotides were detected at the SNP site because the molecular weights of the dyes differ, which allowed separation of the two nucleotides. The next nucleotide detected was a thymidine, which is

complementary to the nucleotide immediately downstream of the SNP site. The next nucleotide detected was a guanine, which was complementary to the nucleotide two bases downstream of the SNP site. Finally, an adenosine was detected, which was complementary to the third nucleotide downstream of the SNP site. Sequence information was obtained not only for the SNP site but for the nucleotide immediately preceding the SNP site and the next three nucleotides.

**[0450]** None of the loci of interest contained a mutation. However, if one of the loci of interest harbored a mutation including but not limited to a point mutation, insertion, deletion, translocation or any combination of said mutations, it could be identified by comparison to the consensus or published sequence. Comparison of the sequences attributed to each of the loci of interest to the native, non-disease related sequence of the gene at each locus of interest determines the presence or absence of a mutation in that sequence. The finding of a mutation in the sequence is then interpreted as the presence of the indicated disease, or a predisposition to develop the same, as appropriate, in that individual. The relative amounts of the mutated vs. normal or non-mutated sequence can be assessed to determine if the subject has one or two alleles of the mutated sequence, and thus whether the subject is a carrier, or whether the indicated mutation results in a dominant or recessive condition.

### **EXAMPLE 3**

**[0451]** Four loci of interest from chromosome 1 and two loci of interest from chromosome 21 were amplified in separate PCR reactions, pooled together, and analyzed. The primers were designed so that each amplified locus of interest was a different size, which allowed detection of the loci of interest.

#### **Preparation of Template DNA**

**[0452]** The template DNA was prepared from a 5 ml sample of blood obtained by venipuncture from a human volunteer with informed consent. Template DNA was isolated using the QIAmp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183). The template DNA was isolated as per instructions included in the kit. Template DNA was isolated from thirty-six human volunteers, and then pooled into a single sample for further analysis.

#### **Primer Design**

**[0453]** SNP TSC 0087315 was amplified using the following primers:

First primer:

5'TTACAATGCATGAATTCATCTTGGTCTCTCAAAGTGC 3' (SEQ ID NO: 15)

Second primer:

5'TGGACCATAAACGGCCAAAACTGTAAG3' (SEQ ID NO: 16).

**[0454]** SNP TSC0214366 was amplified using the following primers:

First primer:

5'ATGACTAGCTATGAATTCGTTCAAGGTAGAAAATGGAA 3' (SEQ ID NO: 13)

Second primer:

5'GAGAATTAGAACGGCCCAAATCCCACTC 3' (SEQ ID NO: 14)

**[0455]** SNP TSC 0413944 was amplified with the following primers:

First primer:

5' TACCTTTTGATCGAATTCAAGGCCAAAAATATTAAGTT 3' (SEQ ID NO: 23)

Second primer:

5' TCGAACTTTAACGGCCTTAGAGTAGAGA 3' (SEQ ID NO: 24)

**[0456]** SNP TSC0095512 was amplified using the following primers:

First primer:

5'AAGTTTAGATCAGAATTCGTGAAAGCAGAAGTTGTCTG 3' (SEQ ID NO: 11)

Second primer:

5'TCTCCAATAACGGCTCATCGAGTAAAG 3' (SEQ ID NO: 12)

**[0457]** SNP HC21S00131 was amplified with the following primers:

First primer:

5' CGATTTTCGATAAGAATTCAAAAGCAGTTCTTAGTTCAG 3' (SEQ ID NO: 25)

Second primer:

5'TGCGAATCTTACGGCTGCATCACATTCA 3' (SEQ ID NO: 26)

**[0458]** SNP HC21S00027 was amplified with the following primers:

First primer:

5' ATAACCGTATGCGAATTCTATAATTTTCCTGATAAAGG 3' (SEQ ID NO: 17)

Second primer:

5' CTTAAATCAGACGGCTAGGTAAACTTCA 3' (SEQ ID NO: 19)

**[0459]** For each SNP, the first primer contained a recognition site for the restriction enzyme EcoRI and had a biotin tag at the extreme 5' end. The second primer used to amplify each SNP contained a recognition site for the restriction enzyme BceA I.



### **PCR Reaction**

[0460] The PCR reactions were performed as described in Example 2 except that the following annealing temperatures were used: the annealing temperature for the first cycle of PCR was 37°C for 30 seconds, the annealing temperature for the second cycle of PCR was 57°C for 30 seconds, and the annealing temperature for the third cycle of PCR was 64°C for 30 seconds. All subsequent cycles had an annealing temperature of 64°C for 30 seconds. Thirty seven (37) cycles of PCR were performed. After PCR, 1/4 of the volume was removed from each reaction, and combined into a single tube.

### **Purification of Fragment of Interest**

[0461] The PCR products (now combined into one sample, and referred to as “the sample”) were separated from the genomic template DNA as described in Example 2 except that the sample was bound to a single well of a Streptawell microtiter plate.

### **Restriction Enzyme Digestion of Isolated Fragments**

[0462] The sample was digested with the restriction enzyme BceA I, which bound the recognition site in the second primer. The restriction enzyme digestions were performed following the instructions supplied with the enzyme. After the restriction enzyme digest, the wells were washed three times with 1X PBS.

### **Incorporation of Nucleotides**

[0463] The restriction enzyme digest described above yielded DNA molecules with a 5' overhang, which contained the SNP site or locus of interest and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide in the presence of a DNA polymerase.

[0464] The following components were used for the fill in reaction: 1 µl of fluorescently labeled ddATP; 1 µl of fluorescently labeled ddTTP; 1 µl of fluorescently labeled ddGTP; 1 µl of fluorescently labeled ddCTP; 2 µl of 10X sequenase buffer, 0.25 µl of Sequenase, and water as needed for a 20 µl reaction. The fill in reaction was performed at 40°C for 10 min. All labeling reagents were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit (US 79565); the concentration of the ddNTPS provided in the kit is proprietary and not published by Amersham). In the presence of fluorescently labeled ddNTPs, the 3' recessed end was filled in by one base, which corresponds to the SNP or locus of interest.

[0465] After the incorporation of nucleotide, the Streptawell was rinsed with 1X PBS (100 µl) three times. The “filled in” DNA fragments were then released from the Streptawell by digestion with the

restriction enzyme EcoRI following the manufacturer's instructions. Digestion was performed for 1 hour at 37 °C with shaking at 120 rpm.

### **Detection of the Locus of Interest**

**[0466]** After release from the streptavidin matrix, 2-3 µl of the 10 µl sample was loaded in a 48 well membrane tray (The Gel Company, catalog number TAM48-01). The sample in the tray was absorbed with a 48 Flow Membrane Comb (The Gel Company, catalog number AM48), and inserted into a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691).

**[0467]** The sample was electrophoresed into the gel at 3000 volts for 3 min. The membrane comb was removed, and the gel was run for 3 hours on an ABI 377 Automated Sequencing Machine. The incorporated nucleotide was detected by fluorescence.

**[0468]** The primers were designed so that each amplified locus of interest differed in size. As shown in FIG. 10, each amplified loci of interest differed by about 5-10 nucleotides, which allowed the loci of interest to be separated from one another by gel electrophoresis. Two nucleotides were detected for SNP TSC0087315, which were guanine and cytosine. These are the two nucleotides reported to exist at SNP TSC0087315 (<http://snp.cshl.org/snpsearch.shtml>). The sample comprised template DNA from 36 individuals and because the DNA molecules that incorporated a guanine differed in molecular weight from those that incorporated a cytosine, distinct bands were seen for each nucleotide.

**[0469]** Two nucleotides were detected at SNP HC21S00027, which were guanine and adenosine (FIG. 10). The two nucleotides reported for this SNP site are guanine and adenosine (<http://snp.cshl.org/snpsearch.shtml>). As discussed above, the sample contained template DNA from thirty-six individuals, and one would expect both nucleotides to be represented in the sample. The molecular weight of the DNA fragments that incorporated a guanine was distinct from the DNA fragments that incorporated an adenosine, which allowed both nucleotides to be detected.

**[0470]** The nucleotide cytosine was detected at SNP TSC0214366 (FIG. 10). The two nucleotides reported to exist at this SNP position are thymidine and cytosine.

**[0471]** The nucleotide guanine was detected at SNP TSC0413944 (FIG. 10). The two nucleotides reported for this SNP are guanine and cytosine (<http://snp.cshl.org/snpsearch.shtml>).

**[0472]** The nucleotide cytosine was detected at SNP TS00095512 (FIG. 10). The two nucleotides reported for this SNP site are guanine and cytosine (<http://snp.cshl.org/snpsearch.shtml>).

[0473] The nucleotide detected at SNP HC21S00131 was guanine. The two nucleotides reported for this SNP site are guanine and adenosine (<http://snp.cshl.org/snpsearch.shtml>).

[0474] As discussed above, the sample was comprised of DNA templates from thirty-six individuals and one would expect both nucleotides at the SNP sites to be represented. For SNP TSC0413944, TSC0095512, TSC0214366 and HC21S00131, one of the two nucleotides was detected. It is likely that both nucleotides reported for these SNP sites are present in the sample but that one fluorescent dye overwhelms the other. The molecular weight of the DNA molecules that incorporated one nucleotide did not allow efficient separation of the DNA molecules that incorporated the other nucleotide. However, the SNPs were readily separated from one another, and for each SNP, a proper nucleotide was incorporated. The sequences of multiple loci of interest from multiple chromosomes, which were treated as a single sample after PCR, were determined.

[0475] A single reaction containing fluorescently labeled ddNTPs was performed with the sample that contained multiple loci of interest. Alternatively, four separate fill in reactions can be performed where each reaction contains one fluorescently labeled nucleotide (ddATP, ddTTP, ddGTP, or ddCTP) and unlabeled ddNTPs (see Example 2, FIGS. 7A-7D and FIGS. 9A-C). Four separate “fill in” reactions will allow detection of any nucleotide that is present at the loci of interest. For example, if analyzing a sample that contains multiple loci of interest from a single individual, and said individual is heterozygous at one or more than one loci of interest, four separate “fill in” reactions can be used to determine the nucleotides at the heterozygous loci of interest.

[0476] Also, when analyzing a sample that contains templates from multiple individuals, four separate “fill in” reactions will allow detection of nucleotides present in the sample, independent of how frequent the nucleotide is found at the locus of interest. For example, if a sample contains DNA templates from 50 individuals, and 49 of the individuals have a thymidine at the locus of interest, and one individual has a guanine, the performance of four separate “fill in” reactions, wherein each “fill in” reaction is run in a separate lane of a gel, such as in FIGS. 9A-9C, will allow detection of the guanine. When analyzing a sample comprised of multiple DNA templates, multiple “fill in” reactions will alleviate the need to distinguish multiple nucleotides at a single site of interest by differences in mass.

[0477] In this example, multiple single nucleotide polymorphisms were analyzed. It is also possible to determine the presence or absence of mutations, including but not limited to point mutations, transitions, transversions, translocations, insertions, and deletions from multiple loci of interest. The multiple loci of interest can be from a single chromosome or from multiple chromosomes. The multiple loci of interest can be from a single gene or from multiple genes.

[0478] The sequence of multiple loci of interest that cause or predispose to a disease phenotype can be determined. For example, one could amplify one to tens to hundreds to thousands of genes implicated in cancer or any other disease. The primers can be designed so that each amplified loci of interest differs in size. After PCR, the amplified loci of interest can be combined and treated as a single sample. Alternatively, the multiple loci of interest can be amplified in one PCR reaction or the total number of loci of interest, for example 100, can be divided into samples, for example 10 loci of interest per PCR reaction, and then later pooled. As demonstrated herein, the sequence of multiple loci of interest can be determined. Thus, in one reaction, the sequence of one to ten to hundreds to thousands of genes that predispose or cause a disease phenotype can be determined.

#### EXAMPLE 4

[0479] The ability to determine the sequence or detect chromosomal abnormalities of a fetus using free fetal DNA in a sample from a pregnant female has been hindered by the low percentage of free fetal DNA. Increasing the percentage of free fetal DNA would enhance the detection of mutation, insertion, deletion, translocation, transversion, monosomy, trisomy, trisomy 21, trisomy 18, trisomy 13, XXY, XXX, other ~~aneuploidies~~aneuploidies, deletion, addition, amplification, translocation and rearrangement. The percent of fetal DNA in plasma obtained from a pregnant female was determined both in the absence and presence of inhibitors of cell lysis. A genetic marker on the Y chromosome was used to calculate the percent of fetal DNA.

#### Preparation of Template DNA

[0480] The DNA template was prepared from a 5 ml sample of blood obtained by venipuncture from a human volunteer with informed consent. The blood was aliquoted into two tubes (Fischer Scientific, 9 ml EDTA Vacuette tubes, catalog number NC9897284). Formaldehyde (25 µl/ml of blood) was added to one of the tubes. The sample in the other tube remained untreated, except for the presence of the EDTA. The tubes were spun at 1000 rpm for ten minutes. Two milliliters of the supernatant (the plasma) of each sample was transferred to a new tube and spun at 3000 rpm for ten minutes. 800 µl of each sample was used for DNA purification. DNA was isolated using the Qiagen Midi Kit for purification of DNA from blood cells (QIAmp DNA Blood Midi Kit, Catalog number 51183). DNA was eluted in 100 µl of distilled water. Two DNA templates were obtained: one from the blood sample treated with EDTA, and one from the blood sample treated with EDTA and formaldehyde.

## Primer Design

[0481] Two different sets of primers were used: one primer set was specific for the Y chromosome, and thus specific for fetal DNA, and the other primer set was designed to amplify the cystic fibrosis gene, which is present on both maternal template DNA and fetal template DNA.

[0482] In this example, the first and second primers were designed so that the entire 5' and 3' sequence of each primer annealed to the template DNA. In this example, the fetus had an XY genotype, and the Y chromosome was used as a marker for the presence of fetal DNA. The following primers were designed to amplify the SRY gene on the Y chromosome.

First primer:

5' TGGCGATTAAGTCAAATTCGC 3' (SEQ ID NO: 263)

Second primer:

5' CCCCTAGTACCCTGACAATGTATT 3' (SEQ ID NO: 264)

[0483] Primers designed to amplify any gene, or region of a region, or any part of any chromosome could be used to detect maternal and fetal DNA. In this example, the following primers were designed to amplify the cystic fibrosis gene:

First primer:

5' CTGTTCTGTGATATTATGTGTGGT 3' (SEQ ID NO: 265)

Second primer:

5' AATTGTTGGCATTCCAGCATTG 3' (SEQ ID NO: 266)

## PCR Reaction

[0484] The SRY gene and the cystic fibrosis gene were amplified from the template genomic DNA using PCR (U.S. Patent Nos. 4,683,195 and 4,683,202). For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by Qiagen (Catalog No. 203443). For amplification of the SRY gene, the DNA eluted from the Qiagen purification column was diluted serially 1:2. For amplification of the cystic fibrosis gene, the DNA from the Qiagen purification column was diluted 1:4, and then serially diluted 1:2. The following components were used for each PCR reaction: 8 µl of template DNA (diluted or undiluted), 1 µl of each primer (5 µM), 10 µl of HotStar Taq mix. The following PCR conditions were used:

- (1) 95°C for 15'
- (2) 94°C for 1'
- (3) 54°C for 15"

- (4) 72°C for 30''
- (5) Repeat steps 2-4 for 45 cycles.
- (6) 10' at 72°C

### **Quantification of Fetal DNA**

**[0485]** The DNA templates that were eluted from the Qiagen columns were serially diluted to the following concentrations: 1:2, 1:4, 1:8, 1:16, 1:32, 1:64, 1:128, 1:256, 1:512, 1:1024, 1:2048, and 1:4096. Amplification of the SRY gene was performed using the templates that were undiluted, 1:2, 1:4, 1:8, 1:16, 1:32, 1:64, 1:128, 1:256, 1:512. Amplification of the cystic fibrosis gene was performed using the DNA templates that were diluted 1:4, 1:8, 1:16, 1:32, 1:64, 1:128, 1:256, 1:512, 1:1024, 1:2048, and 1:4096. The same dilution series was performed with the DNA templates that were purified from the plasma sample treated with EDTA alone and the plasma sample treated with EDTA and formaldehyde.

**[0486]** The results of the PCR reactions using the DNA template that was isolated from the plasma sample treated with EDTA are shown in FIG. 11A. The SRY gene was amplified from the undiluted DNA template, and also in the sample that was diluted 1:2 (FIG. 11A). The SRY gene was not amplified in the next seven serial dilutions. On the other hand, the cystic fibrosis gene was detected in the serial dilutions up to 1:256. A greater presence of the cystic fibrosis gene was expected because of the higher percentage of maternal DNA present in the plasma. The last dilution sample that provided for amplification of the gene product was assumed to have one copy of the cystic fibrosis gene or the SRY gene.

**[0487]** The results of the PCR reactions using the DNA template that was isolated from the plasma sample treated with formaldehyde and EDTA are shown in FIG. 11B. The SRY gene was amplified from the undiluted DNA template, and also in the sample that was diluted 1:2 (FIG. 11B). The SRY gene was not amplified in the next six dilutions. However, in the 1:256 dilution, the SRY gene was detected. It is unlikely that the amplification in the 1:256 sample represents a real signal because the prior six dilution series were all negative for amplification of SRY. Amplification of the SRY gene in this sample was likely an experimental artifact resulting from the high number of PCR cycles used. Thus, the 1:256 sample was not used in calculating the amount of fetal DNA present in the sample.

**[0488]** Amplification of the cystic fibrosis gene was detected in the sample that was diluted 1:16 (FIG. 11B). The presence of the formalin prevents maternal cell lysis, and thus, there is a lower percentage of maternal DNA in the sample. This is in strong contrast to the sample that was treated with only EDTA, which supported amplification up to a dilution of 1:256.

[0489] The percent of fetal DNA present in the maternal plasma was calculated using the following formula:

$$\% \text{ fetal DNA} = (\text{amount of SRY gene} / \text{amount of cystic fibrosis gene}) * 2 * 100.$$

The amount of SRY gene was represented by the highest dilution value in which the gene was amplified. Likewise, the amount of cystic fibrosis gene was represented by the highest dilution value in which it was amplified. The formula contains a multiplication factor of two (2), which is used to normalize for the fact that there is only one copy of the SRY gene (located on the Y chromosome), while there are two copies of the cystic fibrosis gene.

[0490] For the above example, the percentage of fetal DNA present in the sample that was treated with only EDTA was 1.56 % ( $2/256 * 2 * 100$ ). The reported percentage of fetal DNA present in the plasma is between 0.39-11.9 % (Pertl and Bianchi, *Obstetrics and Gynecology*, Vol. 98, No. 3, 483-490 (2001)). The percentage of fetal DNA present in the sample treated with formalin and EDTA was 25% ( $2/16 * 2 * 100$ ). The experiment was repeated numerous times, and each time the presence of formalin increased the overall percentage of fetal DNA.

[0491] The percent fetal DNA from eighteen blood samples with and without formalin was calculated as described above with the exception that serial dilutions of 1:5 were performed. As 1:5 dilutions were performed, the last serial dilution that allowed detection of either the SRY gene or the cystic fibrosis gene may have had one copy of the gene or it may have had 4 copies of the gene. The results from the eighteen samples with and without formalin are summarized in Table V. The low range assumes that the last dilution sample had one copy of the genes and the high range assumes that the last dilution had four copies of the genes.

**Table V. Mean Percentage Fetal DNA with and without formalin.**

Sample	Lower Range	Upper Range
Formalin	19.47	43.69
Without Formalin	7.71	22.1

[0492] An overall increase in fetal DNA was achieved by reducing the maternal cell lysis, and thus, reducing the amount of maternal DNA present in the sample. In this example, formaldehyde was used to prevent lysis of the cells, however any agent that prevents the lysis of cells or increases the structural integrity of the cells can be used. Two or more than two cell lysis inhibitors can be used. The increase in fetal DNA in the maternal plasma allows the sequence of the fetal DNA to be determined, and provides for the rapid detection of abnormal DNA sequences or chromosomal abnormalities including but not

limited to point mutation, reading frame shift, transition, transversion, addition, insertion, deletion, addition-deletion, frame-shift, missense, reverse mutation, and microsatellite alteration, trisomy, monosomy, other aneuploidies, amplification, rearrangement, translocation, transversion, deletion, addition, amplification, fragment, translocation, and rearrangement.

### EXAMPLE 5

**[0493]** A DNA template from an individual with a genotype of trisomy 21 was analyzed. Three loci of interest were analyzed on chromosome 13 and two loci of interest were analyzed on chromosome 21.

#### Preparation of Template DNA

**[0494]** The template DNA was prepared from a 5 ml sample of blood obtained by venipuncture from a human volunteer with informed consent. The human volunteer had previously been genotyped to have an additional chromosome 21 (trisomy 21). Template DNA was isolated using QIAamp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183).

#### Primer Design

**[0495]** The following five single nucleotide polymorphisms were analyzed: SNP TSC 0115603 located on chromosome 21; SNP TSC 03209610 located on chromosome 21; SNP TSC 0198557 located on chromosome 13; and SNP TSC 0200347 located on chromosome 13. The DNA template from another individual was used as an internal control. The SNP TSC 0200347, which was previously identified as being homozygous for guanine, was used as the internal control. The SNP Consortium Ltd database can be accessed at <http://snp.cshl.org/>, website address effective as of April 1, 2002.

**[0496]** SNP TSC 0115603 was amplified using the following primers:

First Primer:

5' GTGCACTTACGTGAATTCAGATGAACGTGATGTAGTAG 3' (SEQ ID NO: 267)

Second Primer:

5' TCCTCGTACTCAACGGCTTTCTCTGAAT 3' (SEQ ID NO: 268)

**[0497]** The first primer was biotinylated at the 5' end, and contained the restriction enzyme recognition site for EcoR I. The second primer contained the restriction enzyme recognition site for the restriction enzyme BceA I.

**[0498]** SNP TSC 0309610 was amplified using the following primers:

First primer:

5' TCCGGAACACTAGAATTCTTATTTACATACACACTTGT 3' (SEQ ID NO: 269)



Second primer:

5' CGAATAAGGTAGACGGCAACAATGAGAA 3' (SEQ ID NO: 270)

[0499] The first primer contained a biotin group at the 5' end, and a restriction enzyme recognition site for the restriction enzyme EcoR I. The second primer contained the restriction enzyme recognition site for BceA I.

[0500] Submitted SNP (ss) 813773 (accession number assigned by the NCBI Submitted SNP (ss) Database) was amplified with the following primers:

First primer:

5' CGGTAAATCGGAGAATTCAGAGGATTTAGAGGAGCTAA 3' (SEQ ID NO: 271)

Second primer:

5' CTCACGTTTCGTTACGGCCATTGTGATAGC 3' (SEQ ID NO: 272)

[0501] The first primer contains a biotin group at the 5' end, and a recognition site for the restriction enzyme EcoR I. The second primer contained the restriction enzyme recognition site for BceA I.

[0502] SNP TSC 0198557 was amplified with the following primers:

First primer:

5' GGGGAAACAGTAGAATTCCATATGGACAGAGCTGTACT 3' (SEQ ID NO: 273)

Second primer:

5' TGAAGCTGTCGGACGGCCTTTGCCCTCTC 3' (SEQ ID NO: 274)

[0503] The first primer contains a biotin group at the 5' end, and a recognition site for the restriction enzyme EcoR I. The second primer contained the restriction enzyme recognition site for BceA I.

[0504] SNP TSC 0197279 was amplified with the following primers:

First primer:

5' ATGGGCAGTTATGAATTCACACTCCCTGTAGCTTGTT 3' (SEQ ID NO: 275)

Second primer:

5' TGATTGGCGCGAACGGCACTCAGAGAAGA 3' (SEQ ID NO: 276)

[0505] The first primer contained a biotin group at the 5' end, and a recognition site for the restriction enzyme for EcoR I. The second primer contained the recognition site for the restriction enzyme BceA I.

[0506] SNP TSC 0200347 was amplified with the following primers:

First primer:

5' CTCAAGGGGACCGAATTCGCTGGGGTCTTCTGTGGGTC 3' (SEQ ID NO: 277)

Second primer:

5' TAGGGCGGCGTGACGGCCAGCCAGTGGT 3' (SEQ ID NO: 278)

**[0507]** The first primer contained a biotin group at the 5' end, and the recognition site for the restriction enzyme EcoR I. The second primer contained the restriction enzyme recognition site for BceA I.

#### **PCR Reaction**

**[0508]** All five loci of interest were amplified from the template genomic DNA using PCR (U.S. Patent Nos. 4,683,195 and 4,683,202). For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction can be optimized for each locus of interest; in this example, 40 ng of template human genomic DNA and 5  $\mu$ M of each primer were used. Thirty-eight cycles of PCR were performed. The following PCR conditions were used for SNP TSC 0115603, SNP TSC 0309610, and SNP TSC 02003437:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 42°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 60°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 69°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty nine (37) times;
- (9) 72°C for 5 minutes.

**[0509]** The following PCR conditions were used for SNP ss813773, SNP TSC 0198557, and SNP TSC 0197279:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 64°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty nine (37) times; and

(9) 72°C for 5 minutes.

**[0510]** In the first cycle of each PCR, the annealing temperature was about the melting temperature of the 3' annealing region of the second primer. The annealing temperature in the second cycle of PCR was about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer. The annealing temperature in the third cycle of PCR was about the melting temperature of the entire sequence of the second primer. Escalating the annealing temperature from TM1 to TM2 to TM3 in the first three cycles of PCR greatly improves specificity. These annealing temperatures are representative, and the skilled artisan will understand the annealing temperatures for each cycle are dependent on the specific primers used. The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results.

#### **Purification of Fragment of Interest**

**[0511]** PCR products were separated from the components of the PCR reaction using Qiagen's MinElute PCR Purification Kit following manufacturer's instructions (Catalog number 28006). The PCR products were eluted in 20 µl of distilled water. For each amplified SNP, one microliter of PCR product, 1 µl of amplified internal control DNA (SNP TSC 0200347), and 8 µl of distilled water were mixed. Five microliters of each sample was placed into two separate reaction wells of a Pierce StreptaWell Microtiter plate (catalog number 15501). The first primers contained a 5' biotin tag so the PCR products bound to the Streptavidin coated wells while the genomic template DNA did not. The streptavidin binding reaction was performed using a Thermomixer (Eppendorf) at 150 rpm for 1 hour at 45°C. Each well was aspirated to remove unbound material, and washed three times with 1X PBS, with gentle mixing (Kandpal et al., Nucl. Acids Res. 18:1789-1795 (1990); Kaneoka et al., Biotechniques 10:30-34 (1991); Green et al., Nucl. Acids Res. 18:6163-6164 (1990)).

#### **Restriction Enzyme Digestion of Isolated Fragments**

**[0512]** The purified PCR products were digested with the restriction enzyme that bound the recognition site that was incorporated into the PCR products from the second primer. The purified PCR products were digested with the restriction enzyme BceA I (New England Biolabs, catalog number R0623S). The digests were performed in the wells of the microtiter plate following the instructions supplied with the restriction enzyme. After digestion with the appropriate restriction enzyme, the wells were washed three times with PBS to remove the cleaved fragments.

### **Incorporation of Labeled Nucleotide**

**[0513]** The restriction enzyme digest described above yielded a DNA fragment with a 5' overhang, which contained the SNP and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase.

**[0514]** For each SNP, two fill in reactions were performed; each reaction contained a different fluorescently labeled dideoxynucleotide (ddATP, ddCTP, ddGTP, or ddTTP, depending on the reported nucleotides to exist at a particular SNP). For example, the nucleotides adenine and thymidine have been reported at SNP TSC 0115603. Therefore, the digested PCR product for SNP TSC 0115603 was mixed with either fluorescently labeled ddATP or fluorescently labeled ddTTP. Each reaction contained fluorescently labeled ddGTP for the internal control. The following components were added to each fill in reaction: 2  $\mu$ l of a ROX-conjugated dideoxynucleotide (depending on the nucleotides reported for each SNP), 2  $\mu$ l of ROX-conjugated ddGTP (internal control), 2.5  $\mu$ l of 10X sequenase buffer, 2  $\mu$ l of Sequenase, and water as needed for a 25  $\mu$ l reaction. All of the fill in reactions were performed at 45°C for 45 min. However, shorter time periods of incorporation can be used. Non-fluorescently labeled ddNTPs were purchased from Fermentas Inc. (Hanover, MD). The ROX-conjugated ddNTPs were obtained from Perkin Elmer. In the presence of fluorescently labeled ddNTPs, the 3' recessed end was extended by one base, which corresponds to the SNP or locus of interest.

**[0515]** After labeling, each Streptawell was rinsed with 1X PBS (100  $\mu$ l) three times. The “filled in” DNA fragments were then released from the Streptawells by digestion with the restriction enzyme EcoR I following manufacturer’s recommendations. Digestion was performed for 1 hour at 37°C with shaking at 120 rpm.

### **Detection of the Locus of Interest**

**[0516]** After release from the streptavidin matrix, 3  $\mu$ l of the 10  $\mu$ l sample was loaded in a 48 well membrane tray (The Gel Company, catalog number TAM48-01). The sample in the tray was absorbed with a 48 Flow Membrane Comb (The Gel Company, catalog number AM48), and inserted into a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691).

**[0517]** The sample was electrophoresed into the gel at 3000 volts for 3 min. The membrane comb was removed, and the gel was run for 3 hours on an ABI 377 Automated Sequencing Machine. The incorporated labeled nucleotide was detected by fluorescence.

**[0518]** As seen in FIG. 12, SNP TSC 0115603 was “filled in” with labeled ddTTP (lane 1) and in a separate reaction with labeled ddATP (lane 3). The calculated ratio between the nucleotides, using the raw data, was 66:34, which is consistent with the theoretical ratio of 66:33 for a SNP on chromosome 21 in an individual with trisomy 21. Both the ddTTP and ddATP were labeled with the same fluorescent dye to minimize variability in incorporation efficiencies of the dyes. However, nucleotides with different fluorescent labels or any detectable label can be used. It is preferable to calculate the coefficients of incorporation when different labels are used.

**[0519]** Each fill in reaction was performed in a separate well so it was possible that there could be variability in DNA binding between the wells of the microtiter plate. To account for the potential variability of DNA binding to the streptavidin-coated plates, an internal control was used. The internal control (SNP TSC 0200347), which is homozygous for guanine, was added to the sample prior to splitting the sample into two separate wells, and thus, an equal amount of the internal control should be present in each well. The amount of incorporated ddGTP can be fixed between the two reactions. If the amount of DNA in each well is equal, the amount of incorporated ddGTP should be equal because the reaction is performed under saturating conditions, with saturating conditions being defined as conditions that support incorporation of a nucleotide at each template molecule. Using the internal control, the ratio of incorporated ddATP to ddTTP was 63.4:36.6. This ratio was very similar to the ratio obtained with the raw data, indicating that there are minor differences in the two fill in reactions for a particular SNP.

**Table VI.**

**Allele Frequencies at Multiple SNPs on DNA Template from Individual with Trisomy 21**

SNP	Allele	Peak Area	Allele Ratio	Internal Control	Normalized Peak Area	Allele Ratio (%)
TSC 0115603	A	5599	66	723	5599	63.4
	T	2951	34	661	3227 ((723/661)*2951)	36.6
TSC 0309610	T	4126	64	1424	4126	66.8
	C	2342	36	1631	2045 ((1424/1631)*2342)	33.2
ss813773	A	4199	46	808	4199	41
	C	4870	54	647	6082 ((808/647)*4870)	59
TSC 0198557	T	3385	55	719	3385	49
	C	2741	45	559	3525 (719/559 *2741)	51
TSC 0197279	T	8085	53	2752	8085	50.7
	C	7202	47	2520	7865 (2752/2520 *7202)	49.3

**[0520]** SNP TSC 0309610 was filled in with ddTTP (lane 3) or ddCTP (lane 4) (FIG. 12). The calculated ratio for the nucleotides, using the raw data, was 64:36. Both ddTTP and ddCTP were labeled

with the same fluorescent dye. After normalization to the internal control, as discussed above, the calculated allele ratio of ddTTP to ddCTP was 66.8:33.2 (Table VI). Again, both the calculated ratio from the raw data and the calculated ratio using the internal control are very similar to the theoretical ratio of 66.6:33.4 for a SNP on chromosome 21 in an individual with trisomy.

**[0521]** To demonstrate that the 66:33 ratios for nucleotides at heterozygous SNPS represented loci on chromosomes present in three copies, SNPs on chromosome 13 were analyzed. The individual from whom the blood sample was obtained had previously been genotyped with one maternal chromosome 13, and one paternal chromosome 13.

**[0522]** Submitted SNP (ss) 813773 was filled in with ddATP (lane 5) or ddCTP (lane 6) (FIG. 12). The calculated ratio for the nucleotides at this heterozygous SNP, using the raw data, was 46:54. This ratio is within 10% of the expected ratio of 50:50. Importantly, the ratio does not approach the 66:33 ratio expected when there is an additional copy of a chromosome.

**[0523]** After normalization to the internal control, the calculated ratio was 41:59. Contrary to the expected result, normalization to the internal control increased the discrepancy between the calculated ratio and the theoretical ratio. This result may represent experimental error that occurred in aliquoting the DNA samples.

**[0524]** Also, it is possible that the restriction enzyme used to generate the overhang, which was used as a template for the “fill-in” reaction, preferentially cut one DNA template over the other DNA template. The two templates differ, with respect to the nucleotide at the SNP site, and this may influence the cutting. The primers can be designed such that the nucleotides adjacent to the cut site are the same, independent of the nucleotide at the SNP site (discussed further in the section entitled “Primer Design”).

**[0525]** SNP TSC 0198557, which is on chromosome 13, was filled in with ddTTP (lane 7) in one reaction and ddCTP (lane 8) in another (FIG. 12). The calculated ratio for the nucleotides at this SNP, using the raw data, was 55:45. After normalization to the internal control, the calculated allele ratio of T:C was 49:51. The normalized ratio was closer to the theoretical ratio of 50:50 for an individual with two copies of chromosome 13.

**[0526]** SNP TSC 0197279, which is on chromosome 13, was filled in with ddTTP (lane 9) in one reaction and ddCTP (lane 10) in another (FIG. 12). The calculated ratio for the nucleotides at this SNP, using the raw data, was 53:47. After normalization to the internal control, the calculated allele ratio of T:C was 50.7:49.3. This is consistent with the theoretical ratio of 50:50 for an individual with only two copies of chromosome 13.

**[0527]** The ratio for the nucleotides at two of the analyzed SNPs on chromosome 13 was approximately 50:50. One SNP, ss813773, showed a ratio of 46:54, and when normalized to the internal control, the ratio was 41:59. These ratios deviate from the expected 50:50, but at the same time, the ratios are not indicative of an extra chromosome, which is indicated with a ratio of 66:33. While the data from this particular SNP is inconclusive, it does not represent a false positive. No conclusion could be drawn on the data from this SNP. However, the other two SNPs provided data that indicated a normal number of chromosomes. It is preferable to analyze multiple SNPs on a chromosome including but not limited to 1-5, 5-10, 10-50, 50-100, 100-200, 200-300, 300-400, 400-500, 500-600, 600-700, 700-800, 800-900, 900-1000, 1000-2000, 2000-3000, and greater than 3000. Preferably, the average of the ratios for a particular chromosome will be used to determine the presence or absence of a chromosomal abnormality. However, it is still possible to analyze one locus of interest. In the event that inconclusive data is obtained, another locus of interest can be analyzed.

**[0528]** The individual from whom the DNA template was obtained had previously been genotyped with trisomy 21, and the allele frequencies at SNPs on chromosome 21 indicate the presence of an additional chromosome 21. The additional chromosome contributes an additional nucleotide for each SNP, and thus alters the traditional 50:50 ratio at a heterozygous SNP. These results are consistent for multiple SNPs, and are specific for those found on chromosome 21. The allele frequencies for SNPs on chromosome 13 gave the expected ratios of approximately 50:50. These results demonstrate that this method of SNP detection can be used to detect chromosomal abnormalities including but not limited to translocations, transversions, monosomies, trisomy 21, trisomy 18, trisomy 13, other ~~aneuploidies~~aneuploidies, deletions, additions, amplifications, translocations and rearrangements.

#### **EXAMPLE 6**

**[0529]** Genomic DNA was obtained from four individuals after informed consent was obtained. Six SNPs on chromosome 13 (TSC0837969, TSC0034767, TSC1130902, TSC0597888, TSC0195492, TSC0607185) were analyzed using the template DNA. Information regarding these SNPs can be found at the following website [www.snp.chsl.org/snpsearch.shtml](http://www.snp.chsl.org/snpsearch.shtml); website active as of February 11, 2003).

**[0530]** A single nucleotide labeled with one fluorescent dye was used to genotype the individuals at the six selected SNP sites. The primers were designed to allow the six SNPs to be analyzed in a single reaction.

## Preparation of Template DNA

[0531] The template DNA was prepared from a 9 ml sample of blood obtained by venipuncture from a human volunteer with informed consent. Template DNA was isolated using the QIAmp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183). The template DNA was isolated as per instructions included in the kit.

## Design of Primers

[0532] SNP TSC0837969 was amplified using the following primer set:

First primer:

5' GGGCTAGTCTCCGAATTCACCTATCCTACCAAATGTC 3' (SEQ ID NO: 30)

Second primer:

5' TAGCTGTAGTTAGGGACTGTTCTGAGCAC 3' (SEQ ID NO: 31)

[0533] The first primer had a biotin tag at the 5' end and contained a restriction enzyme recognition site for EcoRI. The first primer was designed to anneal 44 bases from the locus of interest. The second primer contained a restriction enzyme recognition site for BsmF I.

[0534] SNP TSC0034767 was amplified using the following primer set:

First primer:

5' CGAATGCAAGGCGAATTCGTTAGTAATAACACAGTGCA 3' (SEQ ID NO: 32)

Second primer:

5' AAGACTGGATCCGGGACCATGTAGAATAC 3' (SEQ ID NO: 33)

[0535] The first primer had a biotin tag at the 5' end and contained a restriction enzyme recognition site for EcoRI. The first primer was designed to anneal 50 bases from the locus of interest. The second primer contained a restriction enzyme recognition site for BsmF I.

SNP TSC1130902 was amplified using the following primer set:

First primer:

5' TCTAACCATTGCGAATTCAGGGCAAGGGGGGTGAGATC 3' (SEQ ID NO: 34)

Second primer:

5' TGACTTGGATCCGGGACAACGACTCATCC 3' (SEQ ID NO: 35)

[0536] The first primer had a biotin tag at the 5' end and contained a restriction enzyme recognition site for EcoRI. The first primer was designed to anneal 60 bases from the locus of interest. The second primer contained a restriction enzyme recognition site for BsmF I.

[0537] SNP TSC0597888 was amplified using the following primer set:



First primer:

5' ACCCAGGCGCCAGAATTCTTTAGATAAAAGCTGAAGGGA 3' (SEQ ID NO: 36)

Second primer:

5' GTTACGGGATCCGGGACTCCATATTGATC 3' (SEQ ID NO: 37)

[0538] The first primer had a biotin tag at the 5' end and contained a restriction enzyme recognition site for EcoRI. The first primer was designed to anneal 70 bases from the locus of interest. The second primer contained a restriction enzyme recognition site for BsmF I.

[0539] SNP TSC0195492 was amplified using the following primer set:

First primer:

5'CGTTGGCTTGAGGAATTCGACCAAAAGAGCCAAGAGAA (SEQ ID NO: 38)

Second primer:

5' AAAAAGGGATCCGGGACCTTGACTAGGAC 3' (SEQ ID NO: 39)

[0540] The first primer had a biotin tag at the 5' end and contained a restriction enzyme recognition site for EcoRI. The first primer was designed to anneal 80 bases from the locus of interest. The second primer contained a restriction enzyme recognition site for BsmF I.

[0541] SNP TSC0607185 was amplified using the following primer set:

First primer:

5' ACTTGATTCCGTGAATTCGTTATCAATAAATCTTACAT 3' (SEQ ID NO: 40)

Second primer:

5' CAAGTTGGATCCGGGACCCAGGGCTAACC 3' (SEQ ID NO: 41)

[0542] The first primer had a biotin tag at the 5' end and contained a restriction enzyme recognition site for EcoRI. The first primer was designed to anneal 90 bases from the locus of interest. The second primer contained a restriction enzyme recognition site for BsmF I.

[0543] All loci of interest were amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). In this example, the loci of interest were amplified in separate reaction tubes but they could also be amplified together in a single PCR reaction. For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction can be optimized for each locus of interest but in this example, 40 ng of template human genomic DNA and 5  $\mu$ M of each primer were used. Forty cycles of PCR were performed. The following PCR conditions were used:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 64°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty nine (39) times;
- (9) 72°C for 5 minutes.

**[0544]** In the first cycle of PCR, the annealing temperature was about the melting temperature of the 3' annealing region of the second primers, which was 37°C. The annealing temperature in the second cycle of PCR was about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer, which was 57°C. The annealing temperature in the third cycle of PCR was about the melting temperature of the entire sequence of the second primer, which was 64°C. The annealing temperature for the remaining cycles was 64°C. Escalating the annealing temperature from TM1 to TM2 to TM3 in the first three cycles of PCR greatly improves specificity. These annealing temperatures are representative, and the skilled artisan will understand the annealing temperatures for each cycle are dependent on the specific primers used.

**[0545]** The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results. In this example, the first primer was designed to anneal at various distances from the locus of interest. The skilled artisan understands that the annealing location of the first primer can be 5-10, 11-15, 16-20, 21-25, 26-30, 31-35, 36-40, 41-45, 46-50, 51-55, 56-60, 61-65, 66-70, 71-75, 76-80, 81-85, 86-90, 91-95, 96-100, 101-105, 106-110, 111-115, 116-120, 121-125, 126-130, 131-140, 141-160, 161-180, 181-200, 201-220, 221-240, 241-260, 261-280, 281-300, 301-350, 351-400, 401-450, 451-500, or greater than 500 bases from the locus of interest.

#### **Purification of Fragment of Interest**

**[0546]** The PCR products were separated from the genomic template DNA. After the PCR reaction, 1/4 of the volume of each PCR reaction from one individual was mixed together in a well of a Streptawell, transparent, High-Bind plate from Roche Diagnostics GmbH (catalog number 1 645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog). The first primers contained a 5'

biotin tag so the PCR products bound to the Streptavidin coated wells while the genomic template DNA did not. The streptavidin binding reaction was performed using a Thermomixer (Eppendorf) at 1000 rpm for 20 min. at 37°C. Each well was aspirated to remove unbound material, and washed three times with 1X PBS, with gentle mixing (Kandpal et al., Nucl. Acids Res. 18:1789-1795 (1990); Kaneoka et al., Biotechniques 10:30-34 (1991); Green et al., Nucl. Acids Res. 18:6163-6164 (1990)).

### Restriction Enzyme Digestion of Isolated Fragments

**[0547]** The purified PCR products were digested with the restriction enzyme BsmF I, which binds to the recognition site incorporated into the PCR products from the second primer. The digests were performed in the Streptawells following the instructions supplied with the restriction enzyme. After digestion, the wells were washed three times with PBS to remove the cleaved fragments.

### Incorporation of Labeled Nucleotide

**[0548]** The restriction enzyme digest with BsmF I yielded a DNA fragment with a 5' overhang, which contained the SNP site or locus of interest and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase.

**[0549]** Below, a schematic of the 5' overhang for SNP TSC0837969 is shown. The entire DNA sequence is not reproduced, only the portion to demonstrate the overhang (where R indicates the variable site).

	5' TTAA				
	3' AATT	R	A	C	A
Overhang position		1	2	3	4

**[0550]** The observed nucleotides for TSC0837969 on the 5' sense strand (here depicted as the top strand) are adenine and guanine. The third position in the overhang on the antisense strand corresponds to cytosine, which is complementary to guanine. As this variable site can be adenine or guanine, fluorescently labeled ddGTP in the presence of unlabeled dCTP, dTTP, and dATP was used to determine the sequence of both alleles. The fill-in reactions for an individual homozygous for guanine, homozygous for adenine or heterozygous are diagrammed below.

**[0551]** Homozygous for guanine at TSC 0837969:

Allele 1	5' TTAA	G*			
	3' AATT	C	A	C	A
Overhang position		1	2	3	4
Allele 2	5' TTAA	G*			

	3' AATT	C	A	C	A
Overhang position		1	2	3	4

**[0552]** Labeled ddGTP is incorporated into the first position of the overhang. Only one signal is seen, which corresponds to the molecules filled in with labeled ddGTP at the first position of the overhang.

**[0553]** Homozygous for adenine at TSC 0837969:

Allele 1	5' TTAA	A	T	G*	
	3' AATT	T	A	C	A
Overhang position		1	2	3	4
Allele 2	5' TTAA	A	T	G*	
	3' AATT	T	A	C	A
Overhang position		1	2	3	4

**[0554]** Unlabeled dATP is incorporated at position one of the overhang, and unlabeled dTTP is incorporated at position two of the overhang. Labeled ddGTP was incorporated at position three of the overhang. Only one signal will be seen; the molecules filled in with ddGTP at position 3 will have a different molecular weight from molecules filled in at position one, which allows easy identification of individuals homozygous for adenine or guanine.

**[0555]** Heterozygous at TSC0837969:

Allele 1	5' TTAA	G*			
	3' AATT	C	A	C	A
Overhang position		1	2	3	4
Allele 2	5' TTAA	A	T	G*	
	3' AATT	T	A	C	A
Overhang position		1	2	3	4

**[0556]** Two signals will be seen; one signal corresponds to the DNA molecules filled in with ddGTP at position 1, and a second signal corresponding to molecules filled in at position 3 of the overhang. The two signals can be separated using any technique that separates based on molecular weight including but not limited to gel electrophoresis.

**[0557]** Below, a schematic of the 5' overhang for SNP TSC0034767 is shown. The entire DNA sequence is not reproduced, only the portion to demonstrate the overhang (where R indicates the variable site).

A	C	A	R	GTGT 3'
				CACA 5'
4	3	2	1	Overhang Position

**[0558]** The observed nucleotides for TSC0034767 on the 5' sense strand (here depicted as the top strand) are cytosine and guanine. The second position in the overhang corresponds to adenine, which is complementary to thymidine. The third position in the overhang corresponds to cytosine, which is complementary to guanine. Fluorescently labeled ddGTP in the presence of unlabeled dCTP, dTTP, and dATP is used to determine the sequence of both alleles.

**[0559]** In this case, the second primer anneals upstream of the locus of interest, and thus the fill-in reaction occurs on the anti-sense strand (here depicted as the bottom strand). Either the sense strand or the antisense strand can be filled in depending on whether the second primer, which contains the type IIS restriction enzyme recognition site, anneals upstream or downstream of the locus of interest.

**[0560]** Below, a schematic of the 5' overhang for SNP TSC1130902 is shown. The entire DNA sequence is not reproduced, only a portion to demonstrate the overhang (where R indicates the variable site).

	5' TTCAT				
	3' AAGTA	R	T	C	C
Overhang position		1	2	3	4

**[0561]** The observed nucleotides for TSC1130902 on the 5' sense strand (here depicted as the top strand) are adenine and guanine. The second position in the overhang corresponds to a thymidine, and the third position in the overhang corresponds to cytosine, which is complementary to guanine.

**[0562]** Fluorescently labeled ddGTP in the presence of unlabeled dCTP, dTTP, and dATP is used to determine the sequence of both alleles.

**[0563]** Below, a schematic of the 5' overhang for SNP TSC0597888 is shown. The entire DNA sequence is not reproduced, only the portion to demonstrate the overhang (where R indicates the variable site).

T	C	T	R	ATTC 3'
				TAAG 5'
4	3	2	1	Overhang position

**[0564]** The observed nucleotides for TSC0597888 on the 5' sense strand (here depicted as the top strand) are cytosine and guanine. The third position in the overhang corresponds to cytosine, which is

complementary to guanine. Fluorescently labeled ddGTP in the presence of unlabeled dCTP, dTTP, and dATP is used to determine the sequence of both alleles.

**[0565]** Below, a schematic of the 5' overhang for SNP TSC0607185 is shown. The entire DNA sequence is not reproduced, only the portion to demonstrate the overhang (where R indicates the variable site).

C	C	T	R	TGTC 3'
				ACAG 5'
4	3	2	1	Overhang position

**[0566]** The observed nucleotides for TSC0607185 on the 5' sense strand (here depicted as the top strand) are cytosine and thymidine. In this case, the second primer anneals from the locus of interest, which allows the anti-sense strand to be filled in. The anti-sense strand (here depicted as the bottom strand) will be filled in with guanine or adenine.

**[0567]** The second position in the 5' overhang is thymidine, which is complementary to adenine, and the third position in the overhang corresponds to cytosine, which is complementary to guanine. Fluorescently labeled ddGTP in the presence of unlabeled dCTP, dTTP, and dATP is used to determine the sequence of both alleles.

**[0568]** Below, a schematic of the 5' overhang for SNP TSC0195492 is shown. The entire DNA sequence is not reproduced, only the portion to demonstrate the overhang.

	5' ATCT				
	3' TAGA	R	A	C	A
Overhang position		1	2	3	4

**[0569]** The observed nucleotides at this site are cytosine and guanine (here depicted as the top strand). The second position in the 5' overhang is adenine, which is complementary to thymidine, and the third position in the overhang corresponds to cytosine, which is complementary to guanine. Fluorescently labeled ddGTP in the presence of unlabeled dCTP, dTTP, and dATP is used to determine the sequence of both alleles.

**[0570]** As demonstrated above, the sequence of both alleles of the six SNPs can be determined by labeling with ddGTP in the presence of unlabeled dATP, dTTP, and dCTP. The following components were added to each fill in reaction: 1 µl of fluorescently labeled ddGTP, 0.5 µl of unlabeled dNTPs (40 µM), which contained all nucleotides except guanine, 2 µl of 10X sequenase buffer, 0.25 µl of Sequenase, and water as needed for a 20µl reaction. The fill in reaction was performed at 40°C for 10 min. Non-fluorescently labeled dNTP was purchased from Fermentas Inc. (Hanover, MD). All other labeling

reagents were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565).

**[0571]** After labeling, each Streptawell was rinsed with 1X PBS (100 µl) three times. The "filled in" DNA fragments were then released from the Streptawells by digestion with the restriction enzyme EcoRI, according to the manufacturer's instructions that were supplied with the enzyme. Digestion was performed for 1 hour at 37 °C with shaking at 120 rpm.

#### Detection of the Locus of Interest

**[0572]** After release from the streptavidin matrix, the sample was loaded into a lane of a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691). The sample was electrophoresed into the gel at 3000 volts for 3 min. The gel was run for 3 hours on a sequencing apparatus (Hoefer SQ3 Sequencer). The gel was removed from the apparatus and scanned on the Typhoon 9400 Variable Mode Imager. The incorporated labeled nucleotide was detected by fluorescence.

**[0573]** As shown in FIG. 11, the template DNA in lanes 1 and 2 for SNP TSC0837969 is homozygous for adenine. The following fill-in reaction was expected to occur if the individual was homozygous for adenine:

**[0574]** Homozygous for adenine at TSC 0837969:

5' TTAA	A	T	G*	
3' AATT	T	A	C	A
Overhang position	1	2	3	4

**[0575]** Unlabeled dATP was incorporated in the first position complementary to the overhang. Unlabeled dTTP was incorporated in the second position complementary to the overhang. Labeled ddGTP was incorporated in the third position complementary to the overhang. Only one band was seen, which migrated at about position 46 of the acrylamide gel. This indicated that adenine was the nucleotide filled in at position one. If the nucleotide guanine had been filled in, a band would be expected at position 44.

**[0576]** However, the template DNA in lanes 3 and 4 for SNP TSC0837969 was heterozygous. The following fill-in reactions were expected if the individual was heterozygous:

Heterozygous at TSC0837969:

Allele 1 5' TTAA	G*			
3' AATT	C	A	C	A
Overhang position	1	2	3	4

Allele 25' TTAA	A	T	G*	
3' AATT	T	A	C	A
Overhang position	1	2	3	4

**[0577]** Two distinct bands were seen; one band corresponds to the molecules filled in with ddGTP at position 1 complementary to the overhang (the G allele), and the second band corresponds to molecules filled in with ddGTP at position 3 complementary to the overhang (the A allele). The two bands were separated based on the differences in molecular weight using gel electrophoresis. One fluorescently labeled nucleotide ddGTP was used to determine that an individual was heterozygous at a SNP site. This is the first use of a single nucleotide to effectively detect the presence of two different alleles.

**[0578]** For SNP TSC0034767, the template DNA in lanes 1 and 3 is heterozygous for cytosine and guanine, as evidenced by the two distinct bands. The lower band corresponded to ddGTP filled in at position 1 complementary to the overhang. The second band of slightly higher molecular weight corresponded to ddGTP filled in at position 3, indicating that the first position in the overhang was filled in with unlabeled dCTP, which allowed the polymerase to continue to incorporate nucleotides until it incorporated ddGTP at position 3 complementary to the overhang. The template DNA in lanes 2 and 4 was homozygous for guanine, as evidenced by a single band of higher molecular weight than if ddGTP had been filled in at the first position complementary to the overhang.

**[0579]** For SNP TSC1130902, the template DNA in lanes 1, 2, and 4 is homozygous for adenine at the variable site, as evidenced by a single higher molecular weight band migrating at about position 62 on the gel. The template DNA in lane 3 is heterozygous at the variable site, as indicated by the presence of two distinct bands. The lower band corresponds to molecules filled in with ddGTP at position 1 complementary to the overhang (the guanine allele). The higher molecular weight band corresponds to molecules filled in with ddGTP at position 3 complementary to the overhang (the adenine allele).

**[0580]** For SNP TSC0597888, the template DNA in lanes 1 and 4 was homozygous for cytosine at the variable site; the template DNA in lane 2 was heterozygous at the variable site, and the template DNA in lane 3 was homozygous for guanine. The expected fill-in reactions are diagrammed below:

**[0581]** Homozygous for cytosine:

Allele 1	T	C	T	G	ATTC 3'
		G*	A	C	TAAG 5'
	4	3	2	1	Overhang position
Allele 2	T	C	T	G	ATTC 3'



			G*	A	C	TAAG 5'
	4	3	2	1	Overhang position	
<b>[0582]</b>	Homozygous for guanine:					
	Allele 1	T	C	T	C	ATTC 3'
					G*	TAAG 5'
	4	3	2	1	Overhang position	
	Allele 2	T	C	T	C	ATTC 3'
					G*	TAAG 5'
	4	3	2	1	Overhang position	
<b>[0583]</b>	Heterozygous for guanine/cytosine:					
	Allele 1	T	C	T	G	ATTC 3'
			G*	A	C	TAAG 5'
	4	3	2	1	Overhang position	
	Allele 2	T	C	T	C	ATTC 3'
					G*	TAAG 5'
	4	3	2	1	Overhang position	

**[0584]** Template DNA homozygous for guanine at the variable site displayed a single band, which corresponded to the DNA molecules filled in with ddGTP at position 1 complementary to the overhang. These DNA molecules were of lower molecular weight compared to the DNA molecules filled in with ddGTP at position 3 of the overhang (see lane 3 for SNP TSC0597888). The DNA molecules differed by two bases in molecular weight.

**[0585]** Template DNA homozygous for cytosine at the variable site displayed a single band, which corresponds to the DNA molecules filled in with ddGTP at position 3 complementary to the overhang. These DNA molecules migrated at a higher molecular weight than DNA molecules filled in with ddGTP at position 1 (see lanes 1 and 4 for SNP TSC0597888).

**[0586]** Template DNA heterozygous at the variable site displayed two bands; one band corresponded to the DNA molecules filled in with ddGTP at position 1 complementary to the overhang and was of lower molecular weight, and the second band corresponded to DNA molecules filled in with ddGTP at position 3 complementary to the overhang, and was of higher molecular weight (see lane 3 for SNP TSC0597888).

**[0587]** For SNP TSC0195492, the template DNA in lanes 1 and 3 was heterozygous at the variable site, which was demonstrated by the presence of two distinct bands. The template DNA in lane 2 was

homozygous for guanine at the variable site. The template DNA in lane 4 was homozygous for cytosine. Only one band was seen in lane 4 for this SNP, and it had a higher molecular weight than the DNA molecules filled in with ddGTP at position 1 complementary to the overhang (compare lanes 2, 3 and 4).

**[0588]** The observed alleles for SNP TSC0607185 are reported as cytosine or thymidine. For consistency, the SNP consortium denotes the observed alleles as they appear in the sense strand ([www.snp.cshl.org/shpsearch.shtml](http://www.snp.cshl.org/shpsearch.shtml); website active as of February 11, 2003). For this SNP, the second primer annealed upstream of the locus of interest, which allowed the fill-in reaction to occur on the antisense strand after digestion with BsmF I.

**[0589]** The template DNA in lanes 1 and 3 was heterozygous; the template DNA in lane 2 was homozygous for thymidine, and the template DNA in lane 4 was homozygous for cytosine. The antisense strand was filled in with ddGTP, so the nucleotide on the sense strand corresponded to cytosine.

**[0590]** Molecular weight markers can be used to identify the positions of the expected bands. Alternatively, for each SNP analyzed, a known heterozygous sample can be used, which will identify precisely the position of the two expected bands.

**[0591]** As demonstrated in FIG. 11, one nucleotide labeled with one fluorescent dye can be used to determine the identity of a variable site including but not limited to SNPs and single nucleotide mutations. Typically, to determine if an individual is homozygous or heterozygous at a SNP site, multiple reactions are performed using one nucleotide labeled with one dye and a second nucleotide labeled with a second dye. However, this introduces problems in comparing results because the two dyes have different quantum coefficients. Even if different nucleotides are labeled with the same dye, the quantum coefficients are different. The use of a single nucleotide labeled with one dye eliminates any errors from the quantum coefficients of different dyes.

**[0592]** In this example, fluorescently labeled ddGTP was used. However, the method is applicable for a nucleotide tagged with any signal generating moiety including but not limited to radioactive molecule, fluorescent molecule, antibody, antibody fragment, hapten, carbohydrate, biotin, derivative of biotin, phosphorescent moiety, luminescent moiety, electrochemiluminescent moiety, chromatic moiety, and moiety having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity. In addition, labeled ddATP, ddTTP, or ddCTP can be used.

**[0593]** The above example used the third position complementary to the overhang as an indicator of the second allele. However, the second or fourth position of the overhang can be used as well (see Section on Incorporation of Nucleotides). Furthermore, the overhang was generated with the type IIS

enzyme BsmF I; however any enzymes that cuts DNA at a distance from its binding site can be used including but not limited to the enzymes listed in Table I.

**[0594]** Also, in the above example, the nucleotide immediately preceding the SNP site was not a guanine on the strand that was filled in. This eliminated any effects of the alternative cutting properties of the type IIS restriction enzyme to be removed. For example, at SNP TSC0837969, the nucleotide from the SNP site on the sense strand was an adenine. If BsmF I displayed alternate cutting properties, the following overhangs would be generated for the adenine allele and the guanine allele:

G allele – 11/15 Cut	5' TTA					
	3' AAT	T	C	A	C	
Overhang position		0	1	2	3	
G allele after fill-in	5' TTA	A	G*			
	3' AAT	T	C	A	C	
Overhang position			0	1	2	3
A allele 11/15 Cut	5' TTA					
	3' AAT	T	T	A	C	
Overhang position		0	1	2	3	
A allele after fill-in	5' TTA	A	A	T	G*	
	3' AAT	T	T	A	C	
Overhang position		0	1	2	3	

**[0595]** For the guanine allele, the first position in the overhang would be filled in with dATP, which would allow the polymerase to incorporate ddGTP at position 2 complementary to the overhang. There would be no detectable difference between molecules cut at the 10/14 position or molecules cut at the 11/15 position.

**[0596]** For the adenine allele, the first position complementary to the overhang would be filled in with dATP, the second position would be filled in with dATP, the third position would be filled in with dTTP, and the fourth position would be filled in with ddGTP. There would be no difference in the molecular weights between molecules cut at 10/14 or molecules cut at 11/15. The only differences would correspond to whether the DNA molecules contained an adenine at the variable site or a guanine at the variable site.

**[0597]** As seen in FIG. 11, positioning the annealing region of the first primer allows multiple SNPs to be analyzed in a single lane of a gel. Also, when using the same nucleotide with the same dye, a single fill-in reaction can be performed. In this example, 6 SNPs were analyzed in one lane. However, any

number of SNPs including but not limited to 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 30-40, 40-50, 51-60, 61-70, 71-80, 81-100, 101-120, 121-140, 141-160, 161-180, 181-200, and greater than 200 can be analyzed in a single reaction.

**[0598]** Furthermore, one labeled nucleotide used to detect both alleles can be mixed with a second labeled nucleotide used to detect a different set of SNPs provided that neither of the nucleotides that are labeled occur immediately before the variable site (complementary to nucleotide at position 0 of the 11/15 cut). For example, suppose SNP X can be guanine or thymidine at the variable site and has the following 5' overhang generated after digestion with BsmF I:

SNP X 10/14	5' TTGAC				
G allele	3' AACTG	C	A	C	T
	Overhang position	1	2	3	4
SNP X 11/15	5' TTGA				
G allele	3' AACT	G	C	A	C
	Overhang position	0	1	2	3
SNP X 10/14	5' TTGAC				
T allele	3' AACTG	A	A	C	T
	Overhang position	1	2	3	4
SNP X 11/15	5' TTGA				
T allele	3' AACT	G	A	A	C
	Overhang position	0	1	2	3

**[0599]** After the fill-in reaction with labeled ddGTP, unlabeled dATP, dCTP, and dTTP, the following molecules would be generated:

SNP X 10/14	5' TTGAC	G*			
G allele	3' AACTG	C	A	C	T
	Overhang position	1	2	3	4
SNP X 11/15	5' TTGA	C	G*		
G allele	3' AACT	G	C	A	C
	Overhang position	0	1	2	3
SNP X 10/14	5' TTGAC	T	T	G*	
T allele	3' AACTG	A	A	C	T
	Overhang position	1	2	3	4
SNP X 11/15	5' TTGA	C	T	T	G*

T allele	3'AACT	G	A	A	C
Overhang position		0	1	2	3

**[0600]** Now suppose SNP Y can be adenine or thymidine at the variable site, and has the following 5' overhangs generated after digestion with BsmF I.

SNP Y 10/14	5' GTTT				
A allele	3' CAAA	T	G	T	A
Overhang position		1	2	3	4
SNP Y 11/15	5' GTT				
A allele	3' CAA	A	T	G	T
Overhang position		0	1	2	3

SNP Y 10/14	5' GTTT				
T allele	3' CAAA	A	G	T	A
Overhang position		1	2	3	4
SNP Y 11/15	5' GTT				
T allele	3' CAA	A	A	G	T
Overhang position		0	1	2	3

**[0601]** After fill-in with labeled ddATP and unlabeled dCTP, dGTP, and dTTP, the following molecules would be generated:

SNP Y 10/14	5' GTTT	A*			
A allele	3' CAAA	T	G	T	A
Overhang position		1	2	3	4
SNP Y 11/15	5' GTT	T	A*		
A allele	3' CAA	A	T	G	T
Overhang position		0	1	2	3
SNP Y 10/14	5' GTTT	T	C	A*	
T allele	3' CAAA	A	G	T	A
Overhang position		1	2	3	4
SNP Y 11/15	5' GTT	T	T	C	A*
T allele	3' CAA	A	A	G	T

Overhang position

0      1      2      3

**[0602]** In this example, labeled ddGTP and labeled ddATP are used to determine the identity of both alleles of SNP X and SNP Y respectively. The nucleotide immediately preceding (the complementary nucleotide to position 0 of the overhang from the 11/15 cut SNP X is not guanine or adenine on the strand that is filled-in. Likewise, the nucleotide immediately preceding SNPY is not guanine or adenine on the strand that is filled-in. This allows the fill-in reaction for both SNPs to occur in a single reaction with labeled ddGTP, labeled ddATP, and unlabeled dCTP and dTTP. This reduces the number of reactions that need to be performed and increases the number of SNPs that can be analyzed in one reaction.

**[0603]** The first primers for each SNP can be designed to anneal at different distances from the locus of interest, which allows the SNPs to migrate at different positions on the gel. For example, the first primer used to amplify SNP X can anneal at 30 bases from the locus of interest, and the first primer used to amplify SNP Y can anneal at 35 bases from the locus of interest. Also, the nucleotides can be labeled with fluorescent dyes that emit at spectrums that do not overlap. After running the gel, the gel can be scanned at one wavelength specific for one dye. Only those molecules labeled with that dye will emit a signal. The gel then can be scanned at the wavelength for the second dye. Only those molecules labeled with that dye will emit a signal. This method allows maximum compression for the number of SNPs that can be analyzed in a single reaction.

**[0604]** In this example, the nucleotide preceding the variable site on the strand that was filled-in was not adenine or guanine, and the nucleotide following the variable site can not be adenine or guanine on the sense strand. This method can work with any combination of labeled nucleotides, and the skilled artisan would understand which labeling reactions can be mixed and those that can not. For instance, if one SNP is labeled with thymidine and a second SNP is labeled with cytosine, the SNPs can be labeled in a single reaction if the nucleotide immediately preceding each variable site is not thymidine or cytosine on the sense strand and the nucleotide immediately after the variable site is not thymidine or cytosine on the sense strand.

**[0605]** This method allows the signals from one allele to be compared to the signal from a second allele without the added complexity of determining the degree of alternate cutting, or having to correct for the quantum coefficients of the dyes. This method is especially useful when trying to quantitate a ratio for one allele to another. For example, this method is useful for detecting chromosomal abnormalities. The ratio of alleles at a heterozygous site is expected to be about 1:1 (one A allele and one G allele). However, if an extra chromosome is present the ratio is expected to be about 1:2 (one A allele and 2 G

alleles or 2 A alleles and 1 G allele). This method is especially useful when trying to detect fetal DNA in the presence of maternal DNA.

**[0606]** In addition, this method is useful for detecting two genetic signals in one sample. For example, this method can detect mutant cells in the presence of wild type cells (see Example 5). If a mutant cell contains a mutation in the DNA sequence of a particular gene, this method can be used to detect both the mutant signal and the wild type signal. This method can be used to detect the mutant DNA sequence in the presence of the wild type DNA sequence. The ratio of mutant DNA to wild type DNA can be quantitated because a single nucleotide labeled with one signal generating moiety is used.

### EXAMPLE 7

**[0607]** Non-invasive methods for the detection of various types of cancer have the potential to reduce morbidity and mortality from the disease. Several techniques for the early detection of colorectal tumors have been developed including colonoscopy, barium enemas, and sigmoidoscopy; however the techniques are limited in use because they are invasive, which causes a low rate of patient compliance. Non-invasive genetic tests may be useful in identifying early stage colorectal tumors.

**[0608]** In 1991, researchers identified the Adenomatous Polyposis Coli gene (APC), which plays a critical role in the formation of colorectal tumors (Kinzler *et al.*, Science 253:661-665, 1991). The APC gene resides on chromosome 5q21-22 and a total of 15 exons code for an RNA molecule of 8529 nucleotides, which produces a 300 Kd APC protein. The protein is expressed in numerous cell types and is essential for cell adhesion.

**[0609]** Mutations in the APC gene generally initiate colorectal neoplasia (Tsao, J. *et al.*, Am, J. Pathol. 145:531-534, 1994). Approximately 95% of the mutations in the APC gene result in nonsense/frameshift mutations. The most common mutations occur at codons 1061 and 1309; mutations at these codons account for 1/3 of all germline mutations. With regard to somatic mutations, 60% occur within codons 1286-1513, which is about 10% of the coding sequence. This region is termed the mutation Cluster Region (MCR). Numerous types of mutations have been identified in the APC gene including nucleotide substitutions (see Table VII), splicing errors (see Table VIII), small deletions (see Table IX), small insertions (see Table X), small insertions/deletions (see Table XI), gross deletions (see Table XII), gross insertions (see Table XIII), and complex rearrangements (see Table XIV).

**[0610]** Researchers have attempted to identify cells harboring mutations in the APC gene in stool samples (Traverso, G. *et al.*, New England Journal of Medicine, Vol 346:311-320, 2002). While APC mutations are found in nearly all tumors, about 1 in 250 cells in the stool sample has a mutation in the

APC gene; most of the cells are normal cells that have been shed into the feces. Furthermore, human DNA represents about one-billionth of the total DNA found in stool samples; the majority of DNA is bacterial. The technique employed by Traverso *et al.* only detects mutations that result in a truncated protein.

**[0611]** As discussed above, numerous mutations in the APC gene have been implicated in the formation of colorectal tumors. Thus, a need still exists for a highly sensitive, non-invasive technique for the detection of colorectal tumors. Below, methods are described for detection of two mutations in the APC gene. However, any number of mutations can be analyzed using the methods described herein.

### **Preparation of Template DNA**

**[0612]** The template DNA is purified from a sample containing colon cells including but not limited to a stool sample. The template DNA is purified using the procedures described by Ahlquist *et al.* (Gastroenterology, 119:1219-1227, 2000). If stool samples are frozen, the samples are thawed at room temperature, and homogenized with an Exactor stool shaker (Exact Laboratories, Maynard, Mass.) Following homogenization, a 4 gram stool equivalent of each sample is centrifuged at 2536 x g for 5 minutes. The samples are centrifuged a second time at 16,500 x g for 10 minutes. Supernatants are incubated with 20 µl of RNase (0.5 mg per milliliter) for 1 hour at 37°C. DNA is precipitated with 1/10 volume of 3 mol of sodium acetate per liter and an equal volume of isopropanol. The DNA is dissolved in 5 ml of TRIS-EDTA (0.01 mol of Tris per liter (pH 7.4) and 0.001 mole of EDTA per liter).

### **Design of Primers**

**[0613]** To determine if a mutation resides at codon 1370, the following primers are used:

First primer:

5' GTGCAAAGGCCTGAATTCCCAGGCACAAAGCTGTTGAA 3' (SEQ ID NO: 42)

Second primer:

5' TGAAGCGAACTAGGGACTCAGGTGGACTT (SEQ ID NO: 43)

**[0614]** The first primer contains a biotin tag at the extreme 5' end, and the nucleotide sequence for the restriction enzyme EcoRI. The second primer contains the nucleotide sequence for the restriction enzyme BsmF I.

**[0615]** To determine if a small deletion exists at codon 1302, the following primers are used:

First primer:

5' GATTCCGTAAACGAATTCAGTTCATTATCATCTTTGTC 3' (SEQ ID NO: 44)

Second primer:



5' CCATTGTTAAGCGGGACTTCTGCTATTTG 3' (SEQ ID NO: 45)

**[0616]** The first primer has a biotin tag at the 5' end and contains a restriction enzyme recognition site for EcoRI. The second primer contains a restriction enzyme recognition site for BsmF.I.

#### **PCR Reaction**

**[0617]** The loci of interest are amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). The loci of interest are amplified in separate reaction tubes; they can also be amplified together in a single PCR reaction. For increased specificity, a "hot-start" PCR reaction is used, e.g. by using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction are optimized for each locus of interest but in this example, 40 ng of template human genomic DNA and 5  $\mu$ M of each primer are used. Forty cycles of PCR are performed. The following PCR conditions are used:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 64°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty nine (39) times;
- (9) 72°C for 5 minutes.

**[0618]** In the first cycle of PCR, the annealing temperature is about the melting temperature of the 3' annealing region of the second primers, which is 37°C. The annealing temperature in the second cycle of PCR is about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer, which is 57°C. The annealing temperature in the third cycle of PCR is about the melting temperature of the entire sequence of the second primer, which is 64°C. The annealing temperature for the remaining cycles is 64°C. Escalating the annealing temperature from TM1 to TM2 to TM3 in the first three cycles of PCR greatly improves specificity. These annealing temperatures are representative, and the skilled artisan understands that the annealing temperatures for each cycle are dependent on the specific primers used.

**[0619]** The temperatures and times for denaturing, annealing, and extension, are optimized by trying various settings and using the parameters that yield the best results.

#### **Purification of Fragment of Interest**

**[0620]** The PCR products are separated from the genomic template DNA. Each PCR product is divided into four separate reaction wells of a Streptawell, transparent, High-Bind plate from Roche Diagnostics GmbH (catalog number 1 645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog). The first primers contain a 5' biotin tag so the PCR products bound to the Streptavidin coated wells while the genomic template DNA does not. The streptavidin binding reaction is performed using a Thermomixer (Eppendorf) at 1000 rpm for 20 min. at 37°C. Each well is aspirated to remove unbound material, and washed three times with 1X PBS, with gentle mixing (Kandpal et al., Nucl. Acids Res. 18:1789-1795 (1990); Kaneoka et al., Biotechniques 10:30-34 (1991); Green et al., Nucl. Acids Res. 18:6163-6164 (1990)).

**[0621]** Alternatively, the PCR products are placed into a single well of a streptavidin plate to perform the nucleotide incorporation reaction in a single well.

#### **Restriction Enzyme Digestion of Isolated Fragments**

**[0622]** The purified PCR products are digested with the restriction enzyme BsmF I (New England Biolabs catalog number R0572S), which binds to the recognition site incorporated into the PCR products from the second primer. The digests are performed in the Streptawells following the instructions supplied with the restriction enzyme. After digestion with the appropriate restriction enzyme, the wells are washed three times with PBS to remove the cleaved fragments.

#### **Incorporation of Labeled Nucleotide**

**[0623]** The restriction enzyme digest described above yields a DNA fragment with a 5' overhang, which contains the locus of interest and a 3' recessed end. The 5' overhang functions as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase.

**[0624]** For each locus of interest, four separate fill in reactions are performed; each of the four reactions contains a different fluorescently labeled ddNTP (ddATP, ddTTP, ddGTP, or ddCTP). The following components are added to each fill in reaction: 1 µl of a fluorescently labeled ddNTP, 0.5 µl of unlabeled ddNTPs (40 µM), which contains all nucleotides except the nucleotide that is fluorescently labeled, 2 µl of 10X sequenase buffer, 0.25 µl of Sequenase, and water as needed for a 20µl reaction. The fill are performed in reactions at 40°C for 10 min. Non-fluorescently labeled ddNTP are purchased from

Fermentas Inc. (Hanover, MD). All other labeling reagents are obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565). In the presence of fluorescently labeled ddNTPs, the 3' recessed end is extended by one base, which corresponds to the locus of interest.

**[0625]** A mixture of labeled ddNTPs and unlabeled dNTPs also can be used for the fill-in reaction. The "fill in" conditions are as described above except that a mixture containing 40  $\mu$ M unlabeled dNTPs, 1  $\mu$ l fluorescently labeled ddATP, 1  $\mu$ l fluorescently labeled ddTTP, 1  $\mu$ l fluorescently labeled ddCTP, and 1  $\mu$ l ddGTP are used. The fluorescent ddNTPs are obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565; Amersham does not publish the concentrations of the fluorescent nucleotides). The locus of interest is digested with the restriction enzyme BsmF I, which generates a 5' overhang of four bases. If the first nucleotide incorporated is a labeled ddNTP, the 3' recessed end is filled in by one base, allowing detection of the locus of interest. However, if the first nucleotide incorporated is a dNTP, the polymerase continues to incorporate nucleotides until a ddNTP is filled in. For example, the first two nucleotides may be filled in with dNTPs, and the third nucleotide with a ddNTP, allowing detection of the third nucleotide in the overhang. Thus, the sequence of the entire 5' overhang is determined, which increases the information obtained from each SNP or locus of interest. This type of fill in reaction is especially useful when detecting the presence of insertions, deletions, insertions and deletions, rearrangements, and translocations.

**[0626]** Alternatively, one nucleotide labeled with a single dye is used to determine the sequence of the locus of interest. See Example 6. This method eliminates any potential errors when using different dyes, which have different quantum coefficients.

**[0627]** After labeling, each Streptawell is rinsed with 1X PBS (100  $\mu$ l) three times. The "filled in" DNA fragments are released from the Streptawells by digesting with the restriction enzyme EcoRI, according to the manufacturer's instructions that are supplied with the enzyme. The digestion is performed for 1 hour at 37 °C with shaking at 120 rpm.

#### **Detection of the Locus of Interest**

**[0628]** After release from the streptavidin matrix, the sample is loaded into a lane of a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691). The sample is electrophoresed into the gel at 3000 volts for 3 min. The gel is run for 3 hours using a sequencing apparatus (Hoefer SQ3 Sequencer). The incorporated labeled nucleotide is detected by fluorescence.

**[0629]** To determine if any cells contain mutations at codon 1370 of the APC gene when separate fill-in reactions are performed, the lanes of the gel that correspond to the fill-in reaction for ddATP and ddTTP are analyzed. If only normal cells are present, the lane corresponding to the fill in reaction with ddATP is a bright signal. No signal is detected for the “fill-in” reaction with ddTTP. However, if the patient sample contains cells with mutations at codon 1370 of the APC gene, the lane corresponding to the fill in reaction with ddATP is a bright signal, and a signal is detected from the lane corresponding to the fill in reaction with ddTTP. The intensity of the signal from the lane corresponding to the fill in reaction with ddTTP is indicative of the number of mutant cells in the sample.

**[0630]** Alternatively, one labeled nucleotide is used to determine the sequence of the alleles at codon 1370 of the APC gene. At codon 1370, the normal sequence is AAA, which codes for the amino acid lysine. However, a nucleotide substitution has been identified at codon 1370, which is associated with colorectal tumors. Specifically, a change from A to T (AAA-TAA) typically is found at codon 1370, which results in a stop codon. A single fill-in reaction is performed using labeled ddATP, and unlabeled dTTP, dCTP, and dGTP. A single nucleotide labeled with one fluorescent dye is used to determine the presence of both the normal and mutant DNA sequence that codes for codon 1370. The relevant DNA sequence is depicted below with the sequence corresponding to codon 1370 in bold:

5' CCC**AAA**AGTCCACCTGA (SEQ ID NO: 46)

3' GGGTTTTCAGGTGGACT (SEQ ID NO: 47)

**[0631]** After digest with BsmF I, the following overhang is produced:

5' CCC

3' GGG      T      T      T      T

Overhang position	1	2	3	4
-------------------	---	---	---	---

**[0632]** If the patient sample has no cells harboring a mutation at codon 1370, one signal is seen corresponding to incorporation of labeled ddATP.

5' CCC      **A\***

3' GGG      T      T      T      T

Overhang position	1	2	3	4
-------------------	---	---	---	---

**[0633]** However, if the patient sample has cells with mutations at codon 1370 of the APC gene, one signal is seen, which corresponds to the normal sequence at codon 1370, and a second signal is seen, which corresponds to the mutant sequence at codon 1370. The signals clearly are identified as they differ in molecular weight.

Overhang of normal DNA sequence:    CCC

	GGG	T	T	T	T
Overhang position		1	2	3	4
Normal DNA sequence after fill-in:	CCC	<b>A*</b>			
	GGG	T	T	T	T
Overhang position		1	2	3	4
Overhang of mutant DNA sequence:	CCC				
	GGG	A	T	T	T
Overhang position		1	2	3	4
Mutant DNA sequence after fill-in:	CCC	<b>T</b>	<b>A*</b>		
	GGG	A	T	T	T
Overhang position		1	2	3	4

**[0634]** Two signals are seen when the mutant allele is present. The mutant DNA molecules are filled in one base after the wild type DNA molecules. The two signals are separated using any method that discriminates based on molecular weight. One labeled nucleotide (ddATP) is used to detect the presence of both the wild type DNA sequence and the mutant DNA sequence. This method of labeling reduces the number of reactions that need to be performed and allows accurate quantitation for the number of mutant cells in the patient sample. The number of mutant cells in the sample is used to determine patient prognosis, the degree and the severity of the disease. This method of labeling eliminates the complications associated with using different dyes, which have distinct quantum coefficients. This method of labeling also eliminates errors associated with pipetting reactions.

**[0635]** To determine if any cells contain mutations at codon 1302 of the APC gene when separate fill-in reactions are performed, the lanes of the gel that correspond to the fill-in reaction for ddTTP and ddCTP are analyzed. The normal DNA sequence is depicted below with sequence coding for codon 1302 in bold type-face.

Normal Sequence: 5' ACCCT**GC**AAATAGCAGAA (SEQ ID NO: 48)  
3' TGGGACGTTTATCGTCTT (SEQ ID NO: 49)

**[0636]** After digest, the following 5' overhang is produced:

5' ACCC  
3' TGGG A C G T  
Overhang position 1 2 3 4

**[0637]** After the fill-in reaction, labeled ddTTP is incorporated.

5' ACCC **T\***

3' TGGG	A	C	G	T
Overhang position	1	2	3	4

**[0638]** A deletion of a single base of the APC sequence, which typically codes for codon 1302, has been associated with colorectal tumors. The mutant DNA sequence is depicted below with the relevant sequence in bold:

Mutant Sequence: 5' ACCCGCAAATAGCAGAA (SEQ ID NO: 50)  
 3' TGGGCGTTTATCGTCTT (SEQ ID NO: 51)

**[0639]** After digest:

5' ACC				
3' TGG	G	C	G	T
Overhang position	1	2	3	4

**[0640]** After fill-in:

5' ACC	C*			
3' TGG	G	C	G	T
Overhang position	1	2	3	4

**[0641]** If there are no mutations in the APC gene, signal is not detected for the fill in reaction with ddCTP\*, but a bright signal is detected for the fill-in reaction with ddTTP\*. However, if there are cells in the patient sample that have mutations in the APC gene, signals are seen for the fill-in reactions with ddCTP\* and ddTTP\*.

**[0642]** Alternatively, a single fill-in reaction is performed using a mixture containing unlabeled dNTPs, fluorescently labeled ddATP, fluorescently labeled ddTTP, fluorescently labeled ddCTP, and fluorescently labeled ddGTP. If there is no deletion, labeled ddTTP is incorporated.

5' ACCC	T*			
3' TGGG	A	C	G	T
Overhang position	1	2	3	4

**[0643]** However, if the T has been deleted, labeled ddCTP\* is incorporated.

5' ACCC*				
3' TGGG	C	G	T	
Overhang position	1	2	3	4

**[0644]** The two signals are separated by molecular weight because of the deletion of the thymidine nucleotide. If mutant cells are present, two signals are generated in the same lane but are separated by a single base pair (this principle is demonstrated in FIG 9D). The deletion causes a change in the molecular

weight of the DNA fragments, which allows a single fill in reaction to be used to detect the presence of both normal and mutant cells.

**[0645]** In the above example, methods for the detection of a nucleotide substitution and a small deletion are described. However, the methods can be used for the detection of any type of mutation including but not limited to nucleotide substitutions (see Table VII), splicing errors (see Table VIII), small deletions (see Table IX), small insertions (see Table X), small insertions/deletions (see Table XI), gross deletions (see Table XII), gross insertions (see Table XIII), and complex rearrangements (see Table XIV).

**[0646]** In addition, the above-described methods are used for the detection of any type of disease including but not limited to those listed in Table IV. Furthermore, any type of mutant gene is detected using the inventions described herein including but not limited to the genes associated with the diseases listed in Table IV, BRCA1, BRCA2, MSH6, MSH2, MLH1, RET, PTEN, ATM, H-RAS, p53, ELAC2, CDH1, APC, AR, PMS2, MLH3, CYP1A1, GSTP1, GSTM1, AXIN2, CYP19, MET, NAT1, CDKN2A, NQO1, trc8, RAD51, PMS1, TGFBR2, VHL, MC4R, POMC, NROB2, UCP2, PCSK1, PPARG, ADRB2, UCP3, glur1, cart, SORBS1, LEP, LEPR, SIM1, TNF, IL-6, IL-1, IL-2, IL-3, IL1A, TAP2, THPO, THRB, NBS1, RBM15, LIF, MPL, RUNX1, Her-2, glucocorticoid receptor, estrogen receptor, thyroid receptor, p21, p27, K-RAS, N-RAS, retinoblastoma protein, Wiskott-Aldrich (WAS) gene, Factor V Leiden, Factor II (prothrombin), methylene tetrahydrofolate reductase, cystic fibrosis, LDL receptor, HDL receptor, superoxide dismutase gene, SHOX gene, genes involved in nitric oxide regulation, genes involved in cell cycle regulation, tumor suppressor genes, oncogenes, genes associated with neurodegeneration, genes associated with obesity[[, ]]. Abbreviations correspond to the proteins as listed on the Human Gene Mutation Database, which is incorporated herein by reference [www.archive.uwcm.ac.uk/uwcm](http://www.archive.uwcm.ac.uk/uwcm); website address active as of February 12, 2003).

**[0647]** The above-example demonstrates the detection of mutant cells and mutant alleles from a fecal sample. However, the methods described herein are used for detection of mutant cells from any biological sample including but not limited to blood sample, serum sample, plasma sample, urine sample, spinal fluid, lymphatic fluid, semen, vaginal secretion, ascitic fluid, saliva, mucosa secretion, peritoneal fluid, fecal sample, body exudates, breast fluid, lung aspirates, cells, tissues, individual cells or extracts of the such sources that contain the nucleic acid of the same, and subcellular structures such as mitochondria or chloroplasts. In addition, the methods described herein are used for the detection of mutant cells and mutated DNA from any number of nucleic acid containing sources including but not limited to forensic, food, archeological, agricultural or inorganic samples.

**[0648]** The above example is directed to detection of mutations in the APC gene. However, the inventions described herein are used for the detection of mutations in any gene that is associated with or predisposes to disease (see Table XV).

**[0649]** For example, hypermethylation of the glutathione S-transferase P1 (GSTP1) promoter is the most common DNA alteration in prostate cancer. The methylation state of the promoter is determined using sodium bisulfite and the methods described herein.

**[0650]** Treatment with sodium bisulfite converts unmethylated cytosine residues into uracil, and leaving the methylated cytosines unchanged. Using the methods described herein, a first and second primer are designed to amplify the regions of the GSTP1 promoter that are often methylated. Below, a region of the GSTP1 promoter is shown prior to sodium bisulfite treatment:

Before Sodium Bisulfite treatment:

5' ACCGCTACA  
3' TGGCGATCA

**[0651]** Below, a region of the GSTP1 promoter is shown after sodium bisulfite treatment, PCR amplification, and digestion with the type IIS restriction enzyme BsmF I:

Unmethylated

5' ACC				
3' TGG	U	G	A	T
Overhang position	1	2	3	4

Methylated

5' ACC				
3' TGG	C	G	A	T
Overhang position	1	2	3	4

**[0652]** Labeled ddATP, unlabeled dCTP, dGTP, and dTTP are used to fill-in the 5' overhangs. The following molecules are generated:

Unmethylated

5' ACC	A*			
3' TGG	U	G	A	T
Overhang position	1	2	3	4

Methylated

5' ACC	G	C	T	A*
3' TGG	C	G	A	T



Overhang position                      1                      2                      3                      4

**[0653]** Two signals are seen; one corresponds to DNA molecules filled in with ddATP at position one complementary to the overhang (unmethylated), and the other corresponds to the DNA molecules filled in with ddATP at position 4 complementary to the overhang (methylated). The two signals are separated based on molecular weight. Alternatively, the fill-in reactions are performed in separate reactions using labeled ddGTP in one reaction and labeled ddATP in another reaction.

**[0654]** The methods described herein are used to screen for prostate cancer and also to monitor the progression and severity of the disease. The use of a single nucleotide to detect both the methylated and unmethylated sequences allows accurate quantitation and provides a high level of sensitivity for the methylated sequences, which is a useful tool for earlier detection of the disease.

**[0655]** The information contained in Tables VII-XIV was obtained from the Human Gene Mutation Database. With the information provided herein, the skilled artisan will understand how to apply these methods for determining the sequence of the alleles for any gene. A large number of genes and there associated mutations can be found at the following website: [www.archive.uwcm.ac.uk/uwcm](http://www.archive.uwcm.ac.uk/uwcm).

**TABLE VII: NUCLEOTIDE SUBSTITUTIONS**

Codon	Nucleotide	Amino acid	Phenotype
99	CGG-TGG	Arg-Trp	Adenomatous polyposis coli
121	AGA-TGA	Arg-Term	Adenomatous polyposis coli
157	TGG-TAG	Trp-Term	Adenomatous polyposis coli
159	TAC-TAG	Tyr-Term	Adenomatous polyposis coli
163	CAG-TAG	Gln-Term	Adenomatous polyposis coli
168	AGA-TGA	Arg-Term	Adenomatous polyposis coli
171	AGT-ATT	Ser-Ile	Adenomatous polyposis coli
181	CAA-TAA	Gln-Term	Adenomatous polyposis coli
190	GAA-TAA	Glu-Term	Adenomatous polyposis coli
202	GAA-TAA	Glu-Term	Adenomatous polyposis coli
208	CAG-CGG	Gln-Arg	Adenomatous polyposis coli
208	CAG-TAG	Gln-Term	Adenomatous polyposis coli
213	CGA-TGA	Arg-Term	Adenomatous polyposis coli
215	CAG-TAG	Gln-Term	Adenomatous polyposis coli
216	CGA-TGA	Arg-Term	Adenomatous polyposis coli
232	CGA-TGA	Arg-Term	Adenomatous polyposis coli
233	CAG-TAG	Gln-Term	Adenomatous polyposis coli
247	CAG-TAG	Gln-Term	Adenomatous polyposis coli
267	GGA-TGA	Gly-Term	Adenomatous polyposis coli
278	CAG-TAG	Gln-Term	Adenomatous polyposis coli
280	TCA-TGA	Ser-Term	Adenomatous polyposis coli
280	TCA-TAA	Ser-Term	Adenomatous polyposis coli

283	CGA-TGA	Arg-Term	Adenomatous polyposis coli
302	CGA-TGA	Arg-Term	Adenomatous polyposis coli
332	CGA-TGA	Arg-Term	Adenomatous polyposis coli
358	CAG-TAG	Gln-Term	Adenomatous polyposis coli
405	CGA-TGA	Arg-Term	Adenomatous polyposis coli
414	CGC-TGC	Arg-Cys	Adenomatous polyposis coli
422	GAG-TAG	Glu-Term	Adenomatous polyposis coli
423	TGG-TAG	Trp-Term	Adenomatous polyposis coli
424	CAG-TAG	Gln-Term	Adenomatous polyposis coli
433	CAG-TAG	Gln-Term	Adenomatous polyposis coli
443	GAA-TAA	Glu-Term	Adenomatous polyposis coli
457	TCA-TAA	Ser-Term	Adenomatous polyposis coli
473	CAG-TAG	Gln-Term	Adenomatous polyposis coli
486	TAC-TAG	Tyr-Term	Adenomatous polyposis coli
499	CGA-TGA	Arg-Term	Adenomatous polyposis coli
500	TAT-TAG	Tyr-Term	Adenomatous polyposis coli
541	CAG-TAG	Gln-Term	Adenomatous polyposis coli
553	TGG-TAG	Trp-Term	Adenomatous polyposis coli
554	CGA-TGA	Arg-Term	Adenomatous polyposis coli
564	CGA-TGA	Arg-Term	Adenomatous polyposis coli
577	TTA-TAA	Leu-Term	Adenomatous polyposis coli
586	AAA-TAA	Lys-Term	Adenomatous polyposis coli
592	TTA-TGA	Leu-Term	Adenomatous polyposis coli
593	TGG-TAG	Trp-Term	Adenomatous polyposis coli
593	TGG-TGA	Trp-Term	Adenomatous polyposis coli
622	TAC-TAA	Tyr-Term	Adenomatous polyposis coli
625	CAG-TAG	Gln-Term	Adenomatous polyposis coli
629	TTA-TAA	Leu-Term	Adenomatous polyposis coli
650	GAG-TAG	Glu-Term	Adenomatous polyposis coli
684	TTG-TAG	Leu-Term	Adenomatous polyposis coli
685	TGG-TGA	Trp-Term	Adenomatous polyposis coli
695	CAG-TAG	Gln-Term	Adenomatous polyposis coli
699	TGG-TGA	Trp-Term	Adenomatous polyposis coli
699	TGG-TAG	Trp-Term	Adenomatous polyposis coli
713	TCA-TGA	Ser-Term	Adenomatous polyposis coli
722	AGT-GGT	Ser-Gly	Adenomatous polyposis coli
747	TCA-TGA	Ser-Term	Adenomatous polyposis coli
764	TTA-TAA	Leu-Term	Adenomatous polyposis coli
784	TCT-ACT	Ser-Thr	Adenomatous polyposis coli
805	CGA-TGA	Arg-Term	Adenomatous polyposis coli
811	TCA-TGA	Ser-Term	Adenomatous polyposis coli
848	AAA-TAA	Lys-Term	Adenomatous polyposis coli
876	CGA-TGA	Arg-Term	Adenomatous polyposis coli
879	CAG-TAG	Gln-Term	Adenomatous polyposis coli
893	GAA-TAA	Glu-Term	Adenomatous polyposis coli
932	TCA-TAA	Ser-Term	Adenomatous polyposis coli

932	TCA-TGA	Ser-Term	Adenomatous polyposis coli
935	TAC-TAG	Tyr-Term	Adenomatous polyposis coli
935	TAC-TAA	Tyr-Term	Adenomatous polyposis coli
995	TGC-TGA	Cys-Term	Adenomatous polyposis coli
997	TAT-TAG	Tyr-Term	Adenomatous polyposis coli
999	CAA-TAA	Gln-Term	Adenomatous polyposis coli
1000	TAC-TAA	Tyr-Term	Adenomatous polyposis coli
1020	GAA-TAA	Glu-Term	Adenomatous polyposis coli
1032	TCA-TAA	Ser-Term	Adenomatous polyposis coli
1041	CAA-TAA	Gln-Term	Adenomatous polyposis coli
1044	TCA-TAA	Ser-Term	Adenomatous polyposis coli
1045	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1049	TGG-TGA	Trp-Term	Adenomatous polyposis coli
1067	CAA-TAA	Gln-Term	Adenomatous polyposis coli
1071	CAA-TAA	Gln-Term	Adenomatous polyposis coli
1075	TAT-TAA	Tyr-Term	Adenomatous polyposis coli
1075	TAT-TAG	Tyr-Term	Adenomatous polyposis coli
1102	TAC-TAG	Tyr-Term	Adenomatous polyposis coli
1110	TCA-TGA	Ser-Term	Adenomatous polyposis coli
1114	CGA-TGA	Arg-Term	Adenomatous polyposis coli
1123	CAA-TAA	Gln-Term	Adenomatous polyposis coli
1135	TAT-TAG	Tyr-Term	Adenomatous polyposis coli
1152	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1155	GAA-TAA	Glu-Term	Adenomatous polyposis coli
1168	GAA-TAA	Glu-Term	Adenomatous polyposis coli
1175	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1176	CCT-CTT	Pro-Leu	Adenomatous polyposis coli
1184	GCC-CCC	Ala-Pro	Adenomatous polyposis coli
1193	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1194	TCA-TGA	Ser-Term	Adenomatous polyposis coli
1198	TCA-TGA	Ser-Term	Adenomatous polyposis coli
1201	TCA-TGA	Ser-Term	Adenomatous polyposis coli
1228	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1230	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1244	CAA-TAA	Gln-Term	Adenomatous polyposis coli
1249	TGC-TGA	Cys-Term	Adenomatous polyposis coli
1256	CAA-TAA	Gln-Term	Adenomatous polyposis coli
1262	TAT-TAA	Tyr-Term	Adenomatous polyposis coli
1270	TGT-TGA	Cys-Term	Adenomatous polyposis coli
1276	TCA-TGA	Ser-Term	Adenomatous polyposis coli
1278	TCA-TAA	Ser-Term	Adenomatous polyposis coli
1286	GAA-TAA	Glu-Term	Adenomatous polyposis coli
1289	TGT-TGA	Cys-Term	Adenomatous polyposis coli
1294	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1307	ATA-AAA	Ile-Lys	Colorectal cancer, predisposition to, association
1309	GAA-TAA	Glu-Term	Adenomatous polyposis coli

1317	GAA-CAA	Glu-Gln	Colorectal cancer, predisposition to
1328	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1338	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1342	TTA-TAA	Leu-Term	Adenomatous polyposis coli
1342	TTA-TGA	Leu-Term	Adenomatous polyposis coli
1348	AGG-TGG	Arg-Trp	Adenomatous polyposis coli
1357	GGA-TGA	Gly-Term	Adenomatous polyposis coli
1367	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1370	AAA-TAA	Lys-Term	Adenomatous polyposis coli
1392	TCA-TAA	Ser-Term	Adenomatous polyposis coli
1392	TCA-TGA	Ser-Term	Adenomatous polyposis coli
1397	GAG-TAG	Glu-Term	Adenomatous polyposis coli
1449	AAG-TAG	Lys-Term	Adenomatous polyposis coli
1450	CGA-TGA	Arg-Term	Adenomatous polyposis coli
1451	GAA-TAA	Glu-Term	Adenomatous polyposis coli
1503	TCA-TAA	Ser-Term	Adenomatous polyposis coli
1517	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1529	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1539	TCA-TAA	Ser-Term	Adenomatous polyposis coli
1541	CAG-TAG	Gln-Term	Adenomatous polyposis coli
1564	TTA-TAA	Leu-Term	Adenomatous polyposis coli
1567	TCA-TGA	Ser-Term	Adenomatous polyposis coli
1640	CGG-TGG	Arg-Trp	Adenomatous polyposis coli
1693	GAA-TAA	Glu-Term	Adenomatous polyposis coli
1822	GAC-GTC	Asp-Val	Adenomatous polyposis coli, association with ?
2038	CTG-GTG	Leu-Val	Adenomatous polyposis coli
2040	CAG-TAG	Gln-Term	Adenomatous polyposis coli
2566	AGA-AAA	Arg-Lys	Adenomatous polyposis coli
2621	TCT-TGT	Ser-Cys	Adenomatous polyposis coli
2839	CTT-TTT	Leu-Phe	Adenomatous polyposis coli

TABLE VIII: NUCLEOTIDE SUBSTITUTIONS

Donor/Acceptor	Relative location	Substitution	Phenotype
ds	-1	G-C	Adenomatous polyposis coli
as	-1	G-A	Adenomatous polyposis coli
as	-1	G-C	Adenomatous polyposis coli
ds	+2	T-A	Adenomatous polyposis coli
as	-1	G-C	Adenomatous polyposis coli
as	-1	G-T	Adenomatous polyposis coli
as	-1	G-A	Adenomatous polyposis coli
as	-2	A-C	Adenomatous polyposis coli
as	-5	A-G	Adenomatous polyposis coli
ds	+3	A-C	Adenomatous polyposis coli

as	-1	G-A	Adenomatous polyposis coli
ds	+1	G-A	Adenomatous polyposis coli
as	-1	G-T	Adenomatous polyposis coli
ds	+1	G-A	Adenomatous polyposis coli
as	-1	G-A	Adenomatous polyposis coli
ds	+1	G-A	Adenomatous polyposis coli
ds	+3	A-G	Adenomatous polyposis coli
ds	+5	G-T	Adenomatous polyposis coli
as	-1	G-A	Adenomatous polyposis coli
as	-6	A-G	Adenomatous polyposis coli
as	-5	A-G	Adenomatous polyposis coli
as	-2	A-G	Adenomatous polyposis coli
ds	+2	T-C	Adenomatous polyposis coli
as	-2	A-G	Adenomatous polyposis coli
ds	+1	G-A	Adenomatous polyposis coli
ds	+1	G-T	Adenomatous polyposis coli
ds	+2	T-G	Adenomatous polyposis coli

**TABLE IX: APC SMALL DELETIONS**

**[0656]** Bold letters indicate the codon. Undercase letters represent the deletion. Where deletions extend beyond the coding region, other positional information is provided. For example, the abbreviation 5' UTR represents 5' untranslated region, and the abbreviation E6I6 denotes exon 6/intron 6 boundary.

Location/ codon	Deletion	Phenotype	
77	TTAgataGCAGTAATTT	Adenomatous polyposis coli	<u>SEQ ID NO: 52</u>
97	GGAAGccgggaagGATCTGTATC	Adenomatous polyposis coli	<u>SEQ ID NO: 53</u>
138	GAGAAAGAGAG_E3I3_GTAA	Adenomatous polyposis coli	<u>SEQ ID NO: 54</u>
139	AAAGAgag_E3I3_Gtaactttct	Thyroid cancer	<u>SEQ ID NO: 55</u>
139	AAAGagag_E3I3_GTAACTTTTC	Adenomatous polyposis coli	<u>SEQ ID NO: 56</u>
142	TTTTAAAAAAaAAAAATAG_I3E4_GTCA	Adenomatous polyposis coli	<u>SEQ ID NO: 57</u>
144	AAAATAG_I3E4_GTCatTGCTTCTTGC	Adenomatous polyposis coli	<u>SEQ ID NO: 58</u>
149	GACAaaGAAGAAAAGG	Adenomatous polyposis coli	<u>SEQ ID NO: 59</u>
149	GACAAagaaGAAAAGGAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 60</u>
155	AGGAA^AAAGActggtATTACGCTCA	Adenomatous polyposis coli	<u>SEQ ID NO: 61</u>

169	AAAAGA^ATAGatagTCTTCCTTTA	Adenomatous polyposis coli	<u>SEQ ID NO: 62</u>
172	AGATAGT^CTTcCTTTAACTGA	Adenomatous polyposis coli	<u>SEQ ID NO: 63</u>
179	TCCTTacaaACAGATATGA	Adenomatous polyposis coli	<u>SEQ ID NO: 64</u>
185	ACCaGAAGGCAATT	Adenomatous polyposis coli	<u>SEQ ID NO: 65</u>
196	ATCAGagTTGCGATGGA	Adenomatous polyposis coli	<u>SEQ ID NO: 66</u>
213	CGAGCaCAG_E5I5_GTAAGTT	Adenomatous polyposis coli	<u>SEQ ID NO: 67</u>
298	CACtcTGCACCTCGA	Adenomatous polyposis coli	<u>SEQ ID NO: 68</u>
329	GATaTGTCGCGAAC	Adenomatous polyposis coli	<u>SEQ ID NO: 69</u>
365	AAAGActCTGTATTGTT	Adenomatous polyposis coli	<u>SEQ ID NO: 70</u>
397	GACaaGAGAGGCAGG	Adenomatous polyposis coli	<u>SEQ ID NO: 71</u>
427	CATGAacCAGGCATGGA	Adenomatous polyposis coli	<u>SEQ ID NO: 72</u>
428	GAACCaGGCATGGACC	Adenomatous polyposis coli	<u>SEQ ID NO: 73</u>
436	AATCCaa_E9I9_gTATGTTCTCT	Adenomatous polyposis coli	<u>SEQ ID NO: 74</u>
440	GCTCCtGTTGAACATC	Adenomatous polyposis coli	<u>SEQ ID NO: 75</u>
455	AAACTtTCATTTGATG	Adenomatous polyposis coli	<u>SEQ ID NO: 76</u>
455	AAACtttcaTTTGATGAAG	Adenomatous polyposis coli	<u>SEQ ID NO: 77</u>
472	CTAcAGGCCATTGC	Adenomatous polyposis coli	<u>SEQ ID NO: 78</u>
472	TAAATTAG_I10E11_GGgGACTACAGGC	Adenomatous polyposis coli	<u>SEQ ID NO: 79</u>
478	TTATtGCAAGTGGAC	Adenomatous polyposis coli	<u>SEQ ID NO: 80</u>
486	TACGgGCTTACTAAT	Adenomatous polyposis coli	<u>SEQ ID NO: 81</u>
494	AGTATtACACTAAGAC	Adenomatous polyposis coli	<u>SEQ ID NO: 82</u>
495	ATTACacTAAGACGATA	Adenomatous polyposis coli	<u>SEQ ID NO: 83</u>
497	CTAaGACGATATGC	Adenomatous polyposis coli	<u>SEQ ID NO: 84</u>
520	TGCTCtaTGAAAGGCTG	Adenomatous polyposis coli	<u>SEQ ID NO: 85</u>

		polyposis coli	
526	ATGAGagcactgtgGCCCCAACTAA	Adenomatous polyposis coli	<u>SEQ ID NO: 86</u>
539	GACTTaCAGCAG_E12I12_GTAC	Adenomatous polyposis coli	<u>SEQ ID NO: 87</u>
560	AAAAAgaCGTTGCGAGA	Adenomatous polyposis coli	<u>SEQ ID NO: 88</u>
566	GTTGgaagtGTGAAAGCAT	Adenomatous polyposis coli	<u>SEQ ID NO: 89</u>
570	AAAGCaTTGATGGAAT	Adenomatous polyposis coli	<u>SEQ ID NO: 90</u>
577	TTAGaagtTAAAAAG_E13I13_GTA	Adenomatous polyposis coli	<u>SEQ ID NO: 91</u>
584	ACCCTcAAAAGCGTAT	Adenomatous polyposis coli	<u>SEQ ID NO: 92</u>
591	GCCTtATGGAATTTG	Adenomatous polyposis coli	<u>SEQ ID NO: 93</u>
608	GCTgTAGATGGTGC	Adenomatous polyposis coli	<u>SEQ ID NO: 94</u>
617	GTTggcactcttactaccGGAGCCAGAC	Adenomatous polyposis coli	<u>SEQ ID NO: 95</u>
620	CTTACttacCGGAGCCAGA	Adenomatous polyposis coli	<u>SEQ ID NO: 96</u>
621	ACTTaCCGGAGCCAG	Adenomatous polyposis coli	<u>SEQ ID NO: 97</u>
624	AGCcaGACAAACACT	Adenomatous polyposis coli	<u>SEQ ID NO: 98</u>
624	AGCCagacAAACACTTTA	Adenomatous polyposis coli	<u>SEQ ID NO: 99</u>
626	ACAaacaCTTTAGCCAT	Adenomatous polyposis coli	<u>SEQ ID NO: 100</u>
629	TTAGCcATTATTGAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 101</u>
635	GGAGgTGGGATATTA	Adenomatous polyposis coli	<u>SEQ ID NO: 102</u>
638	ATATtACGGAATGTG	Adenomatous polyposis coli	<u>SEQ ID NO: 103</u>
639	TTACGgAATGTGTCCA	Adenomatous polyposis coli	<u>SEQ ID NO: 104</u>
657	AGAgAaGAACAACACTGT	Adenomatous polyposis coli	<u>SEQ ID NO: 105</u>
659	TATTTTCAG_I14E15_GCaaatcctaagagagAACAACTGTC	Adenomatous polyposis coli	<u>SEQ ID NO: 106</u>
660	AACTgtCTACAAACTT	Adenomatous polyposis coli	<u>SEQ ID NO: 107</u>
665	TTAttACAACACTTA	Adenomatous polyposis coli	<u>SEQ ID NO: 108</u>

668	CACttAAAATCTCAT	Adenomatous polyposis coli	SEQ ID NO: 109
673	AGTttgacaatagtCAGTAATGCA	Adenomatous polyposis coli	SEQ ID NO: 110
768	CACTTaTCAGAAACTT	Adenomatous polyposis coli	SEQ ID NO: 111
769	TTATcAGAAACTTTT	Adenomatous polyposis coli	SEQ ID NO: 112
770	TCAGAAACTTTTGACA	Adenomatous polyposis coli	SEQ ID NO: 113
780	AGTCcCAAGGCATCT	Adenomatous polyposis coli	SEQ ID NO: 114
792	AAGCaAAGTCTCTAT	Adenomatous polyposis coli	SEQ ID NO: 115
792	AAGCAaaGTCTCTATGG	Adenomatous polyposis coli	SEQ ID NO: 116
793	CAAAGTCTCTATGGT	Adenomatous polyposis coli	SEQ ID NO: 117
798	GATTatGTTTTTGACA	Adenomatous polyposis coli	SEQ ID NO: 118
802	GACACcaatcgacatGATGATAATA	Adenomatous polyposis coli	SEQ ID NO: 119
805	CGACatGATGATAATA	Adenomatous polyposis coli	SEQ ID NO: 120
811	TCAGacaaTTTTAATACT	Adenomatous polyposis coli	SEQ ID NO: 121
825	TATtTGAATACTAC	Adenomatous polyposis coli	SEQ ID NO: 122
827	AATAcTACAGTGTTA	Adenomatous polyposis coli	SEQ ID NO: 123
830	GTGTTaccagctcctctTCATCAAGAG	Adenomatous polyposis coli	SEQ ID NO: 124
833	AGCTCcTCTTCATCAA	Adenomatous polyposis coli	SEQ ID NO: 125
836	TCATcAAGAGGAAGC	Adenomatous polyposis coli	SEQ ID NO: 126
848	AAAGAtaGAAGTTTGGA	Adenomatous polyposis coli	SEQ ID NO: 127
848	AAAGatagaagTTTGGAGAGA	Adenomatous polyposis coli	SEQ ID NO: 128
855	GAACgCGGAATTGGT	Adenomatous polyposis coli	SEQ ID NO: 129
856	CGCGgaattGGTCTAGGCA	Adenomatous polyposis coli	SEQ ID NO: 130
856	CGCGgAATTGGTCTA	Adenomatous polyposis coli	SEQ ID NO: 131
879	CAGaTCTCCACCAC	Adenomatous polyposis coli	SEQ ID NO: 132



		polyposis coli	
902	GAAGAcagaAGTTCTGGGT	Adenomatous polyposis coli	<u>SEQ ID NO: 133</u>
907	GGGTcTACCACTGAA	Adenomatous polyposis coli	<u>SEQ ID NO: 134</u>
915	GTGACaGATGAGAGAA	Adenomatous polyposis coli	<u>SEQ ID NO: 135</u>
929	CATACacatTCAAACACTT	Adenomatous polyposis coli	<u>SEQ ID NO: 136</u>
930	ACACAttcaAACACTTACA	Adenomatous polyposis coli	<u>SEQ ID NO: 137</u>
931	CATtCAAACACTTA	Adenomatous polyposis coli	<u>SEQ ID NO: 138</u>
931	CATTcAAACACTTAC	Adenomatous polyposis coli	<u>SEQ ID NO: 139</u>
933	AACacttACAATTTTAC	Adenomatous polyposis coli	<u>SEQ ID NO: 140</u>
935	TACAatttactAAGTCGGAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 141</u>
937	TTCActaaGTCGGAAAAT	Adenomatous polyposis coli	<u>SEQ ID NO: 142</u>
939	AAGtcggAAAATTCAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 143</u>
946	ACATgTTCTATGCCT	Adenomatous polyposis coli	<u>SEQ ID NO: 144</u>
954	TTAGaaTACAAGAGAT	Adenomatous polyposis coli	<u>SEQ ID NO: 145</u>
961	AATgATAGTTTAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 146</u>
963	AGTTTaaATAGTGTCa	Adenomatous polyposis coli	<u>SEQ ID NO: 147</u>
964	TTAaataGTGTCAGTAG	Adenomatous polyposis coli	<u>SEQ ID NO: 148</u>
973	TATGgTAAAAGAGGT	Adenomatous polyposis coli	<u>SEQ ID NO: 149</u>
974	GGTAAaAGAGGTCAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 150</u>
975	AAAAgaGGTCAAATGA	Thyroid cancer	<u>SEQ ID NO: 151</u>
992	AGTAAgTTTTGCAGTT	Thyroid cancer	<u>SEQ ID NO: 152</u>
993	AAGttttgcagttaTGGTCAATAC	Adenomatous polyposis coli	<u>SEQ ID NO: 153</u>
999	CAAtaccagCCGACCTAGC	Adenomatous polyposis coli	<u>SEQ ID NO: 154</u>
1023	ACACcAATAAATTAT	Adenomatous polyposis coli	<u>SEQ ID NO: 155</u>
1030	AAAaATTCAGATGA	Adenomatous	<u>SEQ ID NO: 156</u>

		polyposis coli	
1032	TCAGatgagCAGTTGAACT	Adenomatous polyposis coli	<u>SEQ ID NO: 157</u>
1033	GATGaGCAGTTGAAC	Adenomatous polyposis coli	<u>SEQ ID NO: 158</u>
1049	TGGGcAAGACCCAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 159</u>
1054	CACAtaataGAAGATGAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 160</u>
1055	ATAAtagaaGATGAAATAA	Adenomatous polyposis coli	<u>SEQ ID NO: 161</u>
1056	ATAGaAaGATGAAATAA	Adenomatous polyposis coli	<u>SEQ ID NO: 162</u>
1060	ATAAAacaaaGTGAGCAAAG	Adenomatous polyposis coli	<u>SEQ ID NO: 163</u>
1061	AAAcaaaGTGAGCAAAG	Adenomatous polyposis coli	<u>SEQ ID NO: 164</u>
1061	AAACaaAGTGAGCAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 165</u>
1062	CAAAGtgaGCAAAGACAA	Adenomatous polyposis coli	<u>SEQ ID NO: 166</u>
1065	CAAAGacAATCAAGGAA	Adenomatous polyposis coli	<u>SEQ ID NO: 167</u>
1067	CAAtcaaGGAATCAAAG	Adenomatous polyposis coli	<u>SEQ ID NO: 168</u>
1071	CAAAGtACAACCTTATC	Adenomatous polyposis coli	<u>SEQ ID NO: 169</u>
1079	ACTGagAGCACTGATG	Adenomatous polyposis coli	<u>SEQ ID NO: 170</u>
1082	ACTGAtgATAAACACCT	Adenomatous polyposis coli	<u>SEQ ID NO: 171</u>
1084	GATaaacACCTCAAGTT	Adenomatous polyposis coli	<u>SEQ ID NO: 172</u>
1086	CACCTcAAGTTCCAAC	Adenomatous polyposis coli	<u>SEQ ID NO: 173</u>
1093	TTTGgACAGCAGGAA	Adenomatous polyposis coli	<u>SEQ ID NO: 174</u>
1098	TGTgtTTCTCCATAC	Adenomatous polyposis coli	<u>SEQ ID NO: 175</u>
1105	CGGgGAGCCAATGG	Thyroid cancer	<u>SEQ ID NO: 176</u>
1110	TCAGaAaACAAATCGAG	Adenomatous polyposis coli	<u>SEQ ID NO: 177</u>
1121	ATTAAtaaAATGTAAGCC	Adenomatous polyposis coli	<u>SEQ ID NO: 178</u>
1131	CAAgAAGATGACTA	Adenomatous polyposis coli	<u>SEQ ID NO: 179</u>

1134	GACTAtGAAGATGATA	Adenomatous polyposis coli	<u>SEQ ID NO: 180</u>
1137	GATgataaGCCTACCAAT	Adenomatous polyposis coli	<u>SEQ ID NO: 181</u>
1146	CGTTAcTCTGAAGAAG	Adenomatous polyposis coli	<u>SEQ ID NO: 182</u>
1154	GAAGaagaaGAGAGACCAA	Adenomatous polyposis coli	<u>SEQ ID NO: 183</u>
1155	GAAGaagaGAGACCAACA	Adenomatous polyposis coli	<u>SEQ ID NO: 184</u>
1156	GAAgagaGACCAACAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 185</u>
1168	GAAgagaaACGTCATGTG	Adenomatous polyposis coli	<u>SEQ ID NO: 186</u>
1178	GATTAtagtttaAAATATGCCA	Adenomatous polyposis coli	<u>SEQ ID NO: 187</u>
1181	TTAAaATATGCCACA	Adenomatous polyposis coli	<u>SEQ ID NO: 188</u>
1184	GCCacagaTATTCCTTCA	Adenomatous polyposis coli	<u>SEQ ID NO: 189</u>
1185	ACAgATATTCCTTCA	Adenomatous polyposis coli	<u>SEQ ID NO: 190</u>
1190	TCACAgAAACAGTCAT	Adenomatous polyposis coli	<u>SEQ ID NO: 191</u>
1192	AAAcGTCATTTTCA	Adenomatous polyposis coli	<u>SEQ ID NO: 192</u>
1198	TCAaaGAGTTCATCT	Adenomatous polyposis coli	<u>SEQ ID NO: 193</u>
1207	AAAAcCGAACATATG	Adenomatous polyposis coli	<u>SEQ ID NO: 194</u>
1208	ACCgaacATATGTCTTC	Adenomatous polyposis coli	<u>SEQ ID NO: 195</u>
1210	CATatGTCTTCAAGC	Adenomatous polyposis coli	<u>SEQ ID NO: 196</u>
1233	CCAAGtTCTGCACAGA	Adenomatous polyposis coli	<u>SEQ ID NO: 197</u>
1249	TGCAaaGTTTCTTCTA	Adenomatous polyposis coli	<u>SEQ ID NO: 198</u>
1259	ATAcGACTTATTGT	Adenomatous polyposis coli	<u>SEQ ID NO: 199</u>
1260	CAGACtATTGTGTAGA	Adenomatous polyposis coli	<u>SEQ ID NO: 200</u>
1268	CCAAATATGTTTTTC	Adenomatous polyposis coli	<u>SEQ ID NO: 201</u>
1275	AGTtCATTATCATC	Adenomatous polyposis coli	<u>SEQ ID NO: 202</u>
1294	CAGGAaGCAGATTCTG	Adenomatous	<u>SEQ ID NO: 203</u>

		polyposis coli	
1301	ACCCtGCAAATAGCA	Adenomatous polyposis coli	<u>SEQ ID NO: 204</u>
1306	GAAAtaaaAGAAAAGATT	Adenomatous polyposis coli	<u>SEQ ID NO: 205</u>
1307	ATAaAAGAAAAGAT	Adenomatous polyposis coli	<u>SEQ ID NO: 206</u>
1308	AAAgaaaAGATTGGAAC	Adenomatous polyposis coli	<u>SEQ ID NO: 207</u>
1308	AAAGAAaagaTTGGAAGTAG	Adenomatous polyposis coli	<u>SEQ ID NO: 208</u>
1318	GATCcTGTGAGCGAA	Adenomatous polyposis coli	<u>SEQ ID NO: 209</u>
1320	GTGAGcGAAGTTCCAG	Adenomatous polyposis coli	<u>SEQ ID NO: 210</u>
1323	GTTCcAGCAGTGTCA	Adenomatous polyposis coli	<u>SEQ ID NO: 211</u>
1329	CACCctagaaccAAATCCAGCA	Adenomatous polyposis coli	<u>SEQ ID NO: 212</u>
1336	AGACtgCAGGGTTCTA	Adenomatous polyposis coli	<u>SEQ ID NO: 213</u>
1338	CAGgGTTCTAGTTT	Adenomatous polyposis coli	<u>SEQ ID NO: 214</u>
1340	TCTAgTTTATCTTCA	Adenomatous polyposis coli	<u>SEQ ID NO: 215</u>
1342	TTATcTTCAGAATCA	Adenomatous polyposis coli	<u>SEQ ID NO: 216</u>
1352	GTTgAATTTTCTTC	Adenomatous polyposis coli	<u>SEQ ID NO: 217</u>
1361	CCCTcCAAAAGTGGT	Adenomatous polyposis coli	<u>SEQ ID NO: 218</u>
1364	AGTggtgCTCAGACACC	Adenomatous polyposis coli	<u>SEQ ID NO: 219</u>
1371	AGTCCacCTGAACACTA	Adenomatous polyposis coli	<u>SEQ ID NO: 220</u>
1372	CCACCtGAACACTATG	Adenomatous polyposis coli	<u>SEQ ID NO: 221</u>
1376	TATGttCAGGAGACCC	Adenomatous polyposis coli	<u>SEQ ID NO: 222</u>
1394	GATAgTtTTGAGAGTC	Adenomatous polyposis coli	<u>SEQ ID NO: 223</u>
1401	ATTGCcAGCTCCGTTC	Adenomatous polyposis coli	<u>SEQ ID NO: 224</u>
1415	AGTGGcATTATAAGCC	Adenomatous polyposis coli	<u>SEQ ID NO: 225</u>
1426	AGCCcTGGACAAACC	Adenomatous polyposis coli	<u>SEQ ID NO: 226</u>

1427	CCTGGaCAAACCATGC	Adenomatous polyposis coli	<u>SEQ ID NO: 227</u>
1431	ATGCcACCAAGCAGA	Adenomatous polyposis coli	<u>SEQ ID NO: 228</u>
1454	AAAAAtAAAGCACCTA	Adenomatous polyposis coli	<u>SEQ ID NO: 229</u>
1461	GAAaAGAGAGAGAG	Adenomatous polyposis coli	<u>SEQ ID NO: 230</u>
1463	AGAgagaGTGGACCTAA	Adenomatous polyposis coli	<u>SEQ ID NO: 231</u>
1464	GAGAgTGGACCTAAG	Adenomatous polyposis coli	<u>SEQ ID NO: 232</u>
1464	GAGAgTGGACCTAAGC	Adenomatous polyposis coli	<u>SEQ ID NO: 233</u>
1464	GAGAgTGGACCTAAG	Adenomatous polyposis coli	<u>SEQ ID NO: 234</u>
1492	GCCaCGGAAAGTAC	Adenomatous polyposis coli	<u>SEQ ID NO: 235</u>
1493	ACGGaAGTACTCCAG	Adenomatous polyposis coli	<u>SEQ ID NO: 236</u>
1497	CCAgATGGATTTTC	Adenomatous polyposis coli	<u>SEQ ID NO: 237</u>
1503	TCAAtccaGCCTGAGTGC	Adenomatous polyposis coli	<u>SEQ ID NO: 238</u>
1522	TTAagaataaTGCCTCCAGT	Adenomatous polyposis coli	<u>SEQ ID NO: 239</u>
1536	GAAACagAATCAGAGCA	Adenomatous polyposis coli	<u>SEQ ID NO: 240</u>
1545	TCAAAtgaaaACCAAGAGAA	Adenomatous polyposis coli	<u>SEQ ID NO: 241</u>
1547	GAAaACCAAGAGAA	Adenomatous polyposis coli	<u>SEQ ID NO: 242</u>
1550	GAGAAagaGGCAGAAAAA	Adenomatous polyposis coli	<u>SEQ ID NO: 243</u>
1577	GAATgtATTATTTCTG	Adenomatous polyposis coli	<u>SEQ ID NO: 244</u>
1594	CCAGCcCAGACTGCTT	Adenomatous polyposis coli	<u>SEQ ID NO: 245</u>
1596	CAGACtGCTTCAAAAT	Adenomatous polyposis coli	<u>SEQ ID NO: 246</u>
1823	TTCAaTGATAAGCTC	Adenomatous polyposis coli	<u>SEQ ID NO: 247</u>
1859	AATGAttctTTGAGTTCTC	Adenomatous polyposis coli	<u>SEQ ID NO: 248</u>
1941	CCAGAcagaGGGGCAGCAA	Desmoid tumours	<u>SEQ ID NO: 249</u>
1957	GAAaATACTCCAGT	Adenomatous	<u>SEQ ID NO: 250</u>

		polyposis coli	
1980	AACaATAAAGAAAA	Adenomatous polyposis coli	SEQ ID NO: 251
1985	GAACcTATCAAAGAGA	Adenomatous polyposis coli	SEQ ID NO: 252
1986	CCTaTCAAAGAGAC	Adenomatous polyposis coli	SEQ ID NO: 253
1998	GAACcAAGTAAACCT	Adenomatous polyposis coli	SEQ ID NO: 254
2044	AGCTCcGCAATGCCAA	Adenomatous polyposis coli	SEQ ID NO: 255
2556	TCATCccctctcGAGTAAGCAC	Adenomatous polyposis coli	SEQ ID NO: 256
2643	CTAATttatCAAATGGCAC	Adenomatous polyposis coli	SEQ ID NO: 257

**TABLE X: SMALL INSERTIONS**

Codon	Insertion	Phenotype
157	T	Adenomatous polyposis coli
170	AGAT	Adenomatous polyposis coli
172	T	Adenomatous polyposis coli
199	G	Adenomatous polyposis coli
243	AG	Adenomatous polyposis coli
266	T	Adenomatous polyposis coli
357	A	Adenomatous polyposis coli
405	C	Adenomatous polyposis coli
413	T	Adenomatous polyposis coli
416	A	Adenomatous polyposis coli
457	G	Adenomatous polyposis coli
473	A	Adenomatous polyposis coli
503	ATTC	Adenomatous polyposis coli
519	C	Adenomatous polyposis coli
528	A	Adenomatous polyposis coli
561	A	Adenomatous polyposis coli
608	A	Adenomatous polyposis coli
620	CT	Adenomatous polyposis coli
621	A	Adenomatous polyposis coli
623	TTAC	Adenomatous polyposis coli
627	A	Adenomatous polyposis coli
629	A	Adenomatous polyposis coli
636	GT	Adenomatous polyposis coli
639	A	Adenomatous polyposis coli
704	T	Adenomatous polyposis coli
740	ATGC	Adenomatous polyposis coli
764	T	Adenomatous polyposis coli

779	TT	Adenomatous polyposis coli
807	AT	Adenomatous polyposis coli
827	AT	Adenomatous polyposis coli
831	A	Adenomatous polyposis coli
841	CTTA	Adenomatous polyposis coli
865	CT	Adenomatous polyposis coli
865	AT	Adenomatous polyposis coli
900	TG	Adenomatous polyposis coli
921	G	Adenomatous polyposis coli
927	A	Adenomatous polyposis coli
935	A	Adenomatous polyposis coli
936	C	Adenomatous polyposis coli
975	A	Adenomatous polyposis coli
985	T	Adenomatous polyposis coli
997	A	Adenomatous polyposis coli
1010	TA	Adenomatous polyposis coli
1085	C	Adenomatous polyposis coli
1085	AT	Adenomatous polyposis coli
1095	A	Adenomatous polyposis coli
1100	GTTT	Adenomatous polyposis coli
1107	GGAG	Adenomatous polyposis coli
1120	G	Adenomatous polyposis coli
1166	A	Adenomatous polyposis coli
1179	T	Adenomatous polyposis coli
1187	A	Adenomatous polyposis coli
1211	T	Adenomatous polyposis coli
1256	A	Adenomatous polyposis coli
1265	T	Adenomatous polyposis coli
1267	GATA	Adenomatous polyposis coli
1268	T	Adenomatous polyposis coli
1301	A	Adenomatous polyposis coli
1301	C	Adenomatous polyposis coli
1323	A	Adenomatous polyposis coli
1342	T	Adenomatous polyposis coli
1382	T	Adenomatous polyposis coli
1458	GTAG	Adenomatous polyposis coli
1463	AG	Adenomatous polyposis coli
1488	T	Adenomatous polyposis coli
1531	A	Adenomatous polyposis coli
1533	T	Adenomatous polyposis coli
1554	A	Adenomatous polyposis coli
1555	A	Adenomatous polyposis coli
1556	T	Adenomatous polyposis coli
1563	GACCT	Adenomatous polyposis coli
1924	AA	Desmoid tumours

**TABLE XI: SMALL INSERTIONS/DELETIONS**

Location/ codon	Deletion	Insertion	Phenotype	
538	GAAGAcTTACAGCAGG	gaa	Adenomatous polyposis coli	<u>SEQ ID NO: 258</u>
620	CTTACttaCCGGAGCCAG	ct	Adenomatous polyposis coli	<u>SEQ ID NO: 259</u>
728	AATtctcatGGCAAATAGG ( <u>SEQ ID NO: 260</u> )	ttgcagctttaa ( <u>SEQ ID NO: 261</u> )	Adenomatous polyposis coli	
971	GATGgtTATGGTAAAA	taa	Adenomatous polyposis coli	<u>SEQ ID NO: 262</u>

**TABLE XII: GROSS DELETIONS**

2 kb including ex. 11	Adenomatous polyposis coli
3 kb I10E11-1.5 kb to I12E13-170 bp	Adenomatous polyposis coli
335 bp nt. 1409-1743 ex. 11-13	Adenomatous polyposis coli
6 kb incl. ex. 14	Adenomatous polyposis coli
817 bp I13E14-679 to I13E14+138	Adenomatous polyposis coli
ex. 11-15M	Adenomatous polyposis coli
ex. 11-3'UTR	Adenomatous polyposis coli
ex. 15A – ex. 15F	Adenomatous polyposis coli
ex. 4	Adenomatous polyposis coli
ex. 7, 8 and 9	Adenomatous polyposis coli
ex. 8 to beyond ex. 15F	Adenomatous polyposis coli
ex. 8 - ex. 15F	Adenomatous polyposis coli
ex. 9	Adenomatous polyposis coli
>10mb (del 5q22)	Adenomatous polyposis coli

**TABLE XIII: GROSS INSERTIONS AND DUPLICATIONS**

Description	Phenotype
Insertion of 14 bp nt. 3816	Adenomatous polyposis coli
Insertion of 22 bp nt. 4022	Adenomatous polyposis coli
Duplication of 43 bp cd. 1295	Adenomatous polyposis coli
Insertion of 337 bp of Alu I sequence cd. 1526	Desmoid tumours

**TABLE XIV: COMPLEX REARRANGEMENTS (INCLUDING INVERSIONS)**

A-T nt. 4893 Q1625H, Del C nt. 4897 cd. 1627	Adenomatous polyposis coli
Del 1099 bp I13E14-728 to E14I14+156, ins 126 bp	Adenomatous polyposis coli



Del 1601 bp E14I14+27 to E14I14+1627, ins 180 bp	Adenomatous polyposis coli
Del 310 bp, ins. 15 bp nt. 4394, cd 1464	Adenomatous polyposis coli
Del A and T cd. 1395	Adenomatous polyposis coli
Del TC nt. 4145, Del TGT nt. 4148	Adenomatous polyposis coli
Del. T, nt. 983, Del. 70 bp, nt. 985	Adenomatous polyposis coli
Del. nt. 3892-3903, ins ATTT	Adenomatous polyposis coli

**TABLE XV: DIAGNOSTIC APPLICATIONS**

Cancer Type	Marker	Application	Reference
Breast	Her2/Neu Detection – polymorphism at codon 655 (GTC/valine to ATC/isoleucine [Val(655)Ile])	Using methods described herein, design second primer such that after PCR, and digestion with restriction enzyme, a 5' overhang containing DNA sequence for codon 655 of Her2/Neu is generated. Her2/Neu can be detected and quantified as a possible marker for breast cancer. Methods described herein can detect both mutant allele and normal allele, even when mutant allele is small fraction of total DNA. Herceptin therapy for breast cancer is based upon screening for Her2. The earlier the mutant allele can be detected, the faster therapy can be provided.	D. Xie <i>et al.</i> , <i>J. Natl. Cancer Institute</i> , 92, 412 (2000) K.S. Wilson <i>et al.</i> , <i>Am. J. Pathol.</i> , 161, 1171 (2002) L. Newman, <i>Cancer Control</i> , 9, 473 (2002)
Breast/Ovarian	Hypermethylation of BRCA1	Methods described herein can be used to differentiate between tumors resulting from inherited BRCA1 mutations and those from non-inherited abnormal methylation of the gene	M. Esteller <i>et al.</i> , <i>New England Jnl Med.</i> , 344, 539 (2001)
Bladder	Microsatellite analysis of free tumor DNA in Urine, Serum and Plasma	Methods described herein can be applied to microsatellite analysis and FGFR3 mutation analysis for detection of bladder cancer. Methods described herein provide a non-invasive method for detection of bladder cancer.	W.G. Bas <i>et al.</i> , <i>Clinical Cancer Res.</i> , 9, 257 (2003) M. Utting <i>et al.</i> , <i>Clinical Cancer Res.</i> , 8, 35 (2002) L. Mao, D. Sidransky <i>et al.</i> , <i>Science</i> , 271, 669 (1996)
Lung	Microsatellite analysis of DNA from sputum	Methods described herein can be used to detect mutations in sputum samples, and can markedly boost the accuracy of preclinical lung cancer screening	T. Liloglou <i>et al.</i> , <i>Cancer Research</i> , 61, 1624, (2001) M. Tockman <i>et al.</i> , <i>Cancer Control</i> , 7, 19 (2000) Field <i>et al.</i> , <i>Cancer Research</i> , 59, 2690

			(1999)
Cervical	Analysis of HPV genotype	Methods described herein can be used to detect HPV genotype from a cervical smear preparation.	N. Munoz <i>et al.</i> , <i>New England Jnl Med.</i> , 348, 518 (2003)
Head and Neck	Tumor specific alterations in exfoliated oral mucosal cells (microsatellite markers)	Methods described herein can be used to detect any of 23 microsatellite markers, which are associated with Head and Neck Squamous Cell Carcinoma (HNSCC).	M. Spafford <i>et al.</i> , <i>Clinical Cancer Research</i> , 17, 607 (2001) A. El-Naggar <i>et al.</i> , <i>J. Mol. Diag.</i> , 3, 164 (2001)
Colorectal	Screening for mutation in K-ras2 and APC genes.	Methods described herein can be used to detect K-ras 2 mutations, which can be used as a prognostic indicator for colorectal cancer. APC (see Example 5).	B. Ryan <i>et al.</i> , <i>Gut</i> , 52, 101 (2003)
Prostate	GSTP1 Hypermethylation	Methods described herein can be used to detect GSTP1 hypermethylation in urine from patients with prostate cancer; this can be a more accurate indicator than PSA.	P. Cairns <i>et al.</i> , <i>Clin. Can. Res.</i> , 7, 2727 (2001)

## HIV

Antiretroviral resistance	Screening individuals for mutations in HIV virus – e.g. 154V mutation or CCR5 Δ 32 allele.	Methods described herein can be used for detection of mutations in the HIV virus. Treatment outcomes are improved in individuals receiving anti-retroviral therapy based upon resistance screening.	J. Durant <i>et al.</i> , <i>The Lancet</i> , 353, 2195 (1999)
---------------------------	--	---	--

## Cardiology

Congestive Heart Failure	Synergistic polymorphisms of beta1 and alpha2c adrenergic receptors	Methods described herein can be used to genotype these loci and may help identify people who are at a higher risk of heart failure.	K. Small <i>et al.</i> , <i>New Eng. Jnl. Med.</i> , 347, 1135 (2002)
--------------------------	---	---	---

## EXAMPLE 8

**[0657]** Single nucleotide polymorphisms (SNPs) represent the most common form of sequence variation; three million common SNPs with a population frequency of over 5% have been estimated to be present in the human genome. A genetic map using these polymorphisms as a guide is being developed (<http://research.marshfieldclinic.org/genetics/>; internet address as of February 13, 2003).

**[0658]** The allele frequency varies from SNP to SNP; the allele frequency for one SNP may be 50:50, while the allele frequency for another SNP may be 90:10. The closer the allele frequency is to 50:50, the more likely any particular individual will be heterozygous at that SNP. The SNP consortium provides allele frequency information for some SNPs but not for others. [www.snp.chsl.org](http://www.snp.chsl.org). The allele frequency for a particular SNP provides valuable information as to the utility of that SNP for the non-invasive prenatal screening method described in Example 5. While all SNPs can be used, SNPs with allele frequencies closer to 50:50 are preferable.

**[0659]** Briefly, maternal blood contains fetal DNA. Maternal DNA can be distinguished from fetal DNA by examining SNPs wherein the mother is homozygous. For example, at SNP X, the maternal DNA may be homozygous for guanine. If template DNA obtained from the plasma of a pregnant female is heterozygous, as demonstrated by the detection of signals corresponding to an adenine allele and a guanine allele, the adenine allele can be used as a beacon for the fetal DNA (see Example 5). The closer the allele frequency of a SNP is to 50:50, the more likely there will be allele differences at a particular SNP between the maternal DNA and the fetal DNA.

**[0660]** For example, if at SNP X the observed alleles are adenine and guanine, and the SNP has an allele frequency of 90(A):10(G), it is likely that both mother and father will be homozygous for adenine at that particular SNP. Thus, both the maternal DNA and the fetal DNA will be homozygous for adenine, and there is no distinct signal for the fetal DNA. However, if at SNP X the allele frequency is 50:50, and the mother is homozygous for adenine, the probability is higher that the paternal DNA will contain a guanine allele at SNP X.

**[0661]** Below, a method for determining the allele frequency for a SNP is provided. Seven SNPs located on chromosome 13 were analyzed. The method is applicable for any SNP including but not limited to the SNPs on human chromosomes 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X and Y.

### **Preparation of Template DNA**

**[0662]** To determine the allele frequency of a particular SNP, DNA was obtained from two hundred and fifty individuals after informed consent had been granted. From each individual, a 9 ml blood sample was collected into a sterile tube (Fischer Scientific, 9 ml EDTA Vacuette tubes, catalog number NC9897284). The tubes were spun at 1000 rpm for ten minutes. The supernatant (the plasma) of each sample was removed, and one milliliter of the remaining blood sample, which is commonly referred to as the “buffy-coat” was transferred to a new tube. One milliliter of 1X PBS was added to each sample.

**[0663]** Template DNA was isolated using the QIAmp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183). The template DNA was isolated as per instructions included in the kit. From each individual, 0.76 µg of DNA was pooled together, and the pooled DNA was used in all subsequent reactions.

### **Design of Primers**

**[0664]** SNP TSC0903430 was amplified using the following primer set:

First primer:

5' GTCTTGCATGTAGAAATTCTAGGGACGCTGCTTTTCGTC 3' (SEQ ID NO: 279)

Second primer:

5' CTCCTAGACATCGGGACTAGAATGTCCAC 3' (SEQ ID NO: 280)

**[0665]** The first primer contained a recognition site for the restriction enzyme EcoRI, and was designed to anneal eighty-two bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

**[0666]** SNP TSC0337961 was amplified using the following primer set:

First primer:

5' ACACAAGGCAGAGAATTCCAGTCCTGAGGGTGGGGGCC 3' (SEQ ID NO: 281)

Second primer:

5' CCGTGTTTTAACGGGACAAGCTGTTCTTC 3' (SEQ ID NO: 282)

**[0667]** The first primer contained a recognition site for the restriction enzyme EcoRI, and was designed to anneal ninety-two bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

**[0668]** SNP TSC0786441 was amplified using the following primer set:

First primer:

5' GTAGCGGAGGTTGAATTCTATATGTTGTCTTGGACATT 3' (SEQ ID NO: 283).

Second primer: /

5' CATCAGTAGAGTGGGACGAAAGTTCTGGC 3' (SEQ ID NO: 284)

[0669] The first primer contained a recognition site for the restriction enzyme EcoRI, and was designed to anneal one hundred and four bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0670] SNP TSC1168303 was amplified using the following primer set:

First primer:

5' ATCCACGCCGCAGAATTCGTATTCATGGGCATGTCAAA 3' (SEQ ID NO: 285)

Second primer:

5' CTTGGGACTATTGGGACCAGTGTTCATC 3' (SEQ ID NO: 286)

[0671] The first primer contained a recognition site for the restriction enzyme EcoRI, and was designed to anneal sixty-four bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0672] SNP TSC0056188 was amplified using the following primer set:

First primer:

5' CCAGAAAGCCGTGAATTCGTTAAGCCAACTGACTCCA 3' (SEQ ID NO: 287)

Second primer:

5' TCGGGGTTAGTCGGGACATCCAGCAGCCC 3' (SEQ ID NO: 288)

[0673] The first primer contained a recognition site for the restriction enzyme EcoRI, and was designed to anneal eighty-two bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0674] SNP TSC0466177 was amplified using the following primer set:

First primer:

5' CGAAGGTAATGTGAATTCCAAAACCTTAGTGCCACAATT 3' (SEQ ID NO: 289)

Second primer:

5' ATACCGCCCAACGGGACAGATCCATTGAC 3' (SEQ ID NO: 290)

[0675] The first primer contained a recognition site for the restriction enzyme EcoRI, and was designed to anneal ninety-two bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0676] SNP TSC0197424 was amplified using the following primer set:

First primer:

5' AGAAACCTGTAAGAATTCGATTCCAAATTGTTTTTTGG 3' (SEQ ID NO: 291)

Second primer:

5' CGATCATAGGGGGGGACAGGAGAGAGCAC 3' (SEQ ID NO: 292)

[0677] The first primer contained a recognition site for the restriction enzyme EcoRI, and was designed to anneal one hundred and four bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0678] The first primer was designed to anneal at various distances from the locus of interest. The skilled artisan understands that the annealing location of the first primer can be any distance from the locus of interest including but not limited to 5-10, 11-15, 16-20, 21-25, 26-30, 31-35, 36-40, 41-45, 46-50, 51-55, 56-60, 61-65, 66-70, 71-75, 76-80, 81-85, 86-90, 91-95, 96-100, 101-105, 106-110, 111-115, 116-120, 121-125, 126-130, 131-140, 141-160, 161-180, 181-200, 201-220, 221-240, 241-260, 261-280, 281-300, 301-350, 351-400, 401-450, 451-500, 501-1000, 1001-2000, 2001-3000, or greater than 3000.

[0679] All loci of interest were amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). In this example, the loci of interest were amplified in separate reaction tubes but they can also be amplified together in a single PCR reaction. For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction can be optimized for each locus of interest. In this example, 40 ng of template human genomic DNA (a mixture of template DNA from 245 individuals) and 5  $\mu$ M of each primer were used. Forty cycles of PCR were performed. The following PCR conditions were used:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 64°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty nine (39) times;
- (9) 72°C for 5 minutes.

[0680] In the first cycle of PCR, the annealing temperature was about the melting temperature of the 3' annealing region of the second primers, which was 37°C. The annealing temperature in the second cycle of PCR was about the melting temperature of the 3' region, which anneals to the template DNA, of

the first primer, which was 57°C. The annealing temperature in the third cycle of PCR was about the melting temperature of the entire sequence of the second primer, which was 64°C. The annealing temperature for the remaining cycles was 64°C. Escalating the annealing temperature from TM1 to TM2 to TM3 in the first three cycles of PCR greatly improves specificity. These annealing temperatures are representative, and the skilled artisan will understand the annealing temperatures for each cycle are dependent on the specific primers used.

**[0681]** The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results.

### **Purification of Fragment of Interest**

**[0682]** The PCR products were separated from the unused PCR reagents. After the PCR reaction, 1/2 of the reaction volume for SNP TSC0903430, SNP TSC0337961, and SNP TSC0786441 were mixed together in a single reaction tube. One-half the reaction volumes for SNPs TSC1168303, TSC0056188, TSC0466177, and TSC0197424 were pooled together in a single reaction tube. The un-used primers, and nucleotides were removed from the reaction by using Qiagen MinElute PCR purification kits (Qiagen, Catalog Number 28004). The reactions were performed following the manufacturer's instructions supplied with the columns.

### **Restriction Enzyme Digestion of Isolated Fragments**

**[0683]** The purified PCR products were digested with the restriction enzyme BsmF I, which binds to the recognition site incorporated into the PCR products from the second primer. The digests were performed in eppendorf tubes following the instructions supplied with the restriction enzyme.

### **Incorporation of Labeled Nucleotide**

**[0684]** The restriction enzyme digest with BsmF I yielded a DNA fragment with a 5' overhang, which contained the SNP site or locus of interest and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase.

**[0685]** As discussed in detail in Example 6, the sequence of both alleles of a SNP can be determined with one labeled nucleotide in the presence of the other unlabeled nucleotides. The following components were added to each fill in reaction: 1 µl of fluorescently labeled ddGTP, 0.5 µl of unlabeled ddNTPs (40 µM), which contained all nucleotides except guanine, 2 µl of 10X sequenase buffer, 0.25 µl of Sequenase, and water as needed for a 20µl reaction. The fill in reaction was performed at 40°C for 10 min. Sequenase was the DNA polymerase used in this example. However, any DNA polymerase can be

used for a fill-in reaction including but not limited to E. coli DNA polymerase, Klenow fragment of E. coli DNA polymerase I, T7 DNA polymerase, T4 DNA polymerase, Taq polymerase, Pfu DNA polymerase, Vent DNA polymerase, polymerase from bacteriophage 29, and REDTaq™ Genomic DNA polymerase. Non-fluorescently labeled ddNTP was purchased from Fermentas Inc. (Hanover, MD). All other labeling reagents were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565).

### Detection of the Locus of Interest

**[0686]** The sample was loaded into a lane of a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691). The sample was electrophoresed into the gel at 3000 volts for 3 min. The gel was run for 3 hours on a sequencing apparatus (Hoefer SQ3 Sequencer). The gel was removed from the apparatus and scanned on the Typhoon 9400 Variable Mode Imager. The incorporated labeled nucleotide was detected by fluorescence.

**[0687]** Below, a schematic of the 5' overhang for SNP TSC0056188 is reproduced (where R indicates the variable site). The entire sequence is not shown, only a portion of the overhang.

5'CCA				
3'GGT	R	T	C	C
Overhang position	1	2	3	4

**[0688]** As discussed in detail in Example 6, one nucleotide labeled with one chemical moiety can be used to determine the sequence of the alleles of a locus of interest. The observed nucleotides for TSC0056188 on the 5' sense strand (here depicted as the top strand) are adenine and guanine. The third position in the overhang on the antisense strand is cytosine, which is complementary to guanine. As the variable site can be adenine or guanine, fluorescently labeled ddGTP in the presence of unlabeled dCTP, dTTP, and dATP was used to determine the sequence of both alleles. The fill-in reactions for an individual homozygous for guanine, homozygous for adenine or heterozygous are diagrammed below.

**[0689]** Homozygous adenine:

5'CCA	A	A	G*	
3'GGT	T	T	C	C
Overhang position	1	2	3	4

**[0690]** Homozygous guanine:

5'CCA	G*
-------	----



	3'GGT	C	T	C	C
Overhang position		1	2	3	4
<b>[0691] Heterozygous:</b>					
Allele 1	5'CCA	<b>G*</b>			
	3'GGT	C	T	C	C
Overhang position		1	2	3	4
Allele 2	5'CCA	<b>A</b>	<b>A</b>	<b>G*</b>	
	3'GGT	T	T	C	C
Overhang position		1	2	3	4

**[0692]** As seen in FIG. 14, two bands were detected for SNP TSC0056188. The lower band corresponded to DNA molecules filled in with ddGTP at position one complementary to the overhang, which is representative of the guanine allele. The higher band, separated by a single base from the lower band, corresponded to DNA molecules filled in with ddGTP at position 3 complementary to the overhang. This band represented the adenine allele. The intensity of each band was strong, indicating that each allele was well represented in the population. SNP TSC0056188 is representative of a SNP with high allele frequency.

**[0693]** Below, a schematic of the 5' overhang generated after digestion with BsmF I for SNP TSC0337961 is reproduced (where R indicates the variable site). The entire sequence is not shown, only a portion of the overhang.

	5' GCCA				
	3' CGGT	R	G	C	T
Overhang position		1	2	3	4

**[0694]** The observed nucleotides for SNP TSC0337961 on the 5' sense strand (here depicted as the top strand) are adenine and guanine. The third position in the overhang on the antisense strand was cytosine, which is complementary to guanine. As the variable site can be adenine or guanine, fluorescently labeled ddGTP in the presence of unlabeled dCTP, dTTP, and dATP was used to determine the sequence of both alleles. The fill-in reactions for an individual homozygous for guanine, homozygous for adenine or heterozygous are diagrammed below.

**[0695]** Homozygous for guanine:

	5' GCCA	<b>G*</b>			
	3' CGGT	C	G	C	T
Overhang position		1	2	3	4

**[0696]** Homozygous for adenine:

	5' GCCA	A	C	G*	
	3' CGGT	T	G	C	T
Overhang position		1	2	3	4

**[0697]** Heterozygous

Allele 1	5' GCCA	G*			
	3' CGGT	C	G	C	T
Overhang position		1	2	3	4
Allele 2	5' GCCA	A	C	G*	
	3' CGGT	T	G	C	T
Overhang position		1	2	3	4

**[0698]** As seen in FIG. 14, one band migrating at the position of the expected lower molecular weight band was observed. This band represented the DNA molecules filled in with ddGTP at position one complementary to the overhang, which represents the guanine allele. No band corresponding to the DNA molecules filled in with ddGTP at position 3 complementary to the overhang was detected. SNP TSC0337961 is representative of a SNP that is not highly variable within the population.

**[0699]** Of the seven SNPs analyzed, four of the SNPs (TSC1168303, TSC0056188, TSC0466177, and TSC0197424 had high allele frequencies. Two bands of high intensity were seen for each of the four SNPs, indicating that both alleles were well represented in the population.

**[0700]** However, it is not necessary that the SNPs have allele frequencies of 50:50 to be useful. All SNPs provide useful information. The methods described herein provide a rapid technique for determining the allele frequency of a SNP, or any variable site including but not limited to point mutations. Allele frequencies of 50:50, 51:49, 52:48, 53:47, 54:46, 55:45, 56:46, 57:43, 58:42, 59:41, 60:40, 61:39, 62:38, 63:37, 64:36, 65:35, 66:34, 67:33, 68:32, 69:31, 70:30, 71:29, 72:28, 73:27, 74:26, 75:25, 76:24, 77:23, 78:22, 79:21, 80:20, 81:19, 82:18, 83:17, 84:16, 85:15, 86:14, 87:13, 88:12, 89:11, 90:10, 91:9, 92:8, 93:7, 94:6, 95:5, 96:4, 97:3, 98:2, 99:1 and 100:0 can be useful.

**[0701]** Two bands were seen for SNP TSC0903430. One band, the lower molecular weight band represented the DNA molecules filled in with labeled ddGTP. A band of weaker intensity was seen for the molecules filled in with labeled ddGTP at position 3 complementary to the overhang, which represented the cytosine allele. SNP TSC0903430 represents a SNP with low allele frequency variation. In the population, the majority of individuals carry the guanine allele, but the cytosine allele is still present.

**[0702]** One band of high intensity was seen for SNP TSC0337961 and SNP TSC0786441. The band detected for both SNP TSC0337961 and SNP TSC0786441 corresponded to the DNA molecules filled in with ddGTP at position 1 complementary to the overhang. No signal was detected from DNA molecules that would have been filled in at position 3 complementary to the overhang, which would have represented the second allele. SNP TSC0337961 and SNP TSC0786441 represent SNPs with little variability in the population.

**[0703]** As demonstrated in FIG 14., the first primer used to amplify each locus of interest can be designed to anneal at various distances from the locus of interest. This allows multiple SNPs to be analyzed in the same reaction. By designing the first primer to anneal at specified distances from the loci of interest, any number of loci of interest can be analyzed in a single reaction including but not limited to 1-10, 11-20, 21-30, 31-40, 41-50, 51-60, 61-70, 71-80, 81-90, 91-100, 101-110, 111-120, 121-130, 131-140, 141-150, 151-160, 161-170, 171-180, 181-190, 191-200, 201-300, 301-400, 401-500, and greater than 500.

**[0704]** As discussed in Example 6, some type IIS restriction enzymes display alternate cutting patterns. For example, the type IIS restriction enzyme BsmF I typically cuts 10/14 from its binding site; however, the enzyme also can cut 11/15 from the binding site. To eliminate the effect of the alternate cut, the labeled nucleotide used for the fill-in reaction should be chosen such that it is not complementary to position 0 of the overhang generated by the 11/15 cut (discussed in detail in Example 6). For instance, if you label with ddGTP, the nucleotide preceding the variable site on the strand that is filled in should not be a guanine.

**[0705]** The 11/15 overhang generated by BsmF I for SNP TSC0056188 is depicted below, with the variable site in bold-typeface:

11/15 Overhang for TSC0056188

Allele 1	5'CC				
	3'GG	T	C	T	C
Overhang position		0	1	2	3
Allele 2	5'CC				
	3'GG	T	<b>T</b>	T	C
Overhang position		0	1	2	3

**[0706]** After the fill-in reaction with labeled ddGTP, unlabeled dATP, dTTP, and dCTP, the following molecules were generated:

11/15 Allele 1	5'CC	A	<b>G*</b>
----------------	------	---	-----------

	3'GG	T	C	T	C
	Overhang position	0	1	2	3
11/15 Allele 2	5'CC	A	A	A	G*
	3'GG	T	T	T	C
	Overhang position	0	1	2	3

**[0707]** Two signals were seen; one band corresponded to molecules filled in with ddGTP at position one of the overhang, and the other band corresponded to the molecules filled in with ddGTP at position 3 complementary to the overhang. These are the same DNA molecules generated after the fill-in reaction of the 10/14 overhang. Thus, the two bands can be compared without any ambiguity from the alternate cut. This method of labeling with a single nucleotide eliminates any errors generated from the alternate cutting properties of the enzymes.

**[0708]** The methods described herein is applicable to determining the allele frequency of any SNP including but not limited to SNPs on human chromosomes 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X and Y.

#### EXAMPLE 9

**[0709]** Heterozygous SNPs, by definition, differ by one nucleotide. At a heterozygous SNP, allele 1 and allele 2 may be present at a ratio of 1:1. However, it is possible that DNA polymerases can incorporate one nucleotide at a faster rate than other nucleotides, and thus the observed ratio of a heterozygous SNP may differ from the theoretically expected 1:1 ratio.

**[0710]** Below, methods are described that allow efficient and accurate quantitation for the expected ratio of allele 1 to allele 2 at a heterozygous SNP.

#### Preparation of Template DNA

**[0711]** Template DNA was obtained from twenty-four individuals after informed consent had been granted. From each individual, a 9 ml blood sample was collected into a sterile tube (Fischer Scientific, 9 ml EDTA Vacuette tubes, catalog number NC9897284). The tubes were spun at 1000 rpm for ten minutes without brake. The supernatant (the plasma) of each sample was removed, and one milliliter of the remaining blood sample, which is commonly referred to as the "buffy-coat" was transferred to a new tube. One milliliter of 1X PBS was added to each sample.

**[0712]** Template DNA was isolated using the QIAmp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183). The template DNA was isolated as per instructions included in the kit.

## Design of Primers

[0713] SNP TSC0607185 was amplified using the following primer set:

First primer:

5' ACTTGATTCCGTGAATTCGTTATCAATAAATCTTACAT 3' (SEQ ID NO: 40)

Second primer:

5' CAAGTTGGATCCGGGACCCAGGGCTAACC 3' (SEQ ID NO: 41)

[0714] SNP TSC1130902 was amplified using the following primer set:

First primer:

5' TCTAACCATTGCGAATTCAGGGCAAGGGGGGTGAGATC 3' (SEQ ID NO: 34)

Second primer:

5' TGAAGTTGGATCCGGGACAACGACTCATCC 3' (SEQ ID NO: 35)

[0715] The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI. The second primer contained the recognition site for the restriction enzyme BsmF I. The first primer was designed to anneal at various distances from the locus of interest.

[0716] The first primer for SNP TSC0607185 was designed to anneal ninety bases from the locus of interest. The first primer for SNP TSC1130902 was designed to anneal sixty bases from the locus of interest.

[0717] All loci of interest were amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). In this example, the loci of interest were amplified in separate reaction tubes but they could also be amplified together in a single PCR reaction. For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction can be optimized for each locus of interest but in this example, 40 ng of template human genomic DNA and 5  $\mu$ M of each primer were used. Forty cycles of PCR were performed. The following PCR conditions were used:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 64°C for 30 seconds;
- (7) 95°C for 30 seconds;

(8) Repeat steps 6 and 7 thirty nine (39) times;

(9) 72°C for 5 minutes.

**[0718]** In the first cycle of PCR, the annealing temperature was about the melting temperature of the 3' annealing region of the second primers, which was 37°C. The annealing temperature in the second cycle of PCR was about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer, which was 57°C. The annealing temperature in the third cycle of PCR was about the melting temperature of the entire sequence of the second primer, which was 64°C. The annealing temperature for the remaining cycles was 64°C. Escalating the annealing temperature from TM1 to TM2 to TM3 in the first three cycles of PCR greatly improves specificity. These annealing temperatures are representative, and the skilled artisan will understand the annealing temperatures for each cycle are dependent on the specific primers used.

**[0719]** The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results.

#### **Purification of Fragment of Interest**

**[0720]** The PCR products were separated from the genomic template DNA. One half of the PCR reaction was transferred to a well of a Streptawell, transparent, High-Bind plate from Roche Diagnostics GmbH (catalog number 1 645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog). The first primers contained a 5' biotin tag so the PCR products bound to the Streptavidin coated wells while the genomic template DNA did not. The streptavidin binding reaction was performed using a Thermomixer (Eppendorf) at 1000 rpm for 20 min. at 37°C. Each well was aspirated to remove unbound material, and washed three times with 1X PBS, with gentle mixing (Kandpal et al., Nucl. Acids Res. 18:1789-1795 (1990); Kaneoka et al., Biotechniques 10:30-34 (1991); Green et al., Nucl. Acids Res. 18:6163-6164 (1990)).

#### **Restriction Enzyme Digestion of Isolated Fragments**

**[0721]** The purified PCR products were digested with the restriction enzyme BsmF I, which binds to the recognition site incorporated into the PCR products from the second primer. The digests were performed in the Streptawells following the instructions supplied with the restriction enzyme. After digestion, the wells were washed three times with PBS to remove the cleaved fragments.

### **Incorporation of Labeled Nucleotide**

**[0722]** The restriction enzyme digest with BsmF I yielded a DNA fragment with a 5' overhang, which contained the SNP site or locus of interest and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase.

**[0723]** As discussed in detail in Example 6, the sequence of both alleles of a SNP can be determined by using one labeled nucleotide in the presence of the other unlabeled nucleotides. The following components were added to each fill in reaction: 1 µl of fluorescently labeled ddGTP, 0.5 µl of unlabeled ddNTPs (40 µM), which contained all nucleotides except guanine, 2 µl of 10X sequenase buffer, 0.25 µl of Sequenase, and water as needed for a 20µl reaction. The fill in reaction was performed at 40°C for 10 min. Non-fluorescently labeled ddNTP was purchased from Fermentas Inc. (Hanover, MD). All other labeling reagents were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565).

**[0724]** After labeling, each Streptawell was rinsed with 1X PBS (100 µl) three times. The "filled in" DNA fragments were then released from the Streptawells by digestion with the restriction enzyme EcoRI, according to the manufacturer's instructions that were supplied with the enzyme. Digestion was performed for 1 hour at 37 °C with shaking at 120 rpm.

### **Detection of the Locus of Interest**

**[0725]** The samples were loaded into a lane of a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691). The samples were electrophoresed into the gel at 3000 volts for 3 min. The gel was run for 3 hours on a sequencing apparatus (Hoefer SQ3 Sequencer). The gel was removed from the apparatus and scanned on the Typhoon 9400 Variable Mode Imager. The incorporated labeled nucleotide was detected by fluorescence. A box was drawn around each band and the intensity of the band was calculated using the Typhoon 9400 Variable Mode Imager software.

**[0726]** Below, a schematic of the 5' overhang for SNP TSC0607185 is shown. The entire DNA sequence is not reproduced, only the portion to demonstrate the overhang (where R indicates the variable site).

C	C	T	R	TGTC 3'
				ACAG 5'
4	3	2	1	Overhang position

**[0727]** The observed nucleotides at the variable site for TSC0607185 on the 5' sense strand (here depicted as the top strand) are cytosine and thymidine (depicted here as R). In this case, the second primer anneals from the locus of interest, which allows the fill-in reaction to occur on the anti-sense strand (depicted here as the bottom strand). The antisense strand will be filled in with guanine or adenine.

**[0728]** The second position in the 5' overhang is thymidine, which is complementary to adenine, and the third position in the overhang corresponds to cytosine, which is complementary to guanine.

Fluorescently labeled ddGTP in the presence of unlabeled dCTP, dTTP, and dATP was used to determine the sequence of both alleles. After the fill-in reaction, the following DNA molecules were generated:

C	C	T	C	TGTC 3'	Allele 1
			G*	ACAG 5'	
4	3	2	1	Overhang position	
C	C	T	T	TGTC 3'	Allele 1
	G*	A	A	ACAG 5'	
4	3	2	1	Overhang position	

**[0729]** The overhang generated by BsmF I cutting at 11/15 from the recognition site at TSC0607185 is depicted below:

C	T	R	T	GTC 3' 11/15
				CAG 5'
3	2	1	0	Overhang position

**[0730]** As labeled ddGTP is used for the fill-in reaction, no new signal will be generated from the molecules cut 11/15 from the recognition site. Position 0 complementary to the overhang was filled in with unlabeled dATP. Only signals generated from molecules filled in with labeled ddGTP at position 1 complementary to the overhang or molecules filled in with labeled ddGTP at position 3 complementary to the overhang were seen.

**[0731]** Five of the twenty-four individuals were heterozygous for SNP TSC0607185. As shown in FIG. 15, two bands were detected. The lower molecular weight band corresponded to DNA molecules filled in with ddGTP at position 1 complementary to the overhang. The higher molecular weight band corresponded to DNA molecules filled in with ddGTP at position 3 complementary to the overhang.

**[0732]** The ratio of the two alleles was calculated for each of the five heterozygous samples (see Table XVI). The average ratio of allele 2 to allele 1 was 1.000 with a standard deviation of 0.044. Thus, the allele ratio at SNP TSC0607185 was highly consistent. The experimentally calculated allele ratio for a particular SNP is hereinafter referred to as the "p" value of the SNP. Analysis of SNP TSC0607185



consistently will provide an allele ratio of 1:1, provided that the number of genomes analyzed is of sufficient quantity that no error is generated from statistical sampling.

**[0733]** If the sample contained a low number of genomes, it is statistically possible that the primers will anneal to one chromosome over another chromosome. For example, if the sample contains 40 genomes, which corresponds to a total of 40 chromosomes of allele 1 and 40 chromosomes of allele 2, the primers may anneal to 40 chromosomes of allele 1 but only 35 chromosome of allele 2. This would cause allele 1 to be amplified preferentially to allele 2, which would alter the ratio of allele 1 to allele 2. This problem is eliminated by having a sufficient number of genomes in the sample.

**[0734]** SNP TSC0607185 represents a SNP where the difference in the nucleotide at the variable site does not affect the PCR reaction, or digestion with the restriction enzyme or the fill-in reaction. The use of one nucleotide labeled with one fluorescent dye assures that the bands for one allele can be accurately compared to the bands for the second allele. There is no added complication of having to compare between two different lanes, or having to correct for the quantum coefficients of the dyes. Additionally, any effect from the alternate cutting properties of the type IIS restriction enzymes has been removed.

**TABLE XVI. Ratio of allele 2 to allele 1 at SNPs TSC0607185 and TSC1130902.**

SNP TSC0607185				SNP TSC1130902		
Sample	Allele 1	Allele 2	Allele2/Allele 1	Allele 1	Allele 2	Allele2/Allele 1
1	2382	2313	0.971033	5877	4433	0.754296
2	1581	1533	0.969639	3652	2695	0.737952
3	1795	1879	1.046797	5416	3964	0.730059
4	1921	1855	0.965643	3493	2663	0.762382
5	1618	1701	1.051298	3894	2808	0.721109
Average			1.000882			0.74116
STD			0.044042			0.017018

**[0735]** Below, a schematic of the 5' overhang for SNP TSC1130902 is shown. The entire DNA sequence is not reproduced, only the portion to demonstrate the overhang (where R indicates the variable site).

	5' TTCAT				
	3' AAGTA	R	T	C	C
Overhang position		1	2	3	4

**[0736]** The observed nucleotides for TSC1130902 on the 5' sense strand (here depicted as the top strand) are adenine and guanine. The second position in the overhang corresponds to a thymidine, and the

third position in the overhang corresponds to cytosine, which is complementary to guanine. Fluorescently labeled ddGTP in the presence of unlabeled dCTP, dTTP, and dATP was used to determine the sequence of both alleles. After the fill-in reaction, the following DNA molecules were generated:

Allele 1	5' TTCAT	G*				
	3' AAGTA	C	T	C	C	
Overhang position		1	2	3	4	
Allele 2	5' TTCAT	A	A	G*		
	3' AAGTA	T	T	C	C	
Overhang position		1	2	3	4	

**[0737]** As shown in FIG. 15, two bands were detected. The lower molecular weight band corresponded to DNA molecules filled in with labeled ddGTP at position 1 complementary to the overhang (the G allele). The higher molecular weight band, separated by a single base from the lower band, corresponded to DNA molecules filled in with ddGTP at position 3 complementary to the overhang (the A allele).

**[0738]** Five of the twenty-four individuals were heterozygous for SNP TSC1130902. As seen in FIG. 15, the band corresponding to allele 1 was more intense than the band corresponding to allele 2. This was seen for each of the five individuals. The actual intensity of the band corresponding to allele 1 varied from individual to individual but it was always more intense than the band corresponding to allele 2. For the five individuals, the average ratio of allele 2 to allele 1 was 0.74116, with a standard deviation of 0.017018.

**[0739]** Template DNA was prepared from five different individuals. Separate PCR reactions, separate restriction enzyme digestions, and separate fill-in reactions were performed. However, for each template DNA, the ratio of allele 2 to allele 1 was about 0.75. The “p” value for this SNP was highly consistent.

**[0740]** For example, for SNP TSC1130902, the “p” value was 0.75. Any deviation from this value, provided the sample contains an adequate number of genomes to remove statistical sampling errors, will indicate that there is an abnormal copy number of chromosome 13. If there is an additional copy of allele 2, the “p” value will be higher than the expected 0.75. However, if there is an addition copy of allele 1, the “p” value will be lower than the expected 0.75. With the “p” value quantitated for a particular SNP, that SNP can be used to determine the presence or absence of a chromosomal abnormality. An accurate “p” value measured for a single SNP will be sufficient to detect the presence of a chromosomal abnormality.

**[0741]** There are several possible explanations for why the ratio of one allele to the other allele at some SNPs varies from the theoretically expected ratio of 1:1. First, it is possible that the DNA polymerase incorporates one nucleotide faster than the other nucleotide. As the alleles are being amplified by PCR, even a slight preference for one nucleotide over the other may cause variation from the expected 1:1 ratio. This potential preference for one nucleotide over the other is not seen during the fill-in reaction because a single nucleotide labeled with one dye is used.

**[0742]** It is also possible that the variable nucleotide at the SNP site influences the rate of denaturation of the two alleles. If allele 1 contains a guanine and allele 2 contains an adenine, the difference between the strength of the bonds for these nucleotides may affect the rate at which the DNA strands separate. Again, it is important to mention that the alleles are being amplified by PCR so very subtle differences can make a large impact on the final result. It is also possible that the variable nucleotide at the SNP site influences the rate at which the two strands anneal after separation.

**[0743]** Alternatively, it is possible that the type IIS restriction enzyme cuts one allele preferentially to the other allele. As discussed in detail above, type IIS restriction enzymes cut at a distance from the recognition site. It is possible that the variable nucleotide at the SNP site influences the efficiency of the restriction enzyme digestion. It is possible that at some SNPs the restriction enzyme cuts one allele with an efficiency of 100%, while it cuts the other allele with an efficiency of 90%.

**[0744]** However, the fact that the ratio of allele 1 to allele 2 deviates from the theoretically expected ratio of 1:1, does not influence or reduce the utility of that SNP. As demonstrated above, the “p” value for each SNP is consistent among different individuals.

**[0745]** The “p” value for any SNP can be calculated by analyzing the template DNA of any number of heterozygous individuals including but not limited to 1-10, 11-20, 21-30, 31-40, 41-50, 51-60, 61-70, 71-80, 81-90, 91-100, 101-110, 111-120, 121-130, 131-140, 141-150, 151-160, 161-170, 171-180, 181-190, 191-200, 201-210, 211-220, 221-230, 231-240, 241-250, 251-260, 261-270, 271-280, 281-290, 291-300, and greater than 300.

**[0746]** The methods described herein allow the “p” value for any SNP to be determined. It is possible that some SNPs will behave more consistently than other SNPs. In the human genome, there are over 3 million SNPs; it is not possible to speculate on how each SNP will behave. The “p” value for each SNP will have to be experimentally determined. The methods described herein allow identification of SNPs that have highly consistent, and reproducible “p” values.

## EXAMPLE 10

[0747] As discussed in Example 9, the ratio of one allele to the other allele at a particular SNP may vary from the theoretically expected ratio of 50:50. These SNPs can be used to detect the presence of additional chromosomes provided that the ratio of one allele to the other allele remains linear in individuals with chromosomal disorders. For example, at SNP X if the percentage of allele 1 to allele 2 is 75:25, the expected percentage of allele 1 to allele 2 for an individual with Down's syndrome must be properly adjusted to reflect the variation from the expected percentage at this SNP.

[0748] The percentage of allele 1 to allele 2 for SNP TSC0108992 on chromosome 21 was calculated using template DNA from four normal individuals and template DNA from an individual with Down's syndrome. As demonstrated below, the percentage of one allele to the other allele was consistent and remained linear in an individual with Down's syndrome.

### Preparation of Template DNA

[0749] DNA was obtained from four individuals with a normal genetic karyotype and an individual identified as having an extra copy of chromosome 21 (Down's syndrome). Informed consent was obtained from all individuals. Informed consent also was obtained from the parents of the individual with Down's syndrome.

[0750] From each individual, a 9 ml blood sample was collected into a sterile tube (Fischer Scientific, 9 ml EDTA Vacuette tubes, catalog number NC9897284). Template DNA was isolated using the QIAmp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183). The template DNA was isolated as per instructions included in the kit.

### Design of Primers

[0751] SNP TSC0108992 was amplified using the following primer set:

First primer:

5' CTACTGAGGGCTCGTAGATCCCAATTCCTTCCCAAGCT 3' (SEQ ID NO: 293)

Second primer:

5' AATCCTGCTTTAGGGACCATGCTGGTGGGA 3' (SEQ ID NO: 294)

[0752] The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0753] SNP TSC0108992 was amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq

Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction can be optimized for each locus of interest. In this example, 50 ng of template human genomic DNA and 5  $\mu$ M of each primer were used. Thirty-eight cycles of PCR were performed. The following PCR conditions were used:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 64°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty-seven (37) times;
- (9) 72°C for 5 minutes.

**[0754]** In the first cycle of PCR, the annealing temperature was about the melting temperature of the 3' annealing region of the second primers, which was 37°C. The annealing temperature in the second cycle of PCR was about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer, which was 57°C. The annealing temperature in the third cycle of PCR was about the melting temperature of the entire sequence of the second primer, which was 64°C. The annealing temperature for the remaining cycles was 64°C. Escalating the annealing temperature from TM1 to TM2 to TM3 in the first three cycles of PCR greatly improves specificity. These annealing temperatures are representative, and the skilled artisan will understand the annealing temperatures for each cycle are dependent on the specific primers used.

**[0755]** The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results.

#### **Purification of Fragment of Interest**

**[0756]** The PCR products were separated from the genomic template DNA. Each PCR reaction was split into two samples and transferred to two separate wells of a Streptawell, transparent, High-Bind plate from Roche Diagnostics GmbH (catalog number 1 645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog). For each PCR reaction, there were two replicates; each in a separate well of a microtiter plate. The first primer contained a 5' biotin tag so the PCR products bound to the Streptavidin coated wells while the genomic template DNA did not. The streptavidin binding reaction

was performed using a Thermomixer (Eppendorf) at 1000 rpm for 20 min. at 37°C. Each well was aspirated to remove unbound material, and washed three times with 1X PBS, with gentle mixing (Kandpal et al., Nucl. Acids Res. 18:1789-1795 (1990); Kaneoka et al., Biotechniques 10:30-34 (1991); Green et al., Nucl. Acids Res. 18:6163-6164 (1990)).

#### **Restriction Enzyme Digestion of Isolated Fragments**

**[0757]** The purified PCR products were digested with the restriction enzyme BsmF I, which binds to the recognition site incorporated into the PCR products from the second primer. The digests were performed in the Streptawells following the instructions supplied with the restriction enzyme. After digestion, the wells were washed three times with 1X PBS to remove the cleaved fragments.

#### **Incorporation of Labeled Nucleotide**

**[0758]** The restriction enzyme digest with BsmF I yielded a DNA fragment with a 5' overhang, which contained the SNP site or locus of interest and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase.

**[0759]** As discussed in detail in Example 6, the sequence of both alleles of a SNP can be determined with one labeled nucleotide in the presence of the other unlabeled nucleotides. The following components were added to each fill in reaction: 1 µl of fluorescently labeled ddTTP, 0.5 µl of unlabeled ddNTPs (40 µM), which contained all nucleotides except thymidine, 2 µl of 10X sequenase buffer, 0.25 µl of Sequenase, and water as needed for a 20µl reaction. The fill in reaction was performed at 40°C for 10 min. Non-fluorescently labeled ddNTP was purchased from Fermentas Inc. (Hanover, MD). All other labeling reagents were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565).

**[0760]** After labeling, each Streptawell was rinsed with 1X PBS (100 µl) three times. The "filled in" DNA fragments were then released from the Streptawells by digestion with the restriction enzyme EcoRI, according to the manufacturer's instructions that were supplied with the enzyme. Digestion was performed for 1 hour at 37 °C with shaking at 120 rpm.

#### **Detection of the Locus of Interest**

**[0761]** The samples were loaded into the lanes of a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691). The samples were electrophoresed into the gel at 3000 volts for 3 min. The gel was run for 3 hours on a sequencing apparatus (Hoefer SQ3 Sequencer). The gel was removed from the apparatus and scanned on the

Typhoon 9400 Variable Mode Imager. The incorporated labeled nucleotide was detected by fluorescence. A box was drawn around each band and the intensity of the band was calculated using the Typhoon 9400 Variable Mode Imager software.

**[0762]** Below, a schematic of the 5' overhang for SNP TSC0108992 is shown. The entire DNA sequence is not reproduced, only the portion to demonstrate the overhang (where R indicates the variable site).

				GTCC 3'
G	A	C	R	CAGG 5'
4	3	2	1	Overhang Position

**[0763]** The observed nucleotides for SNP TSC0108992 are adenine and thymidine on the sense strand (here depicted as the top strand). Position 3 of the overhang corresponds to adenine, which is complementary to thymidine. Labeled ddTTP was used in the presence of unlabeled dATP, dCTP, and dGTP. After the fill-in reaction with labeled ddTTP, the following DNA molecules were generated:

	T*	G	A	GTCC 3'	Allele 1
G	A	C	T	CAGG 5'	
4	3	2	1	Overhang Position	
			T*	GTCC 3'	Allele 2
G	A	C	A	CAGG 5'	
4	3	2	1	Overhang Position	

**[0764]** There was no difficulty in comparing the values obtained from allele 1 to allele 2 because one labeled nucleotide was used for the fill-in reaction, and the fill-in reaction for both alleles occurred in a single tube. The alternate cutting properties of BsmF I would not influence this analysis because the 11/15 overhang would be filled in just as the 10/14 overhang. Schematics of the filled-in 11/15 overhangs are depicted below:

T*	G	A	G	TCC 3' 11/15 Allele 1
A	C	T	C	AGG 5'
3	2	1	0	Overhang Position
		T*	G	TCC 3' 11/15 Allele 2
A	C	A	C	AGG 5'
3	2	1	0	Overhang Position

**[0765]** As seen in FIG. 16, two bands were seen for each sample of template DNA. The lower molecular weight band corresponded to the DNA molecules filled in with ddTTP at position one complementary to the overhang, and the higher molecular weight band corresponded to DNA molecules filled in with ddTTP at position 3 complementary to the overhang.

**[0766]** The percentage of allele 2 to allele 1 was highly consistent. (see Table XVII). In addition, for any given individual, the replicates of the PCR reaction showed similar results (see Table XVII). The percentage of allele 2 to allele 1 was calculated by dividing the value of allele 2 by the sum of the values for allele 1 and allele 2 (allele 2/(allele 1+ allele 2)). From four individuals, the average percentage of allele 2 to allele 1 was 0.4773 with a standard deviation of 0.0097. The percentage of allele 2 to allele 1 on template DNA isolated from an individual with Down's syndrome was 0.3086.

**[0767]** The theoretically expected percentage of allele 2 to allele 1 using template DNA from a normal individual is 0.50. However, the experimentally determined percentage was 0.4773. The theoretically expected percentage of allele 2 to allele 1 for an individual with an extra copy of chromosome 21 is 0.33. The experimentally determined percentage of allele 2 to allele 1 for SNP TSC0108992 was 0.3086.

**[0768]** The deviation from the theoretically expected percentage is highly consistent and remains linear. The following formula demonstrates that the percentage of allele 2 to allele 1 at SNP TSC0108992 remains linear even on template DNA obtained from an individual with an extra copy of chromosome 21:

$$\frac{0.47}{0.50} = \frac{X}{0.33}$$

$$X = 0.3102$$

**[0769]** If the percentage of allele 2 to allele 1 using template DNA obtained from a normal individual is determined to be 0.47, then the percentage of allele 2 to allele 1 using template DNA from an individual with Down's syndrome should be 0.3102. The experimentally determined ratio was 0.3086, with a standard deviation of 0.00186. There is no difference between the predicted percentage and the experimentally determined percentage of allele 2 to allele 1 on template DNA from an individual with Down's syndrome.

**[0770]** The percentage of one allele to the other allele at a particular SNP is highly consistent, reproducible, and linear. This demonstrates that any SNP, regardless of the calculated percentage for one allele to another, can be used to determine the presence or absence of a chromosomal disorder.



**TABLE XVII.****Percentage of Allele 2 to Allele 1 at SNP TSC0108992.**

<b>Sample</b>	<b>Allele 2</b>	<b>Allele 1</b>	<b>2/(2+1)</b>
<b>1A</b>	9568886	10578972	0.474933
<b>1B</b>	8330864	9221381	0.474632
<b>2A</b>	9801053	10345444	0.486489
<b>2B</b>	8970942	9603102	0.482983
<b>3A</b>	8676718	9211085	0.485063
<b>3B</b>	10847024	11420943	0.487113
<b>4A</b>	10512420	12227107	0.462297
<b>4B</b>	7883584	9055289	0.465414
		MEAN	0.477366
		STDEV	0.009654
<b>DS</b>	6797400	15138959	0.309869
<b>DS</b>	6025753	13586890	0.307238
		MEAN	0.308554
		STDEV	0.00186

**EXAMPLE 11**

**[0771]** The percentage of allele 2 to allele 1 for a particular SNP is highly consistent. Statistically significant deviation from the experimentally determined ratio indicates the presence of a chromosomal abnormality. Below, the percentage of allele 2 to allele 1 at SNP TSC0108992 on chromosome 21 was calculated using template DNA from a normal individual and template DNA from an individual with Down's syndrome. Mixtures containing various amounts of normal DNA and Down's syndrome DNA were prepared and analyzed in a blind fashion.

**Preparation of Template DNA**

**[0772]** DNA was obtained from an individual with a normal genetic karyotype and an individual identified as having an extra copy of chromosome 21 (Down's syndrome). Informed consent was obtained from both individuals. Informed consent also was obtained from the parents of the individual with Down's syndrome.

**[0773]** From each individual, a 9 ml blood sample was collected into a sterile tube (Fischer Scientific, 9 ml EDTA Vacuette tubes, catalog number NC9897284). Template DNA was isolated using

the QIAmp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183). The template DNA was isolated as per instructions included in the kit.

### **Mixtures of Template DNA**

**[0774]** The template DNA from the individual with the normal karyotype and the template DNA from the individual with an extra copy of chromosome 21 were diluted to a concentration of 10 ng/μl. Four mixtures of normal template DNA and Down's syndrome template DNA were made in the following fashion:

- Mixture 1: 32 μl of Normal DNA + 8 μl of Down's syndrome DNA
- Mixture 2: 28 μl of Normal DNA + 12 μl of Down's syndrome DNA
- Mixture 3: 20 μl of Normal DNA + 20 μl of Down's syndrome DNA
- Mixture 4: 10 μl of Normal DNA + 30 μl of Down's syndrome DNA

**[0775]** Three separate PCR reactions were set up for the normal template DNA and the template DNA from the individual with Down's syndrome. Likewise, for each mixture, three separate PCR reactions were set up.

### **Design of Primers**

**[0776]** SNP TSC0108992 was amplified using the following primer set:

First primer:

5' CTACTGAGGGCTCGTAGATCCCAATTCCTTCCCAAGCT 3' (SEQ ID NO: 293)

Second primer:

5' AATCCTGCTTTAGGGACCATGCTGGTGG 3' (SEQ ID NO: 294)

**[0777]** The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI. The second primer contained the recognition site for the restriction enzyme BsmF I.

**[0778]** SNP TSC0108992 was amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference).

**[0779]** For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction can be optimized for each locus of interest but in this example, 50 ng of template human genomic DNA and 5 μM of each primer were used. Thirty-eight cycles of PCR were performed. The following PCR conditions were used:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 64°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty-seven (37) times;
- (9) 72°C for 5 minutes.

**[0780]** In the first cycle of PCR, the annealing temperature was about the melting temperature of the 3' annealing region of the second primers, which was 37°C. The annealing temperature in the second cycle of PCR was about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer, which was 57°C. The annealing temperature in the third cycle of PCR was about the melting temperature of the entire sequence of the second primer, which was 64°C. The annealing temperature for the remaining cycles was 64°C. Escalating the annealing temperature from TM1 to TM2 to TM3 in the first three cycles of PCR greatly improves specificity. These annealing temperatures are representative, and the skilled artisan will understand the annealing temperatures for each cycle are dependent on the specific primers used.

**[0781]** The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results.

#### **Purification of Fragment of Interest**

**[0782]** The PCR products were separated from the genomic template DNA. Each PCR reaction was split into two samples and transferred to two separate wells of a Streptawell, transparent, High-Bind plate from Roche Diagnostics GmbH (catalog number 1 645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog). For each PCR reaction, there were two replicates, each in a separate well of a microtiter plate. The first primer contained a 5' biotin tag so the PCR products bound to the Streptavidin coated wells while the genomic template DNA did not. The streptavidin binding reaction was performed using a Thermomixer (Eppendorf) at 1000 rpm for 20 min. at 37°C. Each well was aspirated to remove unbound material, and washed three times with 1X PBS, with gentle mixing (Kandpal et al., Nucl. Acids Res. 18:1789-1795 (1990); Kaneoka et al., Biotechniques 10:30-34 (1991); Green et al., Nucl. Acids Res. 18:6163-6164 (1990)).

### **Restriction Enzyme Digestion of Isolated Fragments**

**[0783]** The purified PCR products were digested with the restriction enzyme BsmF I, which binds to the recognition site incorporated into the PCR products from the second primer. The digests were performed in the Streptawells following the instructions supplied with the restriction enzyme. After digestion, the wells were washed three times with 1X PBS to remove the cleaved fragments.

### **Incorporation of Labeled Nucleotide**

**[0784]** The restriction enzyme digest with BsmF I yielded a DNA fragment with a 5' overhang, which contained the SNP site or locus of interest and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase.

**[0785]** As discussed in detail in Example 6, the sequence of both alleles of a SNP can be determined with one labeled nucleotide in the presence of the other unlabeled nucleotides. The following components were added to each fill in reaction: 1  $\mu$ l of fluorescently labeled ddTTP, 0.5  $\mu$ l of unlabeled ddNTPs (40  $\mu$ M), which contained all nucleotides except thymidine, 2  $\mu$ l of 10X sequenase buffer, 0.25  $\mu$ l of Sequenase, and water as needed for a 20 $\mu$ l reaction. The fill in reaction was performed at 40°C for 10 min. Non-fluorescently labeled ddNTP was purchased from Fermentas Inc. (Hanover, MD). All other labeling reagents were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565).

**[0786]** After labeling, each Streptawell was rinsed with 1X PBS (100  $\mu$ l) three times. The "filled in" DNA fragments were then released from the Streptawells by digestion with the restriction enzyme EcoRI, according to the manufacturer's instructions that were supplied with the enzyme. Digestion was performed for 1 hour at 37 °C with shaking at 120 rpm.

### **Detection of the Locus of Interest**

**[0787]** The samples were loaded into the lanes of a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691). The samples were electrophoresed into the gel at 3000 volts for 3 min. The gel was run for 3 hours on a sequencing apparatus (Hoefer SQ3 Sequencer). The gel was removed from the apparatus and scanned on the Typhoon 9400 Variable Mode Imager. The incorporated labeled nucleotide was detected by fluorescence. A box was drawn around each band and the intensity of the band was calculated using the Typhoon 9400 Variable Mode Imager software.

**[0788]** As seen in FIGS. 17 A-F, two bands were seen. The lower molecular weight band corresponded to the DNA molecules filled in with ddTTP at position one complementary to the overhang.

The higher molecular weight band corresponded to DNA molecules filled in with ddTTP at position 3 complementary to the overhang.

**[0789]** The experiment was performed in a blind fashion. The tubes were coded so that it was not known what tube corresponded to what template DNA. After the gels were analyzed, each tube was grouped into the following categories: normal template DNA, Down's syndrome template DNA, 3:1 mixture of Down's syndrome template DNA to normal DNA, 1:1 mixture of normal template DNA to Down's syndrome template DNA, 1:2.3 mixture of Down's syndrome template DNA to normal template DNA, and 1:4 mixture of Down's syndrome template DNA to normal template DNA. Each replicate of each PCR reaction successfully was grouped into the appropriate category, which demonstrates that the method can be used to detect abnormal DNA even if it represents only a small percentage of the total DNA.

**[0790]** The percentage of allele 2 to allele 1 for each replicate of the three PCR reactions from normal template DNA are displayed in Table XVIII (also see FIG. 17A). The average percentage of allele 2 to allele 1 was calculated by dividing the value of allele 2 by the sum of the values for allele 1 and allele 2 ( $\text{allele 2} / (\text{allele 1} + \text{allele 2})$ ), which resulted in an average of 0.50025 with a standard deviation of 0.002897. Thus, allele 1 and allele 2 were present in a ratio of 50:50. While the intensity of the bands varied from one PCR reaction to another (compare reaction 1 with reaction 3), there was no difference in intensity within a PCR reaction. Furthermore, the values obtained for the two replicates of the PCR reactions were very similar. Most of the variation was between PCR reactions and was likely attributable to pipetting errors.

**[0791]** The percentage of allele 2 to allele 1 for each replicate of the three PCR reactions from Down's syndrome template DNA are displayed in Table XVIII (see FIG. 17B). The percentage of allele 2 to allele 1 was calculated by dividing the value of allele 2 by the sum of the values for allele 1 and allele 2 ( $\text{allele 2} / (\text{allele 1} + \text{allele 2})$ ), which resulted in an average of 0.301314 with a standard deviation of 0.012917. It is clear even upon analysis of the gel by the naked eye that allele 1 is present in a higher copy number than allele 2 (see FIG. 17B). Again, most of the variation occurs between PCR reactions and not within the replicate of a PCR reaction. The majority of the statistical variation likely resulted from pipetting errors.

**[0792]** Analysis of a single SNP was sufficient to detect the presence of the chromosomal abnormality. One SNP is sufficient provided that the "p" value of the SNP is known and that there are an adequate number of genomes so that statistical sampling error is not introduced into the analysis. In this experiment, there were approximately 5,000 genomes in each reaction.

**[0793]** The reactions that consisted of a mixture of Down's syndrome template DNA to normal template DNA at a ratio of 3:1 were clearly distinguishable from the normal template DNA, and the other mixtures of DNA (see FIG. 17C). The calculated percentage of allele 2 to allele 1 was 0.319089 with a standard deviation of 0.004346 (see Table XVIII). Likewise, the reactions that consisted of a mixture of Down's syndrome template DNA to normal template DNA at ratios of 1:1, and 1:2.3 were distinguishable (see FIG. 17D and 17E) and the values were statistically significant from all other reactions (see Table XVIII).

**[0794]** As the amount of normal template DNA increased, the percentage of allele 2 to allele 1 increased. With a mixture of Down's syndrome template DNA to normal template DNA of 1:4, the percentage of allele 2 to allele 1 was 0.397642, with a standard deviation of 0.001903 (see FIG 17F). The difference between this value and the value obtained from normal template DNA is statistically significant. Thus, the methods described herein allow the detection of a chromosomal abnormality even when the sample is not a homogeneous sample of abnormal DNA.

**[0795]** As described above, the presence of a small fraction of DNA with an abnormal copy number of chromosomes can be detected even among a large presence of normal DNA. It was clear, even by the naked eye, that as the amount of normal DNA increased and the amount of Down's syndrome DNA decreased, the intensities of the bands that corresponded to alleles 1 and 2 equalized.

**[0796]** The above example analyzed a SNP located on chromosome 21. However, any SNP may be analyzed on any chromosome including but not limited to human chromosomes 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X, and Y and fetal chromosomes 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X, and Y. In addition, chromosomes from non-human organisms can be analyzed using the above methods. Any combination of chromosomes can be analyzed. In the above example, an extra copy of a chromosome was detected. However, the same methods can be used to detect monosomies.

**TABLE XVIII**  
**Percentage of allele 2 to allele 1 at SNP TSC0108992**  
**using normal template DNA and Down's syndrome template DNA.**

	<b>Normal Template DNA</b>		
	<b>Allele 1</b>	<b>Allele 2</b>	<b>2/(2+1)</b>
1A	2602115	2604525	0.500231
1B	2855846	2923860	0.505884
2A	1954765	1941929	0.498353
2B	2084476	2068106	0.498029

3A	2044147	2035719	0.498967
3B	1760291	1760543	0.500036
		Mean	0.50025
		STD	0.002897
		<b>Down's Syndrome</b>	
	<b>Allele 1</b>	<b>Allele 2</b>	<b>2/(2+1)</b>
1A	4046926	1595581	0.282779
1B	4275341	1736260	0.288818
2A	2875698	1299509	0.311244
2B	2453615	1069635	0.303593
3A	3169338	1426643	0.310411
3B	3737440	1687286	0.311036
		Mean	0.301314
		STD	0.012917
		<b>3:1 (Down's: Normal)</b>	
	<b>Allele 1</b>	<b>Allele 2</b>	<b>2/(2+1)</b>
1A	4067623	1980770	0.327487
1B	4058506	1899853	0.318855
2A	2315044	1085860	0.319286
2B	2686984	1243406	0.316357
3A	3880385	1790764	0.315767
3B	3718661	1724189	0.316781
		Mean	0.319089
		STD	0.004346
		<b>1:1 (Down's: Normal)</b>	
	<b>Allele 1</b>	<b>Allele 2</b>	<b>2/(2+1)</b>
1A	3540255	1929840	0.352798
1B	4004085	2161443	0.350569
2A	2358009	1282132	0.35222
2B	2158132	1238377	0.364603
3A	3052330	1648677	0.350707
3B	3852682	2024012	0.344413
		Mean	0.352552
		STD	0.006618
		<b>1:2.3 (Down's: Normal)</b>	
	<b>Allele 1</b>	<b>Allele 2</b>	<b>2/(2+1)</b>
1A	3109326	1942597	0.384526
1B	3392477	2118011	0.38436

2A	2824213	1758428	0.383715
2B	2069889	1249545	0.376433
3A	2335128	1433016	0.380298
3B	2916772	1797965	0.38135
		Mean	0.38178
		STD	0.003128
		<b>1:4 (Down's: Normal)</b>	
	<b>Allele 1</b>	<b>Allele 2</b>	<b>2/(2+1)</b>
1A	3066524	2039636	0.399446
1B	3068284	2038770	0.399207
2A	2325477	1542526	0.398791
2B	2366122	1562218	0.397679
3A	2151205	1403120	0.394764
3B	2397046	1571360	0.395968
		Mean	0.397642
		STD	0.001903

#### EXAMPLE 12

**[0797]** As discussed above in ~~Example~~ Example 9, the ratio for allele 1 to allele 2 at a heterozygous SNP is constant. However, one factor that can influence the ratio of allele 1 to allele 2 at a heterozygous SNP is a low number of genomes. For example, if there are 40 genomes, which means that there are a total of 40 chromosomes of allele 1 and 40 chromosomes of allele 2, it is statistically possible that the primers may anneal to 40 of the chromosomes with allele 1 but only 30 of the chromosomes with allele 2. This will affect the ratio of allele 1 to allele 2, and can erroneously influence the “p” value for a particular SNP.

**[0798]** Typically, whole genomic amplification, which employs degenerate oligonucleotide PCR, is used to increase low quantities of genomic DNA samples. Oligonucleotides of 8, 10, 12, or 14 bases are used to amplify the genome. It is thought that the primers anneal randomly throughout the genome, and will amplify a small genomic DNA sample into hundreds-fold more DNA for genetic analysis.

**[0799]** The methods described herein exploit the fact that typically the whole genome is not of interest. Particular loci of interest located on one chromosome, or on multiple chromosomes or on chromosomes that represent the entire genome are selected for analysis. Even if the loci of interest are located on chromosomes for the entire genome, it is preferential to amplify the region of those chromosomes that contain the loci of interest.



**[0800]** To overcome the limit of a low number of genomes, which is often seen with fetal DNA obtained from the plasma of a pregnant female, a multiplex method can be used to increase the number of genomes. The method described below preferentially amplifies the chromosome or chromosomes that contain the loci of interest.

#### **Preparation of Template DNA**

**[0801]** A 9 ml blood sample was collected into a sterile tube from a human volunteer after informed consent had been granted. (Fischer Scientific, 9 ml EDTA Vacuette tubes, catalog number NC9897284). The tubes were spun at 1000 rpm for ten minutes. The supernatant (the plasma) of each sample was removed, and one milliliter of the remaining blood sample, which is commonly referred to as the “buffy-coat” was transferred to a new tube. One milliliter of 1X PBS was added to each sample. Template DNA was isolated using the QIAmp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183).

#### **Design of Multiplex Primers**

**[0802]** Primers were designed to anneal at various regions on chromosome 21 to increase the copy number of the loci of interest located on chromosome 21. The primers were 12 bases in length. However, primers of any length can be used including but not limited to 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36-45, 46-55, 56-65, 66-75, 76-85, 86-95, 96-105, 106-115, 116-125, and greater than 125 bases. Primers were designed to anneal to both the sense strand and the antisense strand.

**[0803]** Nine SNPs located on chromosome 21 were analyzed: TSC0397235, TSC0470003, TSC1649726, TSC1261039, TSC0310507, TSC1650432, TSC1335008, TSC0128307, and TSC0259757. Any number of SNPs can be analyzed including but not limited to 1-10, 11-20, 21-30, 31-40, 41-50, 51-60, 61-70, 71-80, 81-90, 91-100, 101-200, 201-300, 301-400, 401-500, 501-600, 601-700, 701-800, 801-900, 901-1000, 1001-2000, 2001-3000, 3001-4000, 4001-5000, 5001-6000, 6001-7000, 7001-8000, 8001-9000, 9001-10,000 and greater than 10,000.

**[0804]** For each of the 9 SNPs, a 12 base primer was designed to anneal approximately 130 bases upstream of the loci of interest, and a 12 base primer was designed to anneal approximately 130 bases downstream of the loci of interest (herein referred to as the multiplex primers). The multiplex primers can be designed to anneal at any distance from the loci of interest including but not limited to 10-20, 21-30, 31-40, 41-50, 51-60, 61-70, 71-80, 81-90, 91-100, 101-110, 111-120, 121-130, 131-140, 141-150, 151-160, 161-170, 171-180, 181-190, 191-200, 201-210, 211-220, 221-230, 231-240, 241-250, 251-260, 261-270, 271-280, 281-290, 291-300, 301-310, 311-320, 321-330, 331-340, 341-350, 351-360, 361-370,

371-380, 381-390, 391-400, 401-410, 411-420, 421-430, 431-440, 441-450, 451-460, 461-470, 471-480, 481-490, 491-500, 501-600, 601-700, 701-800, 801-900, 901-1000, 1001-2000, 2001-3000, 3001-4000, 4001-5000, and greater than 5000 bases. In addition, more than one set of multiplex primers can be used for one SNP including but not limited to 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 10-20, 21-30, 31-40, 41-50, and greater than 50.

**[0805]** In addition, 91 sets of forward and reverse primers were used to amplify other regions of chromosome 21, for a total of 100 sets of primers (200 primers in the reaction). These 91 primer sets were used to demonstrate that a large number of primers can be used in a single reaction without producing a large number of non-specific bands. Any number of primers can be used in the reaction including but not limited to 1-10, 11-20, 21-30, 31-40, 41-50, 51-60, 61-70, 71-80, 81-90, 91-100, 101-200, 201-300, 301-400, 401-500, 501-600, 601-700, 701-800, 801-900, 901-1000, 1001-2000, 2001-3000, 3001-4000, 4001-5000, 5001-6000, 6001-7000, 7001-8000, 8001-9000, 9001-10,000, 10,001-20,000, 20,001-30,000 and greater than 30,000.

**[0806]** The multiplex primers were designed to have the same nucleotides at the 3' end of the primer. In this case, the multiplex primers ended in "AA," wherein A indicates adenine. The primers were designed in this manner to minimize primer-dimer formation. However, the primers can terminate in any nucleotides including but not limited to adenine, guanine, cytosine, thymidine, any combination of adenine and guanine, any combination of adenine and cytosine, any combination of adenine and thymidine, any combination of guanine and cytosine, any combination of guanine and thymidine, or any combination of cytosine and thymidine. In addition the multiplex primers can have 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, or more than 10 of the same nucleotides at the 3' end.

**[0807]** The multiplex primers for SNP TSC0397235 were:

Forward Primer:

5' CAAGTGTCTCTAA 3' (SEQ ID NO: 295)

Reverse primer:

5' CAGCTGCTAGAA 3' (SEQ ID NO: 296)

**[0808]** The multiplex primers for SNP TSC0470003 were:

Forward Primer:

5' GGTTGAGGGCAA 3' (SEQ ID NO: 297)

Reverse primer:

5' CACAGCGGGTAA 3' (SEQ ID NO: 298)

**[0809]** The multiplex primers for SNP TSC1649726 were:

Forward Primer:

5' TTGACTTTTAA 3' (SEQ ID NO: 299)

Reverse primer:

5' ACAGAATGGGAA 3' (SEQ ID NO: 300)

**[0810]** The multiplex primers for SNP TSC1261039 were:

Forward Primer:

5' TGCAGGTCACAA 3' (SEQ ID NO: 301)

Reverse primer:

5' TTCTTCTTATAA 3' (SEQ ID NO: 302)

**[0811]** The multiplex primers for SNP TSC0310507 were:

Forward Primer:

5' AGGACAACCTAA 3' (SEQ ID NO: 303)

Reverse primer:

5' TGGTGTTTCAGAA 3' (SEQ ID NO: 304)

**[0812]** The multiplex primers for SNP TSC1650432 were:

Forward Primer:

5' TCAGCATATGAA 3' (SEQ ID NO: 305)

Reverse primer:

5' GTTGCCACACAA 3' (SEQ ID NO: 306)

**[0813]** The multiplex primers for SNP TSC1335008 were:

Forward Primer:

5' CCCAGCTAGCAA 3' (SEQ ID NO: 307)

Reverse primer:

5' GGGTCACTGTAA 3' (SEQ ID NO: 308)

**[0814]** The multiplex primers for SNP TSC0128307 were:

Forward Primer:

5' TTAAATACCCAA 3' (SEQ ID NO: 309)

Reverse primer:

5' TTAGGAGGTAA 3' (SEQ ID NO: 310)

**[0815]** The multiplex primers for SNP TSC0259757 were:

Forward Primer:

5' ACACAGAATCAA 3' (SEQ ID NO: 311)

Reverse primer:

5' CGCTGAGGTCAA 3' (SEQ ID NO: 312)

**[0816]** Ninety-one (91) additional sets of primers, which annealed to various regions along chromosome 21, were included in the reaction:

Set 1:

Forward Primer:

5' AAGTAGAGTCAA 3' (SEQ ID NO: 313)

Reverse primer:

5' CTTCCCATGGAA 3' (SEQ ID NO: 314)

Set 2:

Forward Primer:

5' TTGGTTATTAAA 3' (SEQ ID NO: 315)

Reverse primer:

5' CAACTTACTGAA 3' (SEQ ID NO: 316)

Set 3:

Forward Primer:

5' CACTAAGTGAAA 3' (SEQ ID NO: 317)

Reverse primer:

5' CTCACCTGCCAA 3' (SEQ ID NO: 318)

Set 4:

Forward Primer:

5' ATGCATATATAA 3' (SEQ ID NO: 319)

Reverse primer:

5' AGAGATCAGCAA 3' (SEQ ID NO: 320)

Set 5:

Forward Primer:

5' TATATTTTTCAA 3' (SEQ ID NO: 321)

Reverse primer:

5' CAGAAAGCAGAA 3' (SEQ ID NO: 322)

Set 6:

Forward Primer:

5' GTATTGGGTAA 3' (SEQ ID NO: 323)

Reverse primer:

5' CTGACCCAGGAA 3' (SEQ ID NO: 324)

Set 7:

Forward Primer:

5' CAGTTTTCCCAA 3' (SEQ ID NO: 325)

Reverse primer:

5' AGGGCACAGGAA 3' (SEQ ID NO: 326)

Set 8:

Forward Primer:

5' GTATCAGAGGAA 3' (SEQ ID NO: 327)

Reverse primer:

5' GCATGAAAAGAA 3' (SEQ ID NO: 328)

Set 9:

Forward Primer:

5' GATTGACAGAA 3' (SEQ ID NO: 329)

Reverse primer:

5' TACAGTTTACAA 3' (SEQ ID NO: 330)

Set 10:

Forward Primer:

5' TGTGATTTTAA 3' (SEQ ID NO: 331)

Reverse primer:

5' TTATGTTCTCAA 3' (SEQ ID NO: 332)

Set 11:

Forward Primer:

5' CAAGTACTTGAA 3' (SEQ ID NO: 333)

Reverse primer:

5' CTTGTGTGGCAA 3' (SEQ ID NO: 334)

Set 12:

Forward Primer:

5' AGACTTCTGCAA 3' (SEQ ID NO: 335)

Reverse primer:

5' GTTGTCTTTCAA 3' (SEQ ID NO: 336)

Set 13:

Forward Primer:

5' GGGACACTCCAA 3' (SEQ ID NO: 337)

Reverse primer:

5' ATTATTATTCAA 3' (SEQ ID NO: 338)

Set 14:

Forward Primer:

5' ACATGATGACAA 3' (SEQ ID NO: 339)

Reverse primer:

5' TCAATTATAGAA 3' (SEQ ID NO: 340)

Set 15:

Forward Primer:

5' CTATGGGCTGAA 3' (SEQ ID NO: 341)

Reverse primer:

5' TGTGTGCCTGAA 3' (SEQ ID NO: 342)

Set 16:

Forward Primer:

5' CCATTTGTTGAA 3' (SEQ ID NO: 343)

Reverse primer:

5' TCTCCATCAAAA 3' (SEQ ID NO: 344)

Set 17:

Forward Primer:

5' AATGCTGACAAA 3' (SEQ ID NO: 345)

Reverse primer:

5' TTTCATGTCCAA 3' (SEQ ID NO: 346)

Set 18:

Forward Primer:

5' GGCCTCTTGAA 3' (SEQ ID NO: 347)

Reverse primer:

5' TCATTTTTTGAA 3' (SEQ ID NO: 348)

Set 19:

Forward Primer:

5' GGACTACCATAA 3' (SEQ ID NO: 349)

Reverse primer:

5' AGTCACTCAGAA 3' (SEQ ID NO: 350)

Set 20:

Forward Primer:

5' CCTTGGCAGGAA 3' (SEQ ID NO: 351)

Reverse primer:

5' TTTCTGGTAGAA 3' (SEQ ID NO: 352)

Set 21:

Forward Primer:

5' CCCCCCCCCGAA 3' (SEQ ID NO: 353)

Reverse primer:

5' GCCCAGGCAGAA 3' (SEQ ID NO: 354)

Set 22:

Forward Primer:

5' GAATGCGAAGAA 3' (SEQ ID NO: 355)

Reverse primer:

5' TTAGGTAGAGAA 3' (SEQ ID NO: 356)

Set 23:

Forward Primer:

5' TGCTTTGGTCAA 3' (SEQ ID NO: 357)

Reverse primer:

5' GCCCATTAATAA 3' (SEQ ID NO: 358)

Set 24:

Forward Primer:

5' TGAGATCTTTAA 3' (SEQ ID NO: 359)

Reverse primer:

5' CAGTTTGTTC AA 3' (SEQ ID NO: 360)

Set 25:

Forward Primer:

5' GCTGGGCAAGAA 3' (SEQ ID NO: 361)

Reverse primer:

5' AGTCAAAGTCAA 3' (SEQ ID NO: 362)

Set 26:

Forward Primer:

5' TCTCTGCAGTAA 3' (SEQ ID NO: 363)

Reverse primer:

5' TGAATAACTTAA 3' (SEQ ID NO: 364)

Set 27:

Forward Primer:

5' CGGTTAGAAAAA 3' (SEQ ID NO: 365)

Reverse primer:

5' CATCCCTTTCAA 3' (SEQ ID NO: 366)

Set 28:

Forward Primer:

5' TCTCTTTCTGAA 3' (SEQ ID NO: 367)

Reverse primer:

5' CTCAGATTGTAA 3' (SEQ ID NO: 368)

Set 29:

Forward Primer:

5' TTTGCACCAGAA 3' (SEQ ID NO: 369)

Reverse primer:

5' GGTTAACATGAA 3' (SEQ ID NO: 370)

Set 30:

Forward Primer:

5' ATTATCAACTAA 3' (SEQ ID NO: 371)

Reverse primer:

5' GCCATTTTGTA 3' (SEQ ID NO: 372)

Set 31:

Forward Primer:

5' GATCTAGATGAA 3' (SEQ ID NO: 373)

Reverse primer:

5' TTAATGTATTAA 3' (SEQ ID NO: 374)

Set 32:



Forward Primer:

5' CTAGGGAGACAA 3' (SEQ ID NO: 375)

Reverse primer:

5' TGGAGGAGACAA 3' (SEQ ID NO: 376)

Set 33:

Forward Primer:

5' CATCACATTTAA 3' (SEQ ID NO: 377)

Reverse primer:

5' GGGGTCCTGCAA 3' (SEQ ID NO: 378)

Set 34:

Forward Primer:

5' CAGTTGTGCTAA 3' (SEQ ID NO: 379)

Reverse primer:

5' TCTGCAGCCTAA 3' (SEQ ID NO: 380)

Set 35:

Forward Primer:

5' GAGTCATTTAAA 3' (SEQ ID NO: 381)

Reverse primer:

5' TCTATGGATTAA 3' (SEQ ID NO: 382)

Set 36:

Forward Primer:

5' CAAAAAGTAGAA 3' (SEQ ID NO: 383)

Reverse primer:

5' AATATACTCCAA 3' (SEQ ID NO: 384)

Set 37:

Forward Primer:

5' CGTCCAGCACAA 3' (SEQ ID NO: 385)

Reverse primer:

5' GGATGGTGAGAA 3' (SEQ ID NO: 386)

Set 38:

Forward Primer:

5' TCTCCTTTGTAA 3' (SEQ ID NO: 387)

Reverse primer:

5' TCGTTATTTCAA 3' (SEQ ID NO: 388)

Set 39:

Forward Primer:

5' GATTTTATAGAA 3' (SEQ ID NO: 389)

Reverse primer:

5' AGACATAAGCAA 3' (SEQ ID NO: 390)

Set 40:

Forward Primer:

5' TTCACCTCACAA 3' (SEQ ID NO: 391)

Reverse primer:

5' GGATTGCTTGAA 3' (SEQ ID NO: 392)

Set 41:

Forward Primer:

5' ACTGCATGTGAA 3' (SEQ ID NO: 393)

Reverse primer:

5' TTTATCACAGAA 3' (SEQ ID NO: 394)

Set 42:

Forward Primer:

5' TCAGTAACACAA 3' (SEQ ID NO: 395)

Reverse primer:

5' TACATCTTTGAA 3' (SEQ ID NO: 396)

Set 43:

Forward Primer:

5' TTGTTTCAGTAA 3' (SEQ ID NO: 397)

Reverse primer:

5' TATGAGCATCAA 3' (SEQ ID NO: 398)

Set 44:

Forward Primer:

5' CTCAGCAGGCAA 3' (SEQ ID NO: 399)

Reverse primer:

5' ACCCCTGTATAA 3' (SEQ ID NO: 400)

Set 45:

Forward Primer:

5' TCTGCTCAGCAA 3' (SEQ ID NO: 401)

Reverse primer:

5' GTTCTTTTTTAA 3' (SEQ ID NO: 402)

Set 46:

Forward Primer:

5' GTGATAATCCAA 3' (SEQ ID NO: 403)

Reverse primer:

5' GAGCCCTCAGAA 3' (SEQ ID NO: 404)

Set 47:

Forward Primer:

5' TTTATTGGTTAA 3' (SEQ ID NO: 405)

Reverse primer:

5' GGTACTGGGCAA 3' (SEQ ID NO: 406)

Set 48:

Forward Primer:

5' AGTGTTTTTCAA 3' (SEQ ID NO: 407)

Reverse primer:

5' TGTTATTGGTAA 3' (SEQ ID NO: 408)

Set 49:

Forward Primer:

5' GCGCATTCACAA 3' (SEQ ID NO: 409)

Reverse primer:

5' AAACAAAAGCAA 3' (SEQ ID NO: 410)

Set 50:

Forward Primer:

5' TATATGATAGAA 3' (SEQ ID NO: 411)

Reverse primer:

5' TCCCAGTTCCAA 3' (SEQ ID NO: 412)

Set 51:

Forward Primer:

5' AAAGCCCATAAA 3' (SEQ ID NO: 413)

Reverse primer:

5' TGTCATCCACAA 3' (SEQ ID NO: 414)

Set 52:

Forward Primer:

5' TTGTGAATGCAA 3' (SEQ ID NO: 415)

Reverse primer:

5' GTATTCATACAA 3' (SEQ ID NO: 416)

Set 53:

Forward Primer:

5' TGACATAGGGAA 3' (SEQ ID NO: 417)

Reverse primer:

5' AGCAAATTGCAA 3' (SEQ ID NO: 418)

Set 54:

Forward Primer:

5' AGTAGATGTAA 3' (SEQ ID NO: 419)

Reverse primer:

5' AAAAGATAATAA 3' (SEQ ID NO: 420)

Set 55:

Forward Primer:

5' ACCTCATGGGAA 3' (SEQ ID NO: 421)

Reverse primer:

5' TGGTCGACCTAA 3' (SEQ ID NO: 422)

Set 56:

Forward Primer:

5' TTTGCATGGTAA 3' (SEQ ID NO: 423)

Reverse primer:

5' GCGGCTGCCGAA 3' (SEQ ID NO: 424)

Set 57:

Forward Primer:

5' TCAGGAGTCTAA 3' (SEQ ID NO: 425)

Reverse primer:

5' GCCTACCAGGAA 3' (SEQ ID NO: 426)

Set 58:

Forward Primer:

5' ATCTTCTGTAA 3' (SEQ ID NO: 427)

Reverse primer:

5' AGGTAAGGACAA 3' (SEQ ID NO: 428)

Set 59:

Forward Primer:

5' TGCTTTGAGGAA 3' (SEQ ID NO: 429)

Reverse primer:

5' AACAGTTTAAA 3' (SEQ ID NO: 430)

Set 60:

Forward Primer:

5' TTAAATGTTAA 3' (SEQ ID NO: 431)

Reverse primer:

5' ATAGAAAATCAA 3' (SEQ ID NO: 432)

Set 61:

Forward Primer:

5' GTGTTGTGTTAA 3' (SEQ ID NO: 433)

Reverse primer:

5' GAGGACCTCGAA 3' (SEQ ID NO: 434)

Set 62:

Forward Primer:

5' AGAGGCTGAGAA 3' (SEQ ID NO: 435)

Reverse primer:

5' GGTATTTATTAA 3' (SEQ ID NO: 436)

Set 63:

Forward Primer:

5' ATTTATCTGGAA 3' (SEQ ID NO: 437)

Reverse primer:

5' AGTGCAAACATA 3' (SEQ ID NO: 438)

Set 64:

Forward Primer:

5' TGAACACCTTAA 3' (SEQ ID NO: 439)

Reverse primer:

5' AATTTTTTCTAA 3' (SEQ ID NO: 440)

Set 65:

Forward Primer:

5' TTACTATTATAA 3' (SEQ ID NO: 441)

Reverse primer:

5' TGCTATAGTGAA 3' (SEQ ID NO: 442)

Set 66:

Forward Primer:

5' TGGACTATGGAA 3' (SEQ ID NO: 443)

Reverse primer:

5' CTGCAGTCCGAA 3' (SEQ ID NO: 444)

Set 67:

Forward Primer:

5' GCTACTGCCCAA 3' (SEQ ID NO: 445)

Reverse primer:

5' TCACATGGTGAA 3' (SEQ ID NO: 446)

Set 68:

Forward Primer:

5' GTGGCTCTGGAA 3' (SEQ ID NO: 447)

Reverse primer:

5' GAATTCCATTAA 3' (SEQ ID NO: 448)

Set 69:

Forward Primer:

5' TGGGGTGTCCAA 3' (SEQ ID NO: 449)

Reverse primer:

5' GCAAGCTCCGAA 3' (SEQ ID NO: 450)

Set 70:

Forward Primer:

5' ATGTTTTTTCAA 3' (SEQ ID NO: 451)

Reverse primer:

5' AGATCTGTTGAA 3' (SEQ ID NO: 452)

Set 71:

Forward Primer:

5' AAGTGCTGTGAA 3' (SEQ ID NO: 453)

Reverse primer:

5' ACTTTTTTGGAA 3' (SEQ ID NO: 454)

Set 72:

Forward Primer:

5' AATCGGCAGGAA 3' (SEQ ID NO: 455)

Reverse primer:

5' GGCATGTCACAA 3' (SEQ ID NO: 456)

Set 73:

Forward Primer:

5' AGGAAGAAAGAA 3' (SEQ ID NO: 457)

Reverse primer:

5' CAGTTTCACCAA 3' (SEQ ID NO: 458)

Set 74:

Forward Primer:

5' CACAGAATTAA 3' (SEQ ID NO: 459)

Reverse primer:

5' AAGAATAAGTAA 3' (SEQ ID NO: 460)

Set 75:

Forward Primer:

5' GGGATAGTACAA 3' (SEQ ID NO: 461)

Reverse primer:

5' TTCCCATGATAA 3' (SEQ ID NO: 462)

Set 76:

Forward Primer:

5' TGATTAGTTGAA 3' (SEQ ID NO: 463)

Reverse primer:

5' GCATTCAGTGAA 3' (SEQ ID NO: 464)

Set 77:

Forward Primer:

5' AGGGAATATTAA 3' (SEQ ID NO: 465)

Reverse primer:

5' GACCTTAGGTAA 3' (SEQ ID NO: 466)

Set 78:

Forward Primer:

5' TTCTTTTCACAA 3' (SEQ ID NO: 467)

Reverse primer:

5' CCAAATAAGAA 3' (SEQ ID NO: 468)

Set 79:

Forward Primer:

5' GTGCTCTTAGAA 3' (SEQ ID NO: 469)

Reverse primer:

5' ATGAGTTTAGAA 3' (SEQ ID NO: 470)

Set 80:

Forward Primer:

5' ATGAGCATAGAA 3' (SEQ ID NO: 471)

Reverse primer:

5' GACAAATGAGAA 3' (SEQ ID NO: 472)

Set 81:

Forward Primer:

5' AAACCCAGAGAA 3' (SEQ ID NO: 473)

Reverse primer:

5' CCTCACACAGAA 3' (SEQ ID NO: 474)

Set 82:

Forward Primer:

5' CACACTGTGGAA 3' (SEQ ID NO: 475)

Reverse primer:

5' CACTGTACCCAA 3' (SEQ ID NO: 476)

Set 83:

Forward Primer:



5' GTAGTATTTCAA 3' (SEQ ID NO: 477)

Reverse primer:

5' TGGATACACTAA 3' (SEQ ID NO: 478)

Set 84:

Forward Primer:

5' CCCATGATTCAA 3' (SEQ ID NO: 479)

Reverse primer:

5' TCATAGGAGGAA 3' (SEQ ID NO: 480)

Set 85:

Forward Primer:

5' AGGAAAGAGAAA 3' (SEQ ID NO: 481)

Reverse primer:

5' ATATGGTGATAA 3' (SEQ ID NO: 482)

Set 86:

Forward Primer:

5' GATGCCATCCAA 3' (SEQ ID NO: 483)

Reverse primer:

5' ATACTATTTCAA 3' (SEQ ID NO: 484)

Set 87:

Forward Primer:

5' GTGTGCATGGAA 3' (SEQ ID NO: 485)

Reverse primer:

5' AGGTGTTGAGAA 3' (SEQ ID NO: 486)

Set 88:

Forward Primer:

5' CAGCCTGGGCAA 3' (SEQ ID NO: 487)

Reverse primer:

5' GGAGCTCTACAA 3' (SEQ ID NO: 488)

Set 89:

Forward Primer:

5' AACTAAGGTAA 3' (SEQ ID NO: 489)

Reverse primer:

5' AACTTATGTAA 3' (SEQ ID NO: 490)

Set 90:

Forward Primer:

5' ATCTCAACAGAA 3' (SEQ ID NO: 491)

Reverse primer:

5' TAACAATGTGAA 3' (SEQ ID NO: 492)

Set 91:

Forward Primer:

5' AAGGATCAGGAA 3' (SEQ ID NO: 493)

Reverse primer:

5' CTCAAGTCTTAA 3' (SEQ ID NO: 494)

### **Multiplex PCR**

**[0817]** Regions on chromosome 21 surrounding SNPs TSC0397235, TSC0470003, TSC1649726, TSC1261039, TSC0310507, TSC1650432, TSC1335008, TSC0128307, and TSC0259757 were amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). This PCR reaction used primers that annealed approximately 130 bases upstream and downstream of the loci of interest. It was used to increase the number of copies of the loci of interest to eliminate any errors that may result from a low number of genomes.

**[0818]** For increased specificity, a "hot-start" PCR reaction was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction can be optimized for each locus of interest. In this example, 15 ng of template human genomic DNA and 5  $\mu$ M of each primer were used.

**[0819]** Two microliters of each forward and reverse primer, at concentrations of 5 mM were pooled into a single microcentrifuge tube and mixed. Eight microliters of the primer mix was used in a total PCR reaction volume of 40  $\mu$ l (1.5  $\mu$ l of template DNA, 10.5  $\mu$ l of sterile water, 8  $\mu$ l of primer mix, and 20  $\mu$ l of HotStar Taq). Twenty-five cycles of PCR were performed. The following PCR conditions were used:

- (1) 95°C for 15 minutes;
- (2) 95°C for 30 seconds;
- (3) 4°C for 30 seconds;
- (4) 37°C for 30 seconds;

(5) Repeat steps 2-4 twenty-four (24) times;

(6) 72°C for 10 minutes.

**[0820]** The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results.

**[0821]** In another embodiment, the loci of interest are amplified using 6-base oligonucleotides, 7-base oligonucleotides, 8-base oligonucleotides, 9-base oligonucleotides, 10-base oligonucleotides, 11-base oligonucleotides, 12-base oligonucleotides, 13-base oligonucleotides, 14-base oligonucleotides, or greater than 14-base oligonucleotides. In a preferred embodiment, 6-base oligonucleotides, 7-base oligonucleotides, 8-base oligonucleotides, 9-base oligonucleotides, 10-base oligonucleotides, 11-base oligonucleotides, or 12-base oligonucleotides are used to amplify the loci of interest. In another embodiment, any number of oligonucleotides can be used including but not limited to 1-5, 5-10, 10-15, 15-20, 20-25, 25-30, 30-35, 35-40, 40-45, 45-50, 50-100, 100-500, 500-1000, 1000-2000, 2000-4000, 4000-8000, 8000-10,000 or greater than 10,000. With a small number of random oligos, the concentration of the oligos is large enough to allow efficient amplification, and yet, the number of oligos is small enough that it does not cause interference between the oligos. This allows efficient amplification of the genome.

**[0822]** In another embodiment, the upstream and downstream sequences of the loci of interest are analyzed to identify a 6-base, 7-base, 8-base, 9-base, 10-base, 11-base, or 12-base sequence that is present in the sequence upstream or downstream for each of the loci of interest, which is then used to amplify the loci of interest. In another embodiment, any number of 6-base oligonucleotides can be used to amplify the loci of interest including but not limited to 1-10, 10-50, 50-100, 100-200, 200-500, or greater than 500.

**[0823]** In another embodiment, the number of loci of interest from a small number of genomes can be increased by amplifying a limited number of the loci of interest, followed by removal of the primers, and amplification of the remaining loci of interest. All the loci of interest do not have to be multiplexed in one reaction. Any number of experimentally determined loci of interest can be multiplexed in a single reaction including but not limited to 1-5, 5-10, 10-25, 25-50, 50-100, 100-200, 200-400, or greater than 400. After increasing the number of copies of these loci of interest, the sample can be passed through a column that allows the amplified products to bind and the primers and unused dNTPs to be removed. After eluting the bound products from the column, different loci of interest can be amplified in a single reaction. This reduces the amount of interaction between the primers.

**[0824]** Other methods of genomic amplification can also be used to increase the copy number of the loci of interest including but not limited to primer extension preamplification (PEP) (Zhang *et al.*, PNAS, 89:5847-51, 1992), degenerate oligonucleotide primed PCR (DOP-PCR) (Telenius, *et al.*, Genomics 13:718-25, 1992), strand displacement amplification using DNA polymerase from bacteriophage 29, which undergoes rolling circle replication (Dean *et al.*, Genomic Research 11:1095-99, 2001), multiple displacement amplification (U.S. Patent 6,124,120), REPLI-g™ Whole Genome Amplification kits, and Tagged PCR.

#### **Purification of Fragment of Interest**

**[0825]** The excess primers and nucleotides were removed from the reaction by using Qiagen MinElute PCR purification kits (Qiagen, Catalog Number 28004). The reactions were performed following the manufacturer's instructions supplied with the columns. The DNA was eluted in 100 µl of sterile water.

#### **PCR Reaction Two**

**[0826]** SNP TSC0397235 was amplified using the following primer set:

First Primer:

5' TTAGTCATCGCAGAATTCTACTTCTTTCTGAAGTGGGA 3' (SEQ ID NO: 495)

Second primer:

5' GGACAGCTCGATGGGACTAATGCATACTC 3' (SEQ ID NO: 496)

**[0827]** The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI, and was designed to anneal 103 bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

**[0828]** SNP TSC0470003 was amplified using the following primer set:

First Primer:

5' GTAGCCACTGGTGAATTCGTGCCATCGCAAAGAATAA 3' (SEQ ID NO: 497)

Second primer:

5' ATTAGAATGATGGGGACCCCTGTCTTCCC 3' (SEQ ID NO: 498)

**[0829]** The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI, and was designed to anneal 80 bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

**[0830]** SNP TSC1649726 was amplified using the following primer set:

First Primer:

5' ACGCATAGGAAGGAATTCATTCTGACACGTGTGAGATA 3' (SEQ ID NO: 499)

Second primer:

5' GAAATTGACCACGGGACTGCACACTTTTC 3' (SEQ ID NO: 500)

[0831] The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI, and was designed to anneal 113 bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0832] SNP TSC1261039 was amplified using the following primer set:

First Primer:

5' CGGTAAATCGGAGAATTCAAGTTGAGGCATGCATCCAT 3' (SEQ ID NO: 501)

Second primer:

5' TCGGGGCTCAGCGGGACCACAGCCACTCC 3' (SEQ ID NO: 502)

[0833] The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI, and was designed to anneal 54 bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0834] SNP TSC0310507 was amplified using the following primer set:

First Primer:

5' TCTATGCACCACGAATTCAATATGTGTTCAAGGACATT 3' (SEQ ID NO: 503)

Second primer:

5' TGCTTAATCGGTGGGACTTGTAATTGTAC 3' (SEQ ID NO: 504)

[0835] The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI, and was designed to anneal 93 bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0836] SNP TSC1650432 was amplified using the following primer set:

First Primer:

5' CGCGTTGTATGCGAATTCCCTGGGGTATAAAGATAAGA 3' (SEQ ID NO: 505)

Second primer:

5' CTCACGGGAACCTGGGACACCTGACCCTGC 3' (SEQ ID NO: 506)

[0837] The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI, and was designed to anneal 80 bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0838] SNP TSC1335008 was amplified using the following primer set:

First Primer:

5' GTCTTGCCGCTTGAATTCCCATAGAAGAATGCGCCAAA 3' (SEQ ID NO: 507)

Second primer:

5' TTGAGTAGTACAGGGACACACTAACAGAC 3' (SEQ ID NO: 508)

[0839] The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI, and was designed to anneal 94 bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0840] SNP TSC0128307 was amplified using the following primer set:

First Primer:

5' AATACTGTAGGTGAATTCTTGCCTAAGCATTTTCCCAG 3' (SEQ ID NO: 509)

Second primer:

5' GTGTTGACATTCGGGACTGTAATCTTGAC 3' (SEQ ID NO: 510)

[0841] ~~The first~~ The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI, and was designed to anneal 54 bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0842] SNP TSC0259757 was amplified using the following primer set:

First Primer:

5' TCTGTAGATTCGGAATTCTTTAGAGCCTGTGCGCTGAG 3' (SEQ ID NO: 511)

Second primer:

5' CGTACCAGTACAGGGACGCAAACCTGAGAC 3' (SEQ ID NO: 512)

[0843] The first primer contained a biotin tag at the 5' end and a recognition site for the restriction enzyme EcoRI, and was designed to anneal 100 bases from the locus of interest. The second primer contained the recognition site for the restriction enzyme BsmF I.

[0844] All loci of interest were amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). In this example, the loci of interest were amplified in separate reaction tubes but they can also be amplified together in a single PCR reaction. For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443).

[0845] One microliter of the eluate from the multiplex reaction (PCR product eluted from the MinElute column) was used as template DNA for each PCR reaction. Each SNP was amplified in triplicate when the multiplex sample was used as the template. As a control, each SNP was amplified from 15 ng of the original template DNA (DNA that did not undergo the multiplex reaction). The amount

of template DNA and primer per reaction can be optimized for each locus of interest but in this example, 5  $\mu$ M of each primer was used. Forty cycles of PCR were performed. The following PCR conditions were used:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 64°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty nine (39) times;
- (9) 72°C for 5 minutes.

**[0846]** In the first cycle of PCR, the annealing temperature was about the melting temperature of the 3' annealing region of the second primers, which was 37°C. The annealing temperature in the second cycle of PCR was about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer, which was 57°C. The annealing temperature in the third cycle of PCR was about the melting temperature of the entire sequence of the second primer, which was 64°C. The annealing temperature for the remaining cycles was 64°C. Escalating the annealing temperature from TM1 to TM2 to TM3 in the first three cycles of PCR greatly improves specificity. These annealing temperatures are representative, and the skilled artisan will understand the annealing temperatures for each cycle are dependent on the specific primers used.

**[0847]** The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results.

#### **Agarose Gel Analysis**

**[0848]** Four microliters of a twenty microliter PCR reaction for each SNP from the original template DNA was analyzed by agarose gel electrophoresis (see FIG. 18A). Four microliters of a twenty microliter PCR reaction for each SNP that was amplified from the multiplexed template was analyzed on by agarose gel electrophoresis (see FIG. 18B).

**[0849]** As seen in FIG. 18A, for 8/9 of the SNPs amplified from the original template DNA, a single band of high intensity was seen (lanes 1-3, and 5-9). The band migrated at the correct position for each of the 8 SNPs. Amplification of TSC1261039 from the original template DNA produced a band of high

intensity, which migrated at the correct position, and a faint band of lower molecular weight (lane 4). Only two bands were seen, and the bands could clearly be distinguished based on molecular weight. The PCR method described herein allows clean amplification of the loci of interest from genomic DNA without any concentration or enrichment of the loci of interest.

**[0850]** As seen in FIG. 18B, the primers used to amplify SNPs TSC0397235, TSC0470003, TSC0310507, and TSC0128307 from the multiplexed template DNA produced a single band of high intensity, which migrated at the correct position (lanes 1, 2, 5, and 8). No additional bands were introduced despite the fact that the multiplex reaction contained two hundred primers. While the multiplex primers were 12 bases in length and likely annealed to additional sequences other than those located on chromosome 21, the products were not seen because the bands were not amplified in the second PCR reaction. The second PCR reaction employed primers specific for the loci of interest and used asymmetric oligonucleotides and escalating annealing temperatures, which allows specific amplification from the genome (see Example 1).

**[0851]** Amplification of TSC1649726 from the multiplex template DNA produced one band of high intensity and two weaker bands, which could clearly be distinguished based on molecular weight (see FIG. 18B, lane 3). Amplification of TSC1261039 from the multiplex template DNA produced a high intensity band of the correct molecular weight and a faint band of lower molecular weight (see FIG. 18B, lane 4). The low molecular weight band was the same size as the band seen from the amplification of TSC1261039 from the original template DNA (compare FIG. 18A, lane 4 with FIG. 18B, lane 4). Thus, amplification of TSC1261039 on the multiplex template DNA did not introduce any additional non-specific bands

**[0852]** Amplification of SNPs TSC1650432, TSC1335008, and TSC0259757 from the multiplex template DNA produced one band of high intensity, which migrated at the correct position, and one weaker band (lanes 6, 7, and 9). For SNPs TSC1650432 and TSC0259757, the weaker band was of lower molecular weight, and clearly was distinguishable from the band of interest (see FIG. 18B, lanes 6 and 9). For SNP TSC1335008, the weaker band was of slightly higher molecular weight. However, the correct band can be identified by comparing to the amplification products of TSC1335008 from the original template DNA, (compare FIG. 18A, lane 7 and FIG. 18B, lane 7). The PCR conditions can also be optimized for TSC1335008. All 9 SNPs were amplified under the exact same conditions, which produced clearly distinguishable bands for the amplified SNPs.



### **Purification of Fragment of Interest**

**[0853]** The PCR products were separated from the genomic template DNA. One half of the PCR reaction was transferred to a well of a Streptawell, transparent, High-Bind plate from Roche Diagnostics GmbH (catalog number 1 645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog). The first primers contained a 5' biotin tag so the PCR products bound to the Streptavidin coated wells while the genomic template DNA did not. The streptavidin binding reaction was performed using a Thermomixer (Eppendorf) at 1000 rpm for 20 min. at 37°C. Each well was aspirated to remove unbound material, and washed three times with 1X PBS, with gentle mixing (Kandpal et al., Nucl. Acids Res. 18:1789-1795 (1990); Kaneoka et al., Biotechniques 10:30-34 (1991); Green et al., Nucl. Acids Res. 18:6163-6164 (1990)).

### **Restriction Enzyme Digestion of Isolated Fragments**

**[0854]** The purified PCR products were digested with the restriction enzyme BsmF I, which binds to the recognition site incorporated into the PCR products from the second primer. The digests were performed in the Streptawells following the instructions supplied with the restriction enzyme. After digestion, the wells were washed three times with PBS to remove the cleaved fragments.

### **Incorporation of Labeled Nucleotide**

**[0855]** The restriction enzyme digest with BsmF I yielded a DNA fragment with a 5' overhang, which contained the SNP site or locus of interest and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase.

**[0856]** As discussed in detail in Example 6, the sequence of both alleles of a SNP can be determined by using one labeled nucleotide in the presence of the other unlabeled nucleotides. The following components were added to each fill in reaction: 1 µl of fluorescently labeled ddGTP, 0.5 µl of unlabeled ddNTPs (40 µM), which contained all nucleotides except guanine, 2 µl of 10X sequenase buffer, 0.25 µl of Sequenase, and water as needed for a 20µl reaction. The fill in reaction was performed at 40°C for 10 min. Non-fluorescently labeled ddNTP was purchased from Fermentas Inc. (Hanover, MD). All other labeling reagents were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565).

**[0857]** After labeling, each Streptawell was rinsed with 1X PBS (100 µl) three times. The "filled in" DNA fragments then were released from the Streptawells by digestion with the restriction enzyme EcoRI, according to the manufacturer's instructions that were supplied with the enzyme. Digestion was performed for 1 hour at 37 °C with shaking at 120 rpm.

**Detection of the Locus of Interest**

**[0858]** The samples were loaded into a lane of a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691). The samples were electrophoresed into the gel at 3000 volts for 3 min. The gel was run for 3 hours on a sequencing apparatus (Hoefer SQ3 Sequencer). The gel was removed from the apparatus and scanned on the Typhoon 9400 Variable Mode Imager. The incorporated labeled nucleotide was detected by fluorescence. A box was drawn around each band and the intensity of the band was calculated using the ImageQuant software.

**[0859]** Below, a schematic of the 5' overhang for TSC0470003 after digestion with BsmF I is depicted:

5' CTCT				
3' GAGA	R	A	C	C
Overhang position	1	2	3	4

**[0860]** The observed nucleotides for TSC0470003 are adenine and guanine on the sense strand (herein depicted as the top strand). The third position of the overhang corresponds to cytosine, which is complementary to guanine. Labeled ddGTP was used in the presence of unlabeled dATP, dCTP, and dTTP. Schematics of the DNA molecules after the fill-in reaction are depicted below:

Allele 1	5' CTCT	<b>G*</b>			
	3' GAGA	C	A	C	C
	Overhang position	1	2	3	4
Allele 2	5' CTCT	<b>A</b>	<b>T</b>	<b>G*</b>	
	3' GAGA	T	A	C	C
	Overhang position	1	2	3	4

**[0861]** Two bands were seen; the lower molecular weight band corresponded to the DNA molecules filled in with ddGTP at position 1 complementary to the overhang and the higher molecular weight band corresponded to the DNA molecules filled in with ddGTP at position 3 complementary to the overhang (see FIG. 19).

**[0862]** The percentage of allele 2 to allele 1 at TSC0470003 after amplification from the original template DNA and the multiplexed template DNA was calculated. The use of one fluorescently labeled nucleotide to detect both alleles in a single reaction reduces the amount of error that is introduced through pipetting reactions, and the error that is introduced through the quantum coefficients of different dyes.

**[0863]** For SNP TSC047003, the percentage of allele 2 to allele 1 was calculated by dividing the value of allele 2 by the sum of the values for allele 2 and allele 1. The percentage of allele 2 to allele 1 for TSC047003 on the original template DNA was calculated to be 0.539 (see Table XIX). Three PCR reactions were performed for each SNP on the multiplexed template DNA. The average percentage of allele 2 to allele 1 for TSC047003 on the multiplexed DNA was 0.49 with a standard deviation of 0.0319 (see Table XIX). There was no statistically significant difference between the percentage obtained on the original template DNA and the multiplexed template DNA.

**[0864]** For SNP TSC1261039, the percentage of allele 2 to allele 1 for TSC1261039 on the original template DNA was calculated to be 0.44 (see Table XIX). Three PCR reactions were performed for each SNP on the multiplexed template DNA (see FIG. 19B). The average percentage of allele 2 to allele 1 for TSC1261039 on the multiplexed DNA was 0.468 with a standard deviation of 0.05683 (see Table XIX). There was no statistically significant difference between the percentages of allele 2 to allele 1 obtained on the original template DNA and the multiplexed template DNA.

**[0865]** The variation seen in the percentage of allele 2 to allele 1 for TSC1261039 on the multiplexed template DNA was likely due to pipetting reactions. The variation can be reduced by increasing the number of replicates. With a large number of replicates, a percentage can be obtained with minimum statistical variation.

**[0866]** Likewise, there was no statistical difference between the percentage of allele 2 to allele 1 on the original template DNA and on the multiplexed template DNA for SNPs TSC0310507 and TSC1335008 (see Table XIX, and FIGS. 19C and 19D). Thus, a multiplex reaction can be used to increase the number of chromosomal regions containing the loci of interest without affecting the percentage of one allele to the other at the variable sites.

**TABLE XIX.**

**Percentage of allele 2 to allele 1 at various SNPs with and without multiplexing.**

TSC047003			
	Allele 1	Allele 2	2/(2+1)
IA	5535418	6487873	0.539608748
M1	4804358	4886716	0.504249168
M2	5549389	5958585	0.517778803
M3	8356275	7030245	0.45690936
Mean (M1-M3)			0.49297911
STDEV			0.031961429
TSC1261039			

	Allele 1	Allele 2	2/(2+1)
IA	3488765	2768066	0.442407027
M1	3603388	2573244	0.41660957
M2	4470423	5026872	0.529295131
M3	4306015	36694012	0.46008898
Mean (M1-M3)			0.46866456
STDEV			0.056830136
TSC0310507			
	Allele 1	Allele 2	2/(2+1)
IA	2966511	2688190	0.475390299
M1	4084472	2963451	0.420471535
M2	4509891	4052892	0.47331481
M3	7173191	4642069	0.39288759
Mean (M1-M3)			0.428891312
STDEV			0.040869352
TSC1335008			
	Allele 1	Allele 2	2/(2+1)
IA	2311629	2553016	0.524810341
M1	794790	900879	0.531282343
M2	1261568	1780689	0.5853184
M3	1165156	1427840	0.550653
Mean (M1-M3)			0.555751248
STDEV			0.027376412

**[0867]** The methods described herein used two distinct amplification reactions to amplify the loci of interest. In the first PCR reaction, oligonucleotides were designed to anneal upstream and downstream of the loci of interest. Unlike traditional genomic amplification, these primers were not degenerate and annealed at a specified distance from the loci of interest. However, due to the length of the primers, it is likely that the primers annealed to other regions of the genome. These primers were used to increase the amount of DNA available for genetic analysis.

**[0868]** The second PCR reaction employs the methods described in Examples 1-6. The primers are designed to amplify the loci of interest, and the sequence is determined at the loci of interest. The conditions of the second PCR reaction allowed specific amplification of the loci of interest from the

multiplexed template DNA. If there were any non-specific products from the multiplex reaction, they did not impede amplification of the loci of interest. There was no statistical difference in the percentages of allele 2 to allele 1 at the four SNPs analyzed, regardless of whether the amplification was performed on original template DNA or multiplexed template DNA.

**[0869]** The SNPs analyzed in this example were located on human chromosome 21. However, the methods can be applied to non-human and human DNA including but not limited to chromosomes 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X, and Y. The multiplex methods can also be applied to analysis of genetic mutations including but not limited to nucleotide substitutions, insertions, deletions, and rearrangements.

**[0870]** The above methods can be used to increase the amount of DNA available for genetic analysis whenever the starting template DNA is limiting in quantity. For example, pre-malignant and pre-invasive lesions with malignant cells usually constitute a small fraction of the cells in the specimen, which reduces the number of genetic analyses that can be performed. The methods described herein can be used to increase the amounts of malignant DNA available for genetic analysis. Also, the number of fetal genomes present in the maternal blood is often low; the methods described herein can be used to increase the amount of fetal DNA.

### EXAMPLE 13

**[0871]** Plasma isolated from blood of a pregnant female contains both maternal template DNA and fetal template DNA. As discussed earlier, the percentage of fetal DNA in the maternal plasma varies for each pregnant female. However, the percentage of fetal DNA can be determined by analyzing SNPs wherein the maternal template DNA is homozygous and the template DNA obtained from the plasma displays a heterozygous pattern.

**[0872]** For example, assume SNP X can either be adenine or guanine, and the maternal DNA for SNP X is homozygous for guanine. The labeling method described in Example 6 can be used to determine the sequence of the template DNA in the plasma sample. If the plasma sample contains fetal DNA, which is heterozygous at SNP X, the following DNA molecules are expected after digestion with the type IIS restriction enzyme BsmF I, and the fill-in reaction with labeled ddGTP, unlabeled dATP, dTTP, and dCTP.

Maternal Allele 1	5' GGGT	G*				
	3'CCCA	C	T	C	A	

Maternal Allele 2	5' GGGT	<b>G*</b>				
	3'CCCA	C	T	C	A	
Fetal Allele 1	5' GGGT	<b>G*</b>				
	3'CCCA	C	T	C	A	
Fetal Allele 2	5' GGGT	A	A	<b>G*</b>		
	3'CCCA	T	T	C	A	

**[0873]** Two signals are seen; one signal corresponds to the DNA molecules filled in with ddGTP at position one complementary to the overhang and the second signal corresponds to the DNA molecules filled in with ddGTP at position three complementary to the overhang. However, the maternal DNA is homozygous for guanine, which corresponds to the DNA molecules filled in at position one complementary to the overhang. The signal from the DNA molecules filled in with ddGTP at position three complementary to the overhang corresponds to the adenine allele, which represents the fetal DNA. This signal becomes a beacon for the fetal DNA, and can used to measure the amount of fetal DNA present in the plasma sample.

**[0874]** There is no difference in the amount of fetal DNA from one chromosome to another. For instance, the percentage of fetal DNA in any given individual from chromosome 1 is the same as the percentage of fetal DNA from chromosome 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X and Y. Thus, the allele ratio calculated for SNPs on one chromosome can be compared to the allele ratio for the SNPs on another chromosome.

**[0875]** For example, the allele ratio for the SNPs on chromosome 1 should be equal to the allele ratio for the SNPs on chromosomes 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X, and Y. However, if the fetus has a chromosomal abnormality, including but not limited to a trisomy or monosomy, the ratio for the chromosome that is present in an abnormal copy number will differ from the ratio for the other chromosomes.

**[0876]** Blood from a pregnant female was collected after informed consent had been obtained. The blood sample was used to demonstrate that fetal DNA can be detected in the maternal plasma by analyzing SNPs wherein the maternal DNA was homozygous, and the same SNP displayed a heterozygous pattern from DNA obtained from the plasma of a pregnant woman.

### **Preparation of Plasma from Whole Blood**

[0877] Plasma was isolated from 4 tubes each containing 9 ml of blood (Fischer Scientific, 9 ml EDTA Vacuette tubes, catalog number NC9897284). The blood was obtained by venipuncture from a pregnant female who had given informed consent. After collecting the blood, formaldehyde (25 µl/ml of blood) was added to each of the tubes. The tubes were placed at 4°C until shipment. The tubes were shipped via Federal Express in a foam container containing an ice pack.

[0878] The blood was centrifuged at 1000 rpm for 10 minutes. The brake on the centrifuge was not used. This centrifugation step was repeated. The supernatant was transferred to a new tube and spun at 3,000 rpm for ten minutes. The brake on the centrifuge was not used. The supernatant from each of the four tubes was pooled and aliquoted into two tubes. The plasma was stored at -80°C until the DNA was purified.

[0879] Template DNA was isolated using the QIAmp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183). The template DNA was isolated as per instructions included in the kit. The template DNA from the plasma was eluted in a final volume of 20 microliters.

### **Isolation of Maternal DNA**

[0880] After the plasma was removed from the sample described above, one milliliter of the remaining blood sample, which is commonly referred to as the “buffy-coat,” was transferred to a new tube. One milliliter of 1X PBS was added to the sample. Template DNA was isolated using the QIAmp DNA Blood Midi Kit supplied by QIAGEN (Catalog number 51183).

### **Identification of Homozygous Maternal SNPs**

[0881] Example 8 describes a method for identifying SNPs that are highly variable within the population or for identifying heterozygous SNPs for a given individual. The methods as described in Example 8 were applied to the maternal template DNA to identify SNPs on chromosome 13 wherein the maternal DNA was homozygous. Any number of SNPs can be screened. The number of SNPs to be screened is proportional to the number of heterozygous SNPs in the fetal DNA that need to be analyzed.

[0882] As described in detail in Example 6, one labeled nucleotide can be used to determine the sequence of both alleles at a particular SNP. SNPs for which the sequence can be determined with labeled ddGTP in the presence of unlabeled dATP, dTTP, and dCTP were chosen for this example. However, SNPs for which the sequence can be determined with labeled ddATP, ddCTP or ddTTP can also be used. Additionally, the SNPs to be analyzed can be chosen such that all are labeled with the same

nucleotide or any combination of the four nucleotides. For instance, if 400 SNPs are to be screened, 100 can be chosen such that the sequence is determined with labeled ddATP, 100 can be chosen such that the sequence is determined with labeled ddTTP, 100 can be chosen such that the sequence is determined with labeled ddGTP, and 100 can be chosen such that the sequence is determined with labeled ddCTP, or any combination of the four labeled nucleotides.

**[0883]** Twenty-nine SNPs wherein the maternal DNA was homozygous were identified: TSC0052277, TSC1225391, TSC0289078, TSC1349804, TSC0870209, TSC0194938, TSC0820373, TSC0902859, TSC0501510, TSC1228234, TSC0082910, TSC0838335, TSC0818982, TSC0469204, TSC1084457, TSC0466177, TSC1270598, TSC1002017, TSC1104200, TSC0501389, TSC0039960, TSC0418134, TSC0603688, TSC0129188, TSC1103570, TSC0813449, TSC0701940, TSC0087962, and TSC0660274. Heterozygous SNPs will vary from individual to individual.

### **Design of Multiplex Primers**

**[0884]** A low copy number of fetal genomes typically is present in the maternal plasma. To increase the copy number of the loci of interest located on chromosome 13, primers were designed to anneal at approximately 130 bases upstream and 130 bases downstream of each loci of interest. This was done to reduce statistical sampling error that can occur when working with a low number of genomes, which can influence the ratio of one allele to another (see Example 11). The primers were 12 bases in length. However, primers of any length can be used including but not limited to 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36-45, 46-55, 56-65, 66-75, 76-85, 86-95, 96-105, 106-115, 116-125, and greater than 125 bases. Primers were designed to anneal to both the sense strand and the antisense strand.

**[0885]** The primers were designed to terminate at the 3' end in the dinucleotide "AA" to reduce the formation of primer-dimers. However, the primers can be designed to end in any of the four nucleotides and in any combination of the four nucleotides.

**[0886]** The multiplex primers for SNPTSC0052277 were

Forward primer:

5' GACATGTTGGAA 3' (SEQ ID NO: 513)

Reverse primer:

5' ACTTCCAGTTAA 3' (SEQ ID NO: 514)

**[0887]** The multiplex primers for SNP TSC1225391 were:

Forward primer:



5' GTTTCCTGTAA 3' (SEQ ID NO: 515)

Reverse primer

5' CGATGATGACAA 3' (SEQ ID NO: 516)

**[0888]** The multiplex primers for SNP TSC0289078 were:

Forward primer

5' GAGTAGAGACAA 3' (SEQ ID NO: 517)

Reverse primer

5' TCCCGGATACAA 3' (SEQ ID NO: 518)

**[0889]** The multiplex primers for SNP TSC1349804 were:

Forward primer:

5' CATCCTCTAGAA 3' (SEQ ID NO: 519)

Reverse primer:

5' TATTCCTGAGAA 3' (SEQ ID NO: 520)

**[0890]** The multiplex primers for SNP TSC0870209 were:

Forward primer:

5' AGTTTGTTTAA 3' (SEQ ID NO: 521)

Reverse primer:

5' TATAAACGATAA 3' (SEQ ID NO: 522)

**[0891]** The multiplex primers for SNP TSC0194938 were:

Forward primer:

5' TTTGACCGATAA 3' (SEQ ID NO: 523)

Reverse primer:

5' TGACAGGACCAA 3' (SEQ ID NO: 524)

**[0892]** The multiplex primers for SNP TSC0820373 were:

Forward primer:

5' TTATTCATTCAA 3' (SEQ ID NO: 525)

Reverse primer:

5' AGTTTTTCACAA 3' (SEQ ID NO: 526)

**[0893]** The multiplex primers for SNP TSC0902859 were:

Forward primer:

5' CACCTCCCTGAA 3' (SEQ ID NO: 527)

Reverse primer:

5' CCAGATTGAGAA 3' (SEQ ID NO: 528)

[0894] The multiplex primers for SNP TSC0501510 were:

Forward primer:

5' TGTGTCCACCAA 3' (SEQ ID NO: 529)

Reverse primer:

5' CTTCTATTCCAA 3' (SEQ ID NO: 530)

[0895] The multiplex primers for SNP TSC1228234 were:

Forward primer:

5' TCACAATAGGAA 3' (SEQ ID NO: 531)

Reverse primer:

5' TACAAGTGAGAA 3' (SEQ ID NO: 532)

[0896] The multiplex primers for SNP TSC0082910 were:

Forward primer:

5' GAGTTTTTCGTAA 3' (SEQ ID NO: 533)

Reverse primer:

5' GTGTGCCCCCAA 3' (SEQ ID NO: 534)

[0897] The multiplex primers for SNP TSC0838335 were:

Forward primer:

5' GCACCACTGCAA 3' (SEQ ID NO: 535)

Reverse primer:

5' GAACACAATGAA 3' (SEQ ID NO: 536)

[0898] The multiplex primers for SNP TSC0818982 were:

Forward primer:

5' TATCCTATTCAA 3' (SEQ ID NO: 537)

Reverse primer:

5' CAACCATTATAA 3' (SEQ ID NO: 538)

[0899] The multiplex primers for SNP TSC0469204 were:

Forward primer:

5' TATGCTTTACAA 3' (SEQ ID NO: 539)

Reverse primer:

5' TTTGTTTACCAA 3' (SEQ ID NO: 540)

[0900] The multiplex primers for SNP TSC1084457 were:

Forward primer:

5' AGGAAATTAGAA 3' (SEQ ID NO: 541)

Reverse primer:

5' TGTTAGACTTAA 3' (SEQ ID NO: 542)

[0901] The multiplex primers for SNP TSC0466177 were:

Forward primer:

5' TATTTGGAGGAA 3' (SEQ ID NO: 543)

Reverse primer:

5' GGCATTTGTCAA 3' (SEQ ID NO: 544)

[0902] The multiplex primers for SNP TSC1270598 were:

Forward primer:

5' ATACTCCAGGAA 3' (SEQ ID NO: 545)

Reverse primer:

5' CAGCCTGGACAA 3' (SEQ ID NO: 546)

[0903] The multiplex primers for SNP TSC1002017 were:

Forward primer:

5' CCATTGCAGTAA 3' (SEQ ID NO: 547)

Reverse primer:

5' AGGTTCTCATAA 3' (SEQ ID NO: 548)

[0904] The multiplex primers for SNP TSC1104200 were:

Forward primer:

5' TGTCATCATTA 3' (SEQ ID NO: 549)

Reverse primer:

5' TGGTATTTGCAA 3' (SEQ ID NO: 550)

[0905] The multiplex primers for SNP TSC0501389 were:

Forward primer:

5' TAGGGTTTGTA 3' (SEQ ID NO: 551)

Reverse primer:

5' CCCTAAGTAGAA 3' (SEQ ID NO: 552)

[0906] The multiplex primers for SNP TSC0039960 were:

Forward primer:

5' GTATTTCTTTAA 3' (SEQ ID NO: 553)

Reverse primer:

5' GAGTCTTCCCAA 3' (SEQ ID NO: 554)

[0907] The multiplex primers for SNP TSC0418134 were:

Forward primer:

5' CAGGTAGAGTAA 3' (SEQ ID NO: 555)

Reverse primer:

5' ATAGGATGTGAA 3' (SEQ ID NO: 556)

[0908] The multiplex primers for SNP TSC0603688 were:

Forward primer:

5' CAATGTGTATAA 3' (SEQ ID NO: 557)

Reverse primer:

5' AGAGGGCATCAA 3' (SEQ ID NO: 558)

[0909] The multiplex primers for SNP TSC0129188 were:

Forward primer:

5' CCAGTGGTCTAA 3' (SEQ ID NO: 559)

Reverse primer:

5' TAAACAATAGAA 3' (SEQ ID NO: 560)

[0910] The multiplex primers for SNP TSC1103570 were:

Forward primer:

5' GCACACTTTTAA 3' (SEQ ID NO: 561)

Reverse primer:

5' ATGGCTCTGCAA 3' (SEQ ID NO: 562)

[0911] The multiplex primers for SNP TSC0813449 were:

Forward primer:

5' GTCATCTTGTA 3' (SEQ ID NO: 563)

Reverse primer:

5' TGCTTCATCTAA 3' (SEQ ID NO: 564)

[0912] The multiplex primers for SNP TSC0701940 were:

Forward primer:

5' AGAAAGGGGCAA 3' (SEQ ID NO: 565)

Reverse primer:

5' CTTTTCTTTCAA 3' (SEQ ID NO: 566)

[0913] The multiplex primers for SNP TSC0087962 were:

Forward primer:

5' CTACTCTCTCAA 3' (SEQ ID NO: 567)

Reverse primer:

5' ACAGCATTATAA 3' (SEQ ID NO: 568)

[0914] The multiplex primers for SNP TSC0660274 were:

Forward primer:

5' ACTGCTCTGGAA 3' (SEQ ID NO: 569)

Reverse primer:

5' GCAGAGGCACAA 3' (SEQ ID NO: 570)

Multiplex PCR

[0915] Regions on chromosome 13 surrounding the above-mentioned 29 SNPs were amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). This PCR reaction used primers that annealed approximately 150 bases upstream and downstream of each loci of interest. The fifty-eight primers were mixed together and used in a single reaction to amplify the template DNA. This reaction was done to increase the number of copies of the loci of interest, which eliminates error generated from a low number of genomes.

[0916] For increased specificity, a "hot-start" PCR reaction was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction can be optimized for each locus of interest. In this example, the 20 µl of plasma template DNA was used.

[0917] Two microliters of each forward and reverse primer, at concentrations of 5 mM were pooled into a single microcentrifuge tube and mixed. Four microliters of the primer mix was used in a total PCR reaction volume of 50 µl (20µl of template plasma DNA, 1 µl of sterile water, 4 µl of primer mix, and 25 µl of HotStar Taq. Twenty-five cycles of PCR were performed. The following PCR conditions were used:

- (1) 95°C for 15 minutes;
- (2) 95°C for 30 second;
- (3) 4°C for 30 seconds;
- (4) 37°C for 30 seconds;

(5) Repeat steps 2-4 twenty-four (24) times;

(6) 72°C for 10 minutes.

**[0918]** The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results.

**[0919]** Other methods of genomic amplification can also be used to increase the copy number of the loci of interest including but not limited to primer extension preamplification (PEP) (Zhang *et al.*, PNAS, 89:5847-51, 1992), degenerate oligonucleotide primed PCR (DOP-PCR) (Telenius, *et al.*, Genomics 13:718-25, 1992), strand displacement amplification using DNA polymerase from bacteriophage 29, which undergoes rolling circle replication (Dean *et al.*, Genomic Research 11:1095-99, 2001), multiple displacement amplification (U.S. Patent 6,124,120), REPLI-g™ Whole Genome Amplification kits, and Tagged PCR.

### **Purification of Fragment of Interest**

**[0920]** The unused primers, and nucleotides were removed from the reaction by using Qiagen MinElute PCR purification kits (Qiagen, Catalog Number 28004). The reactions were performed following the manufacturer's instructions supplied with the columns. The DNA was eluted in 100 µl of sterile water.

### **PCR Reaction Two**

#### **Design of Primers**

**[0921]** SNPTSC0052277 was amplified using the following primer set:

First primer:

5' CTCCGTGGTATGGAATTCCACTCAAATCTTCATTCAGA 3' (SEQ ID NO: 571)

Second primer:

5' ACGTCGGGTTACGGGACACCTGATTCCTC 3' (SEQ ID NO: 572)

**[0922]** SNP TSC1225391 was amplified using the following primer set:

First primer:

5' TACCATTGGTTTGAATTCTTGTTTCCTGTAAACCATGC 3' (SEQ ID NO: 573)

Second primer:

5' GCCGAGTTCTACGGGACAGAAAAGGGAGC 3' (SEQ ID NO: 574)

**[0923]** SNP TSC0289078 was amplified using the following primer set:

First primer:

5' TGCAGTGATTTCTGAATTCGAGACAATGCTGCCCAGTCA 3' (SEQ ID NO: 575)

Second primer:

5' TCTAAATTCTCTGGGACCATTCTTCAAC 3' (SEQ ID NO: 576)

**[0924]** SNP TSC1349804 was amplified using the following primer set:

First primer:

5' ACTAACAGCACTGAATTCCATGCTCTTGGACTTTCCAT 3' (SEQ ID NO: 577)

Second primer:

5' TCCCCTAACGTTGGGACACAGAATACTAC 3' (SEQ ID NO: 578)

**[0925]** SNP TSC0870209 was amplified using the following primer set:

First primer:

5' GTCGACGATGGCGAATTCCTGCCACTCATTCAGTTAGC 3' (SEQ ID NO: 579)

Second primer:

5' GAACGGCCCCACAGGGACCTGGCATAACTC 3' (SEQ ID NO: 580)

**[0926]** SNP TSC0194938 was amplified using the following primer set:

First primer:

5' TCATGGTAGCAGGAATTCTGCTTTGACCGATAAGGAGA 3' (SEQ ID NO: 581)

Second primer:

5' ACTGTGGGATTCGGGACTGTCTACTACCC 3' (SEQ ID NO: 582)

**[0927]** SNP TSC0820373 was amplified using the following primer set:

First primer:

5' ACCTCTCGGCCGGAATTCGGAAAAGTGTACAGATCATT 3' (SEQ ID NO: 583)

Second primer:

5' GCCGGATACGAAGGGACGGCTCGTGACTC 3' (SEQ ID NO: 584)

**[0928]** SNP TSC0902859 was amplified using the following primer set:

First primer:

5' CCGTAGACTAAAGAATTCCCTGATGTCAGGCTGTCACC 3' (SEQ ID NO: 585)

Second primer:

5' ATCGGATCAGTCGGGACGGTGTCTTTGCC 3' (SEQ ID NO: 586)

**[0929]** SNP TSC0501510 was amplified using the following primer set:

First primer:

5' GCATAGGCGGGAGAATTCCTGTGTCCACCAAAGTCGG 3' (SEQ ID NO: 587)

Second primer:

5' CCCACATAGGGCGGGACAAAGAGCTGAAC 3' (SEQ ID NO: 588)

[0930] SNP TSC1228234 was amplified using the following primer set:

First primer:

5' GGCTTGCCGAGCGAATTCTAGGAAAGATACGGAATCAA 3' (SEQ ID NO: 589)

Second primer:

5' TAACCCTCATACGGGACTTTCATGGAAGC 3' (SEQ ID NO: 590)

[0931] SNP TSC0082910 was amplified using the following primer set:

First primer:

5' ATGAGCACCCGGGAATTCTGATTGGAGTCTAGGCCAAA 3' (SEQ ID NO: 591)

Second primer:

5' TGCTCACCTTCTGGGACGTGGCTGGTCTC 3' (SEQ ID NO: 592)

[0932] SNP TSC0838335 was amplified using the following primer set:

First primer:

5' ACCGTCTGCCACGAATTCTGGAAAACATGCAGTCTGGT 3' (SEQ ID NO: 593)

Second primer:

5' TACACGGGAGGCGGGACAGGGTGATTAAC 3' (SEQ ID NO: 594)

[0933] SNP TSC0818982 was amplified using the following primer set:

First primer:

5' CTAAAGCTAACGAATTCAGAGCTGTATGAAGATGCTT 3' (SEQ ID NO: 595)

Second primer:

5' AACGCTAAAGGGGGGACAACATAATTGGC 3' (SEQ ID NO: 596)

[0934] SNP TSC0469204 was amplified using the following primer set:

First primer:

5' TTGTAAGAACGAGAATTCTGCAACCTGTCTTTATTGAA 3' (SEQ ID NO: 597)

Second primer:

5' CTTACCACTTTGGGACACTGAAGCCAAC 3' (SEQ ID NO: 598)

[0935] SNP TSC1084457 was amplified using the following primer set:

First primer:

5' AACCATTGATTTGAATTCGAAATGTCCACCAAAGTTCA 3' (SEQ ID NO: 599)

Second primer:

5' TGTCTAGTTCAGGGACGCTGTACTTAC 3' (SEQ ID NO: 600)

[0936] SNP TSC0466177 was amplified using the following primer set:



First primer:

5' CGAAGGTAATGTGAATTCTGCCACAATTAAGACTTGGA 3' (SEQ ID NO: 601)

Second primer:

5' ATACCGGTTTTTCGGGACAGATCCATTGAC 3' (SEQ ID NO: 602)

[0937] SNP TSC1270598 was amplified using the following primer set:

First primer:

5' CCTGAAATCCACGAATTCACCCCTGGCCTCCCAGTGCA 3' (SEQ ID NO: 603)

Second primer:

5' TAGATGGTAGGTGGGACAGGACTGGCTTC 3' (SEQ ID NO: 604)

[0938] SNP TSC1002017 was amplified using the following primer set:

First primer:

5' GCATATCTTAGCGAATTCCTGTGACTAATACAGAGTGC 3' (SEQ ID NO: 605)

Second primer:

5' CCAAATATGGTAGGGACGTGTGAACACTC 3' (SEQ ID NO: 606)

[0939] SNP TSC1104200 was amplified using the following primer set:

First primer:

5' TGCCGCTACAGGGAATTCATATGGCAGATATTCCTGAA 3' (SEQ ID NO: 607)

Second primer:

5' ACGTTGCGGACCGGGACTTCCACAGAGCC 3' (SEQ ID NO: 608)

[0940] SNP TSC0501389 was amplified using the following primer set:

First primer:

5' CTTCGCCCAATGGAATTCGGTACAGGGGTATGCCTTAT 3' (SEQ ID NO: 609)

Second primer:

5' TGCACTTCTGCCGGGACCAGAGGAGAAAC 3' (SEQ ID NO: 610)

[0941] SNP TSC0039960 was amplified using the following primer set:

First primer:

5' TGTGGGTATTCTGAATTCCACAAAATGGACTAACACGC 3' (SEQ ID NO: 611)

Second primer:

5' ACGTCGTTTCAGTGGGACATTAAAAGGCTC 3' (SEQ ID NO: 612)

[0942] SNP TSC0418134 was amplified using the following primer set:

First primer:

5' GGTTATGTGTCAGAATTCTGAACTAGTTTGGAAGTAC 3' (SEQ ID NO: 613)

Second primer:

5' GCCTCAGTTTCGGGGACAGTTCTGAGGAC 3' (SEQ ID NO: 614)

[0943] SNP TSC0603688 was amplified using the following primer set:

First primer:

5' TGTAACACGGCCGAATTCCTCATTTGTATGAAATAGGT 3' (SEQ ID NO: 615)

Second primer:

5' AATCTAACTTGAGGGACCGGCACACACAC 3' (SEQ ID NO: 616)

[0944] SNP TSC0129188 was amplified using the following primer set:

First primer:

5' AGTGTCCCCTTAGAATTCGCAGAGACACCACAGTGTGC 3' (SEQ ID NO: 617)

Second primer:

5' TTTGCTACAGTCGGGACCCTTGTGTGCTC 3' (SEQ ID NO: 618)

[0945] SNP TSC1103570 was amplified using the following primer set:

First primer:

5' AGCACATCACTAGAATTCAATACCATGTGTGAGCTCAA 3' (SEQ ID NO: 619)

Second primer:

5' AATCCTGCTTCCGGGACCTAACTTTGAAC 3' (SEQ ID NO: 620)

[0946] SNP TSC0813449 was amplified using the following primer set:

First primer:

5' TTTCATTTTCTGGAATTCCTCTAATGATTTTCTGGAGC 3' (SEQ ID NO: 621)

Second primer:

5' CGTCGCCGCGTAGGGACTTTTCTTCCAC 3' (SEQ ID NO: 622)

[0947] SNP TSC0701940 was amplified using the following primer set:

First primer:

5' TTACTTAATCCTGAATTCGAGAAAAGCCATGTTGATAA 3' (SEQ ID NO: 623)

Second primer:

5' TCATGGGTCGCTGGGACTTTGCCCTCTGC 3' (SEQ ID NO: 624)

[0948] SNP TSC0087962 was amplified using the following primer set:

First primer:

5' ACTAACAGCACTGAATTCATTTTACTATAATCTGCTAC 3' (SEQ ID NO: 625)

Second primer:

5' GTTAGCCGAGAAGGGACTGTCTGTGAAGC 3' (SEQ ID NO: 626)

**[0949]** SNP TSC0660274 was amplified using the following primer set:

First primer:

5' AAATATGCAGCGGAATTCGTAAGTGACCTATTAATAAC 3' (SEQ ID NO: 627)

Second primer:

5' GCGATGGTTACGGGGACAGCCAGGCAACC 3' (SEQ ID NO: 628)

**[0950]** Each first primer had a biotin tag at the 5' end and contained a restriction enzyme recognition site for EcoRI, and was designed to anneal at a specified distance from the locus of interest. This allows a single reaction to be performed for the loci of interest, as each loci of interest will migrate at a distinct position (based on annealing position of first primer). The second primer contained a restriction enzyme recognition site for BsmF I.

**[0951]** All loci of interest were amplified from the multiplexed template DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). In this example, the loci of interest were amplified in separate reaction tubes but they could also be amplified together in a single PCR reaction. For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443).

**[0952]** The amount of multiplexed template DNA and primer per reaction can be optimized for each locus of interest. One microliter of the multiplexed template DNA eluted from the MinElute column was used in the PCR reaction for each locus of interest, and 5  $\mu$ M of each primer was used. The twenty-nine SNPs described above also were amplified from the maternal DNA (15 ng of DNA was used in the PCR reaction; primer concentrations were as stated above). Forty cycles of PCR were performed. The following PCR conditions were used:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 64°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty nine (39) times;
- (9) 72°C for 5 minutes.

**[0953]** In the first cycle of PCR, the annealing temperature was about the melting temperature of the 3' annealing region of the second primers, which was 37°C. The annealing temperature in the second cycle of PCR was about the melting temperature of the 3' region, which anneals to the template DNA, of the first primer, which was 57°C. The annealing temperature in the third cycle of PCR was about the melting temperature of the entire sequence of the second primer, which was 64°C. The annealing temperature for the remaining cycles was 64°C. Escalating the annealing temperature from TM1 to TM2 to TM3 in the first three cycles of PCR greatly improves specificity. These annealing temperatures are representative, and the skilled artisan will understand the annealing temperatures for each cycle are dependent on the specific primers used.

**[0954]** The temperatures and times for denaturing, annealing, and extension, can be optimized by trying various settings and using the parameters that yield the best results. In this example, the first primer was designed to anneal at various distances from the locus of interest. The skilled artisan understands that the annealing location of the first primer can be 5-10, 11-15, 16-20, 21-25, 26-30, 31-35, 36-40, 41-45, 46-50, 51-55, 56-60, 61-65, 66-70, 71-75, 76-80, 81-85, 86-90, 91-95, 96-100, 101-105, 106-110, 111-115, 116-120, 121-125, 126-130, 131-140, 140-160, 160-180, 180-200, 200-220, 220-240, 240-260, 260-280, 280-300, 300-350, 350-400, 400-450, 450-500, or greater than 500 bases from the locus of interest.

#### **Purification of Fragment of Interest**

**[0955]** The PCR products were separated from the genomic template DNA. Each PCR product was placed into a well of a Streptawell, transparent, High-Bind plate from Roche Diagnostics GmbH (catalog number 1 645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog).

Alternatively, the PCR products can be pooled into a single well because the first primer was designed to allow the loci of interest to separate based on molecular weight. The first primers contained a 5' biotin tag so the PCR products bound to the Streptavidin coated wells while the genomic template DNA did not. The streptavidin binding reaction was performed using a Thermomixer (Eppendorf) at 1000 rpm for 20 min. at 37°C. Each well was aspirated to remove unbound material, and washed three times with 1X PBS, with gentle mixing (Kandpal et al., Nucl. Acids Res. 18:1789-1795 (1990); Kaneoka et al., Biotechniques 10:30-34 (1991); Green et al., Nucl. Acids Res. 18:6163-6164 (1990)).

#### **Restriction Enzyme Digestion of Isolated Fragments**

**[0956]** The purified PCR products were digested with the restriction enzyme BsmF I, which binds to the recognition site incorporated into the PCR products from the second primer. The digests were

performed in the Streptawells following the instructions supplied with the restriction enzyme. After digestion, the wells were washed three times with PBS to remove the cleaved fragments.

### **Incorporation of Labeled Nucleotide**

**[0957]** The restriction enzyme digest with BsmF I yielded a DNA fragment with a 5' overhang, which contained the SNP site or locus of interest and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase.

**[0958]** As demonstrated in Example 6, the sequence of both alleles of a SNP can be determined by filling in the overhang with one labeled nucleotide in the presence of the other unlabeled nucleotides. The following components were added to each fill in reaction: 1 µl of fluorescently labeled ddGTP, 0.5 µl of unlabeled ddNTPs (40 µM), which contained all nucleotides except guanine, 2 µl of 10X sequenase buffer, 0.25 µl of Sequenase, and water as needed for a 20µl reaction. The fill in reaction was performed at 40°C for 10 min. Non-fluorescently labeled ddNTP was purchased from Fermentas Inc. (Hanover, MD). All other labeling reagents were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565).

**[0959]** After labeling, each Streptawell was rinsed with 1X PBS (100 µl) three times. The "filled in" DNA fragments were then released from the Streptawells by digestion with the restriction enzyme EcoRI, according to the manufacturer's instructions that were supplied with the enzyme. Digestion was performed for 1 hour at 37 °C with shaking at 120 rpm.

### **Detection of the Locus of Interest**

**[0960]** After release from the streptavidin matrix, the sample was loaded into a lane of a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691). The sample was electrophoresed into the gel at 3000 volts for 3 min. The gel was run for 3 hours on a sequencing apparatus (Hoefer SQ3 Sequencer). The gel was removed from the apparatus and scanned on the Typhoon 9400 Variable Mode Imager. The incorporated labeled nucleotide was detected by fluorescence.

**[0961]** Below a schematic of the 5' overhang for SNP TSC0838335 is depicted. The entire sequence is not reproduced, only a portion to depict the overhang (where R indicates the variable site).

10/14	5' TAA				
	3' ATT	R	A	C	A
Overhang position		1	2	3	4

**[0962]** The observed nucleotides for TSC0838335 are adenine and guanine on the 5' sense strand (herein depicted as the top strand). The nucleotide in position three of the overhang corresponded to cytosine, which is complementary to guanine. Labeled ddGTP can be used to determine the sequence of both allele in the presence of unlabeled dATP, dCTP, and dTTP.

**[0963]** The restriction enzyme BsmF I was used to create the 5' overhang, which typically cuts 10/14 from the recognition site. At times, BsmF I will cut 11/15 from the recognition site and generate the following overhang:

11/15	5' TA				
	3' AT	T	R	A	C
Overhang position		0	1	2	3

**[0964]** Position 0 in the overhang is thymidine, which is complementary to adenine. Position 0 complementary to the overhang was filled in with unlabeled dATP, and thus after the fill-in reaction, the exact same molecules were generated whether the enzyme cut at 10/14 or 11/15 from the recognition site. The DNA molecules generated after the fill-in reaction are depicted below:

G allele 10/14	5' TAA	G*			
	3' ATT	C	A	C	A
Overhang position		1	2	3	4
G allele 11/15	5' TA	A	G*		
	3' AT	T	C	A	C
Overhang position		0	1	2	3
A allele 10/14	5' TAAA	T	G*		
	3' ATT T	A	C	A	
Overhang position		1	2	3	4
A allele 11/15	5' TA	A	A	T	G*
	3' AT	T	T	A	C
Overhang position		0	1	2	3

**[0965]** The maternal template DNA amplified for TSC0838335 displayed a single band that migrated at the expected position of the higher molecular weight band, which corresponded to the "A" allele (see FIG. 20, lane 1). The maternal template DNA was homozygous for adenine at SNP TSC0838335.

**[0966]** However, in lane 2, amplification of the multiplexed template DNA for TSC0838335 isolated from the plasma of the same individual displayed two bands; a lower molecular weight band, which

corresponded to the “G” allele, and the higher molecular weight band, which corresponded to the “A” allele. The template DNA isolated from the plasma of a pregnant female contains both maternal template DNA and fetal template DNA.

[0967] As seen in FIG. 20, lane 1, the maternal template DNA was homozygous for adenine at this SNP (compare lanes 1 and 2). The “G” allele represented the fetal DNA. Signals from the maternal template DNA and the fetal template DNA clearly have been distinguished. The “G” allele becomes a beacon for the fetal DNA and can be used to measure the amount of fetal DNA present in the sample. Additionally, once the percentage of fetal DNA in the maternal plasma for a given sample has been determined, any deviation from this percentage indicates a chromosomal abnormality. This method provides the first non-invasive method for the detection of fetal chromosomal abnormalities.

[0968] As seen in FIG. 20, lane 3, analysis of the maternal DNA for SNP TSC0418134 generated a single band that migrated at the expected position of the higher molecular weight band, which corresponded to the adenine allele. Likewise, analysis of the multiplexed template DNA isolated from the maternal plasma gave a single band, which migrated at the expected position of the adenine allele (see FIG. 20, lane 4). Both the maternal DNA and the fetal DNA are homozygous for adenine at TSC0418134.

[0969] Below, a schematic of the 5’ overhang for TSC0129188 is depicted, wherein R indicates the variable site:

10/14	5’ TCAT				
	3’ AGTA	R	A	C	T
Overhang position		1	2	3	4

[0970] The nucleotide upstream of the variable site (R) does not correspond to guanine on the sense strand. Thus, the 5’ overhang generated by the 11/15 cutting properties of BsmF I will be filled-in identically to the 5’ overhang generated by the 10/14 cut. Labeled ddGTP in the presence of unlabeled dATP, dTTP, and dCTP was used for the fill-in reaction. The DNA molecules generated after the fill-in reaction are depicted below:

A allele 10/14	5’ TCAT	A	T	G*	
	3’ AGTA	T	A	C	T
Overhang position		1	2	3	4
G allele 10/14	5’ TCAT	G*			
	3’ AGTA	C	A	C	T

Overhang position	1	2	3	4
-------------------	---	---	---	---

**[0971]** Analysis of the maternal DNA for SNP TSC0129188 gave a single band that corresponded to the DNA molecules filled in with ddGTP at position 1 complementary to the overhang, which represented the “G” allele (see FIG. 20, lane 5). No band was detected for adenine allele, indicating that the maternal DNA is homozygous for guanine.

**[0972]** In contrast, analysis of the multiplexed template DNA from the maternal plasma, which contains both maternal DNA, and fetal DNA, gave two distinct bands (see FIG. 20, lane 6). The lower molecular weight band corresponded to the “G” allele, while the higher molecular weight corresponded to the “A” allele. The “A” allele represents the fetal DNA. Thus, a method has been developed that allows separation of maternal DNA and fetal DNA signals without the added complexity of having to isolate fetal cells. In addition, a sample of paternal DNA is not required to detect differences between the maternal DNA and the fetal DNA.

**[0973]** Analysis of the maternal DNA for SNP TSC0501389 gave a single band that migrated at the higher molecular weight position, which corresponded to the “A” allele. No band was detected that corresponded to the “G” allele. Similarly, analysis of the multiplexed template DNA from the maternal plasma for SNP TSC0501389 gave a single band that migrated at the higher molecular weight position, which corresponded to the “A” allele. Both the maternal template DNA and the fetal template DNA were homozygous for adenine at SNP TSC0501389.

**[0974]** The maternal DNA and the template DNA from the plasma originated from the same sample. One sample, which was obtained through a non-invasive procedure, provided a genetic fingerprint for both the mother and the fetus.

**[0975]** Of the twenty-nine SNPs for which the maternal template DNA was homozygous, the fetal template DNA was heterozygous at two of the twenty-nine SNPs. The fetal DNA was homozygous for the same allele as the maternal template DNA at the remaining 27 SNPs (data not shown). Comparing the homozygous allele of the maternal template DNA and the plasma template DNA at a given SNP provides an added level of quality control. It is not possible that the maternal template DNA and the plasma template DNA are homozygous for different alleles at the same SNP. If this is seen, it would indicate that an error in processing had occurred.

**[0976]** The methods described herein demonstrate that the maternal genetic signal can be separated and distinguished from the fetal genetic signal in a maternal plasma sample. The above-example analyzed SNPs located on chromosome 13, however any chromosome can be analyzed including human



chromosome 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X and Y and fetal chromosomes 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X and Y.

**[0977]** In addition, the methods described herein can be used to detect fetal DNA in any biological sample including but not limited to cell, tissue, blood, serum, plasma, saliva, urine, tears, vaginal secretions, umbilical cord blood, chorionic villi, amniotic fluid, embryonic tissues, lymph fluid, cerebrospinal fluid, mucosa secretions, peritoneal fluid, ascitic fluid, fecal matter, or body exudates.

**[0978]** The methods described herein demonstrate that the percentage of fetal DNA in the maternal sample can be determined by analyzing SNPs wherein the maternal DNA is homozygous, and the DNA isolated from the plasma of the pregnant female is heterozygous. The percentage of fetal DNA can be used to determine if the fetal genotype has any chromosomal disorders.

**[0979]** For example, if the percentage of fetal DNA present in the sample is calculated to be 30% by analysis of chromosome 1 (chromosomal abnormalities involving chromosome 1 terminate early in the pregnancy), then any deviation from 30% fetal DNA is indicative of a chromosomal abnormality. For example, if upon analysis of a SNP or multiple SNPs on chromosome 18, the percentage of fetal DNA is higher than 30%, this would indicate that an additional copy of chromosome 18 is present. The calculated percentage of fetal DNA from any chromosome can be compared to any other chromosome. In particular, the percentage of fetal DNA on chromosome 13 can be compared to the percentage of fetal DNA on chromosomes 18 and 21.

**[0980]** This analysis is assisted by knowledge of the expected ratio of one allele to the other allele at each SNP. As discussed in Example 9, not all heterozygous SNPs display ratios of 50:50. Knowledge of the expected ratio of one allele to the other reduces the overall number of variable sites that must be analyzed. However, even without knowledge of the expected ratios for the various SNPs, the percentage of fetal DNA can be calculated by analyzing a large number of SNPs. When the sampling size of SNPs is large enough, the statistical variation arising from the values of the expected ratios will be eliminated.

**[0981]** In addition, heterozygous maternal SNPs also provide valuable information. The analysis is not limited to homozygous maternal SNPs. For example, if at a heterozygous SNP on maternal DNA, the ratio of allele 1 to allele 2 is 1:1, then in the plasma template DNA the ratio should remain 1:1 unless the fetal DNA carries a chromosomal abnormality.

**[0982]** The above methods can also be used to detect mutations in the fetal DNA including but not limited to point mutations, transitions, transversions, translocations, insertions, deletions, and duplications. As seen in FIG. 20, fetal DNA can readily be distinguished from maternal DNA. The above methods can be used to determine the sequence of any locus of interest for any gene.

#### EXAMPLE 14

**[0983]** Plasma isolated from blood of a pregnant female contains both maternal template DNA and fetal template DNA. As discussed above, fetal chromosomal abnormalities can be determined by analyzing SNPs wherein the maternal template DNA is homozygous and the template DNA obtained from the plasma displays a heterozygous pattern.

**[0984]** For example, assume SNP X can either be adenine or guanine, and the maternal DNA for SNP X is homozygous for guanine. The labeling method described in Example 6 can be used to determine the sequence of the DNA in the plasma sample. If the plasma sample contains fetal DNA, which is heterozygous at SNP X, the following DNA molecules are expected after digestion with the type IIS restriction enzyme BsmF I, and the fill-in reaction with labeled ddGTP, unlabeled dATP, dTTP, and dCTP.

Maternal Allele 1	5' GGGT	<b>G*</b>				
	3'CCCA	C	T	C	A	
Maternal Allele 2	5' GGGT	<b>G*</b>				
	3'CCCA	C	T	C	A	
Fetal Allele 1	5' GGGT	<b>G*</b>				
	3'CCCA	C	T	C	A	
Fetal Allele 2	5' GGGT	A	A	<b>G*</b>		
	3'CCCA	T	T	C	A	

**[0985]** Two signals are seen; one signal corresponds to the DNA molecules filled in with ddGTP at position one complementary to the overhang and the second signal corresponds to the DNA molecules filled in with ddGTP at position three complementary to the overhang. However, the maternal DNA is homozygous for guanine, which corresponds to the DNA molecules filled in at position one complementary to the overhang. The signal from the DNA molecules filled in with ddGTP at position three complementary to the overhang corresponds to the adenine allele, which represents the fetal DNA. This signal becomes a beacon for the fetal DNA, and can be used to measure the amount of fetal DNA present in the plasma sample.

**[0986]** There is no difference in the amount of fetal DNA from one chromosome to another. For instance, the percentage of fetal DNA in any given individual from chromosome 1 is the same as the percentage of fetal DNA from chromosome 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X and Y. Thus, the allele ratio calculated for SNPs on one chromosome can be compared to the allele ratio for the SNPs on another chromosome.

[0987] For example, the allele ratio for the SNPs on chromosome 1 should be equal to the allele ratio for the SNPs on chromosomes 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X, and Y. However, if the fetus has a chromosomal abnormality, including but not limited to a trisomy or monosomy, the ratio for the chromosome that is present in an abnormal copy number will differ from the ratio for the other chromosomes.

[0988] To recapitulate the in vivo scenario of blood from a pregnant female, maternal DNA was mixed with DNA isolated from her child, who previously was diagnosed with Trisomy 21, in various ratios to represent varying percentages of fetal DNA. For example, to replicate the in vivo scenario of 50% fetal DNA in maternal blood, equal amounts of maternal DNA were mixed with DNA isolated from her child with Down's syndrome. The maternal DNA was analyzed to identify homozygous SNPs, and these SNPs then were analyzed using the mixture of 50% maternal DNA and 50% Down's syndrome DNA. The ratio of allele 1 to allele 2 at heterozygous SNPs on chromosome 13 was compared to the ratio of allele 1 to allele 2 at heterozygous SNPs on chromosome 21.

[0989] Four different samples were analyzed: a sample with 100% of the DNA from a child with Down syndrome; a sample with 75% DNA from the child with Down syndrome and 25% DNA from the child's mother; a sample with 50% DNA from the child with Down syndrome and 50% DNA from the child's mother; and a sample with 40% DNA from the child with Down syndrome and 60% DNA from the child's mother. The maternal DNA was analyzed to identify homozygous SNPs. The DNA isolated from the child with Down syndrome was genotyped to identify heterozygous SNPs. Then, the samples were genotyped at SNPs where the maternal DNA was homozygous and the DNA from the child was heterozygous. For each sample, these SNPs were analyzed ten times.

### **Collection of Blood Samples**

[0990] An Internal Review Board approved study was designed to allow collection of blood samples from children afflicted with Down's syndrome and their parents. For this study, blood was collected from the mother, the father, and the child with Down's syndrome. Informed consent to collect blood from the child with Down's syndrome was granted by the parents as well as the child. Blood was collected into 9 ml EDTA Vacuette tubes (catalog number NC9897284). The tubes were stored at 4°C until ready for processing.

### **Isolation of Plasma and Maternal Cells**

[0991] ~~he blood~~The blood was stored at 4°C until processing. The tubes were spun at 1000 rpm for ten minutes in a centrifuge with braking power set to zero. The tubes were spun a second time at 1000

rpm for ten minutes. The supernatant (the plasma) of each sample was transferred to a new tube and spun at 3000 rpm for ten minutes with the brake set to zero. The supernatant was transferred to a new tube and stored at -80°C. Approximately two milliliters of the “buffy coat,” which contains maternal cells, was placed into a separate tube and stored at -80°C.

### **Isolation of DNA**

**[0992]** DNA was isolated from the plasma sample using the Qiagen Midi Kit for purification of DNA from blood cells, following the manufacturer’s instructions (QIAmp DNA Blood Midi Kit, Catalog number 51183). DNA was eluted in 100 µl of distilled water. The Qiagen Midi Kit also was used to isolate DNA from the maternal cells contained in the “buffy coat.” Maternal DNA and the plasma DNA were isolated from the same tube of blood.

### **Identification of Maternal Homozygous SNPs**

**[0993]** The maternal DNA was genotyped to identify homozygous SNPs. Seven hundred and sixty-eight SNPs on chromosome 13 and 768 SNPs on chromosome 21 were genotyped using the methods described in Example 6. Any number of SNPs can be analyzed, and the SNPs can be located on human chromosome 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X and Y. Preferably, the SNPs that are genotyped have allele frequencies of 50:50, 60:40, 70:30, 80:20, or 90:10. As described in Example 8, the allele frequency of any given SNP can be determined.

**[0994]** Details regarding the SNPs located on chromosome 13 and 21 can be found at the SNP consortium database, which can be accessed via the internet at <http://www.snp.cshl.org>. The primers were designed following the procedures set forth in the Examples described above, for example, in Examples 1, 2, 3, 5, and 6.

**[0995]** The first primers were designed so that after digestion with a Type II enzyme, the products had different molecular weights as described in Example 6. This allowed the amplified products to be pooled, and run in a single lane of a gel.

**[0996]** For example, the first primer can be designed such that after digestion a 30 base pair product is generated. Likewise, the first primer of a different locus of interest can be designed such that after digestion a 40 base pair product is generated. The first primers can be designed so that in a single reaction, numerous loci can be analyzed in one lane of a gel (30, 40, 50, 60, 70, 80, 90, 100, 110, and 120 base pair products can be run in a single lane). The first primer can be designed to anneal any distance from the locus of interest including but not limited to between 5-10, 10-25, 25-50, 50-75, 75-100,

100-150, 150-200, 200-250, 250-300, 300-350, 350-400, 400-450, 450-500, 500-550, 550-600, 600-650, 650-700, 700-750, 750-800, 800-850, 850-900, 900-950, 950-1000 and greater than 1000 bases.

### **Amplification of the Loci of Interest**

**[0997]** For each SNP that was genotyped, a PCR reaction was used to amplify the loci of interest. The PCR reactions were performed in 96-well plates. The first and second primer (3  $\mu$ l of 1.25  $\mu$ M stock concentration) for each SNP was distributed into a well of a microtiter plate. Eight 96-well PCR plates were set-up for chromosome 21 and eight 96-well plates were set-up for chromosome 13. After the primers had been distributed into the wells of the microtiter plates, a mixture containing the genomic DNA and HotStar PCR reagents was added to each well. Each PCR reaction contained 3  $\mu$ l of each primer, 7.5  $\mu$ l of HotStar Taq Master mix, 0.5  $\mu$ l of water, and 1  $\mu$ l of genomic DNA (10 ng/ $\mu$ l).

**[0998]** The PCR cycling conditions were as follows:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 37°C for 30 seconds;
- (3) 95°C for 30 seconds;
- (4) 52°C for 30 seconds;
- (5) 95°C for 30 seconds;
- (6) 58°C for 30 seconds;
- (7) 95°C for 30 seconds;
- (8) Repeat steps 6 and 7 thirty seven (37) times;
- (9) 72°C for 5 minutes.

### **Purification of Fragment of Interest**

**[0999]** After the PCR reaction, 3  $\mu$ l of a PCR product generated with a first primer designed to produce a 30 base pair product, 3  $\mu$ l of a PCR product generated with a first primer designed to produce a 40 base pair product, 3  $\mu$ l of a PCR product generated with a first primer designed to produce a 50 base pair product, 3  $\mu$ l of a PCR product generated with a first primer designed to produce a 60 base pair product, 3  $\mu$ l of a PCR product generated with a first primer designed to produce a 70 base pair product, 3  $\mu$ l of a PCR product generated with a first primer designed to produce a 80 base pair product, 3  $\mu$ l of a PCR product generated with a first primer designed to produce a 90 base pair product, 3  $\mu$ l of a PCR product generated with a first primer designed to produce a 100 base pair product were mixed together in a well of a Streptawell, transparent, High-Bind plate from Roche Diagnostics GmbH (catalog number 1

645 692, as listed in Roche Molecular Biochemicals, 2001 Biochemicals Catalog). The first primers contained a 5' biotin tag so the PCR products bound to the Streptavidin coated wells while the genomic template DNA did not. The streptavidin binding reaction was performed using a Thermomixer (Eppendorf) at 1000 rpm for 20 min. at 37°C. Each well was aspirated to remove unbound material, and washed three times with 1X PBS, with gentle mixing (Kandpal et al., Nucl. Acids Res. 18:1789-1795 (1990); Kaneoka et al., Biotechniques 10:30-34 (1991); Green et al., Nucl. Acids Res. 18:6163-6164 (1990)).

#### **Restriction Enzyme Digestion of Isolated Fragments**

[1000] The purified PCR products were digested with the restriction enzyme BsmF I, which binds to the recognition site incorporated into the PCR products from the second primer. The digests were performed in the Streptawells following the instructions supplied with the restriction enzyme. After digestion, the wells were washed three times with PBS to remove the cleaved fragments.

#### **Incorporation of Labeled Nucleotide**

[1001] The restriction enzyme digest with BsmF I yielded a DNA fragment with a 5' overhang, which contained the SNP site or locus of interest and a 3' recessed end. The 5' overhang functioned as a template allowing incorporation of a nucleotide or nucleotides in the presence of a DNA polymerase. As discussed in detail in Example 6, a single nucleotide labeled with one chemical moiety can be used to determine the sequence at a SNP.

[1002] The amplified loci of interest were pooled into the streptavidin-well based on size, and on the nucleotide used in the fill-in reaction. The sequence of SNPs that were determined by using a guanine nucleotide were pooled together. Likewise, the sequence of SNPs that were determined by using an adenine nucleotide were pooled together; the sequence of SNPs that were determined by using a thymidine nucleotide were pooled together; and the sequence of SNPs that were determined by using a cytosine nucleotide were pooled together.

[1003] Thus, a typical fill-in reaction contained 8 amplified loci, ranging in size of 30-120 base pair products; the sequence of all eight was determined using a single nucleotide labeled with one chemical moiety. Any number of amplified loci can be pooled together.

[1004] The following components were added to each fill in reaction: 1 µl of fluorescently labeled dideoxynucleotide (ddGTP for G fill-in reactions; ddATP for A fill-in reactions; ddTTP for thymidine fill-in reactions; and ddCTP for cytosine fill-in reactions), 0.5 µl of unlabeled dNTPs (40 µM), which

contained all nucleotides except the labeled nucleotide, 2 µl of 10X sequenase buffer, 0.25 µl of Sequenase, and water as needed for a 20µl reaction.

[1005] The fill in reaction was performed at 40°C for 10 min. Non-fluorescently labeled dNTP was purchased from Fermentas Inc. (Hanover, MD). All other labeling reagents were obtained from Amersham (Thermo Sequenase Dye Terminator Cycle Sequencing Core Kit, US 79565).

[1006] After labeling, each Streptawell was rinsed with 1X PBS (100 µl) three times. The "filled in" DNA fragments were then released from the Streptawells by digestion with the restriction enzyme EcoRI, according to the manufacturer's instructions that were supplied with the enzyme. Digestion was performed for 1 hour at 37 °C with shaking at 120 rpm.

### **Detection of the Locus of Interest**

[1007] After release from the streptavidin matrix, the sample was loaded into a lane of a 36 cm 5% acrylamide (urea) gel (BioWhittaker Molecular Applications, Long Ranger Run Gel Packs, catalog number 50691). The sample was electrophoresed into the gel at 3000 volts for 3 min. The gel was run for 3 hours on a sequencing apparatus (Hoefer SQ3 Sequencer). The gel was removed from the apparatus and scanned on the Typhoon 9400 Variable Mode Imager. The incorporated labeled nucleotide was detected by fluorescence. The homozygous SNPs were identified.

### **Identification of Heterozygous SNPs with the Trisomy 21 template**

[1008] The DNA isolated from the individual with Down syndrome (the child of the mother who was genotyped above) was analyzed to identify heterozygous SNPs. The same seven hundred and sixty-eight SNPs on chromosome 13 and the same 768 SNPs on chromosome 21 that were analyzed with the maternal DNA were genotyped using the methods described in Example 6. Any number of SNPs can be analyzed, and the SNPs can be located on human chromosome 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X or Y. Preferably, the SNPs that are genotyped have allele frequencies of 50:50, 60:40, 70:30, 80:20, or 90:10. As described in Example 8, the allele frequency of any given SNP can be determined.

[1009] The process for genotyping the SNPs with the DNA isolated from the individual with Down syndrome was as described for the maternal DNA. The heterozygous SNPs were identified.

[1010] SNPs that were homozygous for the maternal DNA and heterozygous for the DNA isolated from the individual with Down syndrome were further analyzed using samples that contained mixtures of maternal DNA and Down syndrome DNA.

### **Generation of Samples containing maternal DNA and Down syndrome DNA**

[1011] The DNA and the DNA obtained from her child, who has Down's syndrome, were quantitated using a spectrophotometer. The maternal DNA and the child's DNA were mixed together at various percentages to represent the situation of circulating fetal DNA in the maternal blood. The following percentages were analyzed: 100% Down's syndrome DNA, 75% Down's syndrome DNA, 50% Down's syndrome DNA, and 40% Down's syndrome DNA.

[1012] The ratio at each heterozygous SNP was calculated by dividing the value obtained for allele 1 by the value obtained for allele 2. For example, if SNP X can either be adenine (A) or guanine (G), the ratio at SNP X was calculated by dividing the value obtained for adenine by the value obtained for guanine.

[1013] For the sample containing 100% Down syndrome DNA, sixty-two SNPs on chromosome 13, which were homozygous with the maternal DNA and heterozygous with the DNA isolated from the individual with Down syndrome, were analyzed. For chromosome 21, forty-nine SNPs were analyzed that were homozygous with the maternal DNA and heterozygous with the DNA isolated from the individual with Down syndrome.

[1014] The 62 SNPs on chromosome 13 and 49 SNPs on chromosome 21 were analyzed ten separate times. As shown in Table XX, for each of the ten trials, the ratio of allele 1 to allele 2 on chromosome 13 was approximately 1.0 as expected. For chromosome 13, there is one copy of allele 1 and one copy of allele 2. The average of the ten trials was 1.051 with a standard deviation of 0.085.

[1015] With a Trisomy 21, there are two copies of one allele, which are usually inherited from the mother, and one copy of the other allele. The expected ratio is approximately 0.5 (one copy of allele 1/ two copies of allele 2). As shown in Table XX, the ratio for chromosome 21 varied from a low of 0.462 to a high of 0.634. For every trial, the ratio obtained for chromosome 21 was significantly distinct from the ratio obtained at chromosome 13. The average ratio for the ten trials was 0.531 with a standard deviation of 0.049.

[1016] The experiment was repeated ten times so that a true statistical measurement could be obtained. If ten different genetic samples were used, the SNPs that fit the criteria (maternal homozygous, Down syndrome child heterozygous) would be different, making it difficult to compare from sample to sample.

[1017] Statistical analysis revealed a confidence value of 99.9% that the ratios obtained on chromosome 13 and on chromosome 21 represented true differences, rather than random numerical fluctuations in value. The Ravgen method identified the presence of the chromosomal abnormality.



[1018] For the sample containing 75% Down syndrome DNA and 25% maternal DNA, sixty two SNPS on chromosome 13 and fifty SNPs on chromosome 21 were analyzed, unless stated otherwise. For various trials, not all the SNPS could be quantitated because the bands corresponding to certain SNPs were faint. This may have been caused by poor PCR amplification, poor binding to the streptavidin plate, or a weak fill-in reaction.

[1019] For trial 3, 61 SNPs on chromosome 13 were analyzed. For trial 4, 49 SNPs were analyzed on chromosome 21. With regard to trial 5, 47 SNPs on chromosome 21 were analyzed and 61 SNPs on chromosome 13. For trial 7, 49 SNPs were analyzed on chromosome 21 and 61 SNPs on chromosome 13. For trial 8, 49 chromosomes were analyzed on chromosome 21, and 59 SNPs were analyzed on chromosome 13. For trials 9 and 10, 59 SNPs on chromosome 13 were analyzed.

[1020] The expected ratio on chromosome 13 for a heterozygous SNP is 0.6. If the maternal chromosomes both contain an adenine nucleotide, and the Down syndrome genome is comprised of one chromosome with an adenine nucleotide and one chromosome with a guanine nucleotide, then the ratio of G:A is  $0.75 / (.75 \text{ (Down syndrome A allele)} + 0.25 + 0.25 \text{ (maternal A alleles)})$ , which is 0.6. For the ten trials, the ratios obtained for chromosome 13 varied from 0.567 to 0.645. The average for the ten trials was 0.609 with a standard deviation of 0.032 (see Table XX).

[1021] The expected ratio for chromosome 21 in a Trisomy condition is 0.375. If the maternal chromosomes both contain an adenine nucleotide, and the Down syndrome genome is comprised of two chromosomes with an adenine nucleotide and one chromosome with a guanine nucleotide, then the ratio of G:A is  $0.75 / (0.75 + 0.75 \text{ (Down syndrome A alleles)} + 0.25 + 0.25 \text{ (maternal A alleles)})$ , which is 0.375.

[1022] For the ten trials, the ratios obtained for chromosome 21 varied from 0.350 to 0.4125, with an average of 0.384 and a standard deviation of 0.017 (see Table XX). Statistical analysis revealed a confidence value of 99.9% that the ratios obtained on chromosome 13 and on chromosome 21 represented true differences, rather than random numerical fluctuations in value. The Ravgen method identified the presence of the chromosomal abnormality in the presence of 25% maternal DNA.

[1023] With regard to the sample containing 50% Down syndrome DNA, 46 SNPs on chromosome 13 and 35 SNPs on chromosome 21 were analyzed, unless stated otherwise. For trial 1, 45 SNPs on chromosome 13 were analyzed. For trial 2, 44 SNPs on chromosome 13 were analyzed. For trial 3, 42 SNPs on chromosome 13 were analyzed. For trial 4, 44 SNPs on chromosome 13 and 34 SNPs on chromosome 21 were analyzed. For trial 5, 34 SNPs on chromosome 21 were analyzed. For trials 7 and 8, 44 and 41 SNPs on chromosome 13, respectively, were analyzed. For trial 9, 44 SNPs on chromosome

13 and 34 SNPs on chromosome 21 were analyzed. For trial 10, 44 SNPs on chromosome 13 were analyzed.

**[1024]** The expected ratio at a heterozygous SNP on chromosome 13 for the 50% sample is 0.33. If the maternal chromosomes both contain an adenine nucleotide, and the Down syndrome genome is comprised of one chromosome with an adenine nucleotide and one chromosome with a guanine nucleotide, then the ratio of G:A is  $0.50 / (.50 \text{ (Down syndrome A allele )} + 0.50 + 0.50 \text{ (maternal A alleles)})$ , which is 0.33. For the ten trials, the ratios obtained for chromosome 13 varied from 0.302 to 0.347. The average for the ten trials was 0.324 with a standard deviation of 0.013 (see Table XX).

**[1025]** The expected ratio for chromosome 21 in a Trisomy condition is 0.25. If the maternal chromosomes both contain an adenine nucleotide, and the Down syndrome genome is comprised of two chromosomes with an adenine nucleotide and one chromosome with a guanine nucleotide, then the ratio of G:A is  $0.50 / (0.50 + 0.50 \text{ (Down syndrome A alleles)} + 0.50 + 0.50 \text{ (maternal A alleles)})$ , which is 0.25.

**[1026]** For the ten trials, the ratios obtained for chromosome 21 varied from 0.230 to 0.275, with an average of 0.244 and a standard deviation of 0.015 (see Table XX). Statistical analysis revealed a confidence value of 99.1% that the ratios obtained on chromosome 13 and on chromosome 21 represented true differences, rather than random numerical fluctuations in value. The Ravgen method identified the presence of the chromosomal abnormality in the presence of 50% maternal DNA.

**[1027]** For the sample containing 40% Down syndrome DNA, 60 SNPs on chromosome 13 and 48 SNPs on chromosome 21 were analyzed, unless stated otherwise. For trial 1, 47 SNPs on chromosome 21 were analyzed. For trials 2-4, 59 SNPs on chromosome 13 and 47 SNPs on chromosome 21 were analyzed. For trials 5 and 6, 46 SNPs on chromosome 21 were analyzed. For trial 7, 58 SNPs on chromosome 13 were analyzed. For trial 8, 46 SNPs on chromosome 21 were analyzed and for trials 9 and 10, 47 SNPs on chromosome 21 were analyzed.

**[1028]** The expected ratio at a heterozygous SNP on chromosome 13 for the 40% Down syndrome DNA sample is 0.25. If the maternal chromosomes both contain an adenine nucleotide, and the Down syndrome genome is comprised of one chromosome with an adenine nucleotide and one chromosome with a guanine nucleotide, then the ratio of G:A is  $0.40 / (.40 \text{ (Down syndrome A allele )} + 0.60 + 0.60 \text{ (maternal A alleles)})$ , which is 0.25. For the ten trials, the ratios obtained for chromosome 13 varied from 0.254 to 0.285. The average for the ten trials was 0.269 with a standard deviation of 0.009 (See Table XX).

**[1029]** The expected ratio for chromosome 21 in a Trisomy condition is 0.20. If the maternal chromosomes both contain an adenine nucleotide, and the Down syndrome genome is comprised of two

chromosomes with an adenine nucleotide and one chromosome with a guanine nucleotide, then the ratio of G:A is  $0.40/(0.40 + 0.40 \text{ (Down syndrome A alleles)} + 0.60 + 0.60 \text{ (maternal A alleles)})$ , which is 0.20.

[1030] For the ten trials, the ratios obtained for chromosome 21 varied from 0.216 to 0.249, with an average of 0.23 and a standard deviation of 0.011 (see Table XX). Statistical analysis revealed a confidence value of 94.3% that the ratios obtained on chromosome 13 and on chromosome 21 represented true differences, rather than random numerical fluctuations in value. The Ravgen method identified the presence of the chromosomal abnormality in the presence of 60% maternal DNA.

[1031] The presence of the Trisomy 21 condition was identified with the Ravgen method in numerous samples that contained various percentages of abnormal DNA. Each percentage of abnormal DNA was analyzed ten separate times and each time, the presence of the abnormal condition was identified. The ratio of allele 1 to allele 2 at multiple heterozygous SNPs on chromosome 13 was calculated, and the ratios were averaged. The same was done with the SNPs located on chromosome 21. The ratio obtained for the heterozygous SNPs on chromosome 13 was statistically different from the ratio obtained on chromosome 21. The ratios obtained on both chromosome 13 and 21 were near the mathematically predicted values.

[1032] In this example, the confidence interval for the samples with 100% Down syndrome DNA and 75% Down syndrome DNA was 99.9%, and the confidence interval for the sample with 50% Down syndrome DNA was 99.1%, which is about the accuracy reported for amniocentesis. The confidence interval for the sample containing 40% Down syndrome DNA was 94.3%, which is more accurate than currently marketed non-invasive tests for prenatal diagnostics.

[1033] As discussed above, about 60 SNPs on chromosome 13 and 50 SNPs on chromosome 21 were analyzed. To increase the confidence interval for samples containing 40% fetal DNA or lower, a larger number of SNPs can be analyzed. The Ravgen method provides a highly accurate, cost-effective way to sequence DNA, so sequencing a larger number of SNPs is not difficult. The accuracy of the test is determined by the number of SNPs that are sequenced. For higher accuracy with samples that contain lower percentages of DNA, more SNPs can be analyzed. Alternatively, the methods described in this application can be used to ensure that the samples contain a higher percentage of fetal DNA.

[1034] In this example, a sample containing 40% Down syndrome DNA, which represented the fetal DNA in the maternal blood, was analyzed. Maternal blood samples with any percentage of fetal DNA can be analyzed including but not limited to 0.0001-1%, 1-10%, 10-20%, 20-30%, 30-40%, 40-50%, 50-60%, 60-70%, 70-80%, 80-90%, and 90-100%.

Table XX.

**The Ravgen method identifies chromosomal abnormalities in  
Samples containing 40% Down syndrome DNA**

	100% DS DNA			~ 75% DS DNA		
Ratio at Chrom.	13	21	Expected 21	13	21	Expected 21
Trial 1	0.959	0.5708	0.49	0.637	.3764	0.389
Trial 2	0.916	0.5024	0.48	0.567	.3894	0.362
Trial 3	1029	0.4616	0.51	0.651	.3707	0.394
Trial 4	0.967	0.5123	0.491	0.580	.3901	0.367
Trial 5	1.037	0.6339	0.51	0.645	.4125	0.392
Trial 6	1.111	0.5425	0.53	0.645	.3743	0.392
Trial 7	1.154	0.495	0.54	0.594	.3974	0.373
Trial 8	1.135	0.5276	0.532	0.583	.3901	0.368
Trial 9	1.148	0.5619	0.534	0.579	.3899	0.367
Trial 10	1.057	0.4976	0.52	0.609	.350	0.378
AVG.	1.051	.531	0.512	.609	0.384	0.378
STDEV	.085	.049		.032	.017	
	~ 50% DS DNA			~ 40% DS DNA		
Ratio at Chrom.	13	21	Expected 21	13	21	Expected 21
Trial 1	0.347	0.275	0.258	0.277	0.239	0.217
Trial 2	0.316	0.237	0.24	0.265	0.249	0.21
Trial 3	0.338	0.247	0.253	0.266	0.227	0.21
Trial 4	0.331	0.264	0.249	0.254	0.216	0.202
Trial 5	0.330	0.241	0.248	0.274	0.246	0.215
Trial 6	0.324	0.240	0.244	0.268	0.22	0.211
Trial 7	0.318	0.233	0.241	0.275	0.227	0.216
Trial 8	0.302	0.230	0.231	0.258	0.228	0.21
Trial 9	0.315	0.238	0.240	0.285	0.231	0.222
Trial 10	0.318	0.235	0.241	0.266	0.218	0.21
AVG.	0.324	0.244	0.244	0.269	0.23	0.212
STDEV	0.013	0.015		0.009	0.011	

#### EXAMPLE 15

[1035] As discussed in Example 4 above, the use of cell lysis inhibitors, cell membrane stabilizers, or cross-linking reagents can be used to increase the percentage of fetal DNA in the maternal blood. In this example, methods for the isolation of free fetal DNA are disclosed, which minimize the amount of maternal cell lysis. The effect of formalin on sixty-nine (69) maternal blood samples from twenty-seven clinical practices located in sixteen different states was analyzed. Formalin was added to all samples

collected from the pregnant women, and the percentage of fetal DNA was calculated using serial dilution analysis followed by PCR. A genetic marker on the Y chromosome was used to calculate the percent of fetal DNA.

### **Collection of Blood Samples**

**[1036]** In accordance with an IRB approved study, blood samples were collected from pregnant women after informed consent had been granted. Blood samples were received from 27 different clinical sites operating in 16 different states located throughout the U.S. Blood samples were collected from both women carrying male and female fetuses, however, here, we report results obtained from woman carrying male fetuses, as the Y chromosome is the accepted marker when quantitating percentages of fetal DNA.

**[1037]** Blood is collected by any method or process that results in a substantial increase in the ratio of fetal DNA/maternal DNA in the resulting serum or plasma after appropriate processing. As used herein, a substantial increase in the ratio of fetal DNA/maternal DNA is that which can be detected by the methods as described herein. Such methods or processes typically result in a substantial increase in the ratio of fetal DNA/maternal DNA of about 5%, 10%, 15%, 20%, 30%, 50%, 70%, 80%, 100% or more of the ratio of fetal DNA/maternal DNA found in blood samples collected by standard procedures.

**[1038]** In other embodiments, blood is collected by any method or process that results in a substantial increase in the amount of free fetal DNA compared to the amount of total DNA recovered or detected in the resulting serum or plasma after processing. Such methods or processes typically result in a substantial increase so the fetal DNA recovered or detected is about 10%, 15%, 20%, 25%, 30%, 40%, 50% or more of the total DNA recovered or detected in the processed plasma or serum sample.

**[1039]** All clinical sites were provided with a kit used for the venipuncture procedure, which included 21 gauge needles, 9 ml EDTA Vacuette tubes (catalog number NC9897284) a syringe containing 0.225 ml of 10% neutral buffered solution containing formaldehyde (4% w/v), an icepack, and a shipping container. The clinical sites were instructed to add the formaldehyde immediately after drawing the blood and to gently invert the tubes.

**[1040]** The methods or processes of collecting blood samples may also include other steps that result in lessened or reduced cell lysis. For instance, blood collection devices may be modified to decrease cell lysis due to sheer forces in the collection needle, syringe or tubes used. For instance, needles of large gauge may be employed to reduce cell sheering or vacutainer tubes may be modified to reduce the velocity of blood flow.

### **Isolation of Plasma**

**[1041]** Any method may be used to isolate plasma from the cell components of blood after collection but methods wherein cell lysis is substantially prevented, reduced or inhibited are preferred. The blood was stored at 4°C until processing. Methods for isolation of the plasma were implemented to reduce the amount of maternal cell lysis. The tubes were spun at 1000 rpm for ten minutes in a centrifuge with braking power and acceleration power set at zero to substantially prevent, reduce or inhibit cell lysis and or mixing of blood cell components into the plasma. The tubes were spun a second time at 1000 rpm for ten minutes with braking power (centrifuge stopped by natural deceleration) and acceleration power set to zero. The supernatant (the plasma) of each sample was transferred carefully to a new tube and spun at 3000 rpm for ten minutes with the brake and acceleration power set at zero. The supernatant (the plasma) of each sample was collected via procedures to substantially prevent mixing of cell components into the plasma. Great care was taken to ensure that the buffy-coat was not disturbed. A percentage of the supernatant can be left in the tube including but not limited to 0.001-1%, 1-10%, 10-20%, 20-30%, 30-40%, 40-50%, 50-60%, 60-70%, 70-80% or greater than 80%. In this example, about 0.5 ml of the supernatant was left in the tube to ensure that the buffy-coat was not disturbed. The supernatant was transferred to a new tube and stored at -80°C.

### **Isolation of DNA**

**[1042]** DNA was isolated from the plasma sample using the Qiagen Midi Kit for purification of DNA from blood cells, following the manufacturer's instructions (QIAmp DNA Blood Midi Kit, Catalog number 51183). DNA was eluted in 100 µl of distilled water. However, any method of DNA isolation can be used including cesium chloride gradients, gradients, sucrose gradients, glucose gradients, centrifugation protocols, boiling, Qiagen purification systems, QIA DNA blood purification kit, HiSpeed Plasmid Maxi Kit, QIAfilter plasmid kit, Promega DNA purification systems, MangeSil Paramagnetic Particle based systems, Wizard SV technology, Wizard Genomic DNA purification kit, Amersham purification systems, GFX Genomic Blood DNA purification kit, Invitrogen Life Technologies Purification Systems, CONCERT purification system, Mo Bio Laboratories purification systems, UltraClean BloodSpin Kits, and UltraClean Blood DNA Kit. The skilled artisan understands that the manufacturer's protocols can be modified to increase the yield of DNA. For example, the Qiagen Midi Kit for purification of DNA recommends the use of 1X AL buffer. However, any concentration of AL buffer may be used if the yield of DNA increases including but not limited to 0.1-0.5X AL buffer, 0.5-1X AL buffer, 1X-2X AL buffer, 2-3X AL buffer, 3-4X AL buffer, 4-5X AL buffer, and greater than 5X AL

buffer. The skilled artisan understands that the modifications and manipulations of the reagents are not limited to AL buffer.

### **Quantification of Percentage of Fetal DNA**

[1043] The percentage of fetal DNA present in the maternal plasma sample was calculated using serial dilution analysis followed by PCR. Two different sets of primers were used: one primer set was specific for the Y chromosome, and thus specific for fetal DNA, and the other primer set was designed to amplify the cystic fibrosis gene, which is present on both maternal template DNA and fetal template DNA.

### **Primer Design:**

[1044] The following primers were designed to amplify the SRY gene on the Y chromosome:

Upstream primer:

5' TGGCGATTAAGTCAAATTCGC 3' (SEQ ID NO: 263)

Downstream primer:

5' CCCCTAGTACCCTGACAATGTATT 3' (SEQ ID NO: 264)

[1045] The following primers were designed to amplify the cystic fibrosis gene:

Upstream primer:

5' CTGTTCTGTGATATTATGTGTGGT 3' (SEQ ID NO: 265)

Downstream primer:

5' AATTGTTGGCATTCCAGCATTG 3' (SEQ ID NO: 266)

### **PCR Reaction**

[1046] The SRY gene and the cystic fibrosis gene were amplified from the template genomic DNA using PCR (U.S. Patent Nos. 4,683,195 and 4,683,202). For increased specificity, a "hot-start" PCR was used. PCR reactions were performed using the HotStarTaq Master Mix Kit supplied by Qiagen (Catalog No. 203443). For amplification of the SRY gene, the DNA eluted from the Qiagen purification column was diluted serially 1:2. For amplification of the cystic fibrosis gene, the DNA eluted from the Qiagen purification column was diluted 1:4, and then serially diluted 1:2. The following components were used for each PCR reaction: 8 µl of template DNA (diluted or undiluted), 1 µl of each primer (5 µM), 10 µl of HotStar Taq mix. The following PCR conditions were used:

(1) 95°C for 15'

(2) 94°C for 1'

- (3) 54°C for 15"
- (4) 72°C for 30"
- (5) Repeat steps 2-4 for 45 cycles.
- (6) 10' at 72°C

[1047] Amplification of the SRY gene was performed using the following templates: undiluted, diluted 1:2, diluted 1:4, diluted 1:8, diluted 1:16, diluted 1:32, diluted 1:64, diluted 1:128, diluted 1:256, and diluted 1:512. Amplification of the cystic fibrosis gene was performed using the following templates: diluted 1:4, diluted 1:8, diluted 1:16, diluted 1:32, diluted 1:64, diluted 1:128, diluted 1:256, diluted 1:512, diluted 1:1024, diluted 1:2048, and diluted 1:4096.

[1048] The percent of fetal DNA present in the maternal plasma was calculated using the following formula:

$$\% \text{ fetal DNA} = (\text{amount of SRY gene} / \text{amount of cystic fibrosis gene}) * 2 * 100.$$

[1049] The amount of SRY gene was represented by the highest dilution value in which the gene was amplified. Likewise, the amount of cystic fibrosis gene was represented by the highest dilution value in which it was amplified. The formula contains a multiplication factor of two (2), which is used to normalize for the fact that there is only one copy of the SRY gene (located on the Y chromosome), while there are two copies of the cystic fibrosis gene.

[1050] The effect of formalin on sixty-nine (69) maternal blood samples collected from twenty-seven clinical practices located in sixteen different states, spanning from Washington to Massachusetts is shown in Table XXI. In this study, formalin was added to all samples collected from the pregnant women, and the percentage of fetal DNA was calculated using serial dilution analysis followed by PCR. The serial dilutions and PCR amplifications were performed by four different scientists over a period of five months. The samples were collected from women at gestational ages ranging from 11 weeks to 28 weeks, with the majority of women between 16-19 weeks of gestation. A summary is provided in Table XXIII.

[1051] The average percentage of free fetal DNA for the 69 samples analyzed in the maternal blood was 33.6%. Lo *et al.* reported fetal DNA concentrations of 3.4% in woman in late first to mid-second trimester, which was the gestational age of the majority of women in this study. Thus, the addition of formalin led to approximately a ten-fold increase in the average percentage of fetal DNA.

[1052] While the calculated percentage of fetal DNA in maternal blood is impressive, it is also informative to examine the range of the percentages of fetal DNA observed in this study. About six percent of the women (4/69) had 3.125% of free fetal DNA in the maternal blood, which was the lowest



percentage of fetal DNA observed in this study. Another 10.2% of women had 6.25% fetal DNA, which represents a two-fold increase over the reported average in the literature. The total number of women who had less than 10% fetal DNA in the maternal blood was only 16.0%.

**[1053]** Fifty-eight percent of the women in this study had a percentage of fetal DNA of 25% or greater. Importantly, 26.0% of the women had fifty percent or greater fetal DNA in the maternal blood. Fetal DNA percentages of this magnitude have not been reported, and represent a new tool to the field of prenatal genetics.

**[1054]** There were four samples collected from women at the gestational age of eleven weeks. The percentages of fetal DNA in the maternal blood samples were as follows: two samples at 12.5%; one sample at 25%; and one sample at greater than 50%. Thus, the effect of formalin on the percentages of fetal DNA was observed with samples collected from women in early as well as later gestational periods.

**[1055]** The effect of stabilizing cell membranes and reducing the release of free DNA was not limited to formalin. We have tested several different types of agents, and combinations of agents, that prevent cell lysis and/or stabilize cell membranes, such as glutaraldehyde, and have seen that these agents also reduce the amount of free DNA in the blood sample (data not shown).

**[1056]** The above described methods may also include steps of adding an agent to the blood sample at the time or near to the time of collection to substantially inhibit or impede cell lysis or stabilize cell membranes. Any number of agents that impede cell lysis or stabilize cell membranes or cross-link cell membranes can be added to the maternal blood samples including but not limited to formaldehyde, and derivatives of formaldehyde, formalin, glutaraldehyde, and derivatives of glutaraldehyde, crosslinkers, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA, DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS or the compounds listed in Table XXIII.

Additional cross-linkers that can be used are found at the following website:

[www.piercenet.com/products/](http://www.piercenet.com/products/).

**[1057]** An agent that stabilizes cell membranes may be added to the maternal blood sample to reduce maternal cell lysis including but not limited to aldehydes, urea formaldehyde, phenol formaldehyde, DMAE (dimethylaminoethanol), cholesterol, cholesterol derivatives, high concentrations of magnesium, vitamin E, and vitamin E derivatives, calcium, calcium gluconate, taurine, niacin, hydroxylamine derivatives, bimoclolol, sucrose, astaxanthin, glucose, amitriptyline, isomer A hopane tetral

phenylacetate, isomer B hopane tetral phenylacetate, citicoline, inositol, vitamin B, vitamin B complex, cholesterol hemisuccinate, sorbitol, calcium, coenzyme Q, ubiquinone, vitamin K, vitamin K complex, menaquinone, zonegran, zinc, ginkgo biloba extract, diphenylhydantoin, perftoran, polyvinylpyrrolidone, phosphatidylserine, tegretol, PABA, disodium cromglycate, nedocromil sodium, phenytoin, zinc citrate, mexitil, dilantin, sodium hyaluronate, or polaxamer 188.

**[1058]** Any concentration of agent that stabilizes cell membranes, impedes cell lysis or cross-link cell membranes can be added. In a preferred embodiment, the agent that stabilizes cell membranes, impedes cell lysis, or cross-links cell membranes is added at a concentration that does not impede or hinder subsequent reactions.

**[1059]** While impressive percentages of free fetal DNA in maternal blood samples have been reported, it is thought that higher percentages can be achieved by carefully explaining the importance of the formalin to the physicians. Samples randomly were checked for the presence of formalin and found that about ten percent of the samples did not receive formalin. In addition, aggregates were observed in another ten percent of the samples suggesting that the formalin had not been thoroughly mixed with the collected blood. Thus, while the addition of formalin produced an impressive effect, it is likely that under controlled conditions, the percentage of free fetal DNA may be higher.

**[1060]** In addition, we believe that procedures to minimize hemolysis during the venipuncture procedure and temperature controlled shipping containers (specimens were shipped in a Styrofoam container with ice pack, but there was variation in temperature because samples were shipped from varying distances) may cause a further increase in the percentage of free fetal DNA. Needles designed to reduce hemolysis can be used during the venipuncture procedure.

**[1061]** Also, we hypothesized that procedures for carefully isolating the plasma would help to ensure a minimal amount of maternal DNA in the sample. We implemented procedures, as described above, to reduce cell lysis, such as gentle centrifugation parameters, and allowed the rotors to stop without external force (no brake). Also, we carefully removed the supernatant containing the plasma DNA from the buffy-coat, which contains maternal DNA. These procedures coupled with the addition of formalin to prevent cell lysis resulted in a tremendous increase in the percentage of fetal DNA.

**Table XXI.**

**Formalin increases the percentage of free fetal DNA  
in blood samples collected at numerous clinical sites from women at various stages of gestation.**

Sample	Wks	Fetal Genomes/ ml	% Fetal DNA
--------	-----	-------------------	-------------

	Gestation		
1	16	80	25
2	19	1066	>50
3	17	52	50
4	22	166	25
5	32	457	50
6	19	400	100
7	18	800	100
8	17	100	50
9	16	50	25
10	17	25	12.5
11	16	94.74	12.5
12	16	34.60	50
13	16	22.5	25
14	17	50	12.5
15	17	26.48	12.5
16	17	45.00	25
17	17	94.7	100
18	17	28.13	6.25
19	19	28.13	25
20	20	11.25	12.5
21	15	11.25	12.5
22	11	16.66	12.5
23	18	13.23	25
24	18	12.50	6.25
25	16	112.50	100
26	17	124.13	25
27	14	90.00	50
28	11	100.00	100
29	18	232.00	100
30	19	626.00	100
31	19	112.50	100
32	16	423.50	100
33	16	423.50	25
34	11	105.88	25
35	16	49.60	3.1
36	11	11.84	12.5
37	16	120.00	25
38	18	342.90	100
39	17	51.43	25
40	18	225.00	6.25
41	17	400.00	12.5
42	28	180.00	25
43	17	20.45	12.5
44	18	25.73	25
45	16	68.68	3.1

46	17	218.18	25
47	15	75.00	6.25
48	16	40.58	3.1
49	17	100.00	25
50	17	14.06	12.5
51	22	22.50	12.5
52	15	28.13	12.5
53	17	50.00	3.125
54	18	58.00	50
55	14	100.00	25
56	16	58.08	25
57	16	13.64	12.5
58	16	25.00	6.25
59	20	45.00	25
60	16	23.69	12.5
61	18	5.92	6.25
62	15	28.13	6.25
63	17	50.00	25
63.			
64	16	360.00	50
65	16	25.00	12.5
66	16	48.65	25
67	16	47.38	12.5
68	14	26.45	50
69	17	124.15	25
<b>Average</b>	<b>17</b>	<b>131.15</b>	<b>33.6</b>

**Table XXII.**

**Formalin increases the percentage of free fetal DNA  
in blood samples collected at numerous clinical sites from women at various stages of gestation.**

% Fetal DNA	3.125	6.25	12.5	25	50	Over 50%
Number Women (69)	4	7	18	22	7	11
%	5.8	10.1	26.1	31.9	10.2	15.9

Table XXIII.

A representative list of cross-linkers that can be used to impede maternal cell lysis.

Cross-Linker	Abbreviation
succinimidyl acetylthioacetate	SATA
succinimidyl <i>trans</i> -4-(maleimidylmethyl) cyclohexane-1-carboxylate	SMCC
succinimidyl 3-(2-pyridyldithio)propionate	SPDP
<i>N</i> -((2-pyridyldithio)ethyl)-4-azidosalicylamide	PEAS; AET
4-azido-2,3,5,6-tetrafluorobenzoic acid, succinimidyl ester	ATFB, SE
4-azido-2,3,5,6-tetrafluorobenzoic acid, STP ester, sodium salt	ATFB, STP ester
4-azido-2,3,5,6-tetrafluorobenzyl amine, hydrochloride	
benzophenone-4-isothiocyanate	
benzophenone-4-maleimide	
4-benzoylbenzoic acid, succinimidyl ester	
Disuccinimidylsuberate	DSS
Dithiobis(succinimidylpropionate)	DSP
3,3'-Dithiobis(sulfosuccinimidylpropionate)	DTSSP
Bis[2-(sulfosuccinimidocarbonyloxy)ethyl]sulfone	SULFO BSOECS
Bis[2-(succinimidocarbonyloxy)ethyl]sulfone	BSOECS
Disulfosuccinimidyltartrate	SULFO DST
Disuccinimidyltartrate	DST
Ethylene glycolbis(succinimidylsuccinate)	SULFO EGS
Ethylene glycolbis(sulfosuccinimidylsuccinate)	EGS
1,2-Di[3'-(2'-pyridyldithio)propionamido]butane	DPDPB
Bis(sulfosuccinimidyl)suberate	BSSS
Succinimidyl-4-( <i>p</i> -maleimidophenyl)butyrate	SMPB
Sulfosuccinimidyl-4-( <i>p</i> -maleimidophenyl)butyrate	SULFO SMPB
3-Maleimidobenzoyl- <i>N</i> -hydroxysuccinimide ester	MBS
3-Maleimidobenzoyl- <i>N</i> -hydroxysulfosuccinimide ester	SULFO MBS
<i>N</i> -Succinimidyl(4-iodoacetyl)aminobenzoate	SIAB
<i>N</i> -Sulfosuccinimidyl(4-iodoacetyl)aminobenzoate	SULFO SIAB
Succinimidyl-4-( <i>N</i> -maleimidomethyl)cyclohexane-1-carboxylate	SMCC
Sulfosuccinimidyl-4-( <i>N</i> -maleimidomethyl)cyclohexane-1-carboxylate	SULFO SMCC
Succinimidyl-6-[3-(2-pyridyldithio)propionamido]hexanoate	NHS LC SPDP
Sulfosuccinimidyl-6-[3-(2-pyridyldithio)propionamido]hexanoate	SULFO NHS LC SPDP
<i>N</i> -Succinimidyl-3-(2-pyridyldithio)propionate	SPDP
<i>N</i> -Hydroxysuccinimidylbromoacetate	NHS BROMOACETATE
<i>N</i> -Hydroxysuccinimidyliodoacetate	NHS IODOACETATE
4-( <i>N</i> -Maleimidophenyl)butyric acid hydrazide hydrochloride	MPBH
4-( <i>N</i> -Maleimidomethyl)cyclohexane-1-carboxylic acid hydrazide hydrochloride	MCCH
<i>m</i> -Maleimidobenzoic acid hydrazidehydrochloride	MBH
<i>N</i> -(epsilon-Maleimidocaproyloxy)sulfosuccinimide	SULFO EMCS
<i>N</i> -(epsilon-Maleimidocaproyloxy)succinimide	EMCS
<i>N</i> -( <i>p</i> -Maleimidophenyl)isocyanate	PMPI

N-(kappa-Maleimidoundecanoic acid) hydrazide	KMUH
Succinimidyl-4-(N-maleimidomethyl)-cyclohexane-1-carboxy(6-amidocaproate)	LC SMCC
N-(gamma-Maleimidobutryloxy)sulfosuccinimide ester	SULFO GMBS
Succinimidyl-6-(beta-maleimidopropionamidohexanoate)	SMPH
N-(kappa-Maleimidoundecanoyloxy)sulfosuccinimide ester	SULFO KMUS
N-(gamma-Maleimidobutryloxy)succinimide	GMBS
Dimethyladipimide hydrochloride	DMA
Dimethylpimelimide hydrochloride	DMP
Dimethylsuberimide hydrochloride	DMS
Methyl-p-hydroxybenzimidate hydrochloride, 98%	MHBH(Wood's Reagent)
Amine Reactive	
Bis[sulfosuccinimidyl] suberate	BS3
Bis[2-(succinimidooxycarbonyloxy)ethyl]sulfone	BSOCOES
Disuccinimidyl glutarate	DSG
	DSP (Lomant's Reagent)
1,5-Difluoro-2,4-dinitrobenzene	DFDNB
Dithiobis[succinimidyl]propionate	DTBP
Bis-[b-(4-Azidosalicylamido)ethyl]disulfide	BASED
Sulphydryl Reactive	
BM[PEO] <sub>3</sub> (1,8-bis-Maleimidotriethyleneglycol	BM[PEO] <sub>3</sub>
BM[PEO] <sub>4</sub> (1,11-bis-Maleimidotetraethyleneglycol	BM[PEO] <sub>4</sub>
1,4-bis-Maleimidobutane	BMB
1,4 bis-Maleimidyl-2,3-dihydroxybutane	BMDB
Bis-Maleimidoethane	BMH
Bis-Maleimidoethane	BMOE
1,4-Di-[3'-(2'-pyridyldithio)-propionamido]butane	DPDPB
Dithio-bis-maleimidoethane	DTME
1,6-Hexane-bis-vinylsulfone	HBVS
p-Azidobenzoyl hydrazide	ABH
Amine-Sulphydryl Reactive	
N-[alpha-Maleimidoacetoxy]succinimide ester	AMAS
N-[4-(p-Azidosalicylamido) butyl]-3'-(2'-pyridyldithio)propionamide	APDP
N-[beta-Maleimidopropyl]oxy]succinimide ester	BMPS
N-e-Maleimidocaproic acid	EMCA
N-e-Maleimidocaproyloxy]succinimide ester	EMCS
N-[gamma-Maleimidobutyryloxy]succinimide ester	GMBS
N-k-Maleimidoundecanoic acid	KMUA
Succinimidyl-4-(N-Maleimidomethyl)cyclohexane-1-carboxy-(6-amidocaproate)	LC-SMCC
Succinimidyl 6-(3-[2-pyridyldithio]-propionamido)hexanoate	LC-SPDP
m-Maleimidobenzoyl-N-hydroxysuccinimide ester	MBS
Succinimidyl 3-[bromoacetamido]propionate	SBAP
N-Succinimidyl iodoacetate	SIA

<i>N</i> -Succinimidyl[4-iodoacetyl]aminobenzoate	SIAB
Succinimidyl 4-[ <i>N</i> -maleimidomethyl]cyclohexane-1-carboxylate	SMCC
Succinimidyl 4-[ <i>p</i> -maleimidophenyl]butyrate	SMPB
Succinimidyl-6-[ $\beta$ -maleimidopropionamido]hexanoate	SMPH
4-Succinimidylloxycarbonyl-methyl-a-[2-pyridyldithio]toluene	SMPT
<i>N</i> -Succinimidyl 3-[2-pyridyldithio]-propionamido	SPDP
<i>N</i> -e-Maleimidocaproyloxy]sulfosuccinimide ester	Sulfo-EMCS
<i>N</i> -[ $\gamma$ -Maleimidobutyryloxy]sulfosuccinimide ester	Sulfo-GMBS
<i>N</i> -[ $\kappa$ -Maleimidoundecanoyloxy]sulfosuccinimide ester	Sulfo-KMUS
4-Sulfosuccinimidyl-6-methyl-a-(2-pyridyldithio)toluamido]hexanoate	Sulfo-LC-SMPT
Sulfosuccinimidyl 6-(3'-[2-pyridyldithio]-propionamido)hexanoate	Sulfo-LC-SPDP
m-Maleimidobenzoyl- <i>N</i> -hydroxysulfosuccinimide ester	Sulfo-MBS
<i>N</i> -Sulfosuccinimidyl[4-iodoacetyl]aminobenzoate	Sulfo-SIAB
Sulfosuccinimidyl 4-[ <i>N</i> -maleimidomethyl]cyclohexane-1-carboxylate	Sulfo-SMCC
Sulfosuccinimidyl-4-( <i>P</i> -Maleimidophenyl) Butyrate	Sulfo-SMPB
Amino Groups	
<i>N</i> -5-Azido-2-nitrobenzoyloxysuccinimide	ANB-NOS
Methyl <i>N</i> -succinimidyl adipate	MSA
<i>N</i> -Hydroxysuccinimidyl-4-azidosalicylic acid	NHS-ASA
<i>N</i> -Succinimidyl(4-azidophenyl)-1, 3'-dithiopropionate	SADP
Sulfosuccinimidyl 2-[7-amino-4-methylcoumarin-3-acetamido]ethyl-1,3'-dithiopropionate	SAED
Sulfosuccinimidyl 2[m-azido- <i>o</i> -nitrobenzamido]-ethyl-1,3'-dithiopropionate	SAND
<i>N</i> -Succinimidyl-6-[4'-azido-2'-nitrophenylamino] hexanoate	SANPAH
Sulfosuccinimidyl-2-[ <i>p</i> -azidosalicylamido]ethyl-1,3'-dithiopropionate	SASD
Sulfosuccinimidyl-[perfluoroazidobenzamido]ethyl-1,3'-dithiopropionate	SFAD
<i>N</i> -Hydroxysulfosuccinimidyl-4-azidobenzoate	Sulfo-HSAB
Sulfosuccinimidyl[4-azidosalicylamido]-hexanoate	Sulfo-NHS-LC-ASA
<i>N</i> -Sulfosuccinimidyl(4-azidophenyl)-1, 3'-dithiopropionate	Sulfo-SADP
<i>N</i> -Sulfosuccinimidyl-6-[4'-azido-2'-nitrophenylamino] hexanoate	Sulfo-SANPAH
<i>p</i> -Azidophenyl glyoxal monohydrate	APG
<i>N</i> - $\beta$ -Maleimidopropionic acid	BMPA
Carbohydrate Reactive-Photoreactive	
<i>N</i> -Succinimidyl-S-acetylthiopropionate	SATP
Sulphydryl-Carbohydrate Reactive	
4-(4- <i>N</i> -Maleimidophenyl)butyric acid hydrazide hydrochloride	MPBH
3-(2-Pyridyldithio)propionyl hydrazide	PDPH
Sulphydryl-carbonyl (aldehyde)/carboxyl reactive	
<i>N</i> -[ $\beta$ -Maleimidopropionic acid]hydrazide•TFA	BMPH
<i>N</i> -e-Maleimidocaproic acid]hydrazide	EMCH
<i>N</i> -[ $\kappa$ -Maleimidoundecanoic acid]hydrazide	KMUH
<i>N</i> -[ <i>p</i> -Maleimidophenyl]isocyanate	PMPI
	TFCS

## EXAMPLE 16

**[1062]** Fetal chromosomal abnormalities are determined by analyzing SNPs wherein the maternal template DNA is homozygous and the template DNA obtained from the plasma is heterozygous. Plasma that is isolated from blood of a pregnant female contains both maternal template DNA and fetal template DNA. Any number of SNP detection methods can be used to analyze the maternal and plasma DNA. Any DNA microarray may be used including but not limited to commercially available and non-commercially available arrays.

**[1063]** A DNA microarray can be designed to contain SNPs located on the chromosome or chromosomes of interest including but not limited to a DNA microarray containing SNPs located on chromosomes 13, 18, and 21, a DNA microarray containing SNPS located on chromosomes 13 and 18, a DNA microarray containing SNPS located on chromosomes 13 and 21, a DNA microarray containing SNPS located on chromosomes 18 and 21, a DNA microarray containing SNPS located on chromosomes 13, 18, 21, 15, 22, X, Y, a DNA microarray containing SNPS located on each of the autosomal chromosomes and each of the sex chromosomes, a DNA microarray containing SNPS located on chromosome 13, a DNA microarray containing SNPS located on chromosome 18, a DNA microarray containing SNPS located on chromosome 21, a DNA microarray containing SNPS located on chromosome 15, a DNA microarray containing SNPS located on chromosome 17, a DNA microarray containing SNPS located on chromosome 22, a DNA microarray containing SNPS located on a single chromosome, and a DNA microarray containing SNPS located on multiple chromosomes.

**[1064]** In this example, SNPs are analyzed by GeneChip HuSNP Arrays from Affymetrix, however any number of DNA arrays, including but not limited to GeneChip arrays, GenFlex Tag arrays, Mapping 10K Array, other Affymetrix arrays, and other DNA arrays can be used.

### Collection of Blood Samples

**[1065]** In accordance with an IRB approved study, blood samples are collected from pregnant women after informed consent is granted. Blood is collected into 9 ml EDTA Vacuette tubes (catalog number NC9897284) and 0.225 ml of 10% neutral buffered solution containing formaldehyde (4% w/v), is added to each tube, and each tube gently is inverted. The tubes are stored at 4°C until ready for processing.

**[1066]** Any number of agents that impede cell lysis or stabilize cell membranes or cross-link cell membranes can be added to the tubes including but not limited to formaldehyde, and derivatives of



formaldehyde, formalin, glutaraldehyde, and derivatives of glutaraldehyde, crosslinkers, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA, DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS or compounds listed in Table XXIII. Any concentration of agent that stabilizes cell membranes, impedes cell lysis or cross-link cell membranes can be added. In a preferred embodiment, the agent that stabilizes cell membranes, impedes cell lysis, or cross-links cell membranes is added at a concentration that does not impede or hinder subsequent reactions.

**[1067]** An agent that stabilizes cell membranes may be added to the maternal blood sample to reduce maternal cell lysis including but not limited to aldehydes, urea formaldehyde, phenol formaldehyde, DMAE (dimethylaminoethanol), cholesterol, cholesterol derivatives, high concentrations of magnesium, vitamin E, and vitamin E derivatives, calcium, calcium gluconate, taurine, niacin, hydroxylamine derivatives, bimoclolmol, sucrose, astaxanthin, glucose, amitriptyline, isomer A hopane tetral phenylacetate, isomer B hopane tetral phenylacetate, citicoline, inositol, vitamin B, vitamin B complex, cholesterol hemisuccinate, sorbitol, calcium, coenzyme Q, ubiquinone, vitamin K, vitamin K complex, menaquinone, zonegran, zinc, ginkgo biloba extract, diphenylhydantoin, perftoran, polyvinylpyrrolidone, phosphatidylserine, tegretol, PABA, disodium cromglycate, nedocromil sodium, phenytoin, zinc citrate, mexitol, dilantin, sodium hyaluronate, or polaxamer 188.

### **Isolation of Plasma and Maternal Cells**

**[1068]** The blood is stored at 4°C until processing. The tubes are spun at 1000 rpm for ten minutes in a centrifuge with braking power set at zero. The tubes are spun a second time at 1000 rpm for ten minutes. The supernatant (the plasma) of each sample is transferred to a new tube and spun at 3000 rpm for ten minutes with the brake set at zero. The supernatant is transferred to a new tube and stored at -80°C. Approximately two milliliters of the “buffy coat,” which contains maternal cells, is placed into a separate tube and stored at -80°C.

### **Isolation of DNA**

**[1069]** DNA is isolated from the plasma sample using the Qiagen Midi Kit for purification of DNA from blood cells, following the manufacturer’s instructions (QIAmp DNA Blood Midi Kit, Catalog

number 51183). DNA is eluted in 100 µl of distilled water. The Qiagen Midi Kit also is used to isolate DNA from the maternal cells contained in the “buffy coat.”

## **Identification of Homozygous Maternal SNPs**

### **HuSNP Assay**

**[1070]** The HuSNP assay is done as described by K. Lindblad-Toh et al. (Nature Biotechnology, Vol. 18, 1001-1005). The GeneChip® HuSNP™ Array is thought to enable whole genome surveys by simultaneously tracking nearly 1,500 SNPs dispersed throughout the genome. In this example, HuSNP array is used as a representative Affymetrix array, and is not meant to limit the use of other arrays including but not limited to GeneChip CYP450, and Affymetrix custom arrays that are designed to meet specific user requirements.

### **PCR Amplification**

**[1071]** Maternal DNA is assayed according to the HuSNP protocol supplied by Affymetrix Inc. For each sample, 24 pools of primer pairs (50-100 loci/pool at 50 nM each) are mixed with 5 ng of maternal DNA., 5 mM MgCl<sub>2</sub>, 0.5 mM dNTPs, 1.25 U Amplitaq Gold (PE Biosystems, Foster City, CA), and the supplied buffer in 12.5 µl per pool. Samples are denatured for 5 min at 95°C followed by 30 cycles of 95°C for 30 s, 52°C + 0.2°C/cycle for 55s, and 72°C for 30s; 5 cycles of 95°C for 30 s, 58°C for 55s, and 72°C for 30s and a final extension of 72°C for 7 min. A 1:1000 dilution of each pool is made by adding 1 µl of the amplification product to 999 µl of ddH<sub>2</sub>O. After, 2.5µl of the 1:1000 dilution is transferred to a new plate and amplified with 0.8µM biotinylated T7 and 0.8 µM biotinylated T3 primers, 4mM MgCl<sub>2</sub>, 0.4 mM dNTPs, 2.5 U Taq and the supplied buffer in 25 µl for 8 min at 95°C, followed by 40 cycles of 95°C for 30 s, 55°C for 90s, and 72°C for 30 s, and a final extension of 72°C for 7 min. Then 1.5µl from each pool is tested for amplification on a 3% agarose gel. For each sample, the remainder of each the 24 pools is mixed and loaded on a Microcon-10 spin column (Amicon Bioseparations, Bedford, MA). Samples are concentrated by spinning the column for 20 min at 13,000 g at room temperature and are eluted by inverting the column and centrifuging for 3 min at 3,000 g. Volumes are adjusted to 60 µl.

**[1072]** A custom array can be designed using only the SNPs that are of interest. For example, a custom array may be designed that contains SNPs that are located on chromosomes 1, 13, 21, 18, 15, X, and Y.

**[1073]** Additionally, any number of SNPs can be amplified including SNPs located on any human chromosome including chromosome 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21,

22, X or Y. Two representative SNPs on chromosome 13 and two representative SNPs on chromosome 21 are chosen. The genomic location and sequence of SNPs may be found at the SNP consortium (<http://snp.cshl.org>). If these SNPs are not present on the array, different SNPs can be chosen.

**[1074]** SNP TSC0466917 (C/G), which is located on chromosome 13, is amplified using the following primers:

Upstream Primer:

5' CCAGCTGGTAGAACTT 3' ([SEQ ID NO: 629](#))

Downstream Primer:

5' CCCAATAGACCTATAG 3' ([SEQ ID NO: 630](#))

**[1075]** SNP TSC1172576 (T/A), which is located on chromosome 13, is amplified using the following primers:

Upstream Primer:

5' TAGCAGAATCTCTCAT 3' ([SEQ ID NO: 631](#))

Downstream Primer:

5' AGAGTATCTCATTTGTT 3' ([SEQ ID NO: 632](#))

**[1076]** SNP TSC0271628 (A/G), which is located on chromosome 21, is amplified using the following primers:

Upstream Primer:

5' AGGAAATTGTGAAGTA 3' ([SEQ ID NO: 633](#))

Downstream Primer:

5' TAACTCACTCACTATC 3' ([SEQ ID NO: 634](#))

**[1077]** SNP TSC0069805 (C/T), which is located on chromosome 21, is amplified using the following primers:

Upstream Primer:

5' CTGCTGAGTCATAGTC 3' ([SEQ ID NO: 635](#))

Downstream Primer:

5' TGTTCCTTTGAATCAAC 3' ([SEQ ID NO: 636](#))

### **Hybridization to GeneChip Probe Arrays, Washing and Staining**

**[1078]** 5-30 µl of the sample (depending on the intensity of the chip lot) is diluted in 3 M tetramethylammonium chloride (TMACl), 2mM control oligonucleotide B1 (supplied by Affymetrix), 5X Denhardt's solution, 100 µg/ml herring sperm DNA, 5 mM EDTA pH 8.0, 10 mM Tris pH7.8, and 0.01%

Tween 20 in a volume of 135 µl and is denatured for 10 min at 95°C. After two minutes on ice, the samples are loaded into HuSNP chips and hybridized for 16 h at 44°C and 40 r.p.m.

[1079] Each chip is washed and stained on the Affymetrix fluidics. Chips are washed for two cycles of two mixes with 6X SSPET (Bio Whitaker, Walkersville, MD) (6X SSPE (sodium chloride, sodium phosphate, sodium EDTA) + 0.01% Triton-X-100) at 25°C, and for six cycles of five mixes with 4X SSPET (4X SSPE + 0.01% Triton X-100) at 35°C. Chips are stained for 30 min at 25°C with 50 µg/ml streptavidin-phycoerythrin and 0.25 mg/ml biotinylated anti-streptavidin antibody in 6X SSPE, 1X Denhardt's solution, and 0.01% Tween 20 in a volume of 500 µl. The chip is filled with 6X SSPET following six washes of four mixes with 6X SSPET at 25°C.

[1080] After the hybridization, washing, and staining procedures, the HuSNP probe arrays are scanned using the HP GeneArray Scanner (HuSNP Mapping Assay Manual Affymetrix P/N 700308).

### **Scanning**

[1081] The HuSNP probe arrays are scanned using the HP GeneArray Scanner according to the HuSNP Mapping Assay Manual (Affymetrix P/N 700308). Other scanners may be used including but not limited to the AlphaArray™ Reader. Genotype calls are made automatically from the collected hybridization signal intensities by the Affymetrix Microarray Suite version 5.0 software. Each allele of a SNP is represented by four or five complementary probes with different locations of the SNP base position within the 20-nucleotide probes. Each of these probes, in turn, is paired with a probe of the same sequence except for a central mismatch at or near a SNP position, intended to correct the fluorescence value for non-specific binding to the probe.

[1082] Each SNP is genotyped. SNPs located on chromosomes 13 and 21, wherein the maternal DNA is homozygous, are analyzed with the DNA isolated from the plasma.

### **Analysis of DNA Isolated from Maternal Plasma**

[1083] After the maternal DNA is analyzed and homozygous SNPs are identified, these SNPs are analyzed with the DNA isolated from the plasma. A low copy number of fetal genomes typically is present in the maternal plasma. To increase the copy number of the loci of interest, which are the SNPs at which the maternal DNA is homozygous, primers are designed to anneal at approximately 130 bases upstream and 130 bases downstream of each loci of interest. This is done to reduce statistical sampling error that can occur when working with a low number of genomes, which can influence the ratio of one allele to another (see Example 11).

### **Design of Multiplex Primers**

**[1084]** The primers are 12 bases in length. However, primers of any length can be used including but not limited to 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36-45, 46-55, 56-65, 66-75, 76-85, 86-95, 96-105, 106-115, 116-125, and greater than 125 bases. Primers are designed to anneal to both the sense strand and the antisense strand.

**[1085]** The maternal homozygous SNPs vary from sample to sample so defined sequences are not provided here. Primers are designed to anneal about 130 bases upstream and downstream of the maternal homozygous SNPs. The primers are designed to terminate at the 3' end in the dinucleotide "AA" to reduce the formation of primer-dimers. However, the primers can be designed to end in any of the four nucleotides and in any combination of the four nucleotides.

### **Multiplex PCR**

**[1086]** Regions upstream and downstream of the maternal homozygous SNPs are amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). This PCR reaction uses primers that anneal approximately 130 bases upstream and downstream of each loci of interest. The primers are mixed together and are used in a single reaction to amplify the template DNA. This reaction is done to increase the number of copies of the loci of interest, which eliminates error generated from a low number of genomes.

**[1087]** For increased specificity, a "hot-start" PCR reaction is used. PCR reactions are performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction is optimized for each locus of interest. In this example, the 20 µl of plasma template DNA is used.

**[1088]** Two microliters of each forward and reverse primer, at concentrations of 5 mM are pooled into a single microcentrifuge tube and mixed. Four microliters of the primer mix is used in a total PCR reaction volume of 50 µl (20µl of template plasma DNA, 1 µl of sterile water, 4 µl of primer mix, and 25 µl of HotStar Taq. Twenty-five cycles of PCR are performed. The following PCR conditions are used:

- (1) 95°C for 15 minutes;
- (2) 95°C for 30 second;
- (3) 4°C for 30 seconds;
- (4) 37°C for 30 seconds;
- (5) Repeat steps 2-4 twenty-four (24) times;
- (6) 72°C for 10 minutes.

[1089] The temperatures and times for denaturing, annealing, and extension, are optimized by trying various settings and using the parameters that yield the best results.

[1090] Other methods of genomic amplification can also be used to increase the copy number of the loci of interest including but not limited to primer extension preamplification (PEP) (Zhang *et al.*, PNAS, 89:5847-51, 1992), degenerate oligonucleotide primed PCR (DOP-PCR) (Telenius, *et al.*, Genomics 13:718-25, 1992), strand displacement amplification using DNA polymerase from bacteriophage 29, which undergoes rolling circle replication (Dean *et al.*, Genomic Research 11:1095-99, 2001), multiple displacement amplification (U.S. Patent 6,124,120), REPLI-g™ Whole Genome Amplification kits, and Tagged PCR.

[1091] It is important to ensure that the region amplified (done to increase the copy number of the fetal loci of interest) contains annealing sequences for the primers that are used with the CodeLink assay. Upon purchase of the CodeLink array, each SNP and the primers used to amplify each SNP can be identified. With this knowledge, the multiplex primers are designed to encompass annealing regions for the primers in the HuSNP Array.

#### **Purification of Fragment of Interest**

[1092] The unused primers, and nucleotides are removed from the reaction by using Qiagen MinElute PCR purification kits (Qiagen, Catalog Number 28004). The reactions are performed following the manufacturer's instructions supplied with the columns. The DNA is eluted in 100 µl of sterile water. 5 µl of each amplified loci is mixed together

#### **CodeLink Assay, Washing, Staining and Scanning**

[1093] The pooled DNA is assayed with the CodeLink Array as described above. Washing, staining, and scanning procedures are as described above.

[1094] Each SNP is genotyped. SNPs located on chromosomes 13 and 21, wherein the maternal DNA is homozygous, and DNA isolated from the plasma is heterozygous are quantitated.

#### **Quantification**

[1095] The intensity of the signal for each allele at a heterozygous is SNP is quantitated. As discussed above, the expected ratio of allele 1 to allele 2 can be used to determine the presence or absence of a chromosomal abnormality. If the maternal genome is homozygous at SNP X (A/A), and the plasma DNA is heterozygous at SNP X (A/G), then the G represents the distinct fetal signal. The ratio of G:A depends on the percentage of fetal DNA present in the maternal blood.

[1096] For example, if the sample contains 50% fetal DNA, then the expected ratio is 0.33 (1 fetal G allele/(2 maternal A alleles + 1 fetal A allele)). This ratio should be constant for all chromosomes that are present in two copies. The ratio that is obtained for SNPs on chromosome 13 should be the same as the ratio that is obtained for chromosome 21.

[1097] However, if the fetal genome contains an additional copy of chromosome 21, then the ratio for this chromosome will deviate from the expected ratio. The expected ratio for a Trisomy condition with 50% fetal DNA in the maternal blood is 0.25. Thus, by analyzing SNPs wherein the maternal genome is homozygous, and the DNA that is isolated from the plasma is heterozygous, fetal chromosomal abnormalities can be detected.

[1098] This example explained the use of CodeLink Arrays, but it not intended to limit the use of arrays. Any DNA array may be used including but not limited to the DNA arrays listed in Table XXIII, or DNA arrays available from any of the companies listed in Table XXIV.

#### **EXAMPLE 18**

[1099] Fetal chromosomal abnormalities are determined by analyzing SNPs wherein the maternal template DNA is homozygous and the template DNA obtained from the plasma is heterozygous. Plasma that is isolated from blood of a pregnant female contains both maternal template DNA and fetal template DNA. Any number of SNP detection methods can be used to analyze the maternal and plasma DNA. In this example, SNPs are analyzed using Illumina's BeadArray™ platform, available from Illumina in San Diego, CA.

#### **Collection of Blood Samples**

[1100] In accordance with an IRB approved study, blood samples are collected from pregnant women after informed consent is granted. Blood is collected into 9 ml EDTA Vacuette tubes (catalog number NC9897284) and 0.225 ml of 10% neutral buffered solution containing formaldehyde (4% w/v), is added to each tube, and each tube gently is inverted. The tubes are stored at 4°C until ready for processing.

[1101] Any number of agents that impede cell lysis or stabilize cell membranes can be added to the tubes including but not limited to formaldehyde, and derivatives of formaldehyde, formalin, glutaraldehyde, and derivatives of glutaraldehyde, crosslinkers, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA,

DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS or compounds listed in Table XXIII. Any concentration of agent that stabilizes cell membranes or impedes cell lysis can be added. In a preferred embodiment, the agent that stabilizes cell membranes or impedes cell lysis is added at a concentration that does not impede or hinder subsequent reactions.

**[1102]** An agent that stabilizes cell membranes may be added to the maternal blood sample to reduce maternal cell lysis including but not limited to aldehydes, urea formaldehyde, phenol formaldehyde, DMAE (dimethylaminoethanol), cholesterol, cholesterol derivatives, high concentrations of magnesium, vitamin E, and vitamin E derivatives, calcium, calcium gluconate, taurine, niacin, hydroxylamine derivatives, bimoclomol, sucrose, astaxanthin, glucose, amitriptyline, isomer A hopane tetral phenylacetate, isomer B hopane tetral phenylacetate, citicoline, inositol, vitamin B, vitamin B complex, cholesterol hemisuccinate, sorbitol, calcium, coenzyme Q, ubiquinone, vitamin K, vitamin K complex, menaquinone, zonegran, zinc, ginkgo biloba extract, diphenylhydantoin, perftoran, polyvinylpyrrolidone, phosphatidylserine, tegretol, PABA, disodium cromglycate, nedocromil sodium, phenytoin, zinc citrate, mexitol, dilantin, sodium hyaluronate, or polaxamer 188.

#### **Isolation of Plasma and Maternal Cells**

**[1103]** The blood is stored at 4°C until processing. The tubes are spun at 1000 rpm for ten minutes in a centrifuge with braking power set at zero. The tubes are spun a second time at 1000 rpm for ten minutes. The supernatant (the plasma) of each sample is transferred to a new tube and spun at 3000 rpm for ten minutes with the brake set at zero. The supernatant is transferred to a new tube and stored at -80°C. Approximately two milliliters of the “buffy coat,” which contains maternal cells, is placed into a separate tube and stored at -80°C.

#### **Isolation of DNA**

**[1104]** DNA is isolated from the plasma sample using the Qiagen Midi Kit for purification of DNA from blood cells, following the manufacturer’s instructions (QIAmp DNA Blood Midi Kit, Catalog number 51183). DNA is eluted in 100 µl of distilled water. The Qiagen Midi Kit also is used to isolate DNA from the maternal cells contained in the “buffy coat.”



## **Identification of Homozygous Maternal SNPs**

**[1105]** Illumina's BeadArray™ technology consists of a fiber optic based array system that allegedly allows miniaturized, very-high throughput genetic analysis. Illumina's 96-bundle Sentrix Array™ allegedly enable parallel processing of nearly 150,000 SNPs.

**[1106]** Fiber bundles are manufactured to contain nearly 50,000 individual, light transmitting fiber strands. Each fiber bundle is converted into an array by first chemically etching a microscopic well at the end of each fiber strand within a bundle, which creates up to 50, 000 discrete microscopic wells per bundle.

**[1107]** In a separate process, sensors are created by affixing a specific type of molecule to the beads, each bead approximately 3 microns in diameter. For SNP analysis, a particular DNA sequence is attached to each bead in a batch. Illumina states that hundreds of thousands of molecules of the same type coat each bead. Batches of coated beads are combined to form a pool specific to the type of array desired. For SNP analysis, the array pool allegedly uses DNA sequences that do not cross hybridize with themselves or with known genomic DNA.

**[1108]** Next, the self-assembled array is created. By dipping bundles into a pre-mixed bead pool, the coated beads self-assemble individually, one bead per well, on the end of each fiber in the bundle to create the array. In Illumina's SNP genotyping array, the bead pool consists of up to 1500 sequences, which self assemble in each bundle of 50,000 fibers to create an array with an approximately thirty-fold redundancy.

**[1109]** The BeadArray bundles are assembled into a matrixed device, which is called the Array of Arrays™ platform, where each fiber bundle of the larger array matches a well of a standardized microtiter plate.

**[1110]** Following array assembly, a decoding process is used to determine the bead type that resides in each fiber core. The DNA molecules are synthesized using the Oligator™ custom DNA synthesis technology.

**[1111]** Illumina's SNP genotyping service using the BeadArray technology and other technologies that have sprung from the BeadArray technology are provided at Illumina's facilities or at facilities that have received a license to the BeadArray technology.

**[1112]** Maternal DNA samples are analyzed using BeadArrays that contain oligonucleotide probes to SNPs. The oligonucleotide probes can be for SNPs located on any chromosome including human chromosome 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, X, and Y. The

BeadArrays are analyzed to identify SNPs, wherein the maternal template DNA is homozygous. The identified homozygous SNPs are then analyzed using the DNA isolated from the maternal plasma.

### **Analysis of DNA Isolated from Maternal Plasma**

[1113] After the maternal DNA is analyzed and homozygous SNPs are identified, these SNPs are analyzed with the DNA isolated from the plasma. A low copy number of fetal genomes typically exists in the maternal plasma. To increase the copy number of the loci of interest, which are the SNPs at which the maternal DNA is homozygous, primers are designed to anneal at approximately 130 bases upstream and 130 bases downstream of each loci of interest. This is done to reduce statistical sampling error that can occur when working with a low number of genomes, which can influence the ratio of one allele to another (see Example 11).

### **Design of Multiplex Primers**

[1114] The primers are 12 bases in length. However, primers of any length can be used including but not limited to 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36-45, 46-55, 56-65, 66-75, 76-85, 86-95, 96-105, 106-115, 116-125, and greater than 125 bases. Primers are designed to anneal to both the sense strand and the antisense strand.

[1115] The maternal homozygous SNPs vary from sample to sample so defined sequences are not provided here. Primers are designed to anneal about 130 bases upstream and downstream of the maternal homozygous SNPs. The primers are designed to terminate at the 3' end in the dinucleotide "AA" to reduce the formation of primer-dimers. However, the primers can be designed to end in any of the four nucleotides and in any combination of the four nucleotides.

### **Multiplex PCR**

[1116] Regions upstream and downstream of the maternal homozygous SNPs are amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). This PCR reaction uses primers that anneal approximately 130 bases upstream and downstream of each loci of interest. The primers are mixed together and are used in a single reaction to amplify the template DNA. This reaction is done to increase the number of copies of the loci of interest, which eliminates error generated from a low number of genomes.

[1117] For increased specificity, a "hot-start" PCR reaction is used. PCR reactions are performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of

template DNA and primer per reaction is optimized for each locus of interest. In this example, the 20  $\mu$ l of plasma template DNA is used.

**[1118]** Two microliters of each forward and reverse primer, at concentrations of 5 mM are pooled into a single microcentrifuge tube and mixed. Four microliters of the primer mix is used in a total PCR reaction volume of 50  $\mu$ l (20 $\mu$ l of template plasma DNA, 1  $\mu$ l of sterile water, 4  $\mu$ l of primer mix, and 25  $\mu$ l of HotStar Taq. Twenty-five cycles of PCR are performed. The following PCR conditions are used:

- (1) 95°C for 15 minutes;
- (2) 95°C for 30 second;
- (3) 4°C for 30 seconds;
- (4) 37°C for 30 seconds;
- (5) Repeat steps 2-4 twenty-four (24) times;
- (6) 72°C for 10 minutes.

**[1119]** The temperatures and times for denaturing, annealing, and extension, are optimized by trying various settings and using the parameters that yield the best results.

**[1120]** Other methods of genomic amplification can also be used to increase the copy number of the loci of interest including but not limited to primer extension preamplification (PEP) (Zhang *et al.*, PNAS, 89:5847-51, 1992), degenerate oligonucleotide primed PCR (DOP-PCR) (Telenius, *et al.*, Genomics 13:718-25, 1992), strand displacement amplification using DNA polymerase from bacteriophage 29, which undergoes rolling circle replication (Dean *et al.*, Genomic Research 11:1095-99, 2001), multiple displacement amplification (U.S. Patent 6,124,120), REPLI-g™ Whole Genome Amplification kits, and Tagged PCR.

**[1121]** It is important to ensure that the region amplified contains annealing sequences for the oligonucleotide probes in the BeadArray. Upon purchase of the BeadArray service, each SNP and the primers used to analyze each SNP are identified. With this knowledge, the multiplex primers are designed to encompass annealing regions for the primers in the BeadArray.

### **Purification of Fragment of Interest**

**[1122]** The unused primers, and nucleotides are removed from the reaction by using Qiagen MinElute PCR purification kits (Qiagen, Catalog Number 28004). The reactions are performed following the manufacturer's instructions supplied with the columns. The DNA is eluted in 100  $\mu$ l of sterile water. 5  $\mu$ l of each amplified loci is mixed together

### **BeadArray Technology**

[1123] The pooled DNA is assayed with the BeadArray as described above. Each SNP is genotyped. SNPs located on chromosomes 13 and 21, wherein the maternal DNA is homozygous, and DNA isolated from the plasma is heterozygous are quantitated.

### **Quantification**

[1124] The intensity of the signal for each allele at a heterozygous is SNP is quantitated. As discussed above, the expected ratio of allele 1 to allele 2 can be used to determine the presence or absence of a chromosomal abnormality. If the maternal genome is homozygous at SNP X (A/A), and the plasma DNA is heterozygous at SNP X (A/G), then the G represents the distinct fetal signal. The ratio of G:A depends on the percentage of fetal DNA present in the maternal blood.

[1125] For example, if the sample contains 50% fetal DNA, then the expected ratio is 0.33 (1 fetal G allele/(2 maternal A alleles + 1 fetal A allele)). This ratio should be constant for all chromosomes that are present in two copies. The ratio that is obtained for SNPs on chromosome 13 should be the same as the ratio that is obtained for chromosome 21.

[1126] However, if the fetal genome contains an additional copy of chromosome 21, then the ratio for this chromosome will deviate from the expected ratio. The expected ratio for a Trisomy condition with 50% fetal DNA in the maternal blood is 0.25. Thus, by analyzing SNPs wherein the maternal genome is homozygous, and the DNA that is isolated from the plasma is heterozygous, fetal chromosomal abnormalities can be detected.

[1127] This example explained the use of Illumina's BeadArray Technology, but it not intended to limit the use of arrays. Any DNA array may be used including but not limited to the DNA arrays listed in Table XXIII, or DNA arrays available from any of the companies listed in Table XXIV.

### **EXAMPLE 19**

[1128] Fetal chromosomal abnormalities are determined by analyzing SNPs wherein the maternal template DNA is homozygous and the template DNA obtained from the plasma is heterozygous. Plasma that is isolated from blood of a pregnant female contains both maternal template DNA and fetal template DNA. Any number of SNP detection methods can be used to analyze the maternal and plasma DNA. In this example, SNPs are analyzed using Sequenom's MassArray™ System, which uses Sequenom's homogenous MassCleave™ (hMC) method.

## Collection of Blood Samples

[1129] In accordance with an IRB approved study, blood samples are collected from pregnant women after informed consent is granted. Blood is collected into 9 ml EDTA Vacuette tubes (catalog number NC9897284) and 0.225 ml of 10% neutral buffered solution containing formaldehyde (4% w/v), is added to each tube, and each tube gently is inverted. The tubes are stored at 4°C until ready for processing.

[1130] Any number of agents that impede cell lysis or stabilize cell membranes can be added to the tubes including but not limited to formaldehyde, and derivatives of formaldehyde, formalin, glutaraldehyde, and derivatives of glutaraldehyde, crosslinkers, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA, DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS or compounds listed in Table XXIII. Any concentration of agent that stabilizes cell membranes or impedes cell lysis can be added. In a preferred embodiment, the agent that stabilizes cell membranes or impedes cell lysis is added at a concentration that does not impede or hinder subsequent reactions.

[1131] An agent that stabilizes cell membranes may be added to the maternal blood sample to reduce maternal cell lysis including but not limited to aldehydes, urea formaldehyde, phenol formaldehyde, DMAE (dimethylaminoethanol), cholesterol, cholesterol derivatives, high concentrations of magnesium, vitamin E, and vitamin E derivatives, calcium, calcium gluconate, taurine, niacin, hydroxylamine derivatives, bimoclomol, sucrose, astaxanthin, glucose, amitriptyline, isomer A hopane tetral phenylacetate, isomer B hopane tetral phenylacetate, citicoline, inositol, vitamin B, vitamin B complex, cholesterol hemisuccinate, sorbitol, calcium, coenzyme Q, ubiquinone, vitamin K, vitamin K complex, menaquinone, zonegran, zinc, ginkgo biloba extract, diphenylhydantoin, perftoran, polyvinylpyrrolidone, phosphatidylserine, tegretol, PABA, disodium cromglycate, nedocromil sodium, phenytoin, zinc citrate, mexitil, dilantin, sodium hyaluronate, or polaxamer 188.

## Isolation of Plasma and Maternal Cells

[1132] The blood is stored at 4°C until processing. The tubes are spun at 1000 rpm for ten minutes in a centrifuge with braking power set at zero. The tubes are spun a second time at 1000 rpm for ten minutes. The supernatant (the plasma) of each sample is transferred to a new tube and spun at 3000 rpm

for ten minutes with the brake set at zero. The supernatant is transferred to a new tube and stored at -80°C. Approximately two milliliters of the “buffy coat,” which contains maternal cells, is placed into a separate tube and stored at -80°C.

### **Isolation of DNA**

[1133] DNA is isolated from the plasma sample using the Qiagen Midi Kit for purification of DNA from blood cells, following the manufacturer’s instructions (QIAmp DNA Blood Midi Kit, Catalog number 51183). DNA is eluted in 100 µl of distilled water. The Qiagen Midi Kit also is used to isolate DNA from the maternal cells contained in the “buffy coat.”

### **Identification of Homozygous Maternal SNPs**

#### **Targeted SNP Discovery: hMC Method**

[1134] Sequenom’s hMC method uses nucleotide base-specific cleavage for genotyping. The cleaved fragments are measured using MALDI-TOF to generate a characteristic peak signal, based on the mass of each fragment, for any particular sequence.

### **Primer Design**

[1135] Four primers are needed for the two PCR reactions (one forward reaction and one reverse reaction). The recommended size range for PCR ~~amplicons~~<sup>amplicons</sup> is 300-700 base pairs. The primers contain a T-7 promoter tagged forward or reverse primer to obtain an appropriate product for in vitro transcription. An 8 base insert is included to prevent abortive cycling. The primer that lacks the T-7 promoter contains a 10-mer tag in order to balance the primers.

[1136] The primers for one SNP are provided below. SNP TSC1172576 (T/A), which is located on chromosome 13, is amplified using the following primers for the forward reaction:

Forward reaction:

Upstream Primer:

5' *CAGTAATACGACTCACTATAGGGGTCAGGATTAGCAGAATCTCTCAT* 3' (SEQ ID NO: 637)

Downstream Primer:

5' *GCATTCTATGAGAGTATCTCATTTGTT* 3' (SEQ ID NO: 638)

Reverse Reaction:

Upstream Primer

5' *CAGTAATACGACTCACTATAGGGGTCAGGAAGAGTATCTCATTTGTT* 3' (SEQ ID NO: 639)

#### Downstream Primer

5' GCATTCTATGTAGCAGAATCTCTCAT 3' (SEQ ID NO: 640)

[1137] The sequence of the T-7 promoter are in italics, the 8 base insert is underlined, the 10- base balancing sequence is double underlined, and the gene specific sequences are unmodified.

#### PCR Amplification

[1138] Five nanograms of DNA is amplified in a 5 µl volume using a 384-microtiter format. The following PCR conditions are used:

- 1) 94°C for 15 minutes;
- 2) 94°C for 20 seconds;
- 3) 62°C for 30 seconds;
- 4) 72°C for 1 minute;
- 5) Repeat steps 2-4 44 times; and
- 6) 72°C for 3 minutes.

#### Dephosphorylation

[1139] Shrimp Alkaline Phosphatase (SAP) (2 µl) is added to each 5 µl PCR reaction to dephosphorylate unincorporated dNTPs from the PCR reaction. The plates are incubated at 37°C for 20 minutes. Then, the plates are incubated at 85°C for 5 minutes.

#### In Vitro Transcription

[1140] For each transcription reaction, 2 µl of transcription cocktail and 2 µl of PCR/SAP sample are needed. Add 2 µl of transcription cocktail and 2 µl of PCR/SAP sample to a new microtiter plate. The plates are incubated at 37°C for two hours. For detailed information regarding these protocols see the “Processing homogeneous MassCLEASE Reactions” chapter in the MassARRAY Liquid Handler SNP Discovery User’s Guide for instructions, which is fully incorporated herein by reference.

#### RNase A Cleavage

[1141] RNase A cocktail (2.5 µl) is added to each reaction (T cleavage and C cleavage). The plates are incubated at 37°C for one hour.

[1142] Depending on the nucleotide at the SNP site, various fragments of different weights are generated. For example, the DNA sequence surrounding SNP TSC1172576, which is located on chromosome 13, is as follows:

5' CCGCATA T/A CTCAGCACA 3' (SEQ ID NO: 641)

3'GGCGTAT A/T GAGTCGTGT 5' (SEQ ID NO: 642)

[1143] After PCR, in vitro transcription, and base-specific cleavage, the following fragments for each allele are generated:

T allele	A allele	
Products of forward	1) TATCTCA	1) CTCAGC
Transcription	2) ATC	2) AAC
Products of reverse	1) AGTTA	1) AGATA
transcription	2) GTTATGC	2) ATATGCG

[1144] For the ATC and AAC fragments, the weight difference between T and A is used to determine the genotype at SNP TSC1172576. Likewise, the weight difference between T and A in fragments AGTTA and AGATA is used to determine the genotype at SNP TSC1172576.

### Sample Conditioning

[1145] Double distilled water (20 µl) is added to each sample within the 384-well plate. Clean Resin (6 mg) is added to each well. The plate is rotated for 10 minutes, followed by a centrifugation at 3200 X g. It is recommended that water always be added before the Clean Resin.

### Sample Transfer

[1146] The hMC reaction product (10-15µl) is dispensed onto a 384 element SpectroCHIP®. For further information, see the “Dispensing MassCLEAVE Reaction Products onto SpectroCHIPS” chapter in the MassARRAY Nanodispenser SNP Discovery User’s Guide for instructions.

### Sample Analysis

[1147] Spectra from the four cleavage reactions is acquired using the MassARRAY™ system. For further instructions, see the “Acquiring Spectra” chapter in the MassARRAY Discovery RT Software User’s guide for instructions on acquiring spectra from SpectroCHIPS®.

### SNP Analysis

[1148] The results are analyzed using the SNP Discovery Analysis software. For further instructions, see the “Analyzing SNPs” chapter in the MassARRAY Discovery RT Software User’s Guide for instructions on using the SNP Discovery Analysis software. Components that are useful for the MassARRAY procedure include MassARRAY™ Analyzer (part number 004500), MassARRAY™



Discovery RT Software version 1.2 (part number 11434), MassARRAY™ SNP Discovery Starter Kit (part number 10027), and Liquid Handler SNP Discovery Methods and Macros (part number 11433).

[1149] The SpectroCHIP array is used to genotype the maternal DNA following the manufacturer's recommended protocols and procedures, which are made available after purchase of the SpectroCHIP array. SNPs at which the maternal DNA is homozygous are used to analyze the DNA isolated from the maternal plasma.

#### **Analysis of DNA Isolated from Maternal Plasma**

[1150] After the maternal DNA is analyzed and homozygous SNPs are identified, these SNPs are analyzed with the DNA isolated from the plasma. A low copy number of fetal genomes typically exists in the maternal plasma. To increase the copy number of the loci of interest, which are the SNPs at which the maternal DNA is homozygous, primers are designed to anneal at approximately 130 bases upstream and 130 bases downstream of each loci of interest. This is done to reduce statistical sampling error that can occur when working with a low number of genomes, which can influence the ratio of one allele to another (see Example 11).

#### **Design of Multiplex Primers**

[1151] The primers are 12 bases in length. However, primers of any length can be used including but not limited to 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36-45, 46-55, 56-65, 66-75, 76-85, 86-95, 96-105, 106-115, 116-125, and greater than 125 bases. Primers are designed to anneal to both the sense strand and the antisense strand.

[1152] The maternal homozygous SNPs vary from sample to sample so defined sequences are not provided here. Primers are designed to anneal about 130 bases upstream and downstream of the maternal homozygous SNPs. The primers are designed to terminate at the 3' end in the dinucleotide "AA" to reduce the formation of primer-dimers. However, the primers can be designed to end in any of the four nucleotides and in any combination of the four nucleotides.

#### **Multiplex PCR**

[1153] Regions upstream and downstream of the maternal homozygous SNPs are amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). This PCR reaction uses primers that anneal approximately 130 bases upstream and downstream of each loci of interest. The primers are mixed together and are used

in a single reaction to amplify the template DNA. This reaction is done to increase the number of copies of the loci of interest, which eliminates error generated from a low number of genomes.

[1154] For increased specificity, a "hot-start" PCR reaction is used. PCR reactions are performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction is optimized for each locus of interest. In this example, the 20  $\mu$ l of plasma template DNA is used.

[1155] Two microliters of each forward and reverse primer, at concentrations of 5 mM are pooled into a single microcentrifuge tube and mixed. Four microliters of the primer mix is used in a total PCR reaction volume of 50  $\mu$ l (20 $\mu$ l of template plasma DNA, 1  $\mu$ l of sterile water, 4  $\mu$ l of primer mix, and 25  $\mu$ l of HotStar Taq. Twenty-five cycles of PCR are performed. The following PCR conditions are used:

- (1) 95°C for 15 minutes;
- (2) 95°C for 30 second;
- (3) 4°C for 30 seconds;
- (4) 37°C for 30 seconds;
- (5) Repeat steps 2-4 twenty-four (24) times;
- (6) 72°C for 10 minutes.

[1156] The temperatures and times for denaturing, annealing, and extension, are optimized by trying various settings and using the parameters that yield the best results.

[1157] Other methods of genomic amplification can also be used to increase the copy number of the loci of interest including but not limited to primer extension preamplification (PEP) (Zhang *et al.*, PNAS, 89:5847-51, 1992), degenerate oligonucleotide primed PCR (DOP-PCR) (Telenius, *et al.*, Genomics 13:718-25, 1992), strand displacement amplification using DNA polymerase from bacteriophage 29, which undergoes rolling circle replication (Dean *et al.*, Genomic Research 11:1095-99, 2001), multiple displacement amplification (U.S. Patent 6,124,120), REPLI-g™ Whole Genome Amplification kits, and Tagged PCR.

[1158] It is important to ensure that the region amplified contains annealing sequences for the PCR primers in the targeted SNP discovery, hMC method

### **Purification of Fragment of Interest**

[1159] The unused primers, and nucleotides are removed from the reaction by using Qiagen MinElute PCR purification kits (Qiagen, Catalog Number 28004). The reactions are performed following

the manufacturer's instructions supplied with the columns. The DNA is eluted in 100 µl of sterile water. 5 µl of each amplified loci is mixed together

#### **Targeted SNP Discovery: hMC Method**

[1160] The pooled DNA is assayed with the hMC method as described above. Each SNP is genotyped. SNPs located on chromosomes 13 and 21, wherein the maternal DNA is homozygous, and DNA isolated from the plasma is heterozygous are quantitated. However, SNPs located on other chromosomes can also be quantitated if so desired.

#### **Quantification**

[1161] The intensity of each peak, wherein each peak corresponds to a DNA fragment with a specific molecular weight, is quantitated. As discussed above, the expected ratio of allele 1 to allele 2 is used to determine the presence or absence of a chromosomal abnormality. If the maternal genome is homozygous at SNP X (A/A), and the plasma DNA is heterozygous at SNP X (A/G), then the G represents the distinct fetal signal.

[1162] There will be some fragments that differ in molecular due to the presence of the G nucleotide at SNP X in the fetal genome. The intensity of the peak with the A nucleotide is quantitated and the intensity of the peak that corresponds to fragments with the G nucleotide is quantitated. The ratio of G:A depends on the percentage of fetal DNA present in the maternal blood.

[1163] For example, if the sample contains 50% fetal DNA, then the expected ratio is 0.33 (1 fetal G allele/(2 maternal A alleles + 1 fetal A allele)). This ratio should be constant for all chromosomes that are present in two copies. The ratio that is obtained for SNPs on chromosome 13 should be the same as the ratio that is obtained for chromosome 21.

[1164] However, if the fetal genome contains an additional copy of chromosome 21, then the ratio for this chromosome will deviate from the expected ratio. The expected ratio for a Trisomy condition with 50% fetal DNA in the maternal blood is 0.25. Thus, by analyzing SNPs wherein the maternal genome is homozygous, and the DNA that is isolated from the plasma is heterozygous, fetal chromosomal abnormalities can be detected.

[1165] This example explained the use of Sequenom's hMC method, but it not intended to limit the use of other mass spectrometry techniques.

## EXAMPLE 20

**[1166]** Fetal chromosomal abnormalities are determined by analyzing SNPs wherein the maternal template DNA is homozygous and the template DNA obtained from the plasma is heterozygous. Plasma that is isolated from blood of a pregnant female contains both maternal template DNA and fetal template DNA. Any number of SNP detection methods can be used to analyze the maternal and plasma DNA. In this example, SNPs are analyzed using Sequenom's MassArray™ Homogenous MassEXTEND™ (hME) Assay.

### Collection of Blood Samples

**[1167]** In accordance with an IRB approved study, blood samples are collected from pregnant women after informed consent is granted. Blood is collected into 9 ml EDTA Vacuette tubes (catalog number NC9897284) and 0.225 ml of 10% neutral buffered solution containing formaldehyde (4% w/v), is added to each tube, and each tube gently is inverted. The tubes are stored at 4°C until ready for processing.

**[1168]** Any number of agents that impede cell lysis or stabilize cell membranes can be added to the tubes including but not limited to formaldehyde, and derivatives of formaldehyde, formalin, glutaraldehyde, and derivatives of glutaraldehyde, crosslinkers, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA, DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS or compounds listed in Table XXIII. Any concentration of agent that stabilizes cell membranes or impedes cell lysis can be added. In a preferred embodiment, the agent that stabilizes cell membranes or impedes cell lysis is added at a concentration that does not impede or hinder subsequent reactions.

**[1169]** An agent that stabilizes cell membranes may be added to the maternal blood sample to reduce maternal cell lysis including but not limited to aldehydes, urea formaldehyde, phenol formaldehyde, DMAE (dimethylaminoethanol), cholesterol, cholesterol derivatives, high concentrations of magnesium, vitamin E, and vitamin E derivatives, calcium, calcium gluconate, taurine, niacin, hydroxylamine derivatives, bimoclomol, sucrose, astaxanthin, glucose, amitriptyline, isomer A hopane tetral phenylacetate, isomer B hopane tetral phenylacetate, citicoline, inositol, vitamin B, vitamin B complex, cholesterol hemisuccinate, sorbitol, calcium, coenzyme Q, ubiquinone, vitamin K, vitamin K complex,

menaquinone, zonegran, zinc, ginkgo biloba extract, diphenylhydantoin, perftoran, polyvinylpyrrolidone, phosphatidylserine, tegretol, PABA, disodium cromglycate, nedocromil sodium, phenytoin, zinc citrate, mexilit, dilantin, sodium hyaluronate, or polaxamer 188.

### **Isolation of Plasma and Maternal Cells**

[1170] The blood is stored at 4°C until processing. The tubes are spun at 1000 rpm for ten minutes in a centrifuge with braking power set at zero. The tubes are spun a second time at 1000 rpm for ten minutes. The supernatant (the plasma) of each sample is transferred to a new tube and spun at 3000 rpm for ten minutes with the brake set at zero. The supernatant is transferred to a new tube and stored at -80°C. Approximately two milliliters of the “buffy coat,” which contains maternal cells, is placed into a separate tube and stored at -80°C.

### **Isolation of DNA**

[1171] DNA is isolated from the plasma sample using the Qiagen Midi Kit for purification of DNA from blood cells, following the manufacturer’s instructions (QIAmp DNA Blood Midi Kit, Catalog number 51183). DNA is eluted in 100 µl of distilled water. The Qiagen Midi Kit also is used to isolate DNA from the maternal cells contained in the “buffy coat.”

### **Identification of Homozygous Maternal SNPs**

#### **MassARRAY Homogenous MassEXTEND™ (hME) Assay**

[1172] The Homogenous MassEXTEND™ (hME) Assay uses a beadless, label-free primer extension chemistry for genotyping. Each of the primer products has a unique molecular weight that allows the associated genotype to be precisely identified using mass spectrometry.

### **Template Amplification**

[1173] The isolated maternal DNA is amplified (2.5 ng) in a 5 µl volume using a 384-microtiter plate format. Any number of SNPs can be amplified, either in a single reaction or in multiple reactions. Representative primers that are used to amplify SNP TSC0271628 (A/G), which is located on chromosome 21, are provided below:

Upstream Primer:

5' AGGAAATTGTGAAGTA 3' (SEQ ID NO: 643)

Downstream Primer:

5' TAACTCACTCACTATC 3' (SEQ ID NO: 644)

[1174] The primers can be longer or shorter in nucleotide sequence. PCR conditions recommended by the makers of MassARRAY Homogenous MassEXTEND Assay are followed. Representative PCR conditions are provided below:

- (1) 95°C for 15 minutes and 15 seconds;
- (2) 95°C for 30 seconds;
- (4) 57°C for 30 seconds;
- (5) 72°C for 30 seconds;
- (6) Repeat steps 2-5 thirty two (32) times;
- (7) 72°C for 5 minutes.

### **Dephosphorylation**

[1175] Arctic shrimp alkaline phosphatase is added to the samples, which are then incubated at 37°C for 20 minutes. This step is done to dephosphorylate any remaining nucleotides, which prevents their future incorporation and interference with MassARRAY Homogenous MassEXTEND Assay. Samples are then incubated at 85°C to inactivate the heat-labile SAP.

### **hME Reaction**

[1176] A MassEXTEND primer is designed to anneal close to the polymorphic site, and is designed to identify both alleles of the polymorphic site. For SNP TSC0271628, a representative MassEXTEND primer is:

5' CTTTTTATGCCTTTCCACTCATCCA 3' (SEQ ID NO: 645)

[1177] The length of the MassEXTEND primer is designed according to the instructions provided by the makers of the MassARRAY Homogenous MassEXTEND Assay.

[1178] The ~~MAssEXTEND~~MassEXTEND primer, DNA polymerase, and a cocktail mixture of deoxynucleotides (dNTPs) and dideoxynucleotides (ddNTPs) are added to the initial primer extension reaction. Allele-specific primer products are generated that are generally one to four bases longer than the original MassEXTEND primer.

[1179] A ~~MASSEXTEND~~MassEXTEND primer is hybridized closely adjacent to the polymorphic site following the conditions recommended by the makers of the MassARRAY Homogenous MassEXTEND Assay. Nucleotide mixtures are selected to maximize mass differences for all possible MassEXTEND products. Appropriate dNTPs are incorporated until a single ddNTP is incorporated, and the reaction is terminated. The manufacturer's protocols are followed for all steps of the hME assay.

[1180] Representative reaction products for SNP TSC0271628 are provided below:

#### **A allele Before Primer Extension**

MassEXTEND primer: CT TTTT ATGCCT T TCCACTCATCCA (SEQ ID NO: 646)

Sample DNA: GAAAAATACGGAAAGGTGAGTAGGTTTCC (SEQ ID NO: 647)

[1181] The SNP site is identified in bold. After incubation with DNA polymerase, ddATP, dCTP, dGTP, and dTTP, the following product is generated:

A allele After Primer Extension

MassEXTEND primer: CT TTTT ATGCCT T TCCACTCATCCAA\* (SEQ ID NO: 648)

Sample DNA: GAAAAATACGGAAAGGTGAGTAGGTTTCC (SEQ ID NO: 649)

[1182] ddATP is incorporated into the primer. Either labeled or unlabeled ddNTPs can be used. The asterisk indicates ddATP that is unlabeled. After the incorporation reaction, a 24-mer primer is generated.

G allele Before Primer Extension

MassEXTEND primer: CT TTTT ATGCCT T TCCACTCATCCA (SEQ ID NO: 650)

Sample DNA: GAAAAATACGGAAAGGTGAGTAGGTCTCC (SEQ ID NO: 651)

[1183] The SNP site is identified in bold. After incubation with DNA polymerase, ddATP, dCTP, dGTP, and dTTP, the following product is generated:

G allele After Primer Extension

MassEXTEND primer: CT TTTT ATGCCT T TCCACTCATCCAGA\* (SEQ ID NO: 652)

Sample DNA: GAAAAATACGGAAAGGTGAGTAGGTCTCC (SEQ ID NO: 653)

[1184] After the incorporation reaction, a 25-mer primer is generated. The difference in molecular weight between the reaction product for the A allele (24-mer) and the reaction product for the G allele (25-mer) is used to genotype the locus of interest.

#### **Sample Conditioning**

[1185] SpectroCLEAN™ resin is added to the reaction to remove extraneous salts that interfere with MALDI-TOF analysis.

#### **Sample Transfer**

[1186] 15 nl of sample is transferred from the 384-microtiter plate and spotted onto the pad of the 384 SpectroCHIP™ microarray.

#### **Sample Analysis**

[1187] The SpectroCHIP™ is placed into the MALDI-TOF, which measures the mass of the extension products. Once determined, the genotype is called in real-time with SpectroTYPER™ RT

software. SNPs at which the maternal DNA are homozygous are identified, and analyzed with the DNA that is isolated from the plasma.

### **Analysis of DNA Isolated from Maternal Plasma**

**[1188]** After the maternal DNA is analyzed and homozygous SNPs are identified, these SNPs are analyzed with the DNA isolated from the plasma. A low copy number of fetal genomes typically exists in the maternal plasma. To increase the copy number of the loci of interest, which are the SNPs at which the maternal DNA is homozygous, primers are designed to anneal at approximately 130 bases upstream and 130 bases downstream of each loci of interest. This is done to reduce statistical sampling error that can occur when working with a low number of genomes, which can influence the ratio of one allele to another (see Example 11).

### **Design of Multiplex Primers**

**[1189]** The primers are 12 bases in length. However, primers of any length can be used including but not limited to 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36-45, 46-55, 56-65, 66-75, 76-85, 86-95, 96-105, 106-115, 116-125, and greater than 125 bases. Primers are designed to anneal to both the sense strand and the antisense strand.

**[1190]** The maternal homozygous SNPs vary from sample to sample so defined sequences are not provided here. Primers are designed to anneal about 130 bases upstream and downstream of the maternal homozygous SNPs. The primers are designed to terminate at the 3' end in the dinucleotide "AA" to reduce the formation of primer-dimers. However, the primers can be designed to end in any of the four nucleotides and in any combination of the four nucleotides.

### **Multiplex PCR**

**[1191]** Regions upstream and downstream of the maternal homozygous SNPs are amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). This PCR reaction uses primers that anneal approximately 130 bases upstream and downstream of each loci of interest. The primers are mixed together and are used in a single reaction to amplify the template DNA. This reaction is done to increase the number of copies of the loci of interest, which eliminates error generated from a low number of genomes.

**[1192]** For increased specificity, a "hot-start" PCR reaction is used. PCR reactions are performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of



template DNA and primer per reaction is optimized for each locus of interest. In this example, the 20  $\mu$ l of plasma template DNA is used.

**[1193]** Two microliters of each forward and reverse primer, at concentrations of 5 mM are pooled into a single microcentrifuge tube and mixed. Four microliters of the primer mix is used in a total PCR reaction volume of 50  $\mu$ l (20 $\mu$ l of template plasma DNA, 1  $\mu$ l of sterile water, 4  $\mu$ l of primer mix, and 25  $\mu$ l of HotStar Taq. Twenty-five cycles of PCR are performed. The following PCR conditions are used:

- (1) 95°C for 15 minutes;
- (2) 95°C for 30 second;
- (3) 4°C for 30 seconds;
- (4) 37°C for 30 seconds;
- (5) Repeat steps 2-4 twenty-four (24) times;
- (6) 72°C for 10 minutes.

**[1194]** The temperatures and times for denaturing, annealing, and extension, are optimized by trying various settings and using the parameters that yield the best results.

**[1195]** Other methods of genomic amplification can also be used to increase the copy number of the loci of interest including but not limited to primer extension preamplification (PEP) (Zhang *et al.*, PNAS, 89:5847-51, 1992), degenerate oligonucleotide primed PCR (DOP-PCR) (Telenius, *et al.*, Genomics 13:718-25, 1992), strand displacement amplification using DNA polymerase from bacteriophage 29, which undergoes rolling circle replication (Dean *et al.*, Genomic Research 11:1095-99, 2001), multiple displacement amplification (U.S. Patent 6,124,120), REPLI-g™ Whole Genome Amplification kits, and Tagged PCR.

#### **Purification of Fragment of Interest**

**[1196]** The unused primers, and nucleotides are removed from the reaction by using Qiagen MinElute PCR purification kits (Qiagen, Catalog Number 28004). The reactions are performed following the manufacturer's instructions supplied with the columns. The DNA is eluted in 100  $\mu$ l of sterile water. 5  $\mu$ l of each amplified loci is mixed together

#### **MassARRAY Homogenous MassEXTEND Assay**

**[1197]** The pooled DNA is assayed with the hME assay as described above. Each SNP is genotyped. SNPs located on chromosomes 13 and 21, wherein the maternal DNA is homozygous, and DNA isolated from the plasma is heterozygous are quantitated.

## Quantification

[1198] The intensity of each peak, wherein each peak corresponds to a DNA fragment with a specific molecular weight, is quantitated. As discussed above, the expected ratio of allele 1 to allele 2 is used to determine the presence or absence of a chromosomal abnormality. If the maternal genome is homozygous at SNP X (A/A), and the plasma DNA is heterozygous at SNP X (A/G), then the G represents the distinct fetal signal.

[1199] There will be some fragments that differ in molecular due to the presence of the G nucleotide at SNP X in the fetal genome. The intensity of the peak with the A nucleotide is quantitated and the intensity of the peak that corresponds to fragments with the G nucleotide is quantitated. The ratio of G:A depends on the percentage of fetal DNA present in the maternal blood.

[1200] For example, if the sample contains 50% fetal DNA, then the expected ratio is 0.33 (1 fetal G allele/(2 maternal A alleles + 1 fetal A allele)). This ratio should be constant for all chromosomes that are present in two copies. The ratio that is obtained for SNPs on chromosome 13 should be the same as the ratio that is obtained for chromosome 21.

[1201] However, if the fetal genome contains an additional copy of chromosome 21, then the ratio for this chromosome will deviate from the expected ratio. The expected ratio for a Trisomy condition with 50% fetal DNA in the maternal blood is 0.25. Thus, by analyzing SNPs wherein the maternal genome is homozygous, and the DNA that is isolated from the plasma is heterozygous, fetal chromosomal abnormalities can be detected.

[1202] This example explained the use Sequenom's MassARRAY Homogenous MassEXTEND (hME) assay, but it not intended to limit the use of techniques that differentiate molecules based on molecular weight.

## EXAMPLE 21

[1203] Fetal chromosomal abnormalities are determined by analyzing SNPs wherein the maternal template DNA is homozygous and the template DNA obtained from the plasma is heterozygous. Plasma that is isolated from blood of a pregnant female contains both maternal template DNA and fetal template DNA. Any number of SNP detection methods can be used to analyze the maternal and plasma DNA. In this example, SNPs are analyzed using Orchid's SNP-IT™ Assay. However, other SNP detection methods based on primer extension may also be used.

## Collection of Blood Samples

[1204] In accordance with an IRB approved study, blood samples are collected from pregnant women after informed consent is granted. Blood is collected into 9 ml EDTA Vacuette tubes (catalog number NC9897284) and 0.225 ml of 10% neutral buffered solution containing formaldehyde (4% w/v), is added to each tube, and each tube gently is inverted. The tubes are stored at 4°C until ready for processing.

[1205] Any number of agents that impede cell lysis or stabilize cell membranes can be added to the tubes including but not limited to formaldehyde, and derivatives of formaldehyde, formalin, glutaraldehyde, and derivatives of glutaraldehyde, crosslinkers, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA, DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS or compounds listed in Table XXIII. Any concentration of agent that stabilizes cell membranes or impedes cell lysis can be added. In a preferred embodiment, the agent that stabilizes cell membranes or impedes cell lysis is added at a concentration that does not impede or hinder subsequent reactions.

[1206] An agent that stabilizes cell membranes may be added to the maternal blood sample to reduce maternal cell lysis including but not limited to aldehydes, urea formaldehyde, phenol formaldehyde, DMAE (dimethylaminoethanol), cholesterol, cholesterol derivatives, high concentrations of magnesium, vitamin E, and vitamin E derivatives, calcium, calcium gluconate, taurine, niacin, hydroxylamine derivatives, bimoclomol, sucrose, astaxanthin, glucose, amitriptyline, isomer A hopane tetral phenylacetate, isomer B hopane tetral phenylacetate, citicoline, inositol, vitamin B, vitamin B complex, cholesterol hemisuccinate, sorbitol, calcium, coenzyme Q, ubiquinone, vitamin K, vitamin K complex, menaquinone, zonegran, zinc, ginkgo biloba extract, diphenylhydantoin, perftoran, polyvinylpyrrolidone, phosphatidylserine, tegretol, PABA, disodium cromglycate, nedocromil sodium, phenytoin, zinc citrate, mexitil, dilantin, sodium hyaluronate, or polaxamer 188.

## Isolation of Plasma and Maternal Cells

[1207] The blood is stored at 4°C until processing. The tubes are spun at 1000 rpm for ten minutes in a centrifuge with braking power set at zero. The tubes are spun a second time at 1000 rpm for ten minutes. The supernatant (the plasma) of each sample is transferred to a new tube and spun at 3000 rpm

for ten minutes with the brake set at zero. The supernatant is transferred to a new tube and stored at -80°C. Approximately two milliliters of the “buffy coat,” which contains maternal cells, is placed into a separate tube and stored at -80°C.

### **Isolation of DNA**

[1208] DNA is isolated from the plasma sample using the Qiagen Midi Kit for purification of DNA from blood cells, following the manufacturer’s instructions (QIAmp DNA Blood Midi Kit, Catalog number 51183). DNA is eluted in 100 µl of distilled water. The Qiagen Midi Kit also is used to isolate DNA from the maternal cells contained in the “buffy coat.”

### **Identification of Homozygous Maternal SNPs**

#### **SNP-IT™ Assay**

[1209] The SNP-IT Assay is based on a single base primer extension. Prior to the SNP-IT Assay, a PCR product that includes the SNP of interest is prepared, using one unmodified and one phosphorothiolate modified primer. The PCR product is then rendered single stranded using exonuclease, and the single stranded DNA is annealed to a SNP-IT oligonucleotide immobilized on the surface of a 96-well microtiter plate. Following hybridization, single base extension occurs by the addition of DNA polymerase and the labeled terminators. The incorporated base is detected using antibodies specific to the label followed by colorimetric detection. Data analysis can be done visually or with the use of the absorbance plate reader.

#### **Primer Design**

[1210] For each locus of interest, the SNP-IT™ Assay requires three primers. The primers are designed to produce an amplicon of 100-50 base pairs. The sequence of the SNP-IT primer, which is designed to anneal immediately upstream of the SNP site, is the best sequence available from between the upper and the lower strands. The sequence of the SNP-IT primer is designed to minimize hybridization to self and other sites in the amplicon. In addition, the SNP-IT primer may contain modified bases to prevent self-priming. The length of the primers is designed according to the makers of the SNP-IT Assay.

[1211] Representative primers for the amplification and genotyping of SNP TSC0069085, which is located on chromosome 21, are provided below:

Upstream Primer:

5' ATCACACTGGGGATC 3' (SEQ ID NO: 654)

Downstream Primer:

5' CTAAACCTATGACTC 3' (SEQ ID NO: 655)

SNP-IT primer

5' TTCACAGAGGATATCTTAATA 3' (SEQ ID NO: 656)

[1212] The upstream primer is unmodified and the downstream primer is phosphorothiolate modified.

### **SNP-IT Plate Coating**

[1213] SNP-IT primer is added to coat wells of empty 96-well plates. This reaction typically incubates overnight. The manufacturer's protocols and procedures are followed.

### **PCR**

[1214] Template DNA (15 ng) is amplified either in a reaction vessel including but not limited to an eppendorf tube or a well of a microtiter plate. The manufacturer's protocols and procedures are followed for the PCR reaction.

### **Exonuclease**

[1215] PCR product is treated with exonuclease to degrade the unmodified strand. The protected phosphorothiolate-labeled strand is used in the SNP-IT Assay. The manufacturer's protocols and procedures are followed for the exonuclease reaction.

### **Annealing**

[1216] Single stranded PCR product is transferred to SNP-IT plate and is allowed to form a hybrid with the SNP-IT primer. The annealing reaction typically proceeds for one hour. The manufacturer's protocols and procedures are followed for the annealing reaction.

### **SNP-IT Reaction**

[1217] The extension reagent, which contains DNA polymerase, two terminating nucleotides labeled with either fluorescein or biotin and two unlabeled terminators, is added to the SNP-IT well containing the annealed template and primer complex. For SNP TSC0069085, ddCTP is labeled with fluorescein and ddTTP is labeled with biotin, and the unlabeled terminators are ddATP, and ddGTP. The manufacturer's protocols and procedures are followed for the extension reaction.

**[1218]** The SNP specific base is incorporated by single base extension of the SNP-IT primer. Primers are washed manually or in a plate washer to remove unincorporated material. The manufacturer's protocols and procedures are followed for the washing reaction.

#### **Detection**

**[1219]** Anti-fluorescein labeled with alkaline phosphatase (AP) is added to the plate and allowed to bind to any incorporated fluorescein labeled terminator. The manufacturer's protocols and procedures for the labeling reaction are followed.

**[1220]** The plates are washed, and then color development is performed using pNPP as the detection substrate. The absorbance is read at 405 nm to detect yellow colored pNPP substrate followed by a washing step to remove pNPP detection reagents. The manufacturer's protocols and procedures are followed for color development and washing steps.

**[1221]** Streptavidin labeled with horse radish peroxidase (HRP) is added to the plate and allowed to bind to any incorporated biotin labeled terminator. The manufacturer's protocols and procedures are followed for the labeling reaction.

**[1222]** Following washing, color development is performed using TMB as the detection substrate. The absorbance is read at 620 nm to detect blue colored TMB substrate.

#### **Analysis**

**[1223]** Absorbance is plotted to generate a scatter plot from which genotype calls are made. SNPs at which the maternal DNA is homozygous are identified, and analyzed with the DNA isolated from the maternal plasma.

#### **Analysis of DNA Isolated from Maternal Plasma**

**[1224]** After the maternal DNA is analyzed and homozygous SNPs are identified, these SNPs are analyzed with the DNA isolated from the plasma. A low copy number of fetal genomes typically exists in the maternal plasma. To increase the copy number of the loci of interest, which are the SNPs at which the maternal DNA is homozygous, primers are designed to anneal at approximately 130 bases upstream and 130 bases downstream of each loci of interest. This is done to reduce statistical sampling error that can occur when working with a low number of genomes, which can influence the ratio of one allele to another (see Example 11).

## **Design of Multiplex Primers**

**[1225]** The primers are 12 bases in length. However, primers of any length can be used including but not limited to 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36-45, 46-55, 56-65, 66-75, 76-85, 86-95, 96-105, 106-115, 116-125, and greater than 125 bases. Primers are designed to anneal to both the sense strand and the antisense strand.

**[1226]** The maternal homozygous SNPs vary from sample to sample so defined sequences are not provided here. Primers are designed to anneal about 130 bases upstream and downstream of the maternal homozygous SNPs. The primers are designed to terminate at the 3' end in the dinucleotide "AA" to reduce the formation of primer-dimers. However, the primers can be designed to end in any of the four nucleotides and in any combination of the four nucleotides.

## **Multiplex PCR**

**[1227]** Regions upstream and downstream of the maternal homozygous SNPs are amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). This PCR reaction uses primers that anneal approximately 130 bases upstream and downstream of each loci of interest. The primers are mixed together and are used in a single reaction to amplify the template DNA. This reaction is done to increase the number of copies of the loci of interest, which eliminates error generated from a low number of genomes.

**[1228]** For increased specificity, a "hot-start" PCR reaction is used. PCR reactions are performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction is optimized for each locus of interest. In this example, the 20 µl of plasma template DNA is used.

**[1229]** Two microliters of each forward and reverse primer, at concentrations of 5 mM are pooled into a single microcentrifuge tube and mixed. Four microliters of the primer mix is used in a total PCR reaction volume of 50 µl (20µl of template plasma DNA, 1 µl of sterile water, 4 µl of primer mix, and 25 µl of HotStar Taq. Twenty-five cycles of PCR are performed. The following PCR conditions are used:

- (1) 95°C for 15 minutes;
- (2) 95°C for 30 second;
- (3) 4°C for 30 seconds;
- (4) 37°C for 30 seconds;
- (5) Repeat steps 2-4 twenty-four (24) times;
- (6) 72°C for 10 minutes.

[1230] The temperatures and times for denaturing, annealing, and extension, are optimized by trying various settings and using the parameters that yield the best results.

[1231] Other methods of genomic amplification can also be used to increase the copy number of the loci of interest including but not limited to primer extension preamplification (PEP) (Zhang *et al.*, PNAS, 89:5847-51, 1992), degenerate oligonucleotide primed PCR (DOP-PCR) (Telenius, *et al.*, Genomics 13:718-25, 1992), strand displacement amplification using DNA polymerase from bacteriophage 29, which undergoes rolling circle replication (Dean *et al.*, Genomic Research 11:1095-99, 2001), multiple displacement amplification (U.S. Patent 6,124,120), REPLI-g™ Whole Genome Amplification kits, and Tagged PCR.

#### **Purification of Fragment of Interest**

[1232] The unused primers, and nucleotides are removed from the reaction by using Qiagen MinElute PCR purification kits (Qiagen, Catalog Number 28004). The reactions are performed following the manufacturer's instructions supplied with the columns. The DNA is eluted in 100 µl of sterile water. 5 µl of each amplified loci is mixed together

#### **SNP-IT Assay**

[1233] The pooled DNA is assayed with the SNP-IT assay as described above. Each SNP is genotyped. SNPs located on chromosomes 13 and 21, wherein the maternal DNA is homozygous, and DNA isolated from the plasma is heterozygous are quantitated.

#### **Quantification**

[1234] The fluorescence intensity of each allele is quantitated. As discussed above, the expected ratio of allele 1 to allele 2 is used to determine the presence or absence of a chromosomal abnormality. If the maternal genome is homozygous at SNP X (A/A), and the plasma DNA is heterozygous at SNP X (A/G), then the G represents the distinct fetal signal.

[1235] The intensity of the allele with the A nucleotide is quantitated and the intensity of the allele with the G nucleotide is quantitated. The ratio of G:A depends on the percentage of fetal DNA present in the maternal blood.

[1236] For example, if the sample contains 50% fetal DNA, then the expected ratio is 0.33 (1 fetal G allele/(2 maternal A alleles + 1 fetal A allele)). This ratio should be constant for all chromosomes that are present in two copies. The ratio that is obtained for SNPs on chromosome 13 should be the same as the ratio that is obtained for chromosome 21.



[1237] However, if the fetal genome contains an additional copy of chromosome 21, then the ratio for this chromosome will deviate from the expected ratio. The expected ratio for a Trisomy condition with 50% fetal DNA in the maternal blood is 0.25. Thus, by analyzing SNPs wherein the maternal genome is homozygous, and the DNA that is isolated from the plasma is heterozygous, fetal chromosomal abnormalities can be detected.

[1238] In this example, the terminator nucleotides are labeled with different chemical moieties. However, using the methods described in this application (see Example 6), the SNP-IT assay could be modified to allow detection of both alleles with a single labeled terminator.

[1239] This example explained the use of Orchid's SNP-IT assay, but it not intended to limit the use of other techniques that rely on primer extension. Orchid's SNPstream 25K, as well accompanying software including but not limited to GetGenos™, QCreview™, and ValidGenos™, can also be used to detect the presence of chromosomal abnormalities in the maternal blood. Additional information about these products can be found at:

<http://www.orchidbio.com/products/lsg/products/snpstream.asp>.

## EXAMPLE 22

[1240] Fetal chromosomal abnormalities are determined by analyzing SNPs wherein the maternal template DNA is homozygous and the template DNA obtained from the plasma is heterozygous. Plasma that is isolated from blood of a pregnant female contains both maternal template DNA and fetal template DNA. Any number of SNP detection methods can be used to analyze the maternal and plasma DNA. In this example, SNPs are analyzed using the TaqMan® assay. However, other methods that rely on fluorogenic 5' nuclease assay can be used.

### Collection of Blood Samples

[1241] In accordance with an IRB approved study, blood samples are collected from pregnant women after informed consent is granted. Blood is collected into 9 ml EDTA Vacuette tubes (catalog number NC9897284) and 0.225 ml of 10% neutral buffered solution containing formaldehyde (4% w/v), is added to each tube, and each tube gently is inverted. The tubes are stored at 4°C until ready for processing.

[1242] Any number of agents that impede cell lysis or stabilize cell membranes can be added to the tubes including but not limited to formaldehyde, and derivatives of formaldehyde, formalin, glutaraldehyde, and derivatives of glutaraldehyde, crosslinkers, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive

crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA, DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS or compounds listed in Table XXIII. Any concentration of agent that stabilizes cell membranes or impedes cell lysis can be added. In a preferred embodiment, the agent that stabilizes cell membranes or impedes cell lysis is added at a concentration that does not impede or hinder subsequent reactions.

**[1243]** An agent that stabilizes cell membranes may be added to the maternal blood sample to reduce maternal cell lysis including but not limited to aldehydes, urea formaldehyde, phenol formaldehyde, DMAE (dimethylaminoethanol), cholesterol, cholesterol derivatives, high concentrations of magnesium, vitamin E, and vitamin E derivatives, calcium, calcium gluconate, taurine, niacin, hydroxylamine derivatives, bimoclomol, sucrose, astaxanthin, glucose, amitriptyline, isomer A hopane tetral phenylacetate, isomer B hopane tetral phenylacetate, citicoline, inositol, vitamin B, vitamin B complex, cholesterol hemisuccinate, sorbitol, calcium, coenzyme Q, ubiquinone, vitamin K, vitamin K complex, menaquinone, zonegran, zinc, ginkgo biloba extract, diphenylhydantoin, perftoran, polyvinylpyrrolidone, phosphatidylserine, tegretol, PABA, disodium cromglycate, nedocromil sodium, phenytoin, zinc citrate, mexitil, dilantin, sodium hyaluronate, or polaxamer 188.

#### **Isolation of Plasma and Maternal Cells**

**[1244]** The blood is stored at 4°C until processing. The tubes are spun at 1000 rpm for ten minutes in a centrifuge with braking power set at zero. The tubes are spun a second time at 1000 rpm for ten minutes. The supernatant (the plasma) of each sample is transferred to a new tube and spun at 3000 rpm for ten minutes with the brake set at zero. The supernatant is transferred to a new tube and stored at -80°C. Approximately two milliliters of the “buffy coat,” which contains maternal cells, is placed into a separate tube and stored at -80°C.

#### **Isolation of DNA**

**[1245]** DNA is isolated from the plasma sample using the Qiagen Midi Kit for purification of DNA from blood cells, following the manufacturer’s instructions (QIAmp DNA Blood Midi Kit, Catalog number 51183). DNA is eluted in 100 µl of distilled water. The Qiagen Midi Kit also is used to isolate DNA from the maternal cells contained in the “buffy coat.”

## Identification of Homozygous Maternal SNPs

### TaqMan Assay

[1246] PE Biosystems has two instruments in its Sequence Detection Systems product line, the ABI Prism® 7700 Sequence Detection System and the GeneAmp® 5700 Sequence Detection System. These real-time systems allegedly are capable of detecting PCR products as they accumulate during PCR and so enable the quantitation of DNA in the sample.

[1247] One chemistry available for use on the ABI PRISM® 7700 and GeneAmp® 5700 detection systems is the fluorogenic 5' nuclease assay, or the TaqMan® assay, which uses a fluorogenic probe to enable the detection of a specific PCR product as it accumulates during PCR. PE Biosystems' patented fluorogenic probe design that incorporates the reporter due on the 5' end and the quencher on the 3' end has assisted with the design of TaqMan probes.

[1248] The basis for PCR quantitation in the ABI 7700 instrument is to continuously measure PCR product accumulation using a dual-labeled ~~fluoregenic~~fluorogenic oligonucleotide probe called a TaqMan® probe, which is composed of a short (20-25 bases) oligodeoxynucleotide that is labeled with two different fluorescent dyes. On the 5' terminus is a reporter dye and on the 3' terminus is a quenching dye. This oligonucleotide probe sequence is homologous to an internal sequence present in the PCR amplicon. When the probe is intact, energy transfer occurs between the two fluorophors and emission from the reporter is quenched by the quencher (Livak *et al.*, PCR Methods and Applications, 4:357-362, 1995a; U.S. Patent No. 5,538,848; U.S. Patent No. 5,723,591).

[1249] During the extension phase of PCR, the probe is cleaved by 5' nuclease activity of Taq polymerase thereby releasing the reporter from the oligonucleotide-quencher and producing an increase in reporter emission intensity. The ABI Prism 7700 uses fiber optic systems which connect to each well in a 96-well PCR tray format. The laser light source excites each well and a CCD camera measures the fluorescence spectrum and intensity from each well to generate real-time data during PCR amplification. The ABI 7700 Prism software examines the fluorescence intensity of reporter and quencher dyes and calculates the increase in normalized reporter emission intensity over the course of the amplification. The results are then plotted versus time, represented by cycle number, to produce a continuous measure of PCR amplification. To provide precise quantification of initial target in each PCR reaction, the amplification plot is examined at a point during the early log phase of product accumulation. This is accomplished by assigning a fluorescence threshold above background and determining the time point at which each sample's amplification plot reaches the threshold (defined as the threshold cycle number or

CT). Differences in threshold cycle number are used to quantify the relative amount of PCR target contained within each tube as described previously.

[1250] For SNP analysis, a TaqMan probe can be designed for each allele of the SNP. The reporter emission is used to determine the presence or absence of each allele at the SNP. For example, for a SNP that can either be adenine or guanine, a TaqMan probe will be designed with a complementary nucleotide to the adenine and a separate TaqMan probe will be designed with a complementary nucleotide to the guanine. The two TaqMan probes can be used in separate reaction vessels, which allows the amount of the adenine allele and the amount of the guanine allele to be calculated.

### **Primer and Probe Design**

[1251] Primer and probes can be designed using the Primer Express® software. The probe is designed first, and then the primers are designed as close as possible to the probe without overlapping it. Amplicons of 50-150 base pairs are strongly recommended.

[1252] The primer and probes should be designed following the manufacturer's recommendations. For both the primer and the probes, the G/C is in the range of 20-80%. The primer and probes are designed to avoid runs of an identical nucleotide. This is especially true for guanine, where runs of four or more Gs should be avoided.

[1253] For the probe, the  $T_M$  is about 68-70°C, and is designed so that there is no guanine on the 5' end. Also, the probe is designed so that there are more C than G bases.

[1254] For the primers, the  $T_M$  is about 58-60°C, and the primers are designed so that the five nucleotides at the 3' end have no more than 2 G and/or C bases.

[1255] For example, representative primers and probes for SNP TSC0271628 (A/G), which is located on chromosome 21, are provided below:

Forward Primer ( $T_M$  of 60°C)

5' AGTCTTGTAATACGACAGTCTT 3' (SEQ ID NO: 657)

Reverse Primer ( $T_M$  of 58°C)

5' CCATATCAATCAGTACTCTTG 3' (SEQ ID NO: 658)

TaqMan Probe A allele ( $T_M$  of 68°C; bold indicates variable nucleotide at SNP)

5' CCTTTCCACTCATCCAAAGGTTG 3' (SEQ ID NO: 659)

TaqMan Probe G allele ( $T_M$  of 70°C; bold indicates variable nucleotide at SNP)

5' CCTTTCCACTCATCCAGAGGTTG 3' (SEQ ID NO: 660)

[1256] Information regarding the sequence surrounding SNPs is found at: <http://www.snp.cshl.org>. By independently varying forward and reverse primer concentrations, the concentrations that provide optimal assay conditions can be identified. Primer concentration ranges of 50 nM – 900 nM are tested.

[1257] If the maternal DNA is homozygous for one allele, for example, adenine, then in the sample that contains the TaqMan probe specific for the guanine nucleotide, the reporter is not separated from the quencher because the TaqMan probe does not anneal to the template DNA. However, if the maternal DNA is homozygous, then the reporter will be separated from the quencher in both samples containing the TaqMan probe specific for the guanine allele and samples containing the TaqMan probe specific for the adenine allele.

### **Reagent Solution**

[1258] The polymerase recommended for the TaqMan Assay is AmpliTaq Gold DNA polymerase. It is thought that the use of AmpliTaq Gold DNA polymerase reduces the amount of non-specific product formation. The incorporation of AmpErase® Uracil n-glycosylase (UNG) and dUTP provide protection against PCR carryover contamination. For PCR reactions, the TaqMan Universal PCR Master mix, which is a reagent designed to provide optimal performance for TaqMan assays, is recommended by the manufacturer.

[1259] The TaqMan reaction buffer contains 5.5 mM MgCl<sub>2</sub>, 200 nM each of dATP, dCTP, dGTP, 400 nM dUTP, 0.5 U of uracyl DNA glycosylase, and 1.25 U of AmpliTaq gold.

### **Thermal Cycling Parameters**

[1260] PCR amplification and detection for all primer-probe combinations are performed with the ABI 7700 Sequence Detection System. The recommended cycling parameters for the TaqMan assay are provided below:

- 1) 50°C for 2 min;
- 2) 95°C for 10 min;
- 3) 95°C for 15 sec;
- 4) 60°C for 1 min;
- 5) Repeat steps 3-4 for 40 cycles.

### **TaqMan Quantitation**

[1261] External standards are generated from known quantities of DNA containing an adenine nucleotide at SNP TSC0271628 and a guanine nucleotide at SNP TSC0271628, spanning 6 orders of

magnitude (from  $5 \times 10^0$  to  $5 \times 10^6$  copies). The detection threshold is set at 10 times the standard deviation of the mean baseline emission calculated for PCR cycles 3 to 15 (Stults *et al.*, Applied and Environmental Microbiology, Vol. 67, No. 6, 2781-2789, 2001). Standard curves relating the threshold cycle to DNA concentrations are generated with the ABI Prism 7700 software (available from Perkin Elmer).

[1262] SNPs at which the maternal DNA is homozygous are identified, and are analyzed with the DNA isolated from the plasma sample.

#### **Analysis of DNA Isolated from Maternal Plasma**

[1263] After the maternal DNA is analyzed and homozygous SNPs are identified, these SNPs are analyzed with the DNA isolated from the plasma. A low copy number of fetal genomes typically exists in the maternal plasma. To increase the copy number of the loci of interest, which are the SNPs at which the maternal DNA is homozygous, primers are designed to anneal at approximately 130 bases upstream and 130 bases downstream of each loci of interest. This is done to reduce statistical sampling error that can occur when working with a low number of genomes, which can influence the ratio of one allele to another (see Example 11).

#### **Design of Multiplex Primers**

[1264] The primers are 12 bases in length. However, primers of any length can be used including but not limited to 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36-45, 46-55, 56-65, 66-75, 76-85, 86-95, 96-105, 106-115, 116-125, and greater than 125 bases. Primers are designed to anneal to both the sense strand and the antisense strand.

[1265] The maternal homozygous SNPs vary from sample to sample so defined sequences are not provided here. Primers are designed to anneal about 130 bases upstream and downstream of the maternal homozygous SNPs. The primers are designed to terminate at the 3' end in the dinucleotide "AA" to reduce the formation of primer-dimers. However, the primers can be designed to end in any of the four nucleotides and in any combination of the four nucleotides.

#### **Multiplex PCR**

[1266] Regions upstream and downstream of the maternal homozygous SNPs are amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). This PCR reaction uses primers that anneal approximately 130 bases upstream and downstream of each loci of interest. The primers are mixed together and are used

in a single reaction to amplify the template DNA. This reaction is done to increase the number of copies of the loci of interest, which eliminates error generated from a low number of genomes.

[1267] For increased specificity, a "hot-start" PCR reaction is used. PCR reactions are performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction is optimized for each locus of interest. In this example, the 20  $\mu$ l of plasma template DNA is used.

[1268] Two microliters of each forward and reverse primer, at concentrations of 5 mM are pooled into a single microcentrifuge tube and mixed. Four microliters of the primer mix is used in a total PCR reaction volume of 50  $\mu$ l (20 $\mu$ l of template plasma DNA, 1  $\mu$ l of sterile water, 4  $\mu$ l of primer mix, and 25  $\mu$ l of HotStar Taq. Twenty-five cycles of PCR are performed. The following PCR conditions are used:

- (1) 95°C for 15 minutes;
- (2) 95°C for 30 second;
- (3) 4°C for 30 seconds;
- (4) 37°C for 30 seconds;
- (5) Repeat steps 2-4 twenty-four (24) times;
- (6) 72°C for 10 minutes.

[1269] The temperatures and times for denaturing, annealing, and extension, are optimized by trying various settings and using the parameters that yield the best results.

[1270] Other methods of genomic amplification can also be used to increase the copy number of the loci of interest including but not limited to primer extension preamplification (PEP) (Zhang *et al.*, PNAS, 89:5847-51, 1992), degenerate oligonucleotide primed PCR (DOP-PCR) (Telenius, *et al.*, Genomics 13:718-25, 1992), strand displacement amplification using DNA polymerase from bacteriophage 29, which undergoes rolling circle replication (Dean *et al.*, Genomic Research 11:1095-99, 2001), multiple displacement amplification (U.S. Patent 6,124,120), REPLI-g™ Whole Genome Amplification kits, and Tagged PCR.

[1271] It is important to ensure that the region amplified contains annealing sequences for the oligonucleotide probes in the BeadArray. Upon purchase of the BeadArray service, each SNP and the primers used to analyze each SNP are identified. With this knowledge, the multiplex primers are designed to encompass annealing regions for the primers in the BeadArray.

### **Purification of Fragment of Interest**

[1272] The unused primers, and nucleotides are removed from the reaction by using Qiagen MinElute PCR purification kits (Qiagen, Catalog Number 28004). The reactions are performed following the manufacturer's instructions supplied with the columns. The DNA is eluted in 100 µl of sterile water.

### **TaqMan Assay**

[1273] The amplified DNA is assayed with the TaqMan assay as described above. Each SNP is genotyped. SNPs located on chromosomes 13 and 21, wherein the maternal DNA is homozygous, and DNA isolated from the plasma is heterozygous are quantitated.

### **Quantification**

[1274] The fluorescent intensity of the TaqMan allele specific probe is quantitated. As discussed above, the expected ratio of allele 1 to allele 2 is used to determine the presence or absence of a chromosomal abnormality. If the maternal genome is homozygous at SNP X (A/A), and the plasma DNA is heterozygous at SNP X (A/G), then the G represents the distinct fetal signal.

[1275] The fluorescent intensity of the allele with the A nucleotide is quantitated and the intensity of the allele with the G nucleotide is quantitated. The ratio of G:A depends on the percentage of fetal DNA present in the maternal blood.

[1276] For example, if the sample contains 50% fetal DNA, then the expected ratio is 0.33 (1 fetal G allele/(2 maternal A alleles + 1 fetal A allele)). This ratio should be constant for all chromosomes that are present in two copies. The ratio that is obtained for SNPs on chromosome 13 should be the same as the ratio that is obtained for chromosome 21.

[1277] However, if the fetal genome contains an additional copy of chromosome 21, then the ratio for this chromosome will deviate from the expected ratio. The expected ratio for a Trisomy condition with 50% fetal DNA in the maternal blood is 0.25. Thus, by analyzing SNPs wherein the maternal genome is homozygous, and the DNA that is isolated from the plasma is heterozygous, fetal chromosomal abnormalities can be detected.

[1278] This example explained the use of the TaqMan assay, but it not intended to limit the use of other techniques that employ 5' nuclease activity. For example, the SYBR® Green I double stranded dye can also be used with the ABI PRISM 7700 Sequence Detection System and the GeneAmp® 5700 Sequence Detection system for determining the sequence of maternal and fetal DNA. The SYBR® Green I double stranded dye assay may be used to detect the presence of fetal chromosomal abnormalities in the maternal blood.



[1279] SYBR® Green I double stranded dye is a highly specific double-stranded DNA binding dye that allows the detection of product accumulation during PCR. However, the SYBR® Green I double stranded dye assay detects all double stranded DNA including non-specific reaction products. The advantage of the SYBR® Green I double stranded dye assay is that it does not require a probe.

[1280] The same primers design parameters are recommended for both the TaqMan Assay and the SYBR® Green I double stranded dye assay (see Primer Design section in Example 21). The primer optimization parameters recommended for the TaqMan assay should also be followed for the SYBR® Green I double stranded dye assay. In addition, no template controls should also be run with the various concentrations of primers.

[1281] Furthermore, Applied Biosystems sell other products that may be used to determine the sequence of maternal and fetal DNA including but not limited to Assays-on-Demand™ SNP genotyping products, and Assays-by-Design<sup>SM</sup> Service SNP genotyping products.

[1282] Having now fully described the invention, it will be understood by those of skill in the art that the invention can be performed with a wide and equivalent range of conditions, parameters, and the like, without affecting the spirit or scope of the invention or any embodiment thereof.

[1283] All documents, e.g., scientific publications, patents and patent publications recited herein are hereby incorporated by reference in their entirety to the same extent as if each individual document was specifically and individually indicated to be incorporated by reference in its entirety. Where the document cited only provides the first page of the document, the entire document is intended, including the remaining pages of the document.

### EXAMPLE 23

[1284] Fetal chromosomal abnormalities are determined by analyzing SNPs wherein the maternal template DNA is homozygous and the template DNA obtained from the plasma is heterozygous. Plasma that is isolated from blood of a pregnant female contains both maternal template DNA and fetal template DNA. Any number of SNP detection methods can be used to analyze the maternal and plasma DNA. In this example, SNPs are analyzed by ThirdWave Technologies' Invader™ Assay for Nucleic Acid Detection. However, other techniques that exploit and quantitate biological structures formed in the presence of the correct sequence can be used.

### **Collection of Blood Samples**

[1285] In accordance with an IRB approved study, blood samples are collected from pregnant women after informed consent is granted. Blood is collected into 9 ml EDTA Vacuette tubes (catalog number NC9897284) and 0.225 ml of 10% neutral buffered solution containing formaldehyde (4% w/v), is added to each tube, and each tube gently is inverted. The tubes are stored at 4°C until ready for processing.

[1286] Any number of agents that impede cell lysis or stabilize cell membranes can be added to the tubes including but not limited to formaldehyde, and derivatives of formaldehyde, formalin, glutaraldehyde, and derivatives of glutaraldehyde, crosslinkers, primary amine reactive crosslinkers, sulfhydryl reactive crosslinkers, sulfhydryl addition or disulfide reduction, carbohydrate reactive crosslinkers, carboxyl reactive crosslinkers, photoreactive crosslinkers, cleavable crosslinkers, AEDP, APG, BASED, BM(PEO)<sub>3</sub>, BM(PEO)<sub>4</sub>, BMB, BMDB, BMH, BMOE, BS3, BSOCOES, DFDNB, DMA, DMP, DMS, DPDPB, DSG, DSP, DSS, DST, DTBP, DTME, DTSSP, EGS, HBVS, sulfo-BSOCOES, Sulfo-DST, Sulfo-EGS or compounds listed in Table XXIII. Any concentration of agent that stabilizes cell membranes or impedes cell lysis can be added. In a preferred embodiment, the agent that stabilizes cell membranes or impedes cell lysis is added at a concentration that does not impede or hinder subsequent reactions.

[1287] An agent that stabilizes cell membranes may be added to the maternal blood samples to reduce maternal cell lysis including but not limited to aldehydes, urea formaldehyde, phenol formaldehyde, DMAE (dimethylaminoethanol), cholesterol, cholesterol derivatives, high concentrations of magnesium, vitamin E, and vitamin E derivatives, calcium, calcium gluconate, taurine, niacin, hydroxylamine derivatives, bimoclomol, sucrose, astaxanthin, glucose, amitriptyline, isomer A hopane tetral phenylacetate, isomer B hopane tetral phenylacetate, citicoline, inositol, vitamin B, vitamin B complex, cholesterol hemisuccinate, sorbitol, calcium, coenzyme Q, ubiquinone, vitamin K, vitamin K complex, menaquinone, zonegran, zinc, ginkgo biloba extract, diphenylhydantoin, perftoran, polyvinylpyrrolidone, phosphatidylserine, tegretol, PABA, disodium cromglycate, nedocromil sodium, phenytoin, zinc citrate, mexitil, dilantin, sodium hyaluronate, or polaxamer 188.

### **Isolation of Plasma and Maternal Cells**

[1288] The blood is stored at 4°C until processing. The tubes are spun at 1000 rpm for ten minutes in a centrifuge with braking power set at zero. The tubes are spun a second time at 1000 rpm for ten minutes. The supernatant (the plasma) of each sample is transferred to a new tube and spun at 3000 rpm

for ten minutes with the brake set at zero. The supernatant is transferred to a new tube and stored at -80°C. Approximately two milliliters of the “buffy coat,” which contains maternal cells, is placed into a separate tube and stored at -80°C.

### **Isolation of DNA**

[1289] DNA is isolated from the plasma sample using the Qiagen Midi Kit for purification of DNA from blood cells, following the manufacturer’s instructions (QIAmp DNA Blood Midi Kit, Catalog number 51183). DNA is eluted in 100 µl of distilled water. The Qiagen Midi Kit also is used to isolate DNA from the maternal cells contained in the “buffy coat.”

### **Identification of Homozygous Maternal SNPs**

#### **ThirdWave Technologies Invader™ Assay**

[1290] The Invader™ Assay, which was developed by Third Wave Technologies (Madison, WI), is an isothermal, “PCR-free” approach to the detection and quantitative analysis of DNA. The Invader Assay produces and amplifies an unrelated signal only in the presence of the correct target sequence. The Invader™ Assay relies on a thermostable member of the structure-specific archeobacterial flap endonuclease (FEN) family, which cleaves nucleic acid molecules at specific sites based on structure rather than sequence. When used with structure forming probes for known sequences, the enzymes cleave in a structure and target sequence-specific manner. The nucleases used with Third Wave Technologies’ assays are referred to as “Cleavase®” enzymes.

[1291] The Invader™ Assay uses two target-specific oligonucleotides to create the substrate complex recognized by Cleavase Enzymes (L. DeFrancesco, *The Scientist*, 12(21):16, 1998). The substrate complex is formed when an upstream Invader oligo and a downstream signal probe hybridize in tandem to the nucleic acid. The 3’ end of the Invader oligo must overlap the hybridization site of the signal probe by at least one base (Harrington *et al.*, *Genes and Development*, 8:1344-55, 1994). The 5’ end of the signal probe has additional unpaired bases to form a 5’ flap. Cleavase enzymes cleave the signal probe where it overlaps the Invader oligo, releasing the 5’ arm. Reaction mixtures contain excess signal probe and are carried out near the melting temperature of the probe. Many signal probes can be cleaved for each copy of the target without temperature cycling.

[1292] The overlap between the Invader oligo and the signal probe is important. A mismatch positioned at the site of the overlap will block the cleavage by disrupting the overlap, which may allow discrimination of SNPs and mutations.

[1293] The Invader assay utilizes two sequential cleavage steps. The 5' arm of the signal probe released in the first reaction is not detected directly. Rather, a secondary cleavage product is the actual source of the signal, detected by fluorescence resonance energy transfer (FRET). The primary cleavage product, which is the 5' arm released in the first reaction, is used as an Invader oligo that hybridizes to a supplied FRET probe in the secondary reaction. The FRET probe is labeled with two dyes: a donor fluorophore and a quenching acceptor fluorophore. When the nuclease cleaves the secondary probe, the two fluorophores are separated, quenching is eliminated, and the enhanced fluorescence signal from the donor dye is detected.

### Probe Design

[1294] The following probes are designed to determine the sequence at SNP TSC1172576 (T/A), which is located on chromosome 13:

Invader Oligo for T allele:

5' CATGCAGATATACCGCATAT 3' (SEQ ID NO: 661)

Invader Oligo for A allele:

5' CATGCAGATATACCGCATAA 3' (SEQ ID NO: 662)

[1295] Invader oligonucleotides are designed to be complementary to an 18-22 base region immediately upstream of the signal probe, with an additional one base at the 3' end that "invades" the region hybridized to the signal probe by one base.

Signal probe for T allele:

GGTAGCATCTCTCAGCACAAGAG (SEQ ID NO: 663)

Signal probe for A allele:

GGTAGCATCACTCAGCACAAGAG (SEQ ID NO: 664)

[1296] The signal probes are designed to contain a 3' region that is complementary to the target sequence and a non-complementary 5' arm (the underlined sequence above) that is used for detection. The signal probes are labeled on the 5' end with 6-carboxyfluorescein (TET), hexachloro-6-carboxyfluorescein (HEX), 6-carboxyfluorescein (FAM). However, the 5' end can be labeled with any chemical moiety including but not limited to radioisotope, fluorescent reporter molecule, chemiluminescent reporter molecule, antibody, antibody fragment, hapten, biotin, derivative of biotin, photobiotin, iminobiotin, digoxigenin, avidin, enzyme, acridinium, sugar, enzyme, apoenzyme, homopolymeric oligonucleotide, hormone, ferromagnetic moiety, paramagnetic moiety, diamagnetic moiety, phosphorescent moiety, luminescent moiety, electrochemiluminescent moiety, chromatic moiety,

moiety having a detectable electron spin resonance, electrical capacitance, dielectric constant or electrical conductivity, and combinations thereof

[1297] Both invader and signal probes are complementary to either the sense or the antisense target DNA strand, depending on which results in the formation of the least number of predictable secondary structures.

### PCR Amplification

[1298] A fragment of DNA that surrounds SNP TSC1172576 is amplified by PCR. The sequence of the upstream and downstream primers is provided below:

Upstream Primer:

5' TAGCAGAATCTCTCAT 3' (SEQ ID NO: 665)

Downstream Primer:

5' AGAGTATCTCATTTGTT 3' (SEQ ID NO: 666)

[1299] Amplification reactions are performed in a final volume of 100  $\mu$ l of containing 2  $\mu$ l of genomic DNA, 35 pmol of each primer, 50  $\mu$ M of each deoxynucleotide (Perkin-Elmer Applied Biosystems, Inc., Foster City, CA), 1 X PCR buffer (20 mM Tris-HCl, 50 mM KCl, 1.5 mM MgCl<sub>2</sub>, 0.05% Tween-20, 0.05% NP40), 1 M betaine, 5% dimethylsulfoxide (DMSO), and 2.5 U of Taq polymerase (Roche Boehringer Mannheim, Indianapolis, IN). PCR cycling conditions consist of an initial denaturation step at 95°C for 5 min, 30 cycles of denaturation at 95°C for 1 min, annealing at 68°C for 1 min, and extension at 72°C for 1 min, and a final extension at 72°C for 5 min.

### Invader Reaction

[1300] One microliter of each PCR product is added to 0.5 pmol of the appropriate Invader oligonucleotide, 10 ng human genomic DNA (Promega Corp., Madison, WI) as the carrier, and morpholinepropanesulfonic acid (MOPS) buffer (pH 8.0) at a final concentration of 10 mM in a volume of 7  $\mu$ l. The mixtures are denatured for 5 min at 95°C and then cooled to reaction temperature of 60°C. Invader reactions are ~~initiated~~initiated by the addition of a mixture containing 30 ng of Cleavase VIII (Third Wave Technologies, Inc., Madison, WI) 25 mM MgCl<sub>2</sub>, and 10 pmol of the appropriate signal probe oligonucleotide in a volume of 3  $\mu$ l. Reaction mixtures are incubated for 60 min. The reactions are terminated by the addition of 10  $\mu$ l of 95% formamide -10 mM EDTA (pH 8.0)- 0.05% crystal violet. Following termination, the reactions are diluted 1:10 in reagent-grade water. Samples of 2  $\mu$ l are loaded and electrophoresed in a 24% denaturing polyacrylamide gel (18 cm by 25.5 cm by mm) on an automated

fluorescence sequencing apparatus (model 377, PE-ABI). The data are collected using filter set C and processed with GeneScan software.

[1301] In addition, 5 µl of each sample is electrophoresed in 20% (acrylamide to bisacrylamide, 19:1) denaturing polyacrylamide gels at 20 W. Gel cassettes (20 cm by 20 cm by 0.5 mM) are scanned with a fluorescent scanner (FMBIO-100; Hitachi Corp, San Bruno, CA) by using a 585 nm filter for TET and HEX labeled probes and a 505 nm filter for FAM labeled probes.

[1302] With the Invader assay, it is also possible to perform a second cleavage reaction where the released 5' arm of the signal probe is hybridized to another probe, and fluorescence resonance energy transfer (FRET) is used to detect the presence of a specific nucleic acid. The manufacturer's protocols are followed when using the FRET probe.

[1303] The genotype at each SNP is determined by analyzing the fluorescence intensity of each allele-specific signal probe. For example, for SNP TSC1172576, the presence of the T allele is determined by analyzing the amount of released 5' signal probe using the signal probe from the T allele signal probe (as described above). Likewise, the presence of the A allele is determined by analyzing the amount of released 5' signal probe from the A allele signal probe. The reactions can be performed in a single reaction vessel using two different chemical moieties, which can be analyzed under distinct conditions, or the A and T allele reactions can be performed in two different reaction vessels.

#### **Analysis of DNA Isolated from Maternal Plasma**

[1304] After the maternal DNA is analyzed and homozygous SNPs are identified, these SNPs are analyzed with the DNA isolated from the plasma. A low copy number of fetal genomes typically exist in the maternal plasma. To increase the copy number of the loci of interest, which are the SNPs at which the maternal DNA is homozygous, primers are designed to anneal at approximately 130 bases upstream and 130 bases downstream of each loci of interest. This is done to reduce statistical sampling error that can occur when working with a low number of genomes, which can influence the ratio of one allele to another (see Example 11).

#### **Design of Multiplex Primers**

[1305] The primers are 12 bases in length. However, primers of any length can be used including but not limited to 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36-45, 46-55, 56-65, 66-75, 76-85, 86-95, 96-105, 106-115, 116-125, and greater than 125 bases. Primers are designed to anneal to both the sense strand and the antisense strand.

[1306] The maternal homozygous SNPs vary from sample to sample so defined sequences are not provided here. Primers are designed to anneal about 130 bases upstream and downstream of the maternal homozygous SNPs. The primers are designed to terminate at the 3' end in the dinucleotide "AA" to reduce the formation of primer-dimers. However, the primers can be designed to end in any of the four nucleotides and in any combination of the four nucleotides.

### **Multiplex PCR**

[1307] Regions upstream and downstream of the maternal homozygous SNPs are amplified from the template genomic DNA using the polymerase chain reaction (PCR, U.S. Patent Nos. 4,683,195 and 4,683,202, incorporated herein by reference). This PCR reaction uses primers that anneal approximately 130 bases upstream and downstream of each loci of interest. The primers are mixed together and are used in a single reaction to amplify the template DNA. This reaction is done to increase the number of copies of the loci of interest, which eliminates error generated from a low number of genomes.

[1308] For increased specificity, a "hot-start" PCR reaction is used. PCR reactions are performed using the HotStarTaq Master Mix Kit supplied by QIAGEN (catalog number 203443). The amount of template DNA and primer per reaction is optimized for each locus of interest. In this example, the 20 µl of plasma template DNA is used.

[1309] Two microliters of each forward and reverse primer, at concentrations of 5 mM are pooled into a single microcentrifuge tube and mixed. Four microliters of the primer mix is used in a total PCR reaction volume of 50 µl (20µl of template plasma DNA, 1 µl of sterile water, 4 µl of primer mix, and 25 µl of HotStar Taq. Twenty-five cycles of PCR are performed. The following PCR conditions are used:

- (1) 95°C for 15 minutes;
- (2) 95°C for 30 second;
- (3) 4°C for 30 seconds;
- (4) 37°C for 30 seconds;
- (5) Repeat steps 2-4 twenty-four (24) times;
- (6) 72°C for 10 minutes.

[1310] The temperatures and times for denaturing, annealing, and extension, are optimized by trying various settings and using the parameters that yield the best results.

[1311] Other methods of genomic amplification can also be used to increase the copy number of the loci of interest including but not limited to primer extension preamplification (PEP) (Zhang *et al.*, PNAS, 89:5847-51, 1992), degenerate oligonucleotide primed PCR (DOP-PCR) (Telenius, *et al.*, Genomics

13:718-25, 1992), strand displacement amplification using DNA polymerase from bacteriophage 29, which undergoes rolling circle replication (Dean *et al.*, Genomic Research 11:1095-99, 2001), multiple displacement amplification (U.S. Patent 6,124,120), REPLI-g™ Whole Genome Amplification kits, and Tagged PCR.

[1312] It is important to ensure that the region amplified contains annealing sequences for the oligonucleotide probes in the BeadArray. Upon purchase of the BeadArray service, each SNP and the primers used to analyze each SNP are identified. With this knowledge, the multiplex primers are designed to encompass annealing regions for the primers in the BeadArray.

#### **Purification of Fragment of Interest**

[1313] The unused primers, and nucleotides are removed from the reaction by using Qiagen MinElute PCR purification kits (Qiagen, Catalog Number 28004). The reactions are performed following the manufacturer's instructions supplied with the columns. The DNA is eluted in 100 µl of sterile water.

#### **Invader Assay**

[1314] The amplified DNA is assayed with the Invader assay as described above. Each SNP is genotyped. SNPs located on chromosomes 13 and 21, wherein the maternal DNA is homozygous, and DNA isolated from the plasma is heterozygous are quantitated.

#### **Quantification**

[1315] The fluorescent intensity of the allele specific signal probe is quantitated. As discussed above, the expected ratio of allele 1 to allele 2 is used to determine the presence or absence of a chromosomal abnormality. If the maternal genome is homozygous at SNP X (A/A), and the plasma DNA is heterozygous at SNP X (A/G), then the G represents the distinct fetal signal.

[1316] The fluorescent intensity of the allele with the A nucleotide is quantitated and the intensity of the allele with the G nucleotide is quantitated. The ratio of G:A depends on the percentage of fetal DNA present in the maternal blood.

[1317] For example, if the sample contains 50% fetal DNA, then the expected ratio is 0.33 (1 fetal G allele/ (2 maternal A alleles + 1 fetal A allele)). This ratio should be constant for all chromosomes that are present in two copies. The ratio that is obtained for SNPs on chromosome 13 should be the same as the ratio that is obtained for chromosome 21.

[1318] However, if the fetal genome contains an additional copy of chromosome 21, then the ratio for this chromosome will deviate from the expected ratio. The expected ratio for a Trisomy condition



with 50% fetal DNA in the maternal blood is 0.25. Thus, by analyzing SNPs wherein the maternal genome is homozygous, and the DNA that is isolated from the plasma is heterozygous, fetal chromosomal abnormalities can be detected.

### **ABSTRACT OF THE DISCLOSURE**

The invention provides a method useful for detection of genetic disorders. The method comprises determining the sequence of alleles of a locus of interest, and quantitating a ratio for the alleles at the locus of interest, wherein the ratio indicates the presence or absence of a chromosomal abnormality. The present invention also provides a non-invasive method for the detection of chromosomal abnormalities in a fetus. The invention is especially useful as a non-invasive method for determining the sequence of fetal DNA. The invention further provides methods of isolation of free DNA from a sample.